


LINET

Instructions for Use and Technical Description

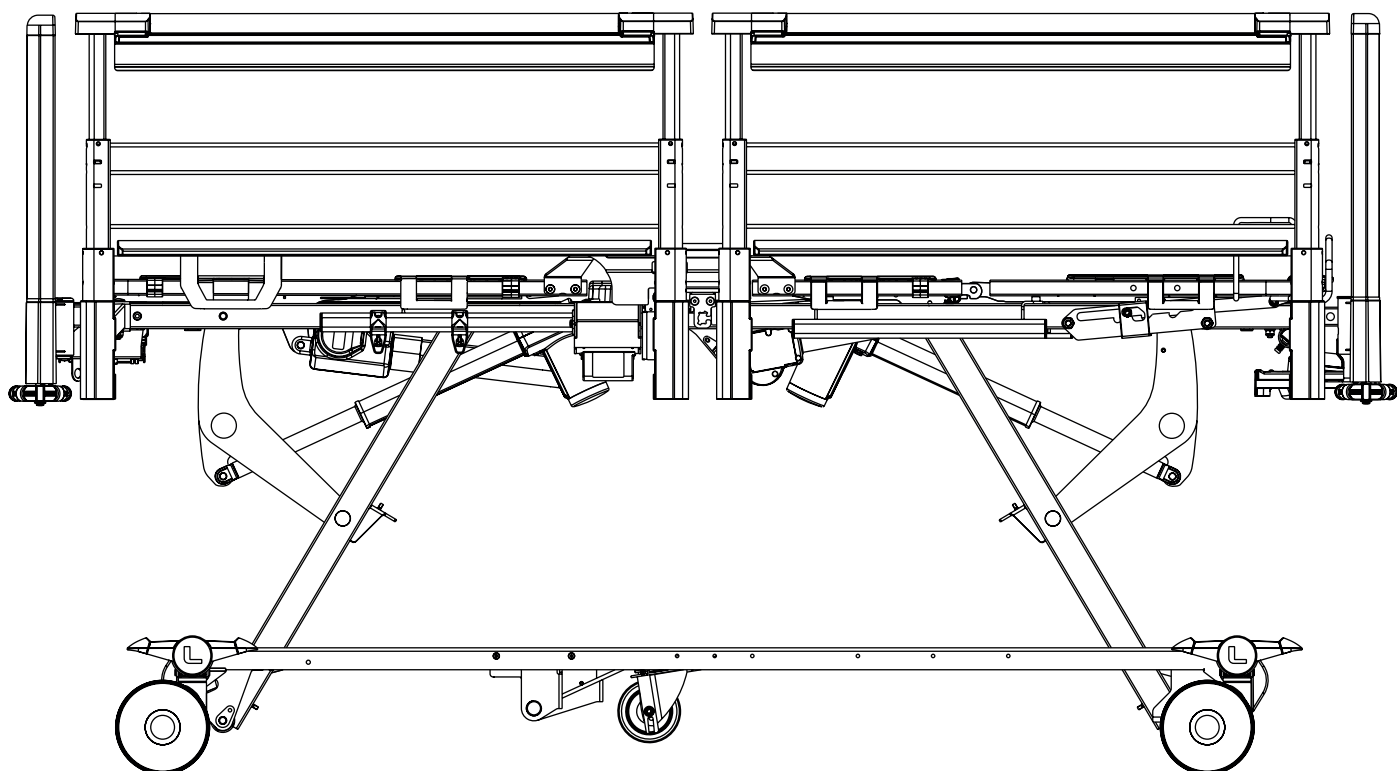

Essenza[®] 300 W

Medical Bed for Acute Care and Long-term Care



D9U001EZW-0101

Version: 02

Publication Date: 2025-07

Manufacturer:

L I N E T spol. s r.o.
Želečnice 5
274 01 Slaný

Tel.: +420 312 576 111
Fax: +420 312 522 668

E-mail: info@linet.cz
<http://www.linet.com>
Service department: service@linetgroup.com

Authorized Representative in Great Britain:

LINET UK Ltd
11 Brunel Way
Segensworth East
Fareham
Hampshire
PO15 5TX
United Kingdom

Authorized Representative in Switzerland:

Bigla Care AG
Bernstrasse 3
CH-3421 Lyssach
Switzerland

Authorized Representative in Malaysia:

Interscience Sdn. Bhd.
2, Jalan Sg. Kayu Ara 32/38
Berjaya Industrial Park
40460 Shah Alam
Selangor, Malaysia

Authorized Representative in Peru:

Drogería Emergo Peru S.R.L.
Calle Las Orquídeas Nro. 585, Int. 1301
San Isidro, Lima
Director Técnico Renato Delgado Rivera
RUC: 205 52 75 65 35
Peru

Authorized Representative in Philippines:

Biomedica Healthcare Inc.
Unit 2103 City land 10, Tower 1
Ayala Corner, Dela Costa Street
Makati City, Metro Manila
Philippines

Essenza® 300 W

Medical Bed for Acute Care and Long-term Care

Author: L I N E T, s.r.o.
Related links: www.linet.com

D9U001EZW-0101

Version: 02
Publication Date: 2025-07

Copyright © L I N E T, s.r.o., 2025
Translation © L I N E T, s.r.o., 2025
All rights reserved.

All trademarks and brands are the property of the appropriate owners. The manufacturer reserves the right to changes in the contents of the instructions for use that relate to the product's technical regulations. It is for this reason that the contents of the instructions for use may indicate differences from the current manufacture of the product. Reproduction, also excerpts, only with prior permission of the publisher. Subject to changes due to technical developments. All technical data are rated data and are subject to construction and manufacturing tolerances.

Table of Contents

1 Symbols and Definitions	7	11.2.1 Description of the Telescopic Siderails without DIN Rails	58
1.1 Warning Notices	7	11.3 Collapsible Siderails	61
1.1.1 Types of Warning Notices	7	11.3.1 Description of the Collapsible Siderails	62
1.1.2 Structure of Warning Notices	7	11.4 Castor Control	64
1.2 Instructions	7	12 Control Elements	66
1.3 Lists	7	12.1 Attendant Control Panel for Essenza 300 W with Lateral Tilt	68
1.4 Symbols on the Package	8	12.2 Attendant Control Panel for Essenza 300 W	69
1.5 Symbols on the Bed	9	12.3 Siderail Control Panel for Essenza 300 W with Lateral Tilt	70
1.6 Serial Label with UDI	14	12.4 Siderail Control Panel for Essenza 300 W	71
1.6.1 Essenza 300 W	14	12.5 Handset Control Panel	72
1.7 Acoustic Signalling	14	12.5.1 Cable of the Handset Control Panel	73
1.8 Visual Signalling	15	12.5.2 Moving the Handset Control Panel to the opposite bed side	74
1.9 Definitions	17	12.6 Patient Control Panel	81
1.10 Abbreviations	18	13 Bed Positioning	82
2 Safety Instructions	19	13.1 Backrest Adjustment	84
3 Intended use of the Essenza 300 W Medical Bed	23	13.2 Thighrest Adjustment	87
3.1 User population	23	13.3 Calfrest Adjustment	89
3.2 Contraindications	23	13.4 Autocontour Adjustment	91
3.3 Operator	23	13.5 Bed Height Adjustment	93
4 Product Description	24	13.6 Emergency Trendelenburg Position	95
4.1 Essenza 300 W with plastic siderails	25	13.7 Trendelenburg Tilt	97
4.2 Essenza 300 W with telescopic siderails	26	13.8 Reverse-Trendelenburg Tilt	99
4.3 Essenza 300 W with collapsible siderails	27	13.9 Bed Extension Adjustment	101
5 Technical Specification	28	13.10 Calfrest Extension Adjustment	109
5.1 Identification of Type B Applied Parts	28	13.10.1 Calfrest Extension Adjustment	110
5.2 Mechanical Specifications	28	13.11 Emergency CPR Position	111
5.2.1 Mechanical Specifications of the Essenza 300 W	29	13.12 Cardiac Chair Position	117
5.2.2 Mechanical Specifications of the Essenza 300 W with Lateral Tilt	30	13.13 Mobilisation Position	119
5.3 Electrical Specifications	31	13.14 Lateral Tilt (only Essenza 300 W with Lateral Tilt)	121
5.4 Environment Conditions	31	13.15 Extra Low Position	124
5.5 Electromagnetic Compatibility	32	13.16 Ergoframe	126
5.5.1 Manufacturer instructions - electromagnetic emissions	32	14 Locking and Unlocking the Bed Positioning Functions	127
5.5.2 Manufacturer instructions - electromagnetic susceptibility	33	14.1 Locking	129
6 Use and Storage Conditions	34	14.2 Unlocking	130
7 Scope of Delivery and Product Variants	34	15 Equipment	131
7.1 Delivery	34	15.1 Backrest stop at 30 degrees	131
7.2 Scope of Delivery	34	15.2 Accessory Rail with plastic hooks	131
7.3 Essenza 300 W Variants	35	15.3 Urinary Bag Holders in Plastic Siderails	133
8 Putting into Service	36	15.4 Bed Illumination and Handset Control Panel Illumination	134
8.1 Accumulator Activation	37	15.5 Safe Height Visual Signalling	135
8.1.1 Removing the Isolating Foil	37	15.6 Fifth Castor	136
8.2 Plastic Foot Board with one Safety Lock	38	15.7 Angle Indicators in Plastic Siderails	137
8.3 Plastic Head Board with one Safety Lock	39	15.8 Plastic Linen Shelf	138
8.4 HPL Foot Board	41	15.9 Telescopic Linen Shelf	140
8.5 HPL Head Board	42	15.10 Brake Signal	142
8.6 Plastic Foot Board with two Safety Locks	44	15.11 SafeSense® 3 ready bed	143
8.7 Plastic Head Board with two Safety Locks	45	16 Mattress	144
8.8 Mattress Support Platform	47	16.1 Passive Mattress	146
8.8.1 Handles and Patient Restraint Points	49	16.2 Active Mattress	146
8.9 Potential Interconnection	49	16.3 Cleaning of Mattress	147
8.10 Before Use	50	16.3.1 General Guidance	147
8.11 Transport	50	16.3.2 Routine Cleaning and Disinfection	148
8.11.1 Transport Position	51	16.3.3 Complete Cleaning and Disinfection	148
8.12 Firmware	51	16.3.4 Mattress Core	148
9 Power Supply Cord	52	17 Accessories	149
10 Accumulator of the bed	53	18 Cleaning and Disinfection	150
10.1 Replacing the accumulator	54	18.1 Cleaning the Essenza 300 W bed	151
10.2 Removing the Bed from Service	54	18.1.1 Cleaning before Changing Patients	151
11 Manipulation	55	18.1.2 Daily Cleaning	151
11.1 Plastic Siderails	55	18.1.3 Complete Cleaning and Disinfection	151
11.1.1 Description of the Plastic Siderails	56	18.2 Machine Washing	152
11.2 Telescopic Siderails	58	18.2.1 Preparing for cleaning	153

18.2.2 Machine washing	155
18.2.3 After machine washing	155
19 Troubleshooting	156
20 Maintenance	157
20.1 Regular maintenance	157
20.2 Spare Parts	157
20.3 Electrical Safety Check and Preventive Maintenance	157
20.3.1 Electrical Safety Check	157
20.3.2 Preventive Maintenance	158
20.4 Instructions related to machine washing	159
20.4.1 Recommended places for greasing	159
20.5 Control Unit Leak Tester (8300-9004) for the Essenza 300 W washable bed	163
21 Disposal	164
21.1 Environment Protection	164
21.2 Disposal	164
21.2.1 Within Europe	164
21.2.2 Outside Europe	164
22 Warranty	165
23 Standards and Regulations	165
23.1 Essenza 300 W	165
23.2 Manufacturer	165

List of Figures

Fig. Accumulator Activation Instructions: Remove isolating foil from the control unit under the Backrest.....	13
Fig. Visual Signalling of Attendant Control Panel for Essenza 300 W with Lateral Tilt.....	16
Fig. Position of Isolating Foil under the Backrest	37
Fig. Shape of Isolating Foil	37
Fig. Standard Plastic Foot Board	38
Fig. Standard Plastic Head Board	39
Fig. Narrow Plastic Head Board	40
Fig. HPL Foot Board	41
Fig. HPL Head Board	42
Fig. Two Head Board Safety Locks	43
Fig. Correctly inserted HPL Head Board.....	43
Fig. Plastic Foot Board with two Safety Locks	44
Fig. Plastic Head Board with two Safety Locks	45
Fig. Two Head Board Safety Locks	46
Fig. 4-part Mattress Support Platform with Plastic Siderails.....	47
Fig. 4-part Mattress Support Platform with Telescopic Siderails	48
Fig. Potential equalisation male connector at head end of the bed.....	49
Fig. Potential equalisation female connector.....	49
Fig. Hook for hanging Power Supply Cord.....	52
Fig. Manipulation with Plastic Siderails.....	56
Fig. Correct positions of the Handset Control Panel on the Plastic Siderails	57
Fig. Correct positions of the Attendant Control Panel on the Plastic Siderails.....	57
Fig. Manipulation with Telescopic Siderails without DIN Rails	58
Fig. Red Buttons of Head Siderail Release Mechanism.....	59
Fig. Telescopic Siderail in highest (1) and lowest (4) positions	59
Fig. Telescopic Siderail in two interpositions (2, 3).....	60
Fig. Correct positions of the Handset Control Panel on the Telescopic Siderails	60
Fig. Correct positions of the Attendant Control Panel on the Telescopic Siderails.....	61
Fig. Manipulation with Collapsible Siderail	62
Fig. Manipulation with Siderail Release Mechanism in version with collapsible siderails.....	63
Fig. Correct position of the Handset Control Panel on the Collapsible Siderail.....	63
Fig. Correct position of the Attendant Control Panel on the Collapsible Siderail	64
Fig. Three Positions of Castor Control Lever	65
Fig. Attendant Control Panel for Essenza 300 W with Lateral Tilt.....	68
Fig. Attendant Control Panel for Essenza 300 W	69
Fig. Siderail Control Panel for Essenza 300 W with Lateral Tilt	70
Fig. Siderail Control Panel for Essenza 300 W.....	71
Fig. Handset Control Panel	72
Fig. Plastic Lock of the Handset Control Panel Cable	73
Fig. Handset Control Panel Cable Holder with unlocked plastic lock	73
Fig. Handset Control Panel Cable Holder with locked plastic lock.....	73
Fig. Position of the Handset Control Panel Cable Holder on the bed with telescopic siderails	75
Fig. Position of the Handset Control Panel Cable Holder on the bed with plastic siderails.....	75
Fig. Moving the Handset Control Panel Cable to the opposite bed side through the bed seat section.....	76
Fig. Position of the Handset Control Panel Cable Holder on the bed with telescopic siderails	77
Fig. Handset Control Panel Cable Routing between the seat section and thighrest on the bed with telescopic siderails ...	78
Fig. Position of the Handset Control Panel Cable Holder on the bed with collapsible siderails	79
Fig. Handset Control Panel Cable Routing between the seat section and thighrest on the bed with collapsible siderails...	80
Fig. Patient Control Panel	81
Fig. Bed adjustment unsuitable for sitting in the middle of the bed.....	86
Fig. Calfrest Adjustment.....	89
Fig. Calfrest Handle	90
Fig. Position of Simple Bed Extension Safety Catch on the bed with plastic siderails	102
Fig. Position of Advanced Bed Extension Handles under the Foot Board on the bed with plastic siderails	104
Fig. Position of Simple Bed Extension Safety Catch on the bed with telescopic siderails	106
Fig. Position of Advanced Bed Extension Handles under the Foot Board on the bed with telescopic siderails	108
Fig. Adjustment of Calfrest Extension Holder.....	109
Fig. Two steps of the Calfrest Extension Adjustment	110
Fig. CPR lever under the Backrest of the bed with Plastic siderails	112
Fig. Manipulation with CPR lever in version with plastic siderails	112
Fig. CPR lever under the Seat Section of the bed with Telescopic siderails	114
Fig. Manipulation with CPR lever in version with telescopic siderails.....	114
Fig. CPR lever under the seat section of the bed with collapsible siderails	116
Fig. Manipulation with CPR lever in version with collapsible siderails	116

Fig. Lock Button and Lock LED Indicators on the Attendant Control Panel.....	127
Fig. Lock Button and Lock LED Indicators on the Siderail Control Panel	128
Fig. Lock LED Indicators on the Handset Control Panel	128
Fig. Accessory Rails with plastic hooks on the bed with Plastic Siderails.....	131
Fig. Accessory Rails with plastic hooks on the bed with Telescopic Siderails	132
Fig. Accessory Rails with plastic hooks on the bed with Collapsible Siderails	132
Fig. Urinary Bag Holders in Plastic Siderails	133
Fig. Bed Illumination source on the Essenza 300 W medical bed	134
Fig. Green light source on the Essenza 300 W medical bed	135
Fig. Activation of the Fifth Castor	136
Fig. Angle Indicators in Plastic Siderails.....	137
Fig. Plastic Linen Shelf at bed foot end.....	138
Fig. Attendant Control Panel in the Plastic Linen Shelf.....	139
Fig. Telescopic Linen Shelf at bed foot end	140
Fig. Telescopic Linen Shelf in spread out position	141
Fig. Telescopic Linen Shelf in retracted position	142
Fig. SafeSense 3 ready Essenza 300 W bed.....	143
Fig. Places for fixation of the mattress on the mattress support platform of the Essenza 300 W bed	145
Fig. Arrow-like symbol indicating bed centre on Head Board and Foot Board	146
Fig. Essenza 300 W bed intended for machine washing.....	152
Fig. Recommended position of bed for machine washing	154

1 Symbols and Definitions

1.1 Warning Notices

1.1.1 Types of Warning Notices

Warning notices are differentiated by the type of danger using the following signal words:

- ▶ **CAUTION** warns about the risk of material damage.
- ▶ **WARNING** warns about the risk of physical injury.
- ▶ **DANGER** warns about the risk of fatal injury.

1.1.2 Structure of Warning Notices



DANGER!
Dangerous electricity!
 ▶ Measures to avoid the risk.



WARNING!
Type and source of risk!
 ▶ Measures to avoid the risk, if necessary.



CAUTION!
Type and source of risk!
 ▶ Measures to avoid the risk, if necessary.

1.2 Instructions

Structure of instructions:



- ▶ Perform this step.
 Results, if necessary.

1.3 Lists


Structure of bulleted lists:


- List level 1
 - List level 2
 - List level 3



1.4 Symbols on the Package

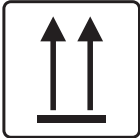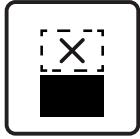 	FRAGILE, HANDLE WITH CARE
	THIS WAY UP
 	KEEP DRY (PROTECT FROM HUMIDITY)
	PAPER RECYCLING SYMBOL
	DO NOT USE HAND TRUCK HERE
	DO NOT STACK DURING STORAGE
	OVERSEAS PACKAGE: STACKING LIMIT BY NUMBER (4 PACKAGES FOR STORAGE)
	OVERSEAS PACKAGE: STACKING LIMIT BY NUMBER (2 PACKAGES FOR TRANSPORT)
	READ INSTRUCTIONS FOR USE
	DO NOT POLLUTE THE ENVIRONMENT

1.5 Symbols on the Bed

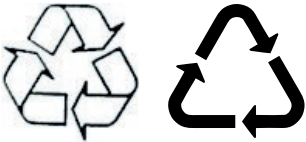	READ INSTRUCTIONS FOR USE
	SAFE WORKING LOAD
	WARNING AGAINST CRUSHING OR TRAPPING
	JACK FOR ATTACHMENT OF CONDUCTOR FOR POTENTIAL EQUALISATION
	GENERAL WARNING SIGN
	CAUTION SIGN
	ONLY SUITABLE FOR INDOOR USE
	MAXIMUM WEIGHT OF PATIENT
	WEIGHT OF BED (depending on configuration)
	RECYCLING SYMBOLS
	DO NOT POLLUTE THE ENVIRONMENT

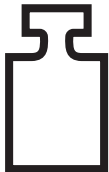	MANUFACTURER
	MANUFACTURING DATE
	REFERENCE NUMBER (PRODUCT TYPE DEPENDING ON CONFIGURATION)
	SERIAL NUMBER
	MEDICAL DEVICE (COMPATIBLE WITH MEDICAL DEVICE REGULATION)
	WEEE SYMBOL (RECYCLE AS ELECTRONIC WASTE, DO NOT PUT INTO THE HOUSEHOLD WASTE)
	BATCH NUMBER (ACCESSORIES)
	UNIQUE DEVICE IDENTIFICATION (FOR MEDICAL DEVICES)
	MAXIMUM MASS OF MOBILE HOSPITAL BED (MAXIMUM MASS OF EMPTY Bed + SAFE WORKING LOAD)
	CE MARKING FOR Essenza 300 W (PRODUCT NORMATIVELY HARMONIZED FOR EUROPEAN ECONOMIC AREA)

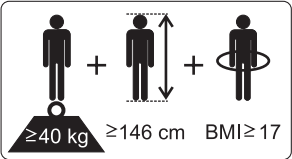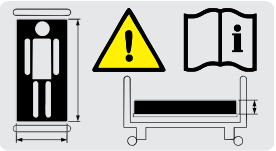	UK CONFORMITY ASSESSED (UKCA) MARKING (PRODUCT NORMATIVELY HARMONIZED FOR GREAT BRITAIN ECONOMIC AREA)
	AUTHORIZED REPRESENTATIVE IN GREAT BRITAIN
	AUTHORIZED REPRESENTATIVE IN SWITZERLAND
	GO BUTTON (PRESS TO ACTIVATE CONTROL ELEMENT)
	STOP BUTTON (PRESS TO INTERRUPT BED POSITIONING)
	DESIGNATION OF MEDICAL BED FOR ADULTS (PHYSICAL DESCRIPTION OF AN ADULT)
	ONLY USE MATTRESS RECOMMENDED BY MANUFACTURER
	DO NOT PUT ANY OBJECTS ON UNDERCARRIAGE
	LEVER INTENDED FOR EMERGENCY CPR POSITION ADJUSTMENT
	THERMAL PROTECTION FOR TRANSFORMER

	GENERAL SAFETY ISOLATING TRANSFORMER
	DO NOT OPEN A CASE MARKED WITH THIS SYMBOL
	CALFREST LOAD LIMIT
	EARTH GROUND
	APPLIED PARTS TYPE B
	WARNING, READ INSTRUCTIONS FOR USE
	EAC MARKING (EURASIAN CONFORMITY)
	BED IS INTENDED FOR CLEANING IN DECONTAMINATION DEVICES, WHICH MEETS CRITERIA STATED IN STANDARD AK:BWA 2022, 9TH EDITION (LABEL IS PLACED ON THE UNDERCARRIAGE UNDER THE BED FOOT BOARD)

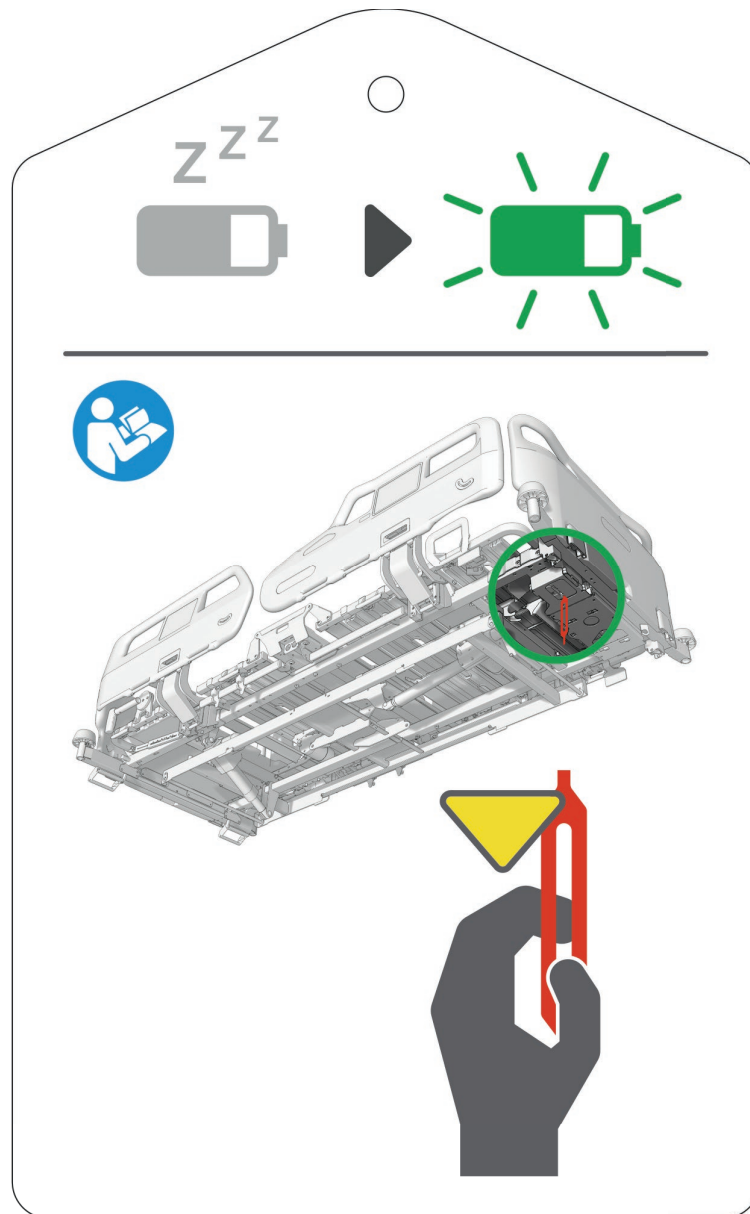


Fig. Accumulator Activation Instructions: Remove isolating foil from the control unit under the Backrest

1.6 Serial Label with UDI

1.6.1 Essenza 300 W

Serial label of the Essenza 300 W is placed on the left head end side of the mattress support platform frame. The serial label contains information about Address of Manufacturer, Manufacturing Date (Year-Month-Day), Product Reference Number, Product Serial Number, Global Trade Item Number (GTIN), Unique Device Identification (UDI), symbols, weight specifications and electrical specifications.

1.7 Acoustic Signalling

SOUND	MEANING
CONTINUOUS SOUND	overheating
	accumulator overcurrent
	actuator overload
REPEATED BEEP: 0,6s sound / 2,6s silence	STOP error (all STOP buttons are disabled)
BEEP lasting 0,3s	confirmation
	stopping or locked function
	transition from tilt (lateral tilt, Trendelenburg, Reverse-Trendelenburg) to horizontal position
4 TIMES REPEATED BEEP lasting 0,3s	disconnected from the mains
	positioning powered by the accumulator
BEEP lasting 0,5s	adjustment to the Extra Low Position started
	start of service mode or end of service mode
	keyboard error (positioning blocked)
BEEP lasting 3s	system error
REPEATED BEEP during 3 minutes: 1,1s sound / 1,1s silence	Brake Signal: acoustic signal indicating the bed is not braked when plugged into the mains (only version with Brake Signal)

1.8 Visual Signalling

MAINS POWER LED INDICATOR (1)	MEANING
lit (ON) when bed is plugged into the mains	bed is mains powered
not lit (OFF)	bed is not mains powered
flashing: 0,5s lit / 0,5s not lit	control unit overheated
flashing: 0,05s lit / 0,05s not lit	calibration mode

ACCUMULATOR CHARGE STATUS LED INDICATOR (2)	MEANING
lit (ON) when bed is plugged into the mains	accumulator fault: accumulator not connected or low voltage
not lit (OFF) when bed is plugged into the mains	accumulator works correctly
flashing (0,2s lit / 1,6s not lit) when bed is plugged into the mains	accumulator is charging
flashing (0,1s lit / 0,1s not lit) when bed is only accumulator powered	accumulator low (Positioning does not work. Beep sounds anytime a positioning button was pressed.)
not lit (OFF) when bed is only accumulator powered	accumulator works correctly
flashing: 0,05s lit / 0,05s not lit	calibration mode

BACKREST ADJUSTMENT LOCK LED INDICATOR (4)	MEANING
lit (ON)	locked function
not lit (OFF)	unlocked function
flashing once per second	unlocked function but in the mode of lock selection
flashing 4 times per second	locked function but in the mode of lock selection

THIGHREST ADJUSTMENT LOCK LED INDICATOR (5)	MEANING
lit (ON)	locked function
not lit (OFF)	unlocked function
flashing once per second	unlocked function but in the mode of lock selection
flashing 4 times per second	locked function but in the mode of lock selection

BED HEIGHT ADJUSTMENT LOCK LED INDICATOR / LATERAL TILT LOCK LED INDICATOR / TRENDELENBURG TILT LOCK LED INDICATOR / REVERSE-TRENDELENBURG TILT LOCK LED INDICATOR (6)	MEANING
lit (ON)	locked function
not lit (OFF)	unlocked function
flashing once per second	unlocked function but in the mode of lock selection
flashing 4 times per second	locked function but in the mode of lock selection

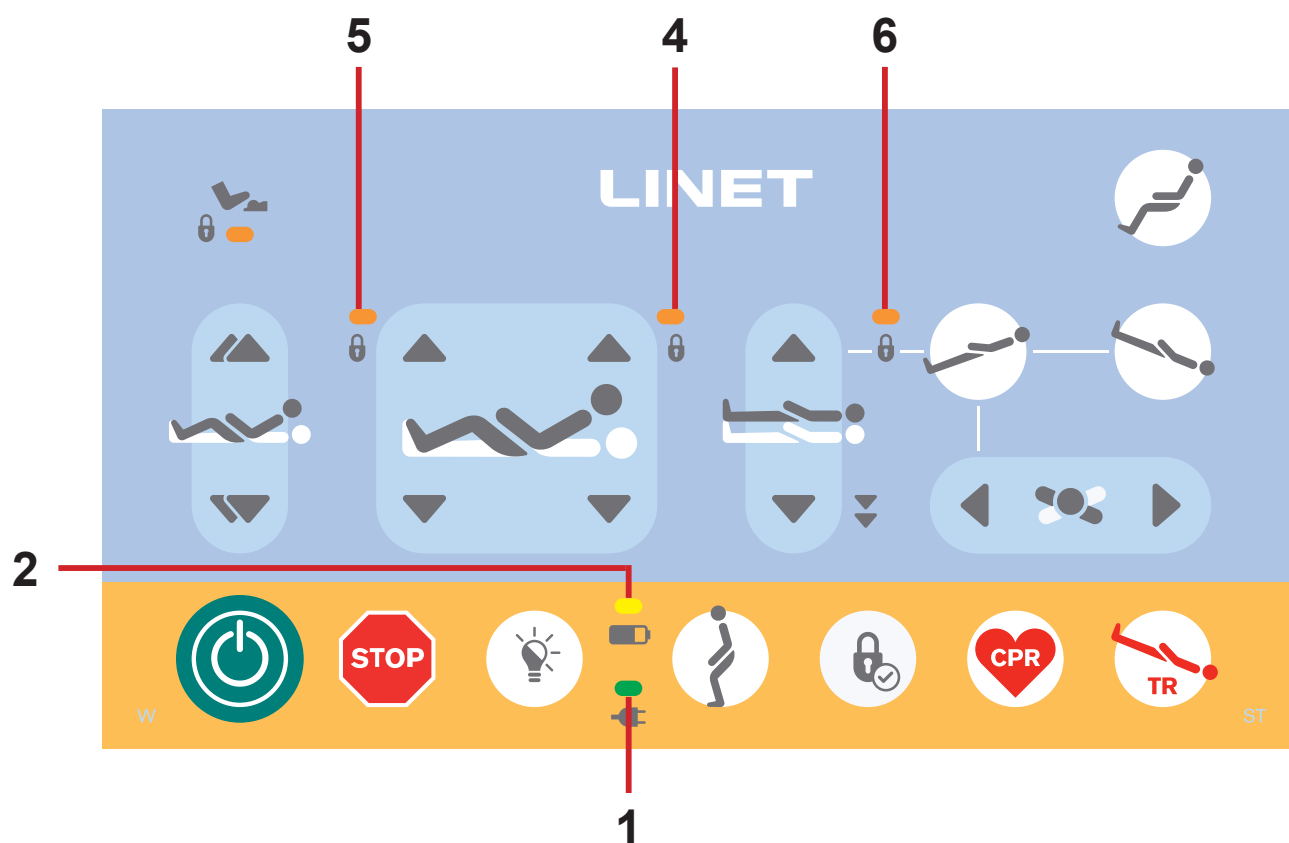


Fig. Visual Signalling of Attendant Control Panel for Essenza 300 W with Lateral Tilt

1.9 Definitions

Adult	Patient having a physical size equal to or more than 146 cm, a mass equal to or more than 40 kg and a body mass index (BMI) equal to or more than 17 (according to IEC 60601-2-52).
Basic Bed Configuration	The pricelist model configuration, not including a mattress.
Clearance of Undercarriage	The height from the floor to the lowest point of the undercarriage between the castors, for the manipulation of accessories under a braked bed in the standard position.
Duty Cycle	Cycle of operation of the motor: time of activity/time of rest.
Ergoframe	Ergoframe is the kinematic system of Mattress support platform Adjustment whose effect is the reduction of pressure on the patient's abdomen and pelvic area and frictional forces on the patient's back and legs.
Maximum Mass of Mobile Hospital Bed	Sum of Empty Bed Maximum Mass and Safe Working Load.
Reference number	Reference number depends on configuration.
Safe Working Load	The highest allowable load on the bed. The Safe Working Load is the sum of: – the patient – the mattress – the accessories of the medical bed (only if they are supported by the support system of the medical bed) and – the load supported by those accessories (excluding patient weight).
Siderail Height	The height of the upper crossbar or the edges of the siderails (not the highest point of the siderail controls) from the patient surface.
Standard Bed Position	– The mattress support platform with regard to the floor is in the middle position – The Mattress Support Platform, including its individual parts, has to be in a horizontal position (level - 0°). – The siderails are always locked in the upper position.
Bed Weight	The value depends on the product configuration, accessories or customer adjustments.
Type B Applied Parts	The degree of protection against electric shock regarding the product parts in contact with patient.
Machine Washing	This bed can be cleaned in the compatible automatic washing tunnel.

1.10 Abbreviations

AC (~)	Alternating Current
ALT	Automatic Lateral Therapy
CE	European Conformity
CPR	Cardiopulmonary Resuscitation
dB	Sound Intensity Unit
DC (—)	Direct Current
CUC	Configuration number
EMC	Electromagnetic Compatibility
FET	Field-effect transistor
HF	High Frequency
HPL	High Pressure Laminate
ICU	Intensive Care Unit
INT.	Duty Cycle
IP	Ingress Protection
IV	Intravenous
LCD	Liquid Crystal Display
LED	Light Emitting Diodes
ME	Medical Electrical (Equipment)
ON	Activation
OFF	Deactivation
ppm	Parts per million, millionth (1000 ppm = 0,1%)
REF	Reference Number (product type depending on configuration)
SCU	System Control Unit (integrated mattress)
SN	Serial Number
SWL	Safe Working Load
UDI	Unique Device Identification (for medical devices)
USB	Universal Serial Bus
WEEE	Waste Electrical and Electronic Equipment

2 Safety Instructions



WARNING!

Risk of injury due to incorrect use!

- ▶ **Essenza 300 W bed should be left in its lowest position when the patient is unattended in order to reduce risk of injury due to falls!**



WARNING!

Risk of injury due to incorrect use!

- ▶ **Siderails of Essenza 300 W should be located in the „up“ position to reduce the risk of the patient accidentally slipping or rolling off the mattress!**



WARNING!

Risk of injury due to incorrect use!

- ▶ **Lateral Tilt of the Essenza 300 W with Lateral Tilt is a positioning function intended for trained medical personnel only. The use of the Lateral Tilt, including its angle of adjustment, must be preceded by an evaluation of the patient's condition by the staff. This evaluation must take into account all risks resulting from the physical and mental state of the patient, especially the risk of fall, and loss of bed stability.**



WARNING!

Risk of injury due to incorrect use!

- ▶ **Staff expert assessment is needed to consider all individual cases of contraindications!**



WARNING!

Risk of trapping or squeezing because of patient's body constitution disproportionate to the size of mattress support platform!



WARNING!

Risk of injury due to incorrect use!

- ▶ **Certain bed positions are not suitable for specific diagnosis/medical conditions. Fowler position is not suitable for spinal cord injuries! Trendelenburg position is not suitable for patients with higher intracranial pressure!**



WARNING!

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established!



WARNING!

This medical device is not intended for oxygen enriched environment!



WARNING!

This medical device is not intended for use with flammable substances!



WARNING!

This medical device is not portable medical electrical equipment!



WARNING!

Patient is allowed to use selected control elements only if hospital personnel had assessed that the patient's physical and psychological state is in accordance with use of them and only if the hospital personnel had trained the patient in accordance with the instructions for use!



WARNING!

An additional multiple socket-outlet or extension cord shall not be connected to the medical electrical system of Essenza 300 W.



WARNING!

Inappropriate handling of the power supply cord of Essenza 300 W, e. g. by kinking, shearing or other mechanical damages is hazardous!



WARNING!

During specific investigations or treatments, the significant risks of reciprocal interference posed by Essenza 300 W (ME equipment) may occur.



WARNING!

To avoid the risk of electric shock, this equipment must only be connected to a supply mains with protective earth.



WARNING!

Incompatible siderails and mattresses can cause an entrapment hazard!



WARNING!

Incompatible mattresses can create hazards!



WARNING!

The bed is intended for one adult.

► Follow the Intended use chapter.



WARNING!

When routing cables from other equipment in the Essenza 300 W bed avoid squeezing those between parts of the Essenza 300 W bed!

**WARNING!**

No modification of this equipment is allowed.

Do not modify this equipment without authorization of the manufacturer.

**WARNING!**

If this Essenza 300 W medical bed (ME equipment) is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.

**WARNING!**

Only authorised and trained person using an appropriate tool is allowed to change fuses, accumulators, control units and power supplies!

**WARNING!**

Make sure the duty cycle (2 min ON/18 min OFF) is not exceeded during bed positioning!

**WARNING!**

No part of the Essenza 300 W ME equipment shall be serviced or maintained while in use with a patient!

**WARNING!**

Patient is allowed to use selected control elements of the ME equipment unless servicing and maintenance of the ME equipment is performed!

**CAUTION!**

For cleaning purposes this bed can be used with automatic washing systems.

Additional Instructions for correct use:

- ▶ Follow the instructions for use carefully.
- ▶ Use the bed exclusively if it is in perfect working order.
- ▶ If necessary, check the bed functions daily or at each shift change.
- ▶ Ensure any user has read and understood this manual completely before operating the product.
- ▶ Use the bed exclusively with the correct mains supply.
- ▶ Ensure that the bed is operated exclusively by qualified personnel who have been trained according to the instructions for use.
- ▶ Ensure the bed is operated by qualified personnel trained in using the bed in accordance with the instructions for use by manufacturer or by person approved by the manufacturer, eventually by authorized representative or by person approved by the authorized representative.
- ▶ Ensure that the patient (health permitting) has been informed about the operation of the bed and all applicable safety instructions.
- ▶ Contact service department of the manufacturer immediately to let it replace any damaged parts by the original spare parts.
- ▶ Move the bed exclusively on even, hard-surfaced floors.
- ▶ Replace any damaged parts immediately with original spare parts.
- ▶ Ensure that maintenance and installation are performed exclusively by qualified personnel who have been trained by the manufacturer.
- ▶ Before peak loads or unavoidable excess loads (CPR), place Mattress Support Platform to the lowest position.
- ▶ Ensure that only one adult patient lies on the bed at any time.
- ▶ Take care to avoid injuries or squeezing when operating moving parts.
- ▶ When using lifting pole or infusion stand, ensure that nothing will be damaged when you move or adjust the bed.
- ▶ Brake the castors when the bed is occupied.
- ▶ Keep the Mattress Support Platform in the lowest position when the patient is unattended by healthcare personnel in order to minimize risk of patient falls.
- ▶ Ensure that siderails are operated exclusively by healthcare personnel.
- ▶ Never use the bed in areas where there is a hazard of explosion.
- ▶ Enable or disable functions on control elements as appropriate for the patient's physical and mental state. Afterwards ensure the lock function was used as required.
- ▶ Ensure that parts of the bed intended for movement are not blocked.
- ▶ Never handle the mains plug with wet hands.
- ▶ Disconnect the product from the mains exclusively by pulling the mains plug.
- ▶ When pulling the mains plug, always hold the plug, not the cable.
- ▶ Position the mains cable so that there are no loops or kinks in the cable. Protect the cable from mechanical wear and tear.
- ▶ Ensure that the stipulated duty cycle of motor is not exceeded.
- ▶ To change fuses or cables contact service organisation authorized by manufacturer.
- ▶ To prevent failures, use exclusively the manufacturer's original accessories and mattresses.
- ▶ Ensure that the stipulated safe working load is not exceeded.
- ▶ If the patient's condition could lead to an entrapment, leave the mattress support platform in the flat position whilst unattended.
- ▶ Adjust bed height when transporting the bed in order to facilitate overcoming possible obstacles.
- ▶ Do not exceed maximum load of 80 kg for mattress support platform extension.
- ▶ Do not hang anything on any cable.
- ▶ Do not modify bed and its components without the manufacturer's approval.
- ▶ Use the mattress exclusively as specified in this manual and in perfect working order.
- ▶ Use the mattress exclusively in its original state and do not modify it in any way.
- ▶ Have the mattress used exclusively by or under supervision of trained and qualified nursing personnel.
- ▶ Have the mattress serviced and installed exclusively by qualified personnel trained and authorised by the manufacturer.
- ▶ Do not exceed the maximum patient weight limit (see Mechanical Specifications).
- ▶ Do not use the bed in the case its parts have been removed unless these parts are designed to be removed.
- ▶ To avoid injury or crushing, take extra caution when operating any moving parts of the bed.
- ▶ Never use the mattress replacement system compatible with the Essenza 300 W bed near radiators or other heat sources.
- ▶ Select a suitable location for the placement of bed accessories and other objects to prevent involuntary activation of buttons or controls which may result in the adjustment of bed positioning.
- ▶ Never place any accessories or Handset Control Panel on the siderails where keyboards are located.
- ▶ After each emergency situation always check if any of the control elements is not involuntarily pressed by the bed accessories or by the mattress.
- ▶ To avoid unintended activation of moving parts during any use of the bed always check that none of the control elements of the bed is pressed by persons, mattress or other objects.
- ▶ When adjusting the Lateral Tilt in the case of Essenza 300 W with Lateral Tilt, ensure that the process does not pose any risk of the patient falling or getting injured, or of any lines or tracheal tubes getting pulled out.

3 Intended use of the Essenza 300 W Medical Bed

The intended use is the hospitalization of the patient in the acute and long-term care units, which includes above all the following aspects:

- ▶ Adjustment of the specific positions needed for the preventive reasons, routine nursing, treatments, mobilization, physiotherapy, examinations, sleeping, and relaxation.
- ▶ Providing the safe environment for the patient during all relevant procedures. The particular requirements on patient safety are the subject of the clinical evaluation, including evaluation of the risk/benefit ratio. The relevant safety issues are the part of the risk management file.
- ▶ Patient in-bed indoor transport out of the patient room.
- ▶ Providing the suitable working conditions for the caregivers to perform the routine and specific tasks during the patient hospitalization.

3.1 User population

- ▶ Adult patients (weight ≥ 40 kg, height ≥ 146 cm, BMI ≥ 17) in the acute and long-term care units (Application Environment 2 and 3, as in IEC 60601-2-52), alternatively in Application Environment 1 and 5
- ▶ Caregivers (nurses, doctors, technical personnel, transport personnel, cleaning personnel)

3.2 Contraindications

- ▶ The medical device is not intended for the pediatric patients use.
- ▶ Limited patient handling leads to development of pressure injury.
- ▶ Certain positions are not suitable for specific diagnoses/medical conditions (e.g. spinal cord injuries vs. Fowler position, higher ICP patients vs. Trendelenburg). Staff expert assessment / nursing consideration is needed in all individual case of contraindication.


3.3 Operator

- ▶ Caregiver
- ▶ Patient (based on individual patient status assessment by caregiver the patient can utilize dedicated device functions)

4 Product Description


Following pictures show basic visible features of the Essenza 300 W medical bed together with 2 types of the siderails.

4.1 Essenza 300 W with plastic siderails

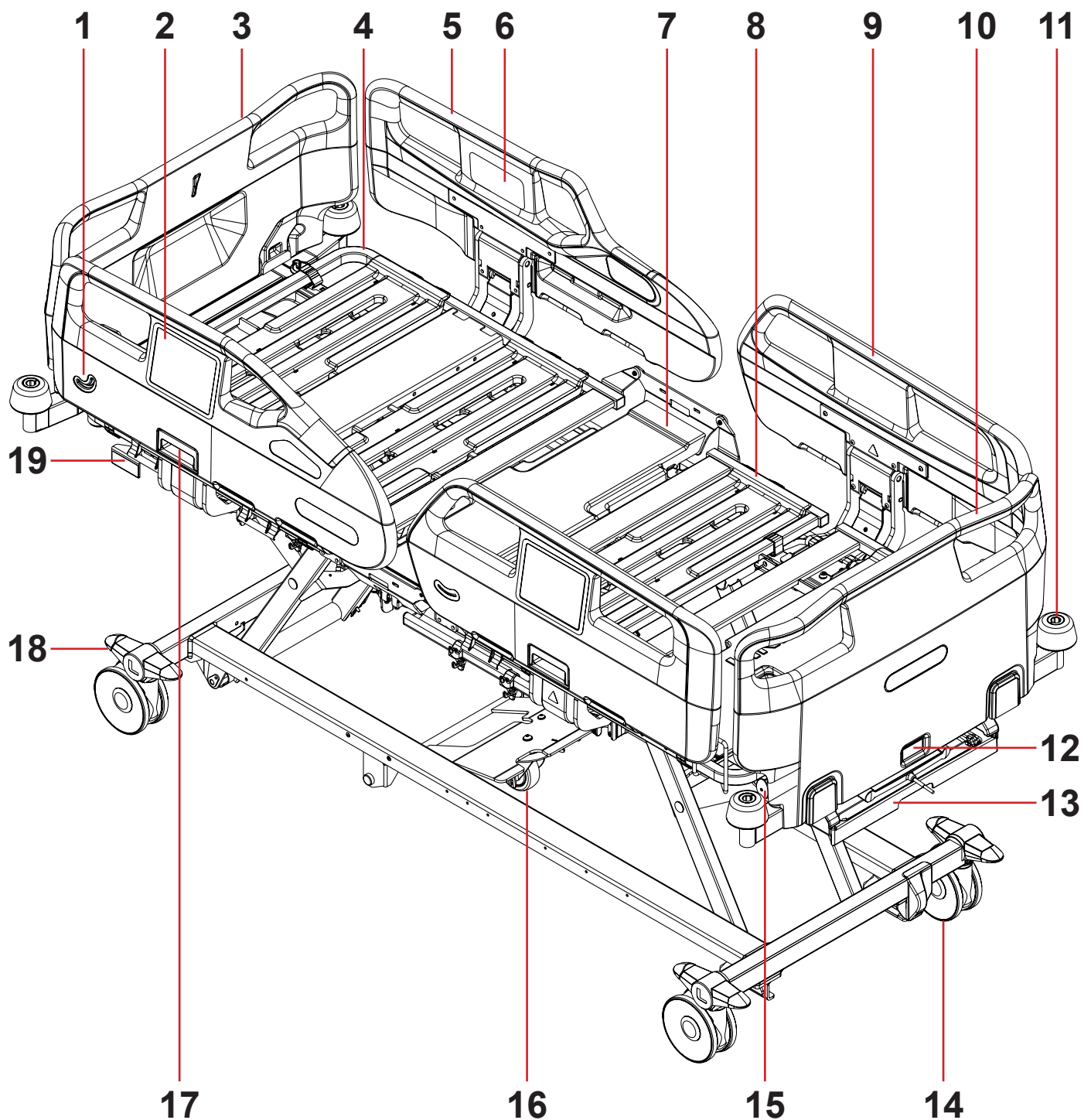

- | | | | |
|-----|--------------------------------------|-----|--|
| 1. | Angle Indicator | 12. | Foot Board Safety Lock |
| 2. | Siderail Control Panel | 13. | Plastic Linen Shelf (also for Attendant Control Panel) |
| 3. | Head Board | 14. | Castor |
| 4. | Backrest | 15. | Bed Extension Safety Catch (Simple Bed Extension) |
| 5. | Head Plastic Siderail | 16. | Fifth Castor |
| 6. | Patient Control Panel | 17. | Siderail Release Mechanism |
| 7. | Seat Section | 18. | Castor Control Lever |
| 8. | Thighrest | 19. | CPR Lever for Emergency Backrest Release |
| 9. | Foot Plastic Siderail | | |
| 10. | Foot Board | | |
| 11. | Corner Bumper with accessory bushing | | |


4.2 Essenza 300 W with telescopic siderails

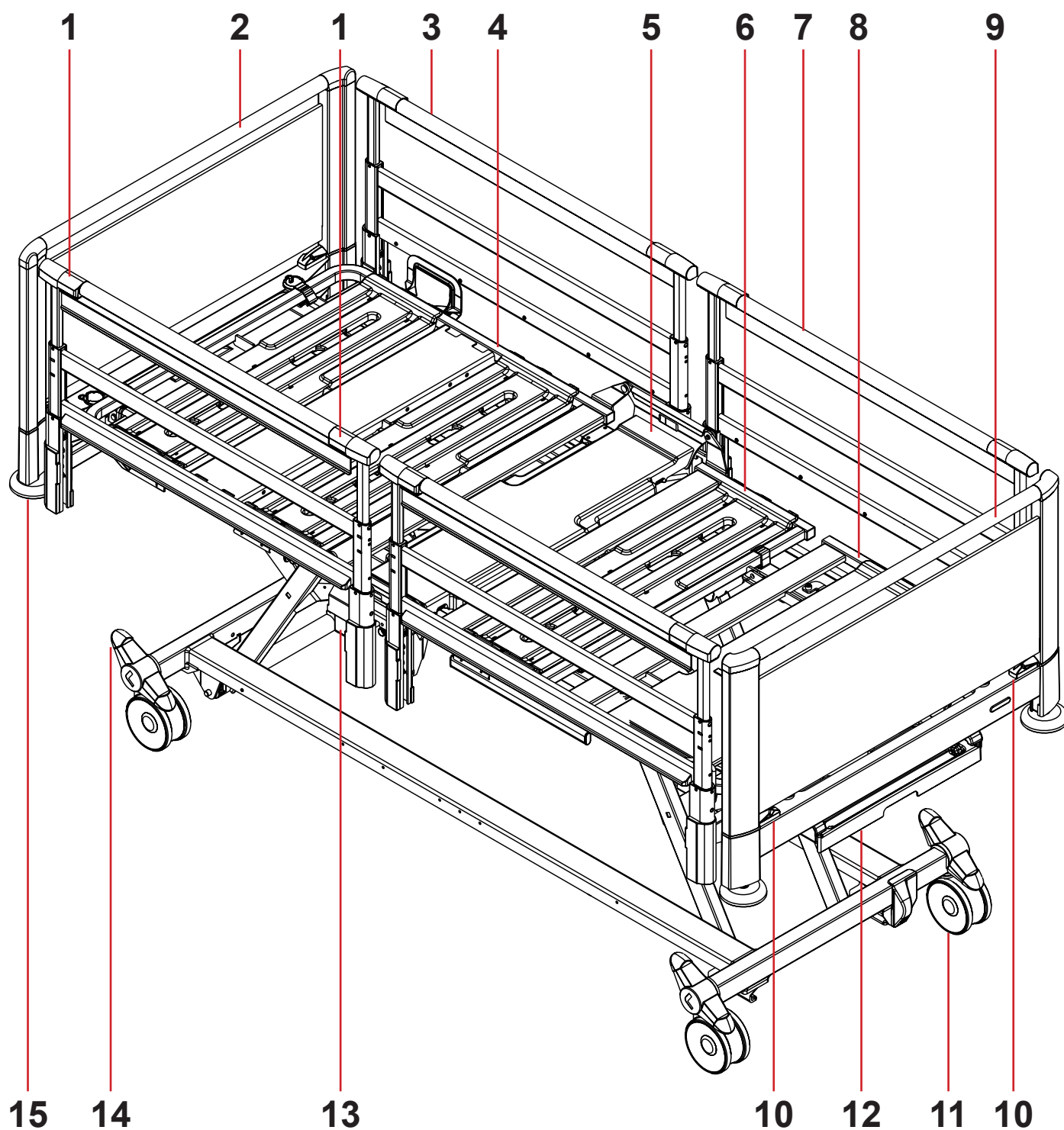

- 1. Siderail Release Mechanism
- 2. HPL Head Board
- 3. Head Telescopic Siderail
- 4. Backrest
- 5. Seat Section
- 6. Thighrest
- 7. Foot Telescopic Siderail
- 8. Calfrest
- 9. HPL Foot Board
- 10. Foot Board Safety Lock

- 11. Castor
- 12. Plastic Linen Shelf (also for Attendant Control Panel)
- 13. CPR Lever for Emergency Backrest Release
- 14. Castor Control Lever
- 15. Corner Bumper without accessory bushing

4.3 Essenza 300 W with collapsible siderails

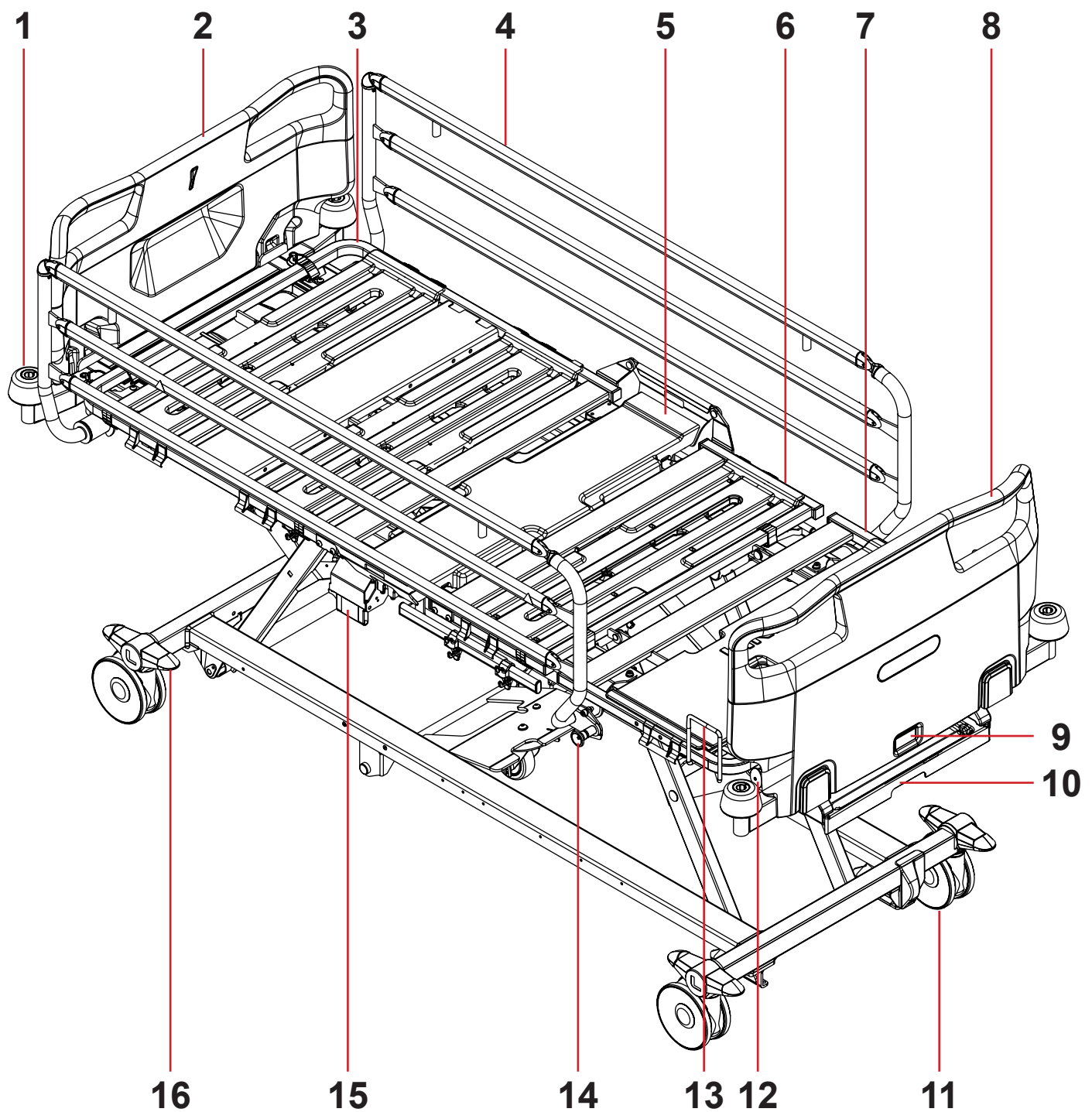

1. Corner Bumper with accessory bushing
2. Head Board
3. Backrest
4. Single Collapsible Siderail
5. Seat Section
6. Thighrest
7. Calfrest
8. Foot Board
9. Foot Board Safety Lock

10. Plastic Linen Shelf (also for Attendant Control Panel)
11. Castor
12. Bed Extension Safety Catch (Simple Bed Extension)
13. Calfrest Handle
14. Siderail Release Mechanism
15. CPR Lever for Emergency Backrest Release
16. Castor Control Lever

5 Technical Specification

All technical data are rated data and are subject to construction and manufacturing tolerances.

5.1 Identification of Type B Applied Parts

All part of the bed (and accessories) the patient can reach are type B Applied Parts.

- Mattress Support Platform Frame, Parts of Mattress Support Platform
- Siderails
- Head Board and Foot Board
- Mattress

5.2 Mechanical Specifications



WARNING!

Risk of injury due to mattress incompatibility!

- ▶ Value of the *Maximum Siderail Height above Mattress Support Platform* must be checked by medical personnel to confirm the compatibility of the medical bed with a mattress prepared for patient!

5.2.1 Mechanical Specifications of the Essenza 300 W

Parameter	Value
External Bed Length	219 cm (221 cm for bed with Accessory Adaptor with 2 Bushings for Lifting Pole at bed head end)
External Bed Width	99,5 cm
Maximum Siderail Height above Mattress Support Platform: with plastic siderails with telescopic siderails with collapsible siderails	42 cm 43 cm 38 cm
Siderail Dimensions (length x height): plastic siderails (head siderail; foot siderail) telescopic siderails (head siderail; foot siderail) collapsible siderails	101,3 cm; 88,5 cm x 43,5 cm 100 cm; 100 cm x 43 cm 162 cm x 38 cm
Maximum Mattress Dimensions (length x width)	200 cm x 90 cm
Maximum Mattress Height: with plastic siderails with telescopic siderails with collapsible siderails	20 cm 20 cm 16 cm
Bed Extension	0 cm / 20 cm / 30 cm
Calfrest Extension	0 cm – 20 cm
Castor Diameter: 125 Tente Linea double castor 150 Tente Linea double castor	12,5 cm 15 cm
Minimum — Middle — Maximum Mattress Support Platform Height above floor (without mattress): Essenza 300 W with 125 Tente Linea double castors Essenza 300 W with 150 Tente Linea double castors	25,5 cm – 37,4 cm – 75,9 cm 28 cm – 39,9 cm – 78,4 cm
Maximum Clearance under the Undercarriage in Standard Position (with 125 Tente Linea double castors)	14,1 cm
Maximum Clearance under the Undercarriage in Standard Position (with 150 Tente Linea double castors)	16,6 cm
Maximum Backrest Angle	70°
Maximum Thighrest Angle	35°
Maximum Calfrest Angle	25°
Maximum Angle between Calfrest and Thighrest	300°
Ergoframe Distances (Backrest Distance / Thighrest Distance)	5 cm / 3 cm
Maximum Trendelenburg Tilt Angle	14°
Maximum Reverse-Trendelenburg Tilt Angle	14°
Average Bed Weight	165 kg
Maximum Lifting Pole Load	75 kg
Safe Working Load (SWL)	300 kg
Maximum Mass of Mobile Hospital Bed (Maximum Mass of Empty Bed + Safe Working Load)	485 kg
Maximum Patient Weight Application environment 1, 2 Application environment 3, 5	235 kg 265 kg
Maximum Sound Pressure Level	54 dBA

5.2.2 Mechanical Specifications of the Essenza 300 W with Lateral Tilt

Parameter	Value
External Bed Length	219 cm (221 cm for bed with Accessory Adaptor with 2 Bushings for Lifting Pole at bed head end)
External Bed Width	99,5 cm
Maximum Siderail Height above Mattress Support Platform: with plastic siderails with telescopic siderails with collapsible siderails	42 cm 43 cm 38 cm
Siderail Dimensions (length x height): plastic siderails (head siderail; foot siderail) telescopic siderails (head siderail; foot siderail) collapsible siderails	101,3 cm; 88,5 cm x 43,5 cm 100 cm; 100 cm x 43 cm 162 cm x 38 cm
Maximum Mattress Dimensions (length x width)	200 cm x 90 cm
Maximum Mattress Height: with plastic siderails with telescopic siderails with collapsible siderails	20 cm 20 cm 16 cm
Bed Extension	0 cm / 20 cm / 30 cm
Calfrest Extension	0 cm – 20 cm
Castor Diameter: 125 Tente Linea double castor 150 Tente Linea double castor	12,5 cm 15 cm
Minimum — Middle — Maximum Mattress Support Platform Height above floor (without mattress): Essenza 300 W with 125 Tente Linea double castors Essenza 300 W with 150 Tente Linea double castors	39,4 cm – 51,3 cm – 89,8 cm 41,9 cm – 53,8 cm – 92,3 cm
Maximum Clearance under the Undercarriage in Standard Position (with 125 Tente Linea double castors)	14,1 cm
Maximum Clearance under the Undercarriage in Standard Position (with 150 Tente Linea double castors)	16,6 cm
Maximum Backrest Angle	70°
Maximum Thighrest Angle	35°
Maximum Calfrest Angle	25°
Maximum Lateral Tilt Angle	15°
Maximum Angle between Calfrest and Thighrest	300°
Ergoframe Distances (Backrest Distance / Thighrest Distance)	5 cm / 3 cm
Maximum Trendelenburg Tilt Angle	14°
Maximum Reverse-Trendelenburg Tilt Angle	14°
Average Bed Weight	180 kg
Maximum Lifting Pole Load	75 kg
Safe Working Load (SWL)	270 kg
Maximum Mass of Mobile Hospital Bed (Maximum Mass of Empty Bed + Safe Working Load)	470 kg
Maximum Patient Weight Application environment 1, 2 Application environment 3, 5	205 kg 235 kg
Maximum Sound Pressure Level	54 dBA

5.3 Electrical Specifications

Parameter	Value
Input Voltage, Frequency	100 — 240 V AC, 50/60 Hz
Maximum Power Input	max. 480 VA
Ingress Protection according to EN 60529	IPX6
Electric Protection Class	Class I (with type B Applied Parts)
Duty Cycle	max. 2 minutes ON / 18 minutes OFF
Bed Accumulator	AKU Pb VRLA 2 x 12 V / 1,2 Ah / Fuse 15A
Control Unit Fuses	2x T4A L 250 V

5.4 Environment Conditions


WARNING!

Risk of damaging the product due to incorrect environment conditions!

► Do not use the Essenza 300 W bed under the environmental conditions outside of those specified in the Environment Conditions chapter!


CAUTION!

Risk of damaging the product if its packaging is exposed to environmental conditions outside of those specified in the Environment Conditions chapter!

Parameter	Value
Use Conditions	
Ambient Temperature	10°C — 40°C
Relative Humidity	30% — 75 %
Atmospheric Pressure	795 — 1060 hPa
Storage and Transport Conditions	
Ambient Temperature	-20°C — 50°C
Relative Humidity	20% — 90 % (non-condensing environment)
Atmospheric Pressure	795 — 1060 hPa

5.5 Electromagnetic Compatibility

Bed is intended for hospitals except for near active HF surgical equipment and the RF shielded room of a medical system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Bed has defined no essential performance based on the bed electric system.



WARNING!

It is recommended to avoid the use of this device next to or in block with other device, because it could lead to improper operation. If such use is needed, this device and the other equipment should be under surveillance to verify proper operation.

List of used cables:

- Mains cable, maximum length 5 m



WARNING!

Use of the accessories, converters and other cables, than specified and provided by manufacturer of this bed could lead to increase of electromagnetic emission or lower the electromagnetic immunity of this bed and lead to improper operation.



WARNING!

Mobile RF communication device (including end use devices like antenna cables and external antenna) should not be used closer than 30 cm (12 inches) from any part of this Essenza 300 W bed, including cables specified by manufacturer. Otherwise this could lead to deterioration of functionality of this bed.



WARNING!

Do not overload the bed (SWL), respect the duty cycle (INT. on the serial label) and consider Maintenance chapter (see „20 Maintenance“ on the page 157) in order to maintain the basic safety with regard to electromagnetic disturbances for the expected service life.

5.5.1 Manufacturer instructions - electromagnetic emissions

Emission Test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies

5.5.2 Manufacturer instructions - electromagnetic susceptibility

Immunity Tests	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV for contact discharge ± 15 kV for air discharge
Radiated RF IEC 61000-4-3 Proximity fields from RF wireless communications equipment IEC 61000-4-3	3 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz See Table 1
Fast electrical transients / burst IEC 61000-4-4	±2 kV for power line repetition frequency 100 kHz
Surge IEC 61000-4-5	± 1 kV Line-to-line ± 2 kV Line-to-ground
Conducted RF IEC 61000-4-6	3 V (0,15 MHz – 80 MHz) 6 V in ISM bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m
Voltage dips, short interruptions on power supply input lines IEC 61000-4-11	0 % UT; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycle Single phase: at 0° 0 % UT; 250/300 cycle

Table 1 - IMMUNITY to RF wireless communications equipment

Test frequency (MHz)	Band (MHz)	Service	Modulation	Immunity Test Level V/m
385	380 - 390	TETRA 400	Pulse modulation 18 Hz	27
450	430 - 470	GMRS 460, FRS 460	FM ± 5 kHz deviation 1 kHz sine	28
710 745 780	704 - 787	LTE band 13, 17	Pulse modulation 217 Hz	9
810 870 930	800 - 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE band 5	Pulse modulation 18 Hz	28
1 720 1 845 1 970	1 700 - 1 990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	28
2 450	2 400 - 2 570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE band 7	Pulse modulation 217 Hz	28
5 240 5 500 5 785	5 100 - 5 800	WLAN 802.11 a/n	Pulse modulation 217 Hz	9

NOTE There are applied no deviations to requirements of IEC 60601-1-2 ed. 4.

NOTE There are no known other measures for keeping the basic safety based on EMC phenomena.

6 Use and Storage Conditions



DANGER!

Danger to life due to electric shock!

To ensure the bed's class I protection against electric shocks:

- ▶ Ground the mains.
- ▶ Use exclusively Hospital Grade or Hospital Only receptacles for grounding.

Essenza 300 W is designed for use in rooms for medical purposes. Electrical installations connected to the Essenza 300 W must therefore meet local norms laying down the necessary conditions for electrical installations.

To increase the safety of electrical equipment:

- ▶ Disconnect the Essenza 300 W from the mains in exceptional cases (i.e. lightnings, earthquake).

Respect values of the parameters connected with environment conditions in the chapter Technical Specification during use and storage of the product.

Essenza 300 W is not suitable for indoor environments containing flammable gases (except oxygen cylinders).

Essenza 300 W is not suitable for continuous operation. Duty cycle equals to 2 minutes in operation during 20 minutes.

7 Scope of Delivery and Product Variants

7.1 Delivery

- ▶ Upon receipt, check that the shipment is complete as specified on the delivery note.
- ▶ Notify the carrier and supplier of any deficiencies or damages immediately as well as in writing or make a note on the delivery note.

7.2 Scope of Delivery

- Essenza 300 W Medical bed
- Instructions for Use

7.3 Essenza 300 W Variants

Variable bed features:

- **Lateral Tilt**
 - Essenza 300 W without Lateral Tilt
 - Essenza 300 W with Lateral Tilt
- **Castors**
 - Central Castor Control System (castors with diameter 125 mm, castors with diameter 150 mm, double castors)
 - Fifth Castor
 - Brake Signal
- **Mattress Support Platform Extension**
 - Simple Bed Extension
 - Advanced Bed Extension
- **Patient Restraint Holders**
- **Backrest Handles**
- **Accessory Rail**
 - at foot end of the bed
 - at foot end of the bed and at head end of the bed
- **Siderails**
 - Plastic Siderails (compatible with Accessory Adaptor with 2 Bushings for Lifting Pole in the undercarriage at bed head end)
 - Telescopic Siderails with DIN Rails (not compatible with Accessory Adaptor with 2 Bushings for Lifting Pole in the undercarriage at bed head end)
 - Telescopic Siderails without DIN Rails (not compatible with Accessory Adaptor with 2 Bushings for Lifting Pole in the undercarriage at bed head end)
 - Collapsible Siderails (not compatible with Accessory Adaptor with 2 Bushings for Lifting Pole in the undercarriage at bed head end)
- **Head Board and Foot Board**
 - Standard Plastic Head Board and Foot Board with One Safety Lock (not compatible with Telescopic Siderails)
 - Narrow Plastic Head Board and Foot Board with One Safety Lock (compatible with Accessory Adaptor with 2 Bushings for Lifting Pole at bed head end)
 - Plastic Head Board and Foot Board with Two Safety Locks (only compatible with Telescopic Siderails)
 - HPL Head Board and Foot Board (only compatible with Telescopic Siderails)
- **Corner Bumpers**
 - with accessory bushings
 - 4 Corner Bumpers in outer corners of the mattress support platform
 - 2 Corner Bumpers in foot end outer corners of the mattress support platform (compatible with Accessory Adaptor with 2 Bushings for Lifting Pole at bed head end)
 - without accessory bushings
 - 4 Corner Bumpers with vertical bumpers in outer corners of the mattress support platform
 - 2 Corner Bumpers with vertical bumpers in head end outer corners of the mattress support platform and 2 Corner Bumpers without vertical bumpers in foot end outer corners of the mattress support platform
- **Bushings for Lifting Pole**
 - 2 Bushings for Lifting Pole in inner corners of the mattress platform at bed head end
 - 2 Bushings for Lifting Pole in inner corners of the mattress platform at bed head end and Accessory Adaptor with 2 external Bushings for Lifting Pole at bed head end (only compatible with Plastic Siderails and with Narrow Plastic Head Board without side corner bumpers at bed head end)
- **Control Elements**
 - Attendant Control Panel for Essenza 300 W with Lateral Tilt
 - Attendant Control Panel for Essenza 300 W
 - Siderail Control Panel for Essenza 300 W with Lateral Tilt
 - Siderail Control Panel for Essenza 300 W
 - Handset Control Panel
 - Patient Control Panel
- **Back-up Accumulator**
- **Bed Illumination**
- **Safe Height Visual Signalling (green light)**
- **Extra Low Position**
- **Plastic Linen Shelf**
- **Telescopic Linen Shelf**

8 Putting into Service



WARNING!

Risk of injury when working on the bed!

- ▶ Ensure that the bed is disconnected from the mains connection prior to putting into service, putting out of service and maintenance.
- ▶ Ensure that the castors are locked prior to putting into service, putting out of service and maintenance.



CAUTION!

Risk of material damage due to incorrect putting into service!

- ▶ Ensure that putting into service is performed exclusively by manufacturer's customer service or trained hospital personnel.



CAUTION!

Risk of material damage due to temperature difference!

- ▶ If there is a considerable temperature difference between the bed and the place of operation (after transport/storage), leave bed unconnected for 24 hours for the difference to balance itself.

NOTE For safe, easy handling, LINET ® recommends having two technicians put the bed into service.

Set up the bed as follows:

- ▶ Unpack the bed.
 - ▶ Check the delivery (see „7 Scope of Delivery and Product Variants“ on the page 34).
 - ▶ Remove isolating foil from the accumulator box (see „8.1 Accumulator Activation“ on the page 37).
 - ▶ Raise siderails up.
 - ▶ Install equipment and accessories.
 - ▶ In case of delivery with dismantled head board and foot board, mount the head board and foot board correctly.
 - ▶ Set up the bed exclusively on a suitable floor surface (see „8.11 Transport“ on the page 50).
 - ▶ Ensure that the Power Supply Cord is not in collision with something and it is not stretched when adjusting the bed.
 - ▶ Check that the plug is inserted correctly into the mains.
 - ▶ Do not leave any extension cords or power strips loose on the floor.
 - ▶ Ensure that all of the required mechanical and electrical prevention mechanisms are available on site.
 - ▶ There is no mains switch on the bed, i.e. the Power Supply Cord is the only means to isolate the bed from the mains.
- Ensure that the Power Supply Cord is always accessible.
- ▶ Have the separable plug of the Power Supply Cord changed and maintained exclusively by qualified and trained service technicians authorised by the manufacturer if the Power Supply Cord is damaged or most likely damaged.

8.1 Accumulator Activation

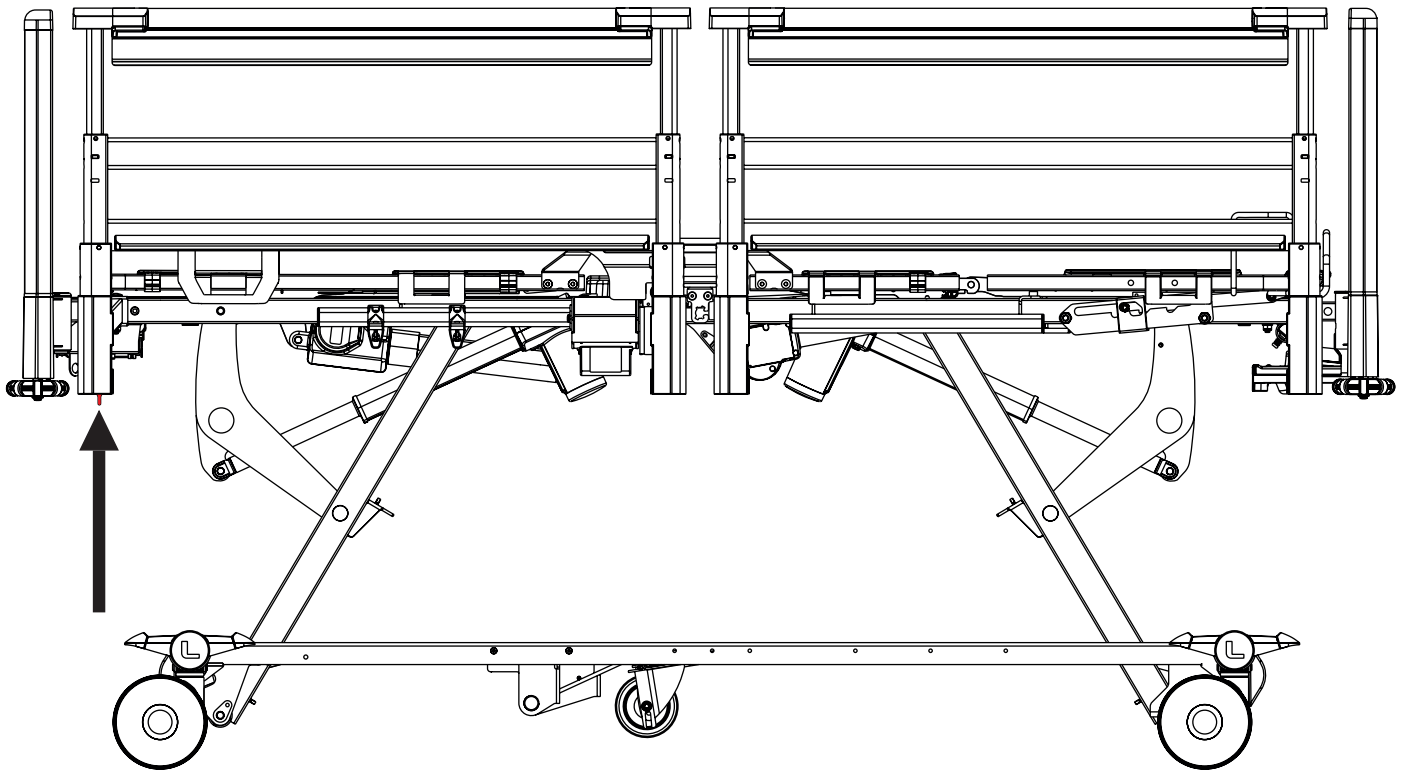


Fig. Position of Isolating Foil under the Backrest

8.1.1 Removing the Isolating Foil

To remove isolating foil:

- ▶ Remove isolating foil from the control unit under the Backrest by pulling the strap.
- ▶ Check if isolating foil is complete and undamaged as shown in figure Shape of Isolating Foil.
- ▶ If isolating foil is damaged, contact the manufacturer's service department immediately.

NOTE Isolating Foil is sharp-edged. Remove it carefully not to cut yourself.


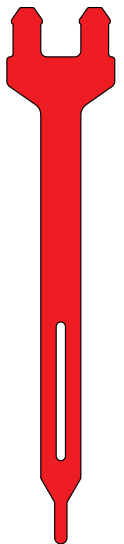


Fig. Shape of Isolating Foil

8.2 Plastic Foot Board with one Safety Lock

There are 2 types of the Plastic Head Board and Foot Board depending on the compatibility of Head Board with additional 2 Bushings for Lifting Pole in the undercarriage at bed head end. Only the Narrow Head Board without side corner bumpers is compatible with the additional 2 Bushings for Lifting Pole. Standard Head Board and Foot Board are interchangeable in the case of the compatible configuration without the additional 2 Bushings for Lifting Pole.

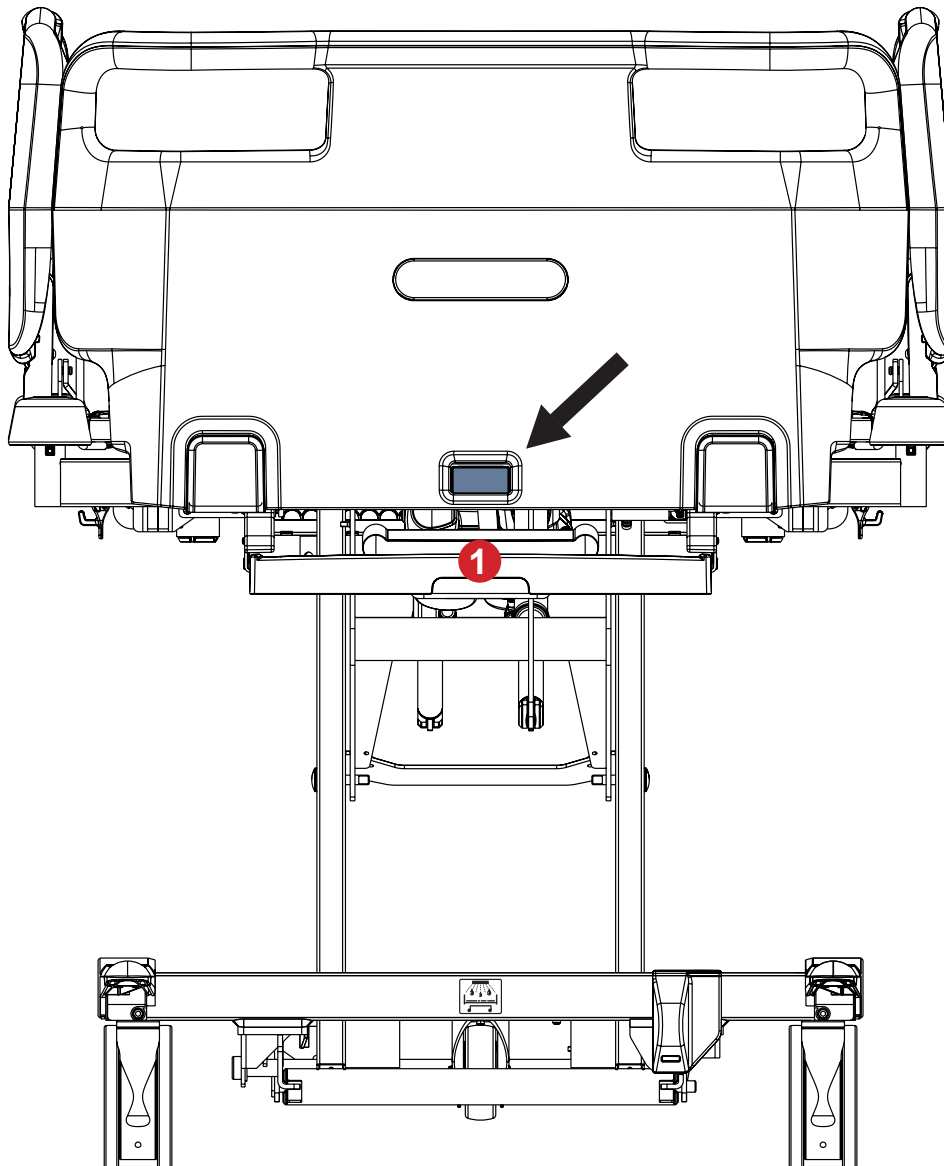

1. Attendant Control Panel in the Plastic Linen Shelf

Fig. Standard Plastic Foot Board

Dismount the Foot Board as follows:

- ▶ Unlock Foot Board Safety Lock.
- ▶ Pull Foot Board from sleeve fittings.

Install the Foot Board as follows:

- ▶ Unlock Foot Board Safety Lock.
- ▶ Slide Foot Board into sleeve fittings.
- ▶ Lock Foot Board Safety Lock.

8.3 Plastic Head Board with one Safety Lock



CAUTION!

Risk of collision due to a Standard Head Board incompatible with the additional 2 Bushings for Lifting Pole in the undercarriage at bed head end!

- Ensure the Head Board you want to fix at bed head end is compatible with your bed configuration if Standard Head Board and Narrow Head Board are available in your ward!

There are 2 types of the Plastic Head Board and Foot Board depending on the compatibility of Head Board with additional 2 Bushings for Lifting Pole in the undercarriage at bed head end. Only the Narrow Head Board without side corner bumpers is compatible with the additional 2 Bushings for Lifting Pole. Standard Head Board and Foot Board are interchangeable in the case of the compatible configuration without the additional 2 Bushings for Lifting Pole.

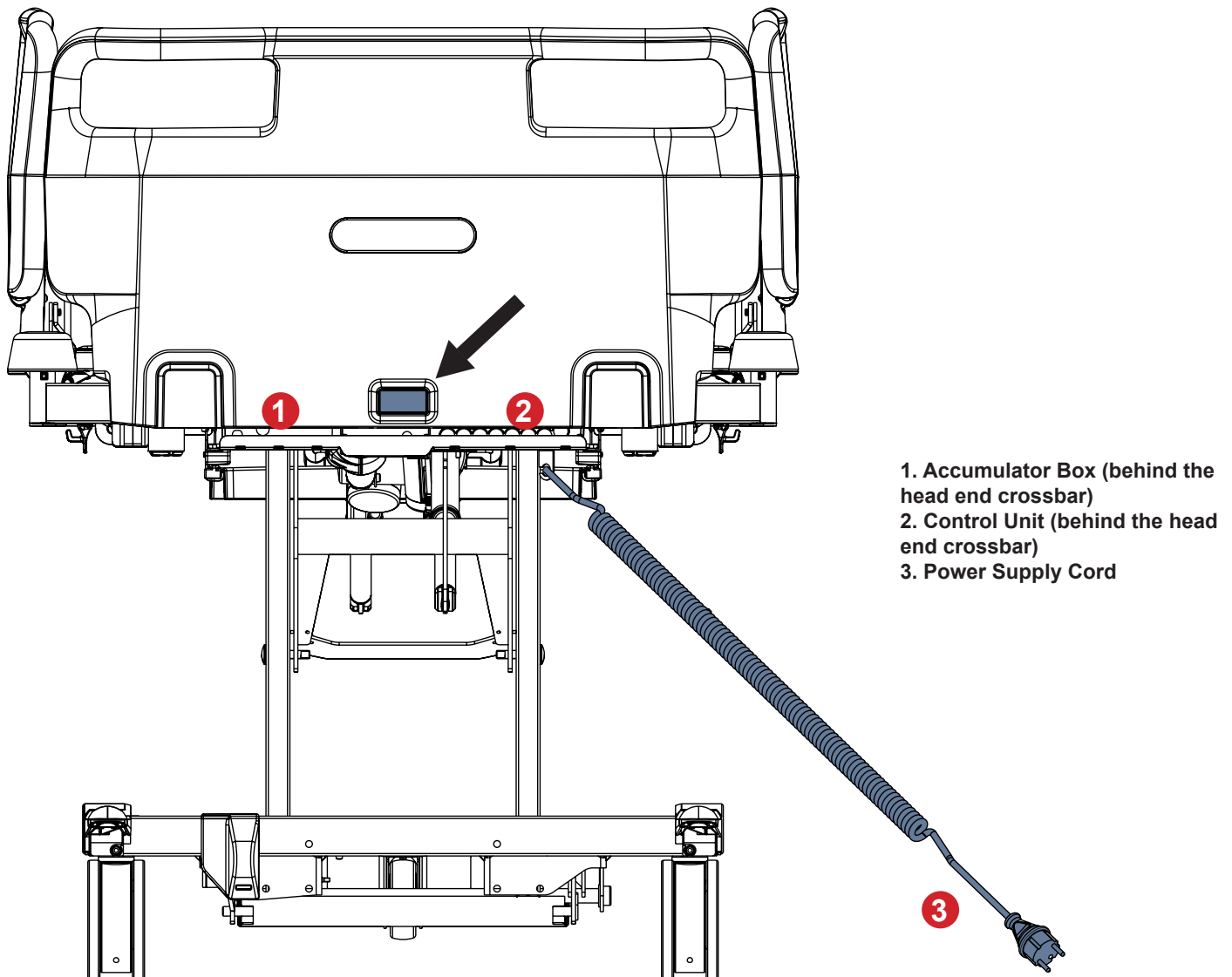


Fig. Standard Plastic Head Board

Dismount the Head Board as follows:

- Unlock Head Board Safety Lock.
- Pull Head Board from sleeve fittings.

Install the Head Board as follows:

- Unlock Head Board Safety Lock.
- Slide Head Board into sleeve fittings.
- Lock Head Board Safety Lock.


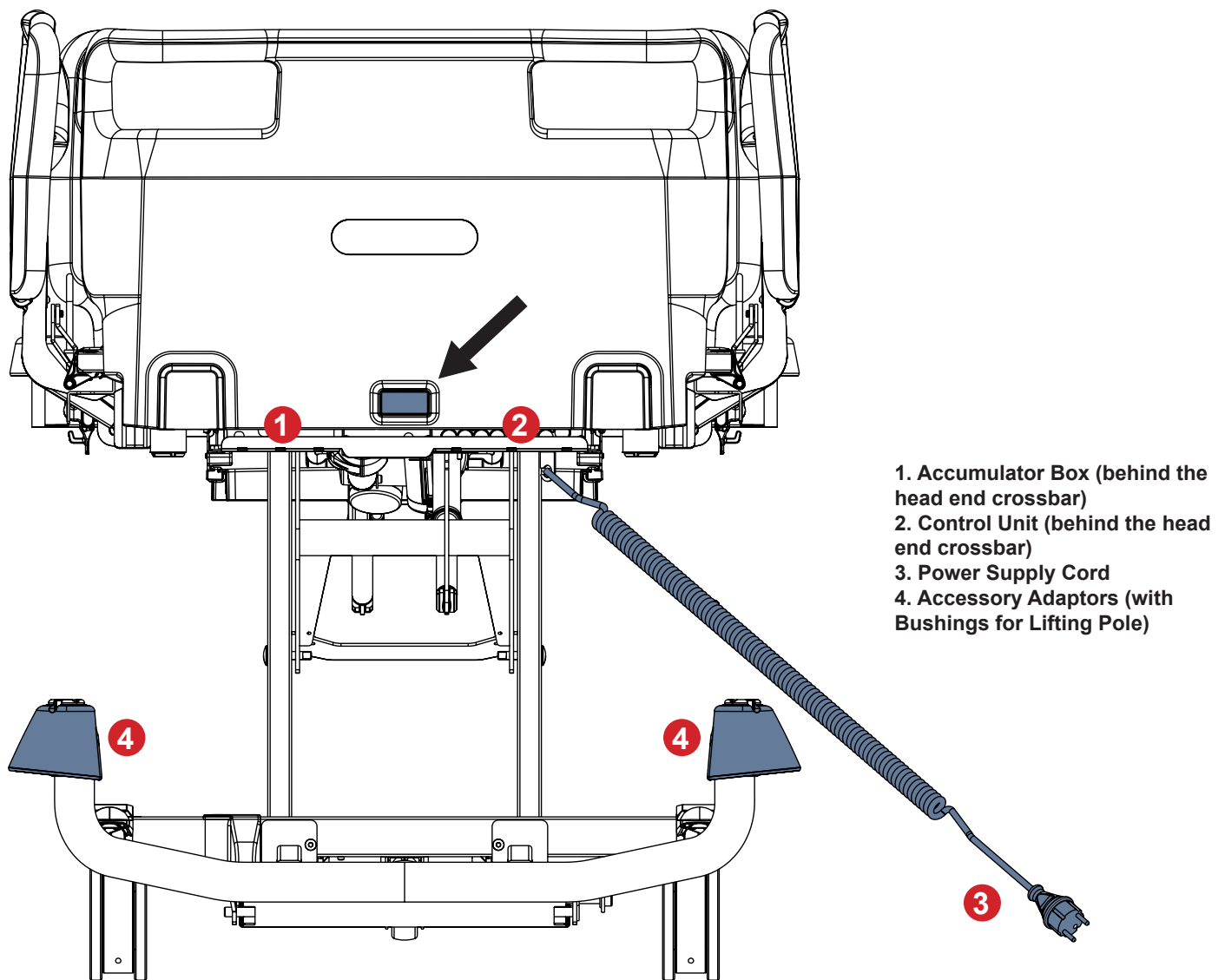


Fig. Narrow Plastic Head Board

8.4 HPL Foot Board

HPL Head Board and HPL Foot Board are only compatible with telescopic siderails. Essenza 300 W bed compatible with HPL Head Board and HPL Foot Board is equipped with two Head Board Locks and two Foot Board Locks. HPL Head Board and HPL Foot Board are interchangeable in the case of the compatible configuration with telescopic siderails.

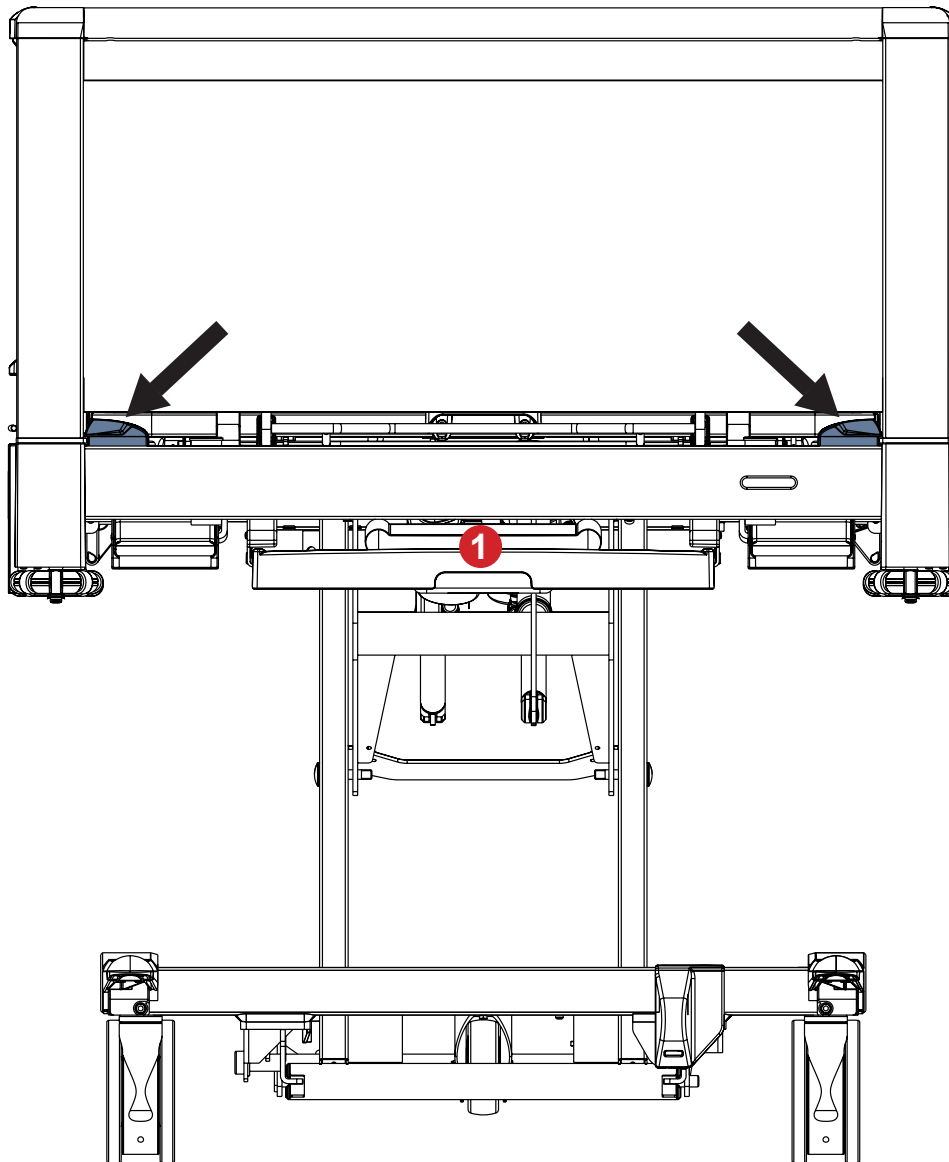

1. Attendant Control Panel in the Plastic Linen Shelf

Fig. HPL Foot Board

Insert the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Slide head board or foot board into slots on corner posts.
- ▶ Lock safety locks at corner posts.

Remove the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Pull head board or foot board upward.

8.5 HPL Head Board

HPL Head Board and HPL Foot Board are only compatible with telescopic siderails. Essenza 300 W bed compatible with HPL Head Board and HPL Foot Board is equipped with two Head Board Locks and two Foot Board Locks. HPL Head Board and HPL Foot Board are interchangeable in the case of the compatible configuration with telescopic siderails.

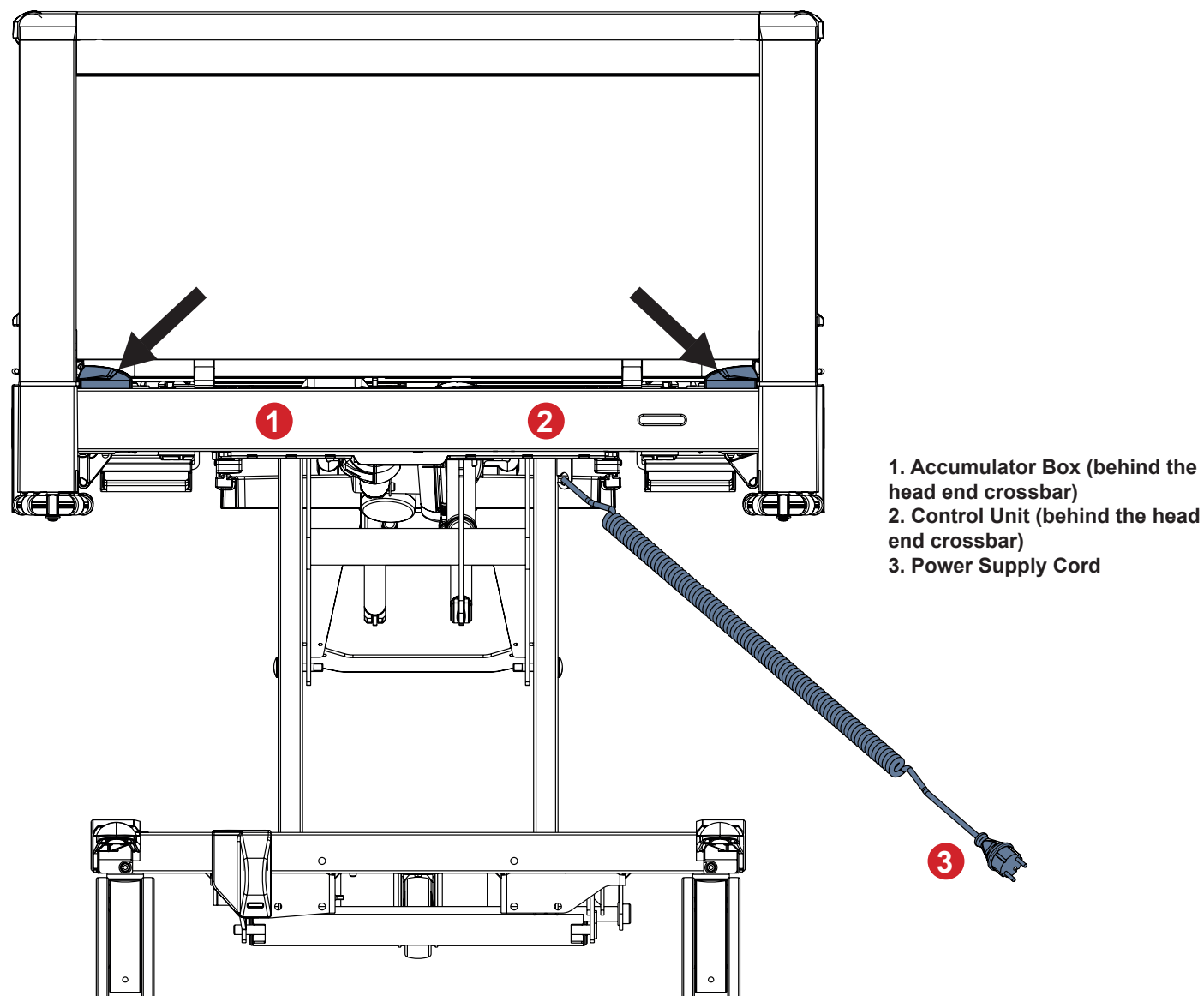


Fig. HPL Head Board

Insert the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Slide head board or foot board into slots on corner posts.
- ▶ Lock safety locks at corner posts.

Remove the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Pull head board or foot board upward.

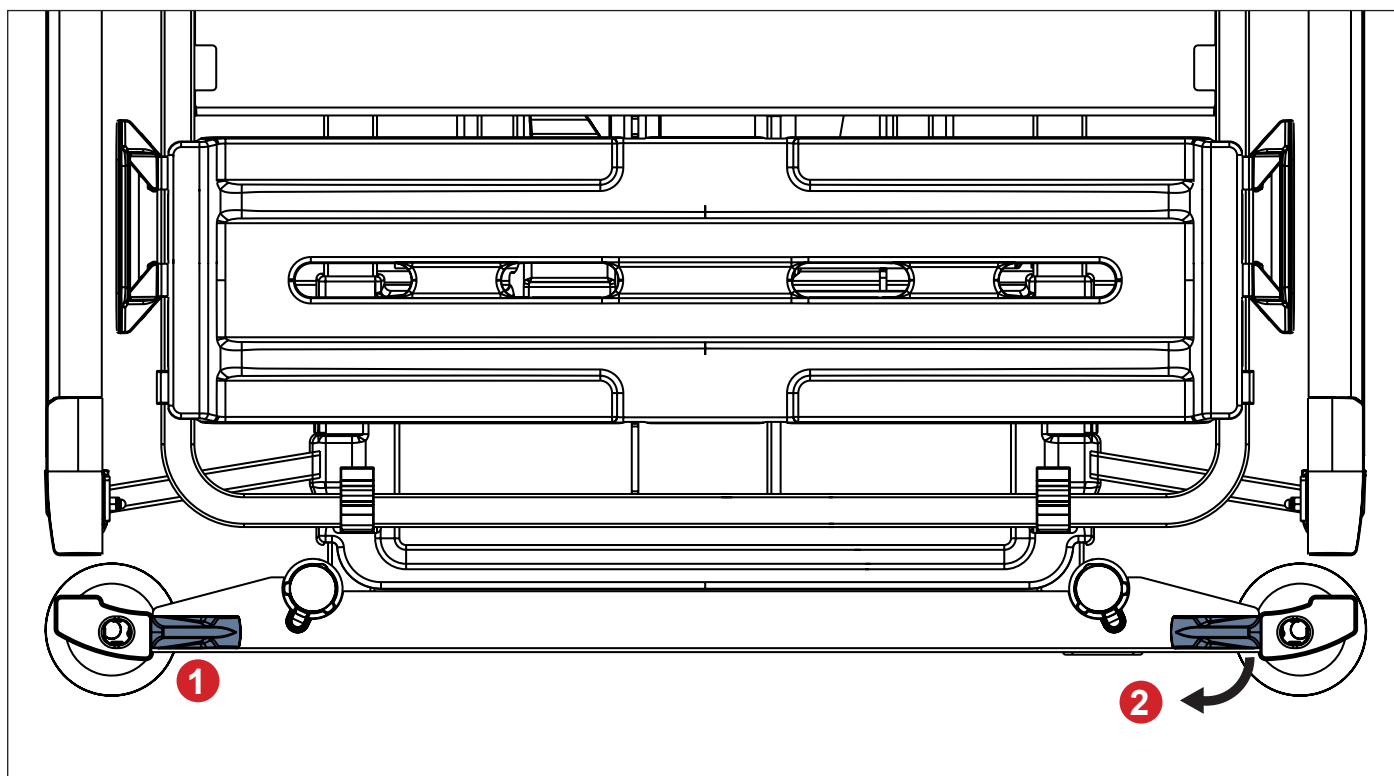


Fig. Two Head Board Safety Locks

- 1. Locked Lock
- 2. Unlocking

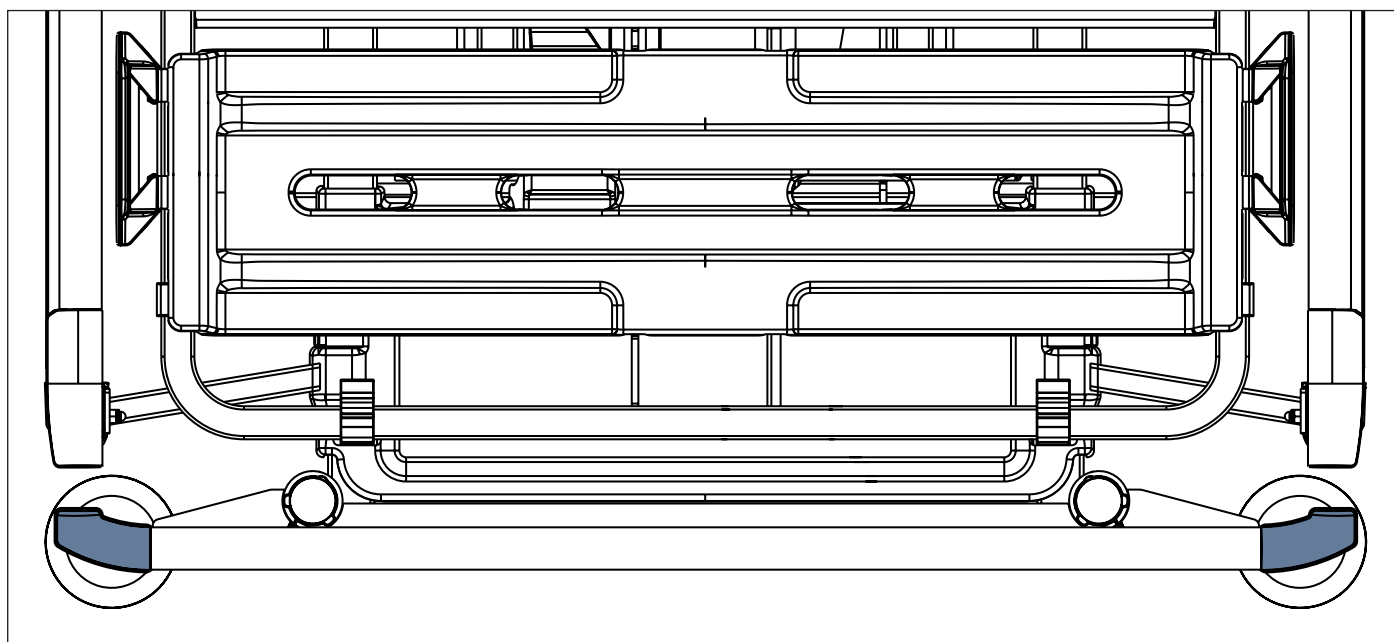


Fig. Correctly inserted HPL Head Board

8.6 Plastic Foot Board with two Safety Locks

Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks are only compatible with bed equipped with corner bumpers without accessory bushings. Essenza 300 W bed compatible with Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks is equipped with two Head Board Locks and two Foot Board Locks. Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks are interchangeable in the case of the compatible bed configuration with corner bumpers without accessory bushings.

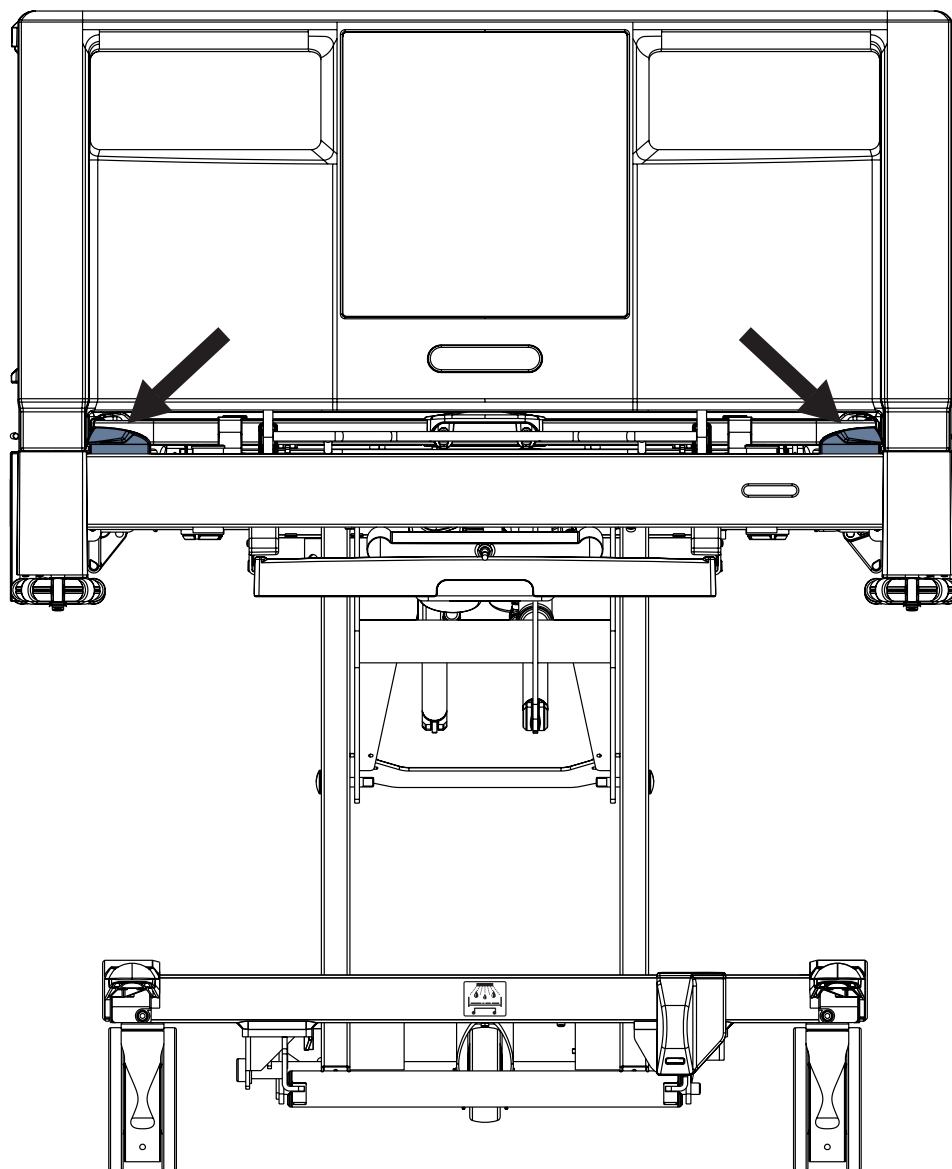


Fig. Plastic Foot Board with two Safety Locks

Insert the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Slide head board or foot board into slots on corner posts.
- ▶ Lock safety locks at corner posts.

Remove the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Pull head board or foot board upward.

8.7 Plastic Head Board with two Safety Locks

Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks are only compatible with bed equipped with corner bumpers without accessory bushings. Essenza 300 W bed compatible with Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks is equipped with two Head Board Locks and two Foot Board Locks. Plastic Head Board with two Safety Locks and Plastic Foot Board with two Safety Locks are interchangeable in the case of the compatible bed configuration with corner bumpers without accessory bushings.

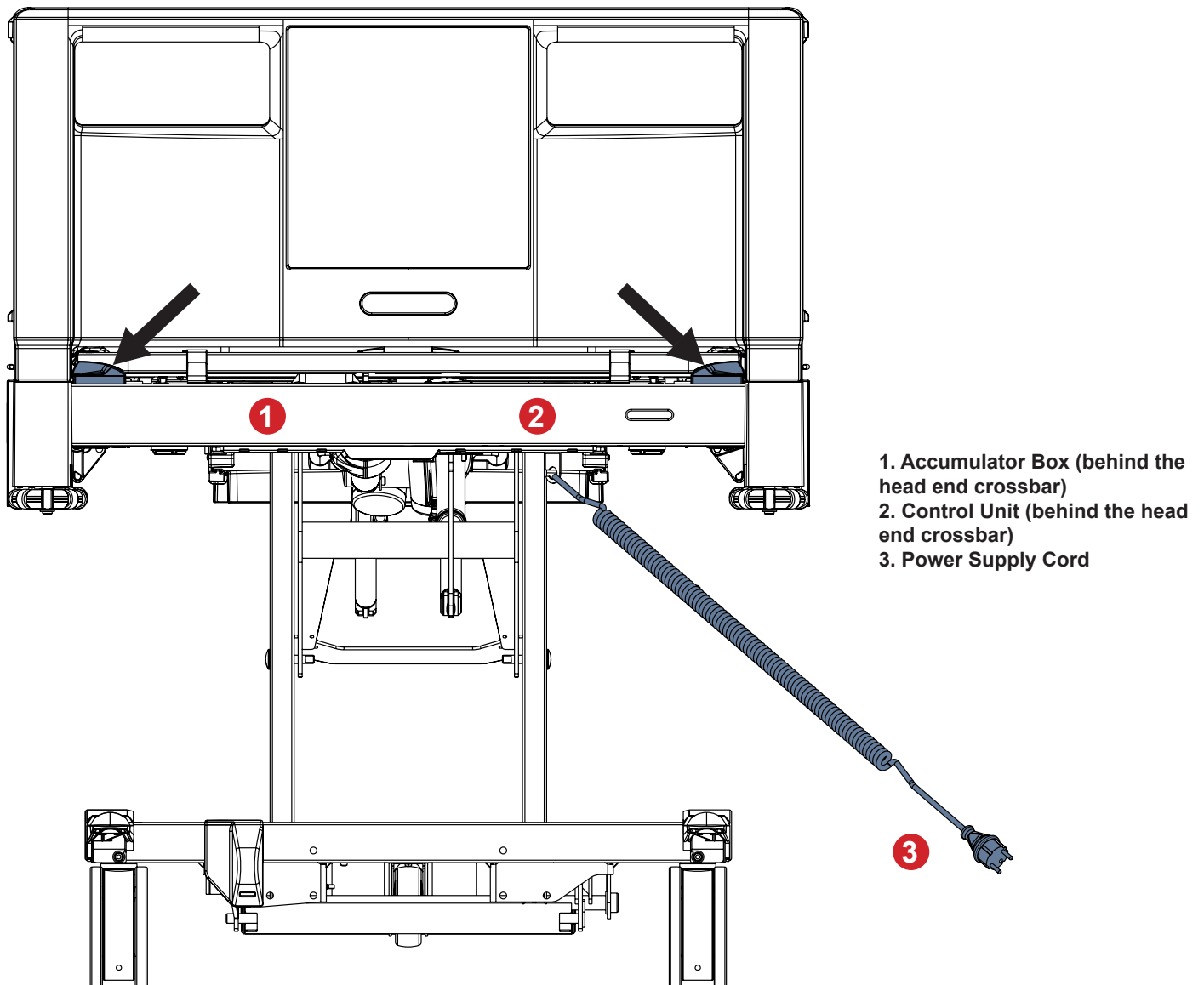


Fig. Plastic Head Board with two Safety Locks

Insert the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Slide head board or foot board into slots on corner posts.
- ▶ Lock safety locks at corner posts.

Remove the head board or foot board as follows:

- ▶ Unlock safety locks at corner posts.
- ▶ Pull head board or foot board upward.

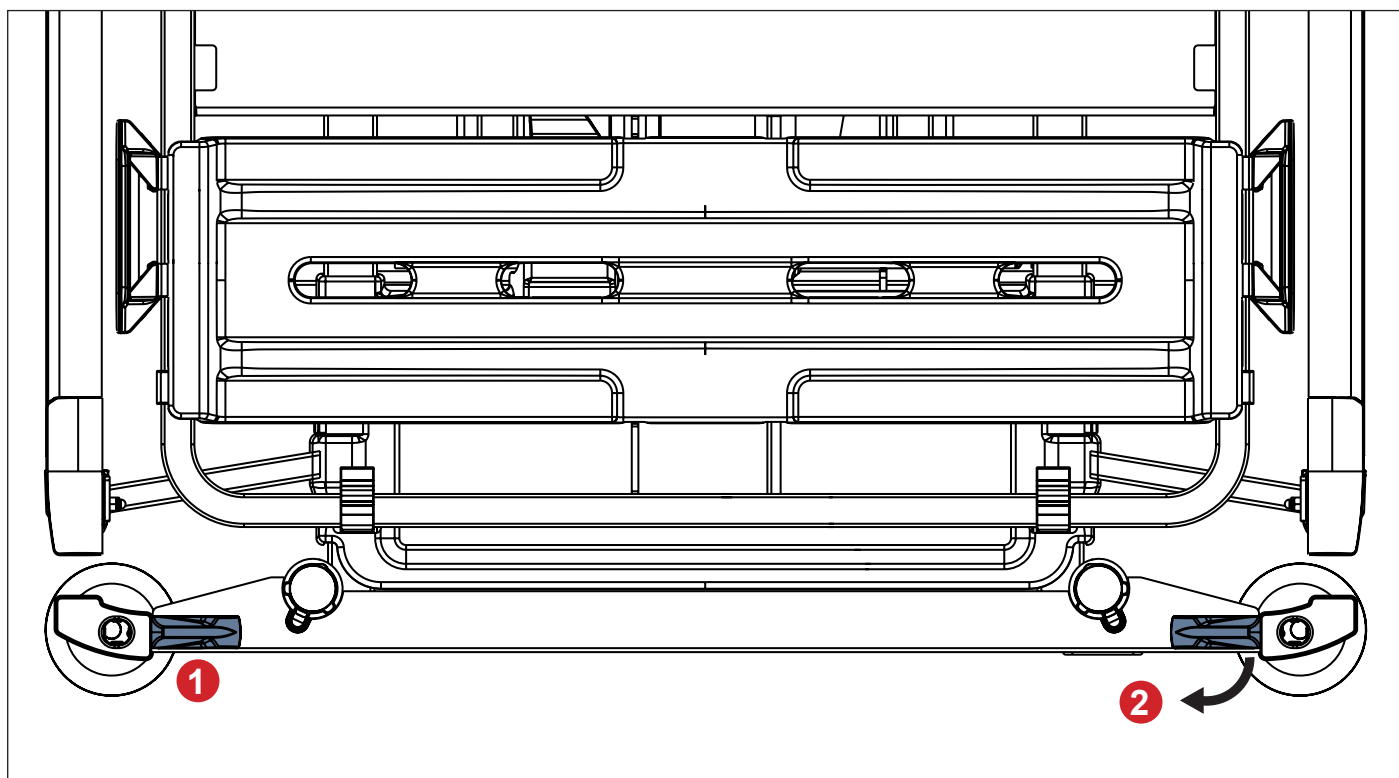


Fig. Two Head Board Safety Locks

- 1. Locked Lock
- 2. Unlocking

8.8 Mattress Support Platform

Essenza 300 W has Mattress Support Platform with four sections.


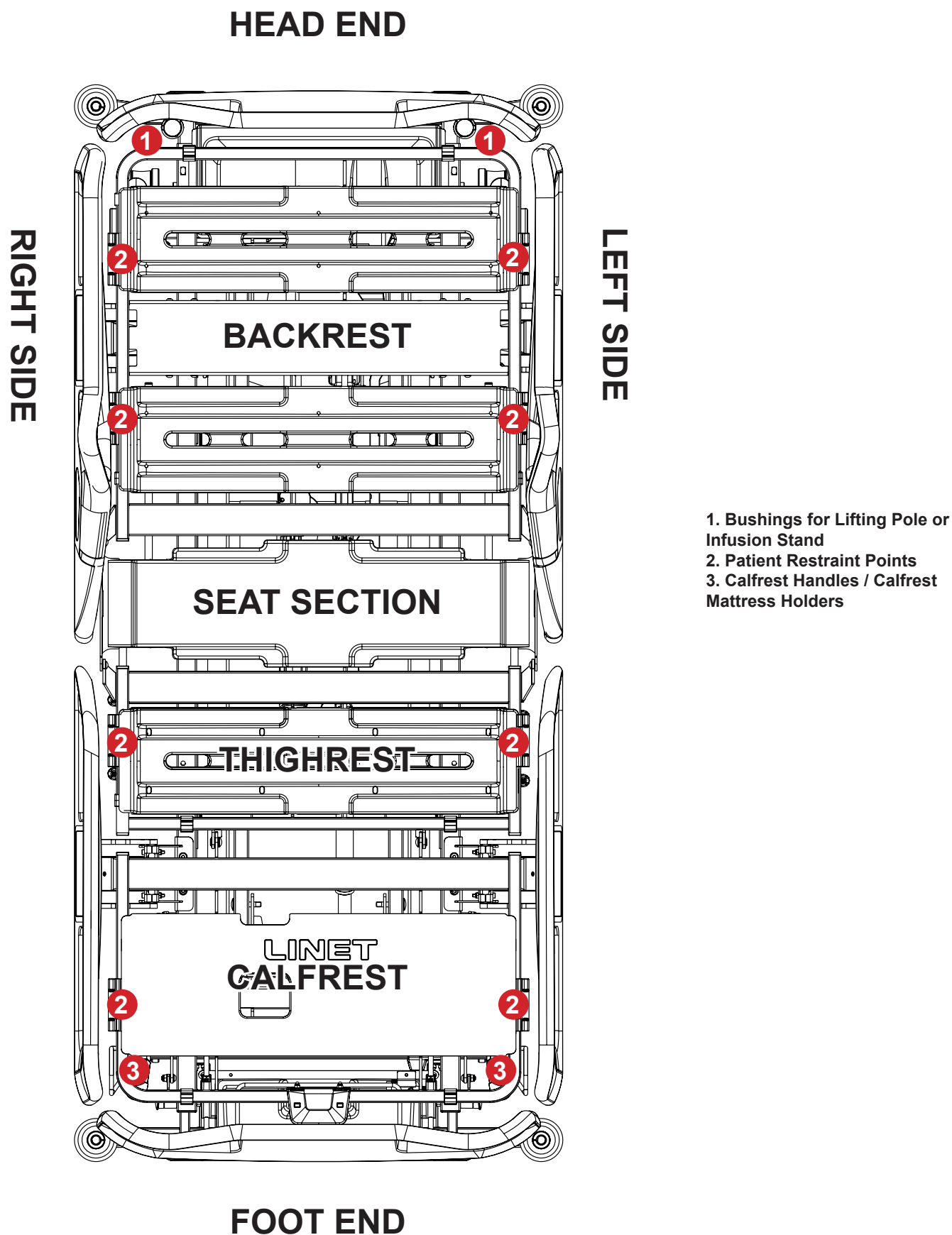


Fig. 4-part Mattress Support Platform with Plastic Siderails

HEAD END

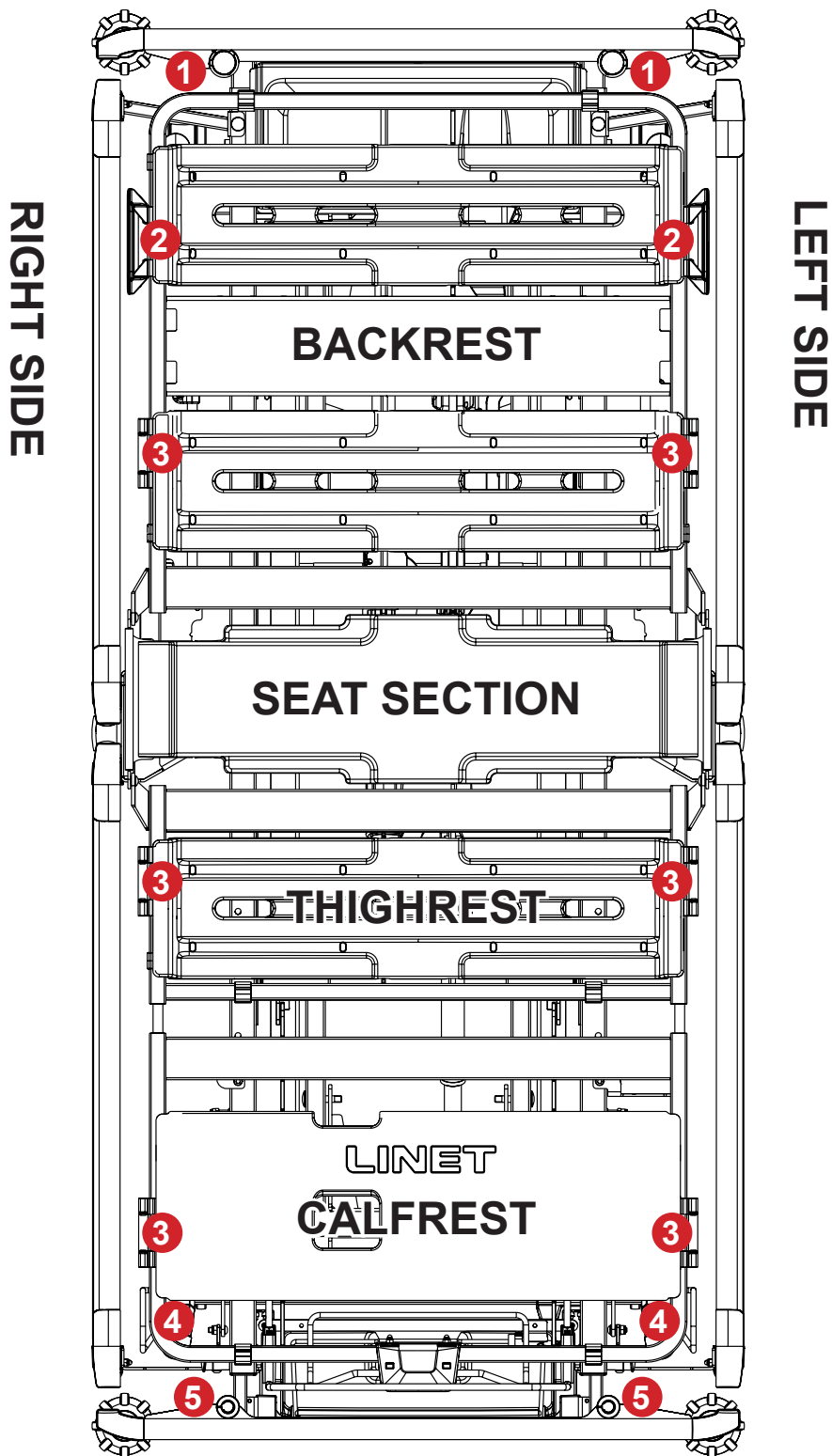

- 1. Bushings for Lifting Pole or Infusion Stand
- 2. Backrest Handles / Backrest Mattress Holders
- 3. Patient Restraint Points
- 4. Calfrest Handles / Calfrest Mattress Holders
- 5. Bushings for accessories

FOOT END

Fig. 4-part Mattress Support Platform with Telescopic Siderails

8.8.1 Handles and Patient Restraint Points

Variants of Handles / Mattress Holders:

- on the sides of the Backrest upper part (version with telescopic siderails only)
- on the sides of the Calfrest lower part

Variants of Patient Restraint Points:

- on the sides of the Backrest lower part
- on the sides of the Thighrest upper part
- on the sides of the Calfrest upper part

8.9 Potential Interconnection



DANGER!

Danger to life due to electric shock!

To protect patient against static electric shock:

- ▶ Connect potential equalisation male connector to the potential equalisation female connector.

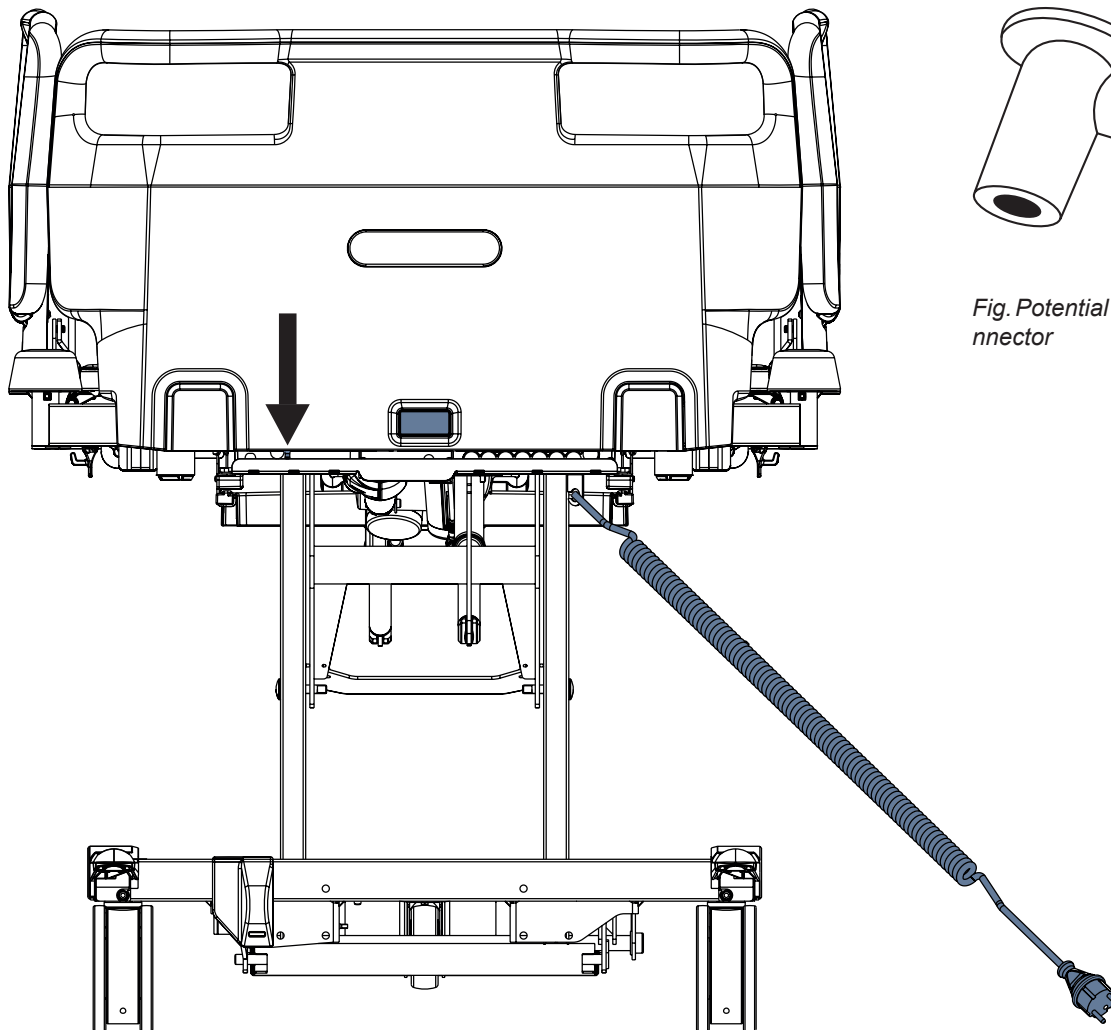


Fig. Potential equalisation male connector at head end of the bed

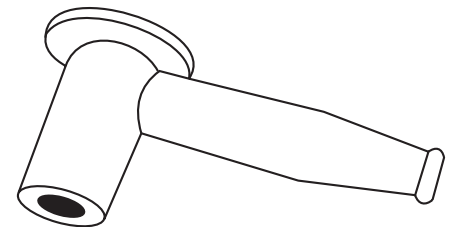


Fig. Potential equalisation female connector

The bed is equipped with a standard protective connector. This connector is used for potential equalisation between the bed and any intravascular or intracardiac device connected to the patient to protect the patient from static electric shocks.

Use equalisation connector if the patient is connected to any intravascular or intracardiac device.

Before connecting the patient to an intravascular/intracardiac device:

- ▶ Connect the ground wire of the device to the potential equalisation connector on the bed on which the patient in question is lying.
- ▶ Use a standard hospital connector.
- ▶ Make sure that the connectors match.
- ▶ Make sure that there is no possibility for inadvertent disconnection.

Before moving the bed:

- ▶ Disconnect the patient from the intravascular or intracardiac device.
- ▶ Disconnect the potential equalisation connector.

8.10 Before Use

Prepare the bed for use as follows:

- ▶ Connect the bed to the mains.
- ▶ Charge the accumulator.
- ▶ Raise and tilt the mattress support platform to the highest position.
- ▶ Lower and tilt the mattress support platform to the lowest position.
- ▶ Check that the castors as well as brake levers work correctly.
- ▶ Check that the bed extension works correctly.
- ▶ Check that it is possible to remove the head board and foot board.
- ▶ Check all of the functions on the control elements.
- ▶ Check that the siderails function properly.
- ▶ Dispose of all packaging (see „21 Disposal“ on the page 164).

8.11 Transport

For a safe transport, observe the following:

- ▶ Ensure that no cables are run over when moving the bed.
- ▶ Ensure that the power supply cord is attached with a hook (at the head end of the bed).
- ▶ Adjust bed height to at least 20 cm below maximum height.
- ▶ Push bed by handles on Head Board or Foot Board.
- ▶ Move the bed exclusively on suitable floor surfaces.
- ▶ Ensure the bed is braked when it should not move (see „11.4 Castor Control“ on the page 64).
- ▶ For longer distances, ensure the castor direction determining function is activated.
- ▶ Ensure that the castors are not braked before moving the bed.

Suitable surfaces:

- Tile
- Hard linoleum
- Poured flooring

Unsuitable surfaces:

- Too soft, unsealed or defective flooring
- Soft wooden flooring
- Soft and porous stone floors
- Carpeted floors with underlay
- Soft linoleum

8.11.1 Transport Position



WARNING!

Hospital staff is responsible for assessing the suitable adjustment of the mattress support platform in accordance with patient's state and needs!



WARNING!

Placement of the accessories and of the compatible medical devices must not cause any hazard for patient or for hospital staff or for the other persons and any damage to the product or to the surroundings during transport of the Essenza 300 W bed!



WARNING!

Siderails in down position can cause patient's fall from the Essenza 300 W bed during transport of the Essenza 300 W bed!

Transport Position of the Essenza 300 W bed:

- ▶ siderails up and locked
- ▶ height of the mattress support platform in accordance with ergonomically suitable drive
- ▶ not tilted mattress support platform
- ▶ parts of the mattress support platform adjusted according to the patient's state and needs
- ▶ Power Supply Cord disconnected from the mains power and correctly inserted to the Hook for hanging Power Supply Cord

8.12 Firmware

The control unit of the bed includes firmware that can be updated only by an authorised service technician.

9 Power Supply Cord



WARNING!

Where the integrity of the external protective conductor in the installation or its arrangement is in doubt
► Use the other power source that is not in doubt.



WARNING!

Inappropriate handling of the power supply cord, e. g. by kinking, shearing or other mechanical damages is hazardous!



WARNING!

When routing Power Supply Cord in a Essenza 300 W avoid squeezing the cable between parts of the Essenza 300 W!



WARNING!

An additional multiple socket-outlet or extension cord shall not be connected to the Essenza 300 W!



WARNING!

Use exclusively the power supply cord in perfect condition!
► Ensure the power supply cord is not damaged before each use!




WARNING!

Power Supply Cord must not be tightened or stretched during use!
► Ensure the Power Supply Cord hangs loosely between the Essenza 300 W and the wall where the Power Supply Cord is connected to the mains.



WARNING!

Plug of the Power Supply Cord must not be wet or damp before use!
► Ensure the plug is completely dry before connecting it to the mains.

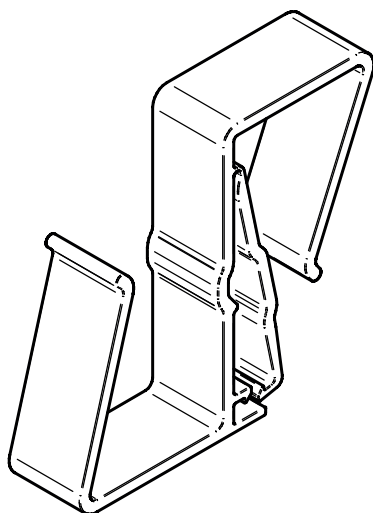

Attachment plug is means of connecting and disconnecting Essenza 300 W from the mains.

Fig. Hook for hanging Power Supply Cord

10 Accumulator of the bed



WARNING!

Risk of material damage to the bed accumulator due to incorrect charging!

- ▶ When the bed is not connected to the mains and accumulator is not sufficiently charged all electrical functions of the bed are blocked! After some long period without charging the battery, the battery will be deeply discharged and lose its full functionality!



WARNING!

Risk of damaging the product due to incorrect maintenance!

- ▶ Only authorised and trained service technician is allowed to replace the bed accumulator! Shape of the accumulator box does not allow different ways of insertion.

Purpose

The accumulator serves as a backup during power failures or for emergency bed positioning.

- ▶ Use only accumulator approved by the manufacturer.
- ▶ Check the functionality of the accumulator at least once a month and have the accumulator changed if necessary.

The manufacturer will assume no responsibility for any damage to the bed or the accumulator caused by:

- non-observance of the manufacturer's instructions in the instructions for use.
- using accumulator not approved by the manufacturer.

Warranty

The manufacturer provides a **6-month** warranty for the full function of the accumulator.

Bed Accumulator lifetime could be up to 5 years if operated under optimum conditions.

Accumulator capacity can be significantly reduced if:

- ▶ too high ambient temperature
- ▶ many accumulator charge/discharge cycles
- ▶ recurrence of deep discharge
- ▶ bed is often powered only by the accumulator

Charging

The accumulator supplied with the bed is delivered insufficiently charged.

Charging of the accumulator before use of the bed takes approximately 4 hours.

Charging of the accumulator runs automatically when the power cable is connected to the mains.

To charge the accumulator:

- ▶ Connect the bed to the mains.

Storage

For declared lifetime period of leaded accumulators is recommended during storage:

- ▶ To prevent accumulator from deep discharging and to keep accumulator at least partly charged by regular recharging.
- ▶ To store accumulator on the dry places with temperature from 10°C to 40°C.
- ▶ To prevent accumulators from being in the sunshine.

Signalling

The corresponding indicators on the Attendant Control Panel and the Siderail Control Panel indicate the accumulator charge status.

10.1 Replacing the accumulator



CAUTION!

Damage to the bed due to incorrect accumulator replacement!

- ▶ Have the accumulator replaced exclusively by qualified personnel.
- ▶ Exclusively use the accumulator approved by the manufacturer.



CAUTION!

Material damage due to overheating!

If the accumulator is faulty, degassing may occur. In rare cases this might cause deformations of the accumulator case, control panel housing or cable.

- ▶ Stop using the bed immediately (see „10.2 Removing the Bed from Service“ on the page 54).
- ▶ Inform the manufacturer's service department.



CAUTION!

Risk of reducing accumulator durability due to incorrect use!

- ▶ Use bed on accumulator only in crisis situations (e.g.: power blackout, patient complications during transport, etc.)
- ▶ After reconnecting bed to the mains charge accumulator to full capacity (see chart Accumulator charge status).
- ▶ Have the accumulator replaced exclusively by a qualified service organisation of the manufacturer.

Faulty accumulator

The accumulator is regarded as faulty if at least one of the following conditions applies:

- ▶ Accumulator charging constantly
- ▶ Low voltage on accumulator
- ▶ Low charging current of accumulator

This status is indicated by the accumulator status indicator being constantly lit. These statuses are summarised and written to black box in the control unit.

To cancel this status:

- ▶ Press **STOP** button.

Discharged accumulator

This status is indicated by the accumulator status indicator flashing quickly.

This status will be cancelled automatically when the bed switches to sleep mode (bed is disconnected from the mains, no button is pressed).

To cancel this status:

- ▶ Press **STOP** button.

10.2 Removing the Bed from Service

Remove the bed from use as follows:

- ▶ Disconnect the bed from the mains.
- ▶ Disconnect the ground wire.
- ▶ Deactivate the accumulator.
- ▶ Remove accessories.

To prevent damage during storage:

- ▶ Pack or cover the bed and accessories.
- ▶ Ensure that storage conditions are the same as the operating conditions.

11 Manipulation



WARNING!

Risk of injury when adjusting the bed position and siderail position!

- ▶ Ensure that there are no body parts between the Mattress Support Platform elements and the Mattress Support Platform frame when adjusting the bed.
- ▶ Ensure that there are no body parts between the siderails and the Mattress Support Platform frame when adjusting the bed position and siderail position.
- ▶ Ensure that there are no body parts below the Mattress Support Platform frame before adjusting the bed.



CAUTION!

Risk of damaging the bed when adjusting the bed position and siderail position!

- ▶ Avoid any collision during bed positioning and siderail adjustment!



CAUTION!

Risk of damaging the bed when overloading the Calfrest!

- ▶ Respect the maximum allowed load of the calfrest!

11.1 Plastic Siderails

The siderails are components of the bed in contact with patient.

The nursing personnel are responsible for the siderails being raised up while the patient is in the bed.



WARNING!

Risk of injury, damaging or unintentional movement of the bed due to incorrect placement of accessories or Handset Control Panel!

- ▶ Never place any accessories or Handset Control Panel on the siderails in the area where keyboards are located.



WARNING!

Risk of injury due to incorrectly latched siderail!

- ▶ Ensure that siderail is secured in the upper or lower position. Audible „click“ indicates the siderail locked in the upper position. Push the siderail forward and backward to ensure the siderail is fixed in the upper position!



WARNING!

Risk of injury due to incorrect position of siderails!

- ▶ Ensure that siderails are folded up while the patient is in bed. Prevent patient from falling off the bed when a siderail is folded down! Hospital personnel is responsible for meeting this requirement.
- ▶ Ensure the CPR lever is accessible when necessary!

11.1.1 Description of the Plastic Siderails

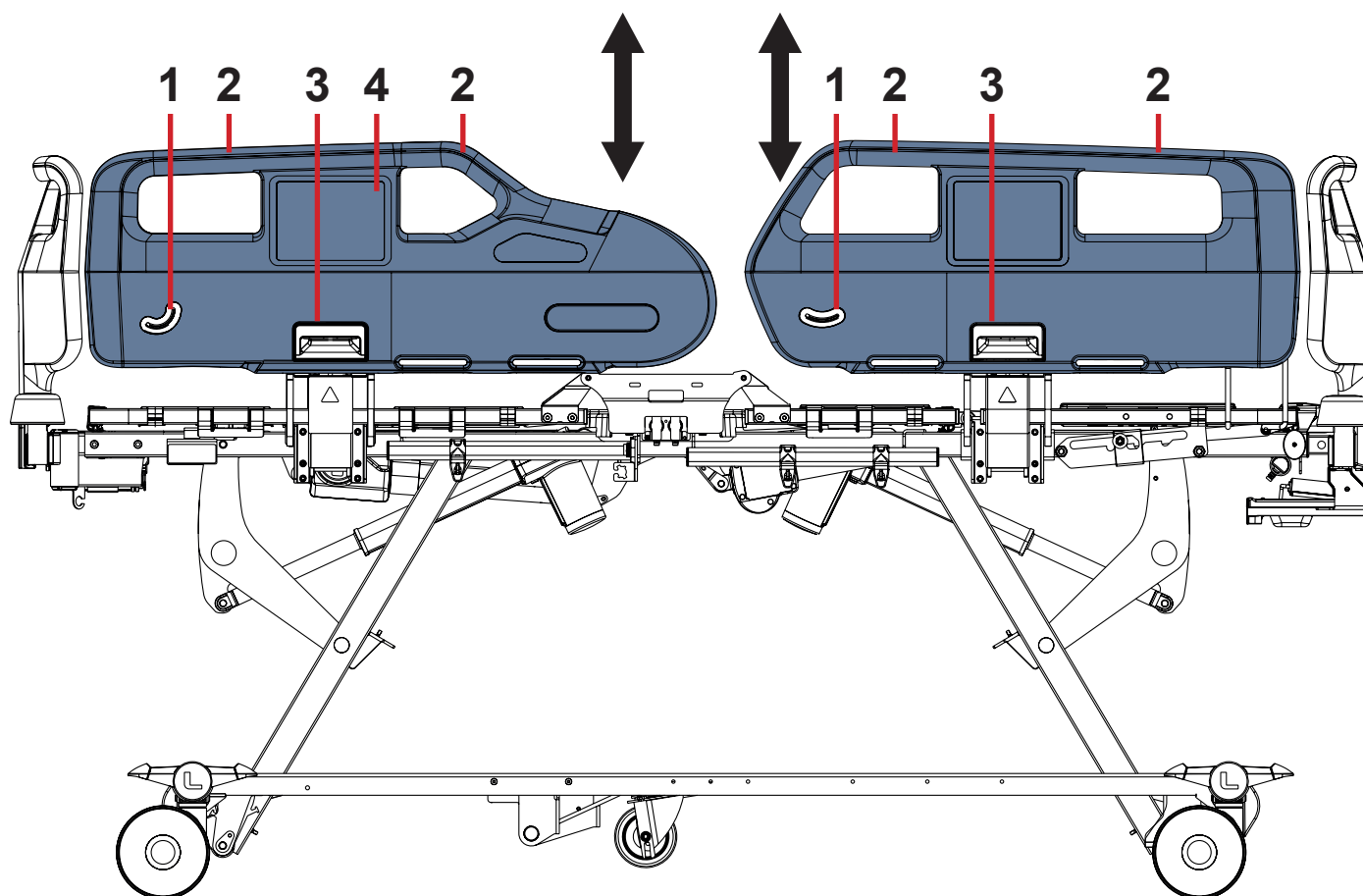


Fig. Manipulation with Plastic Siderails

- 1. Angle Indicator
- 2. Siderail Handle
- 3. Siderail Release Mechanism with Siderail Release Handle
- 4. Siderail Control Panel

MANIPULATION

To raise siderail up:

- ▶ Grab siderail by Siderail Handle (2).
- ▶ Pull siderail up until it latches. You will hear audible „click“.
- ▶ Push the siderail forward and backward to ensure the siderail is fixed in the upper position!

To release siderail down:

- ▶ Grab siderail by Siderail Handle (2).
- ▶ Unlock siderail by pulling Siderail Release Handle (3) to yourself.
- ▶ Fold down siderail slowly.

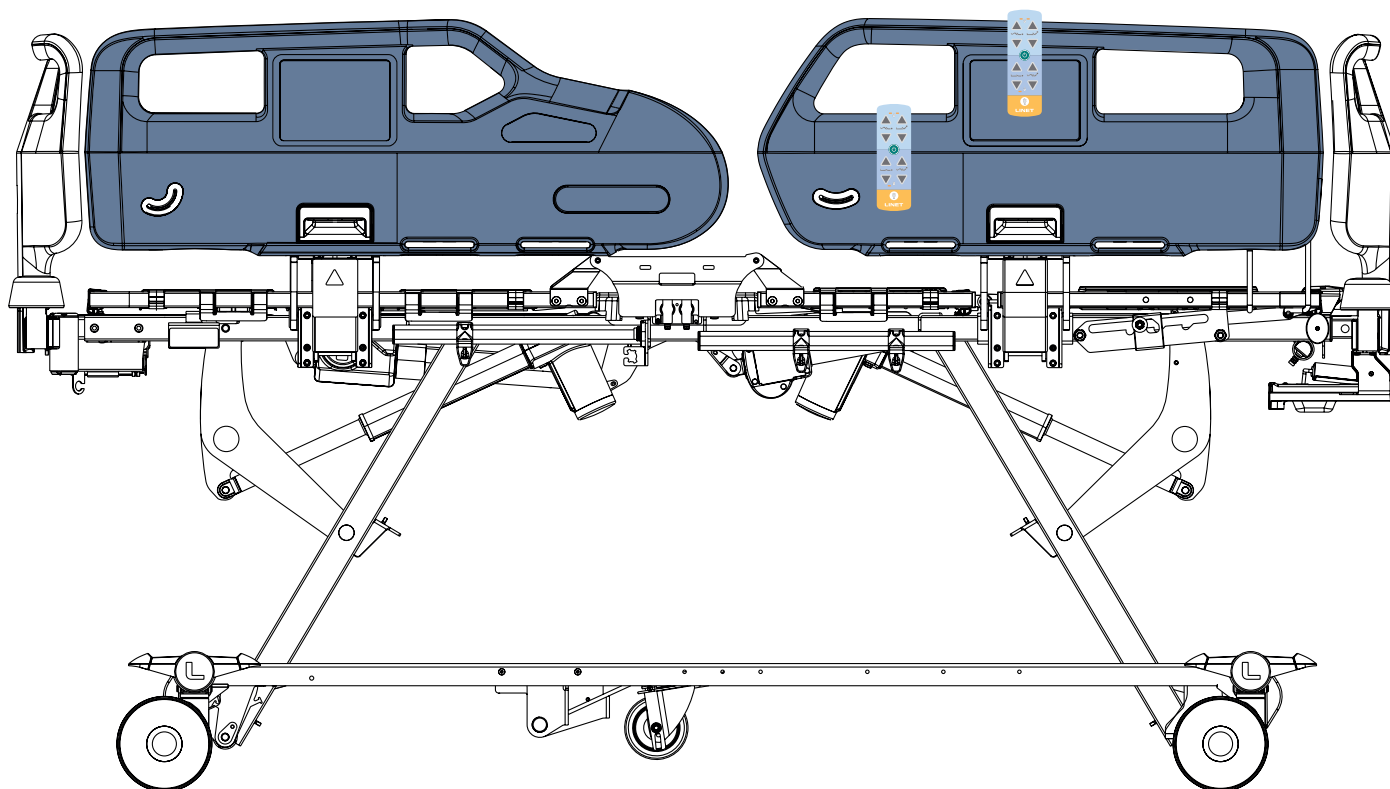


Fig. Correct positions of the Handset Control Panel on the Plastic Siderails

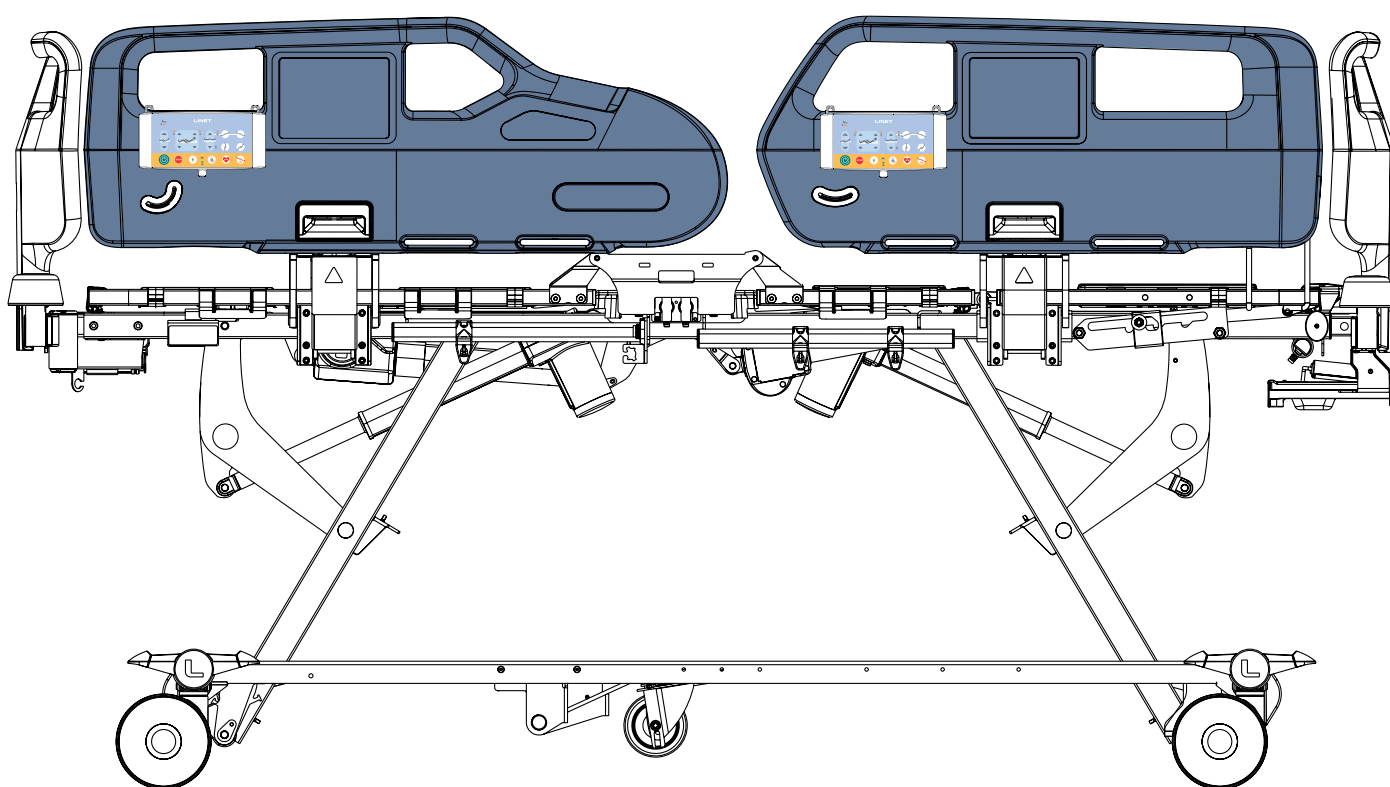


Fig. Correct positions of the Attendant Control Panel on the Plastic Siderails

11.2 Telescopic Siderails

The siderails are components of the bed in contact with patient.

The nursing personnel are responsible for the siderails being raised up while the patient is in the bed.



WARNING!

Risk of injury due to incorrectly latched siderail!

- Ensure that siderail is secured in the upper position or folded down to the lowest position. Audible „click“ indicates the siderail locked in the upper position. Push the siderail up and down to ensure the siderail is fixed in the required position!



WARNING!

Risk of injury due to incorrect position of siderails!

- Ensure that siderails are folded up while the patient is in the bed. Prevent patient from falling off the bed when a siderail is folded down! Hospital personnel is responsible for meeting this requirement.
- Hospital staff is responsible for assessing if the physical and psychological state of the patient is in accordance with use of the telescopic siderails in middle positions!



WARNING!

Risk of injury due to incorrect manipulation with siderails!

- Ensure that there is no body part between siderail bars during folding siderail down.

11.2.1 Description of the Telescopic Siderails without DIN Rails

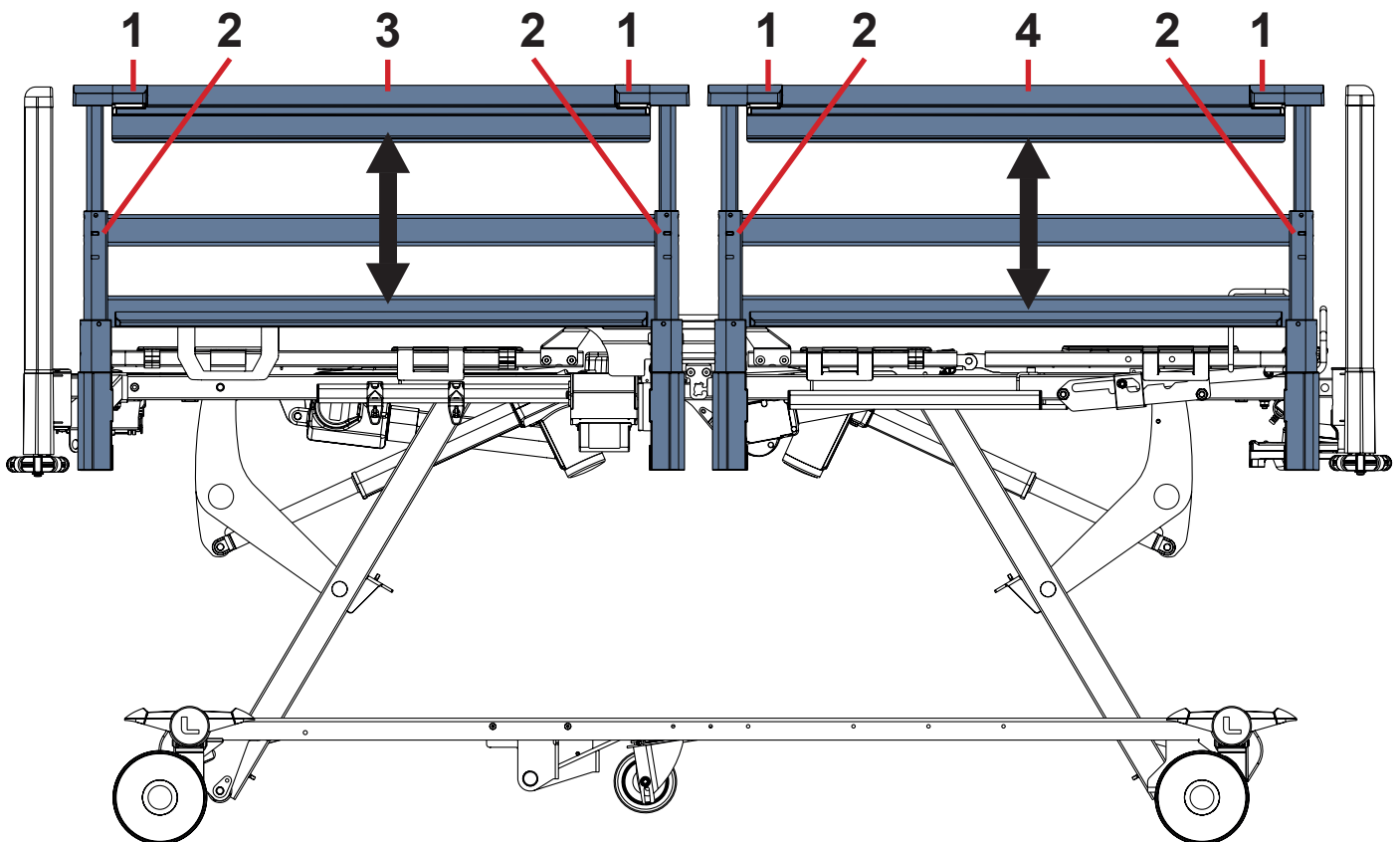


Fig. Manipulation with Telescopic Siderails without DIN Rails

1. Siderail Release Mechanism (with red button)
2. Check Latch (with yellow catch)
3. Head Telescopic Siderail
4. Foot Telescopic Siderail

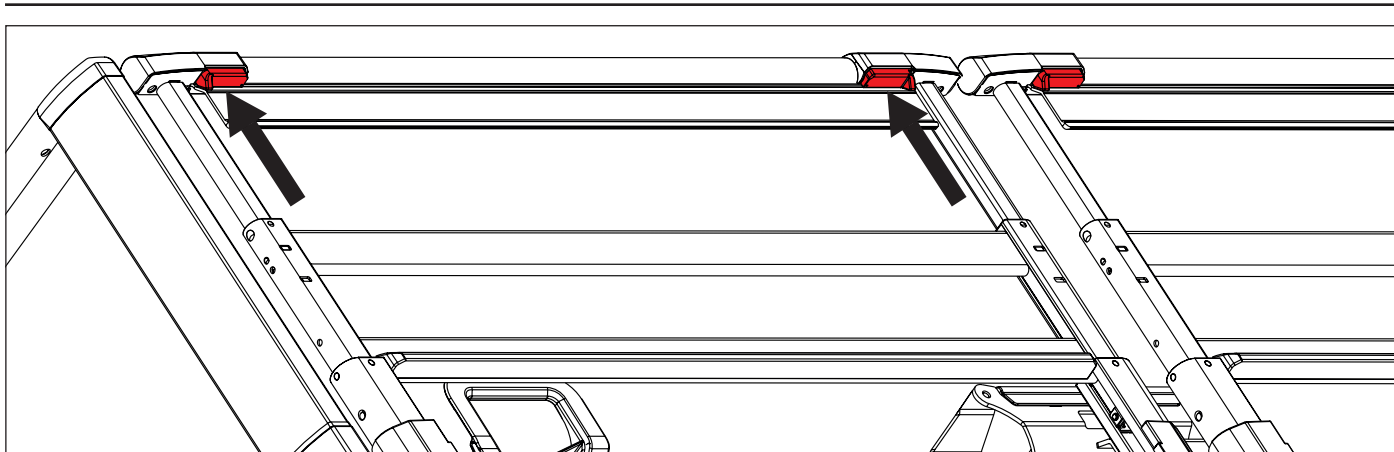


Fig. Red Buttons of Head Siderail Release Mechanism

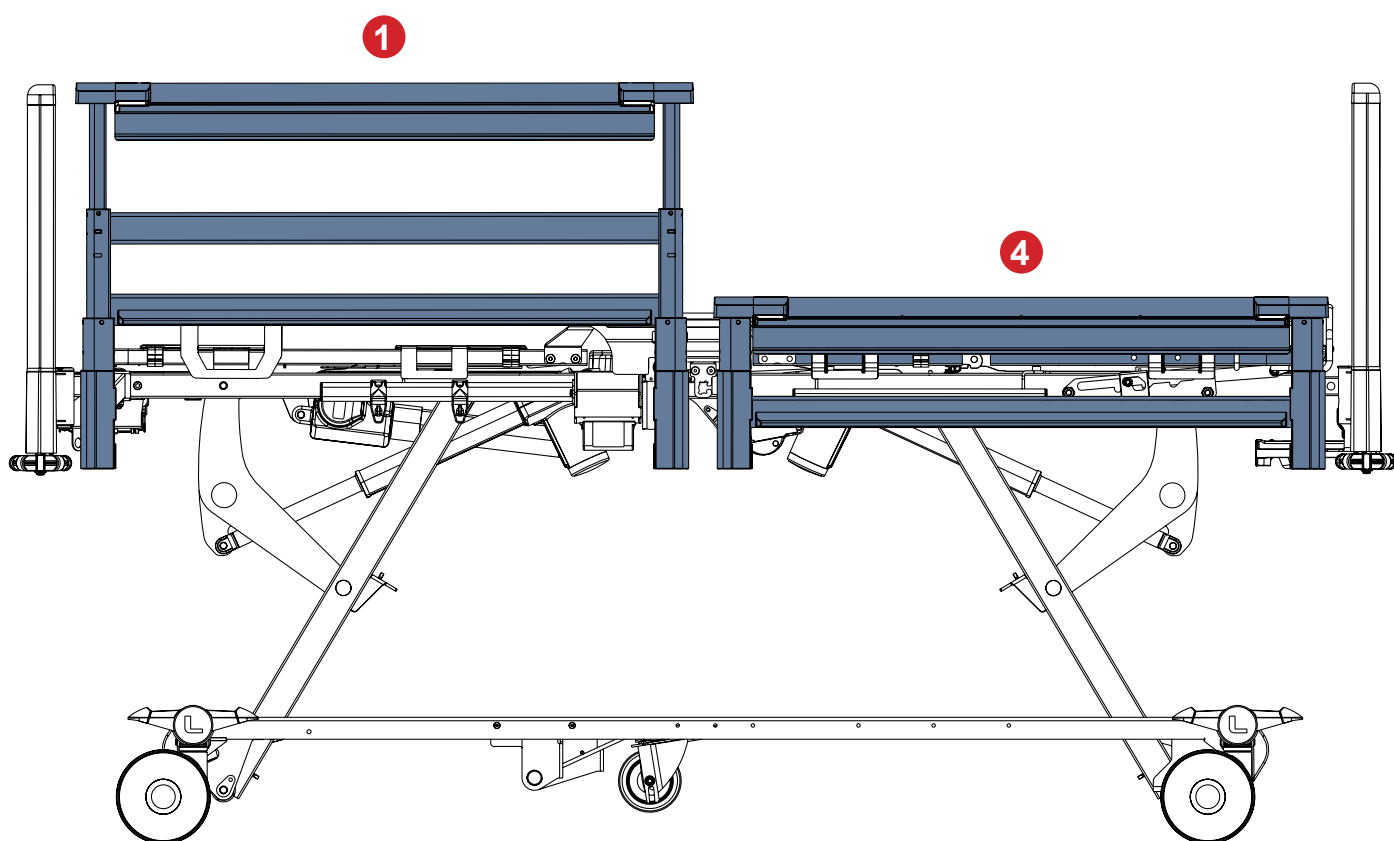


Fig. Telescopic Siderail in highest (1) and lowest (4) positions

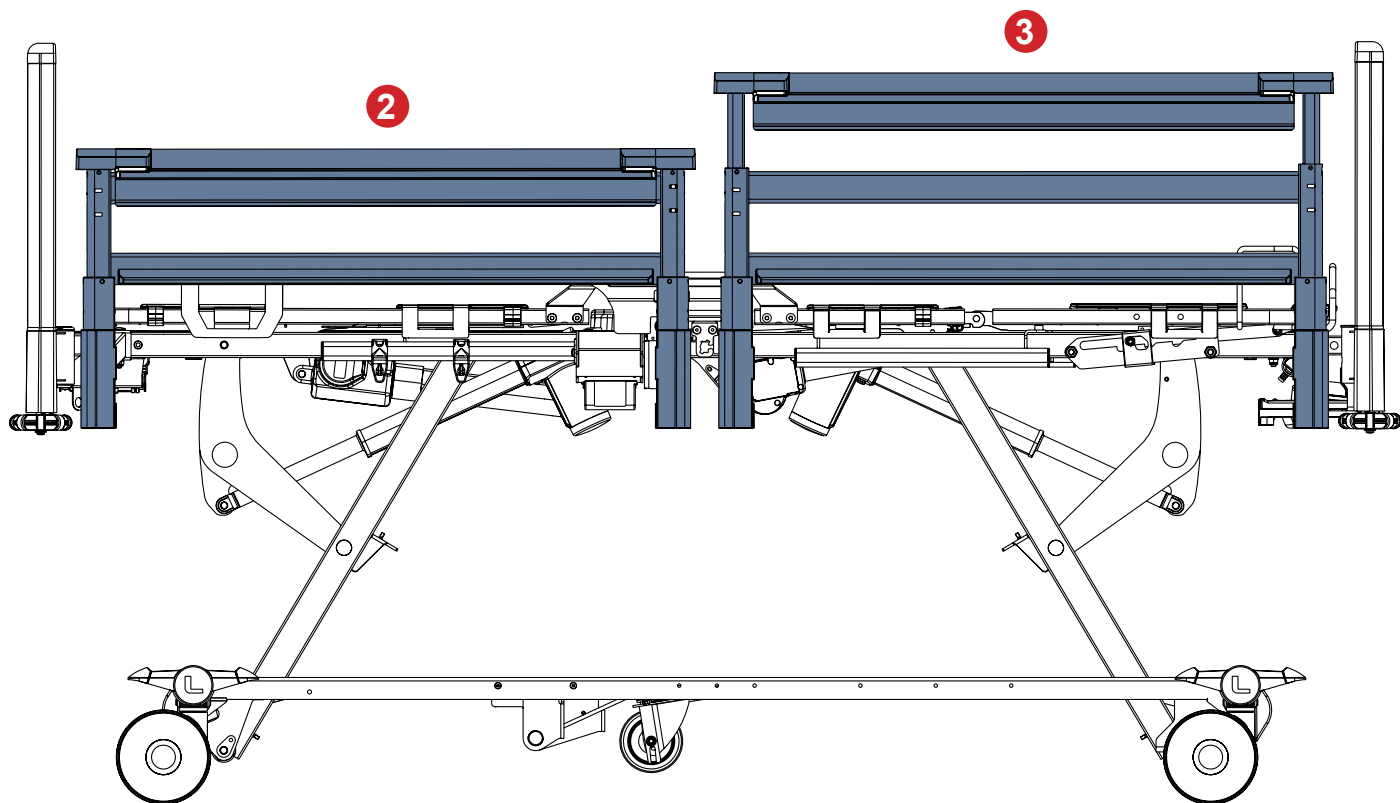


Fig. Telescopic Siderail in two interpositions (2, 3)

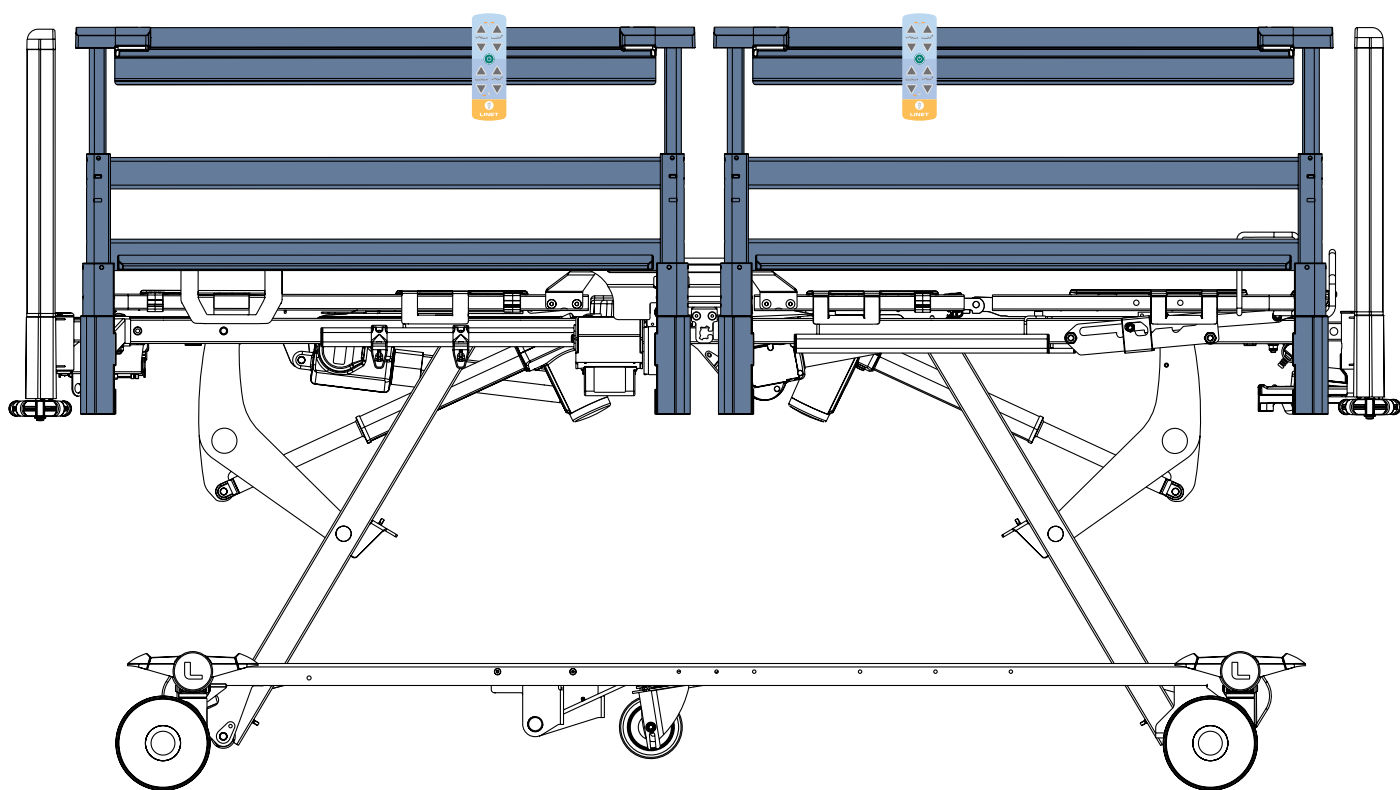


Fig. Correct positions of the Handset Control Panel on the Telescopic Siderails

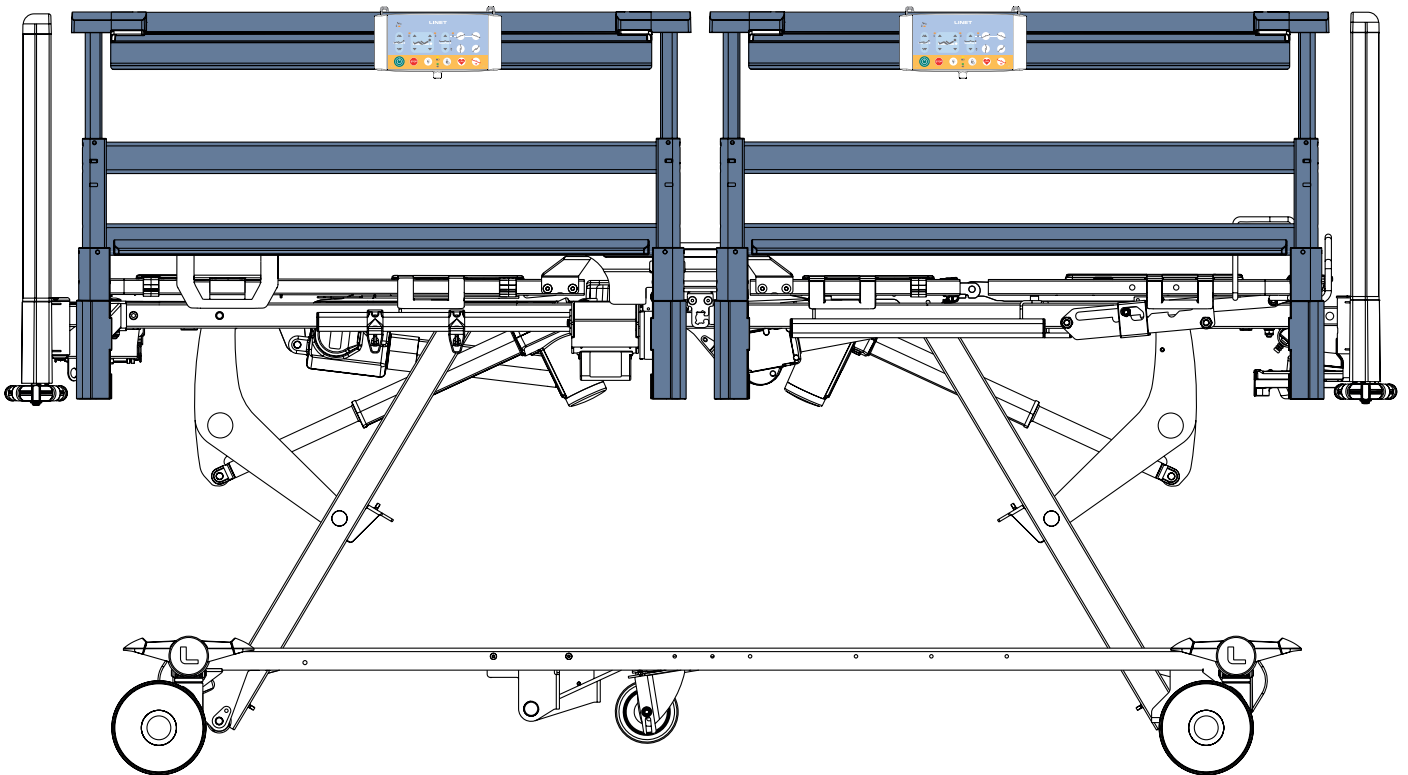


Fig. Correct positions of the Attendant Control Panel on the Telescopic Siderails

11.3 Collapsible Siderails

The siderails are components of the bed in contact with patient.

The nursing personnel are responsible for the siderails being raised up while the patient is in the bed.



WARNING!

Risk of injury due to incorrectly latched siderail!

- Ensure that siderail is secured in the upper position or folded down to the lowest position. Audible „click“ indicates the siderail locked in the upper position. Push the siderail towards head end and towards foot end to ensure the siderail is fixed in the upper position!



WARNING!

Risk of injury due to incorrect position of siderails!

- Ensure that siderails are folded up while the patient is in the bed. Prevent patient from falling off the bed when a siderail is folded down! Hospital personnel is responsible for meeting this requirement.



WARNING!

Risk of injury due to incorrect manipulation with siderails!

- Ensure that there is no body part between siderail bars during folding the siderail down.

11.3.1 Description of the Collapsible Siderails

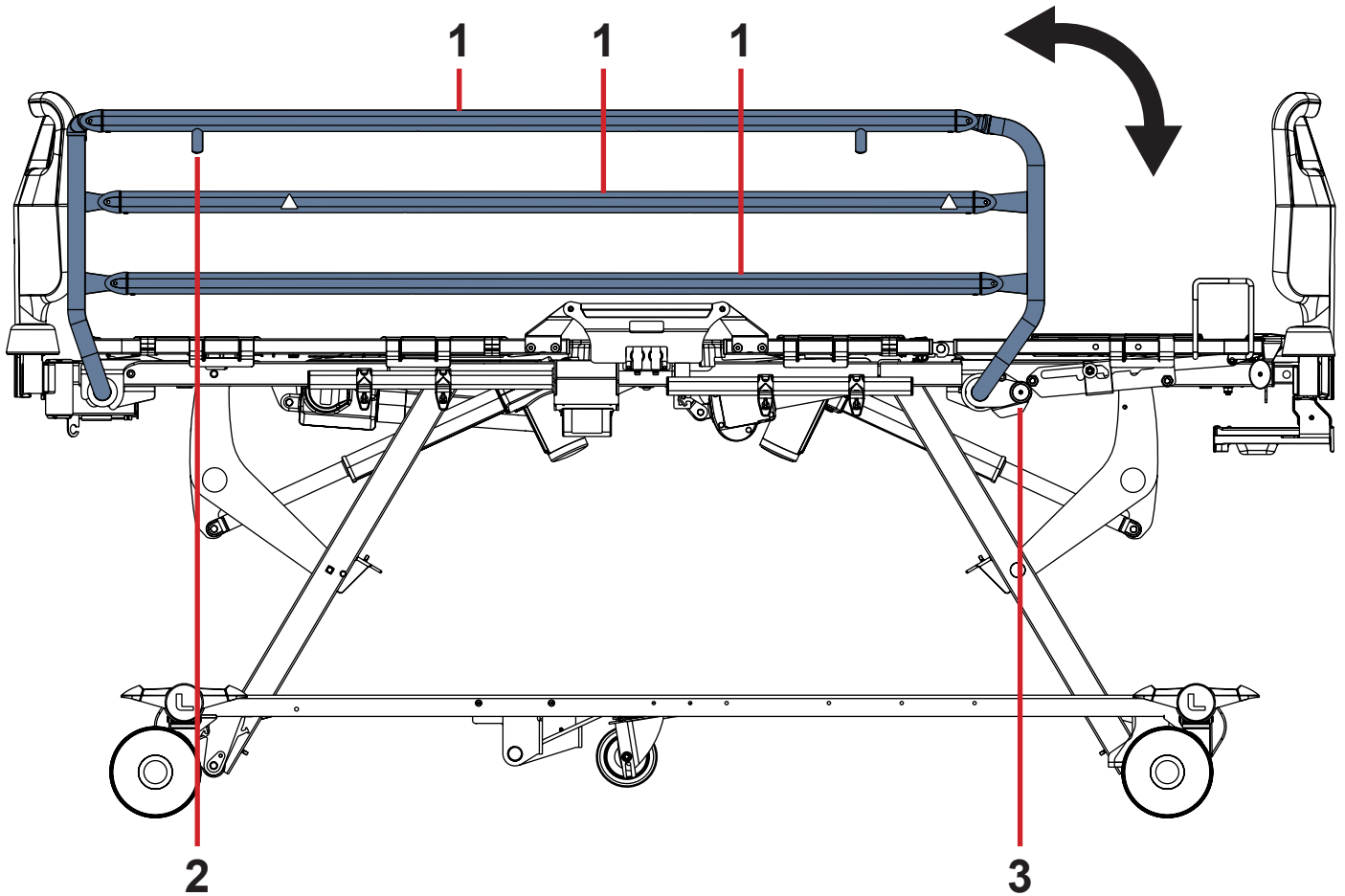


Fig. Manipulation with Collapsible Siderail

- 1. Siderail Bar
- 2. Siderail Stop
- 3. Siderail Release Mechanism

MANIPULATION

To raise siderail up:

- Grab an upper Siderail Bar and push the siderail up. Siderail will click into place and lock automatically. You will hear audible „click“.
- Push the siderail towards head end and towards foot end to ensure the siderail is fixed in place.

To fold siderail down:

- Unlock siderail by pulling out the lever of Siderail Release Mechanism (3).
- Fold down siderail slowly.

Do not let the siderail fall down when it is unblocked!

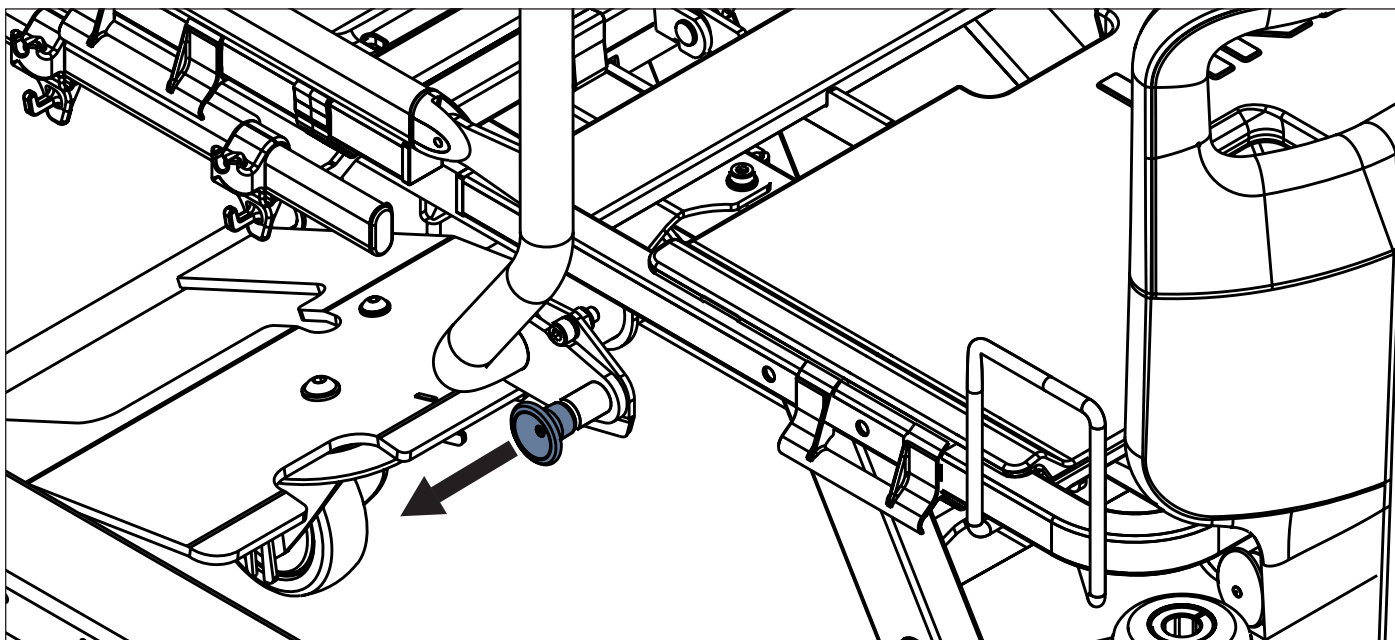


Fig. Manipulation with Siderail Release Mechanism in version with collapsible siderails

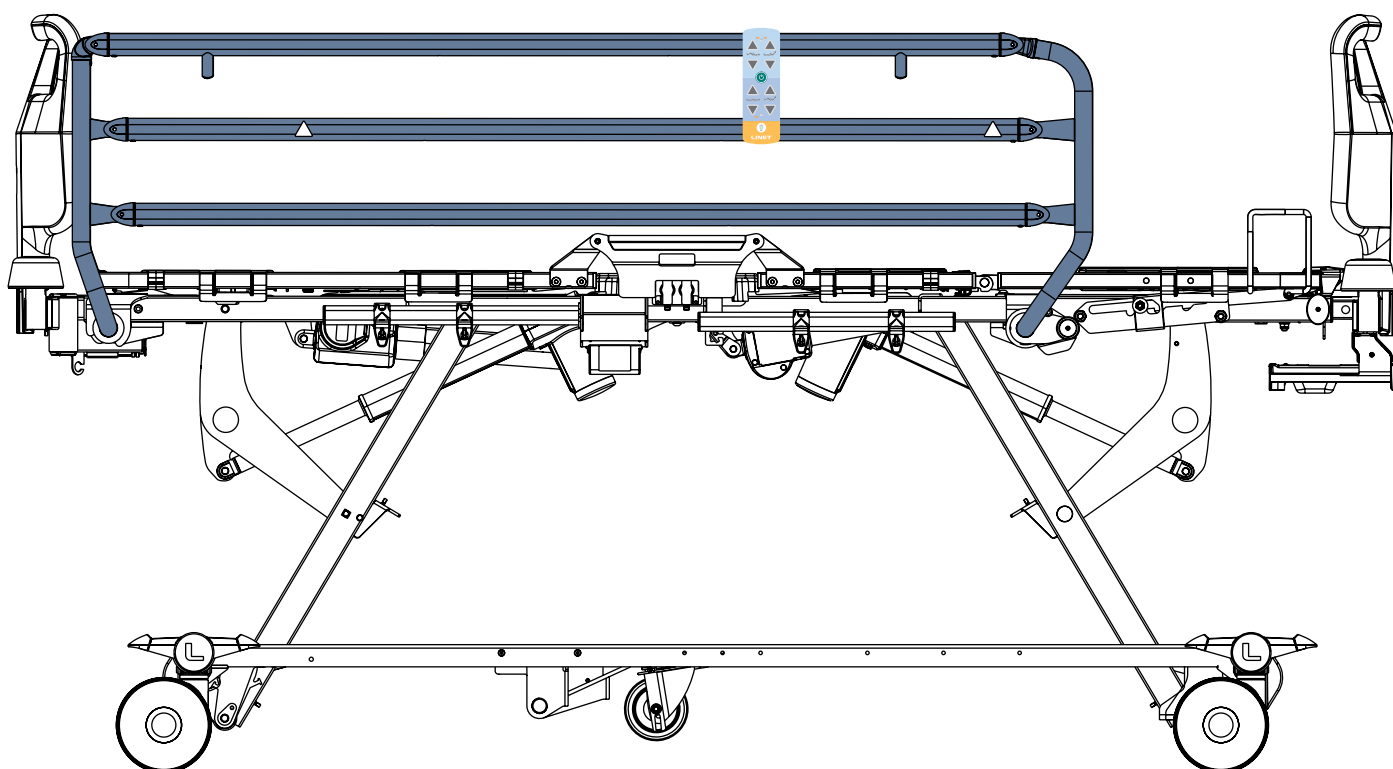


Fig. Correct position of the Handset Control Panel on the Collapsible Siderail

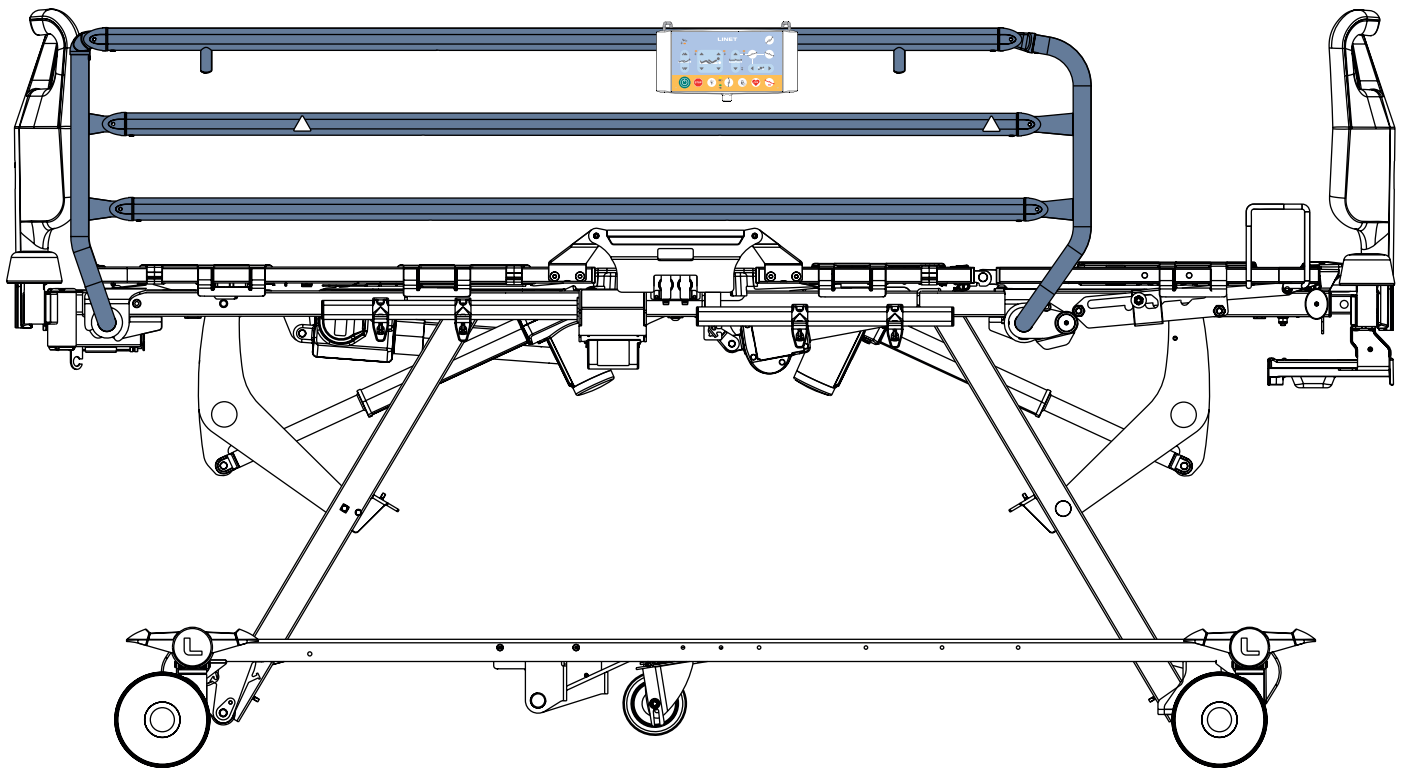


Fig. Correct position of the Attendant Control Panel on the Collapsible Siderail

11.4 Castor Control



CAUTION!

Risk of material damage due to incorrect transport and involuntary movement!

- ▶ Prior to transport, ensure that the bed is disconnected from the mains.
- ▶ Hang the power supply cord on the transport hook on the Head Board before bed transport.
- ▶ Ensure the castors are braked prior to putting into service, removal from service and maintenance.
- ▶ Ensure the castors are braked when the bed is occupied.
- ▶ Ensure the castors are braked when the bed should not move.
- ▶ Ensure the castors are braked before bed positioning.
- ▶ Have the bed transported exclusively by nursing personnel.
- ▶ Observe the path for any obstacles and avoid collisions and possible damages to any bed parts on the undercarriage.

The bed is equipped with central castor control and brake system.

A directional castor can be situated at bed head end or bed foot end depending on bed configuration.

The castor control levers are located in the four corners of the undercarriage.

Castor Configurations of the Essenza 300 W bed:

- ☐ 150 mm Tente Linea double castors + Fifth castor
- ☐ 150 mm Tente Linea double castors
- ☐ 125 mm Tente Linea double castors + Fifth castor
- ☐ 125 mm Tente Linea double castors

Each castor control lever is equipped with green pedal and red pedal.

Red colour refers to braking and green colour refers to steering.

Each castor control lever has 3 control positions.

All castor control levers are interconnected such that all pedal functions belong to each pedal.

In following table the pedal functions are described.

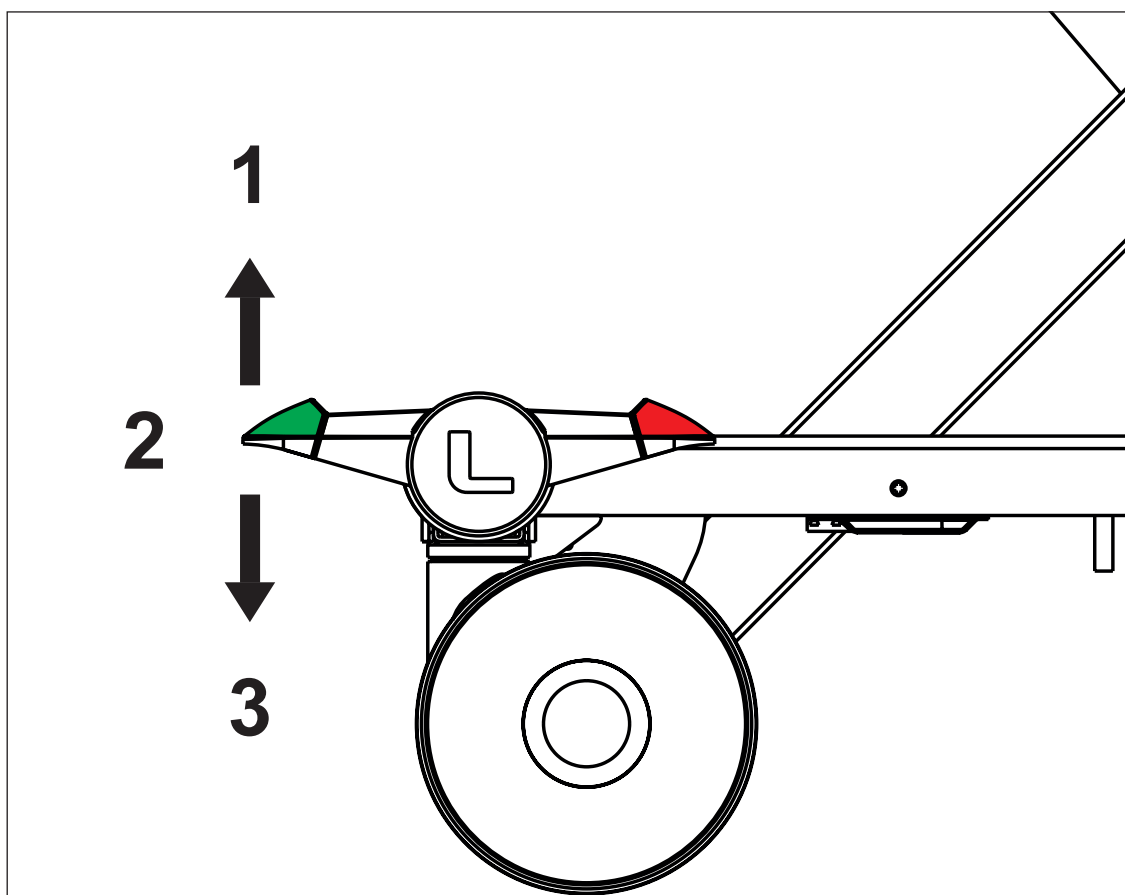


Fig. Three Positions of Castor Control Lever

Pedal Colour	Upper Position (1)	Middle Position (2)	Lower Position (3)
GREEN	BRAKED	UNRESTRICTED MOVEMENT	STEERING / FIFTH CASTOR
RED	STEERING / FIFTH CASTOR	UNRESTRICTED MOVEMENT	BRAKED

12 Control Elements

The Essenza 300 W medical bed is operated by different control elements.

Control elements depending on the Essenza 300 W configuration:

- Attendant Control Panel for Essenza 300 W with Lateral Tilt
- Attendant Control Panel for Essenza 300 W
- Siderail Control Panel for Essenza 300 W with Lateral Tilt
- Siderail Control Panel for Essenza 300 W
- Handset Control Panel
- Patient Control Panel

Control Elements for Essenza 300 W without Lateral Tilt

	Attendant Control Panel for Essenza 300 W	Siderail Control Panel for Essenza 300 W	Handset Control Panel	Patient Control Panel
Backrest	✓	✓	✓	✓
Thighrest	✓	✓	✓	✓
Calfrest				
Autocontour	✓	✓	✓	
Bed Height	✓	✓	✓	✓
Emergency Trendelenburg Position	✓	✓		
Trendelenburg Tilt	✓	✓		
Reverse-Trendelenburg Tilt	✓	✓		
Bed Extension				
Emergency CPR Position	✓	✓		
Cardiac Chair Position	✓	✓		
Mobilization Position	✓	✓		
Extra Low Position	✓	✓		

Control Elements for Essenza 300 W with Lateral Tilt

	Attendant Control Panel for Essenza 300 W with Lateral Tilt	Siderail Control Panel for Essenza 300 W with Lateral Tilt	Handset Control Panel	Patient Control Panel
Backrest	✓	✓	✓	✓
Thighrest	✓	✓	✓	✓
Calfrest				
Autocontour	✓	✓	✓	
Bed Height	✓	✓	✓	✓
Emergency Trendelenburg Position	✓	✓		
Trendelenburg Tilt	✓	✓		
Reverse-Trendelenburg Tilt	✓	✓		
Bed Extension				
Emergency CPR Position	✓	✓		
Cardiac Chair Position	✓	✓		
Mobilization Position	✓	✓		
Lateral Tilt	✓	✓		
Extra Low Position	✓	✓	✓	✓

12.1 Attendant Control Panel for Essenza 300 W with Lateral Tilt

Ensure that exclusively trained nursing staff operates the Attendant Control Panel for Essenza 300 W with Lateral Tilt. The use of the Attendant Control Panel for Essenza 300 W with Lateral Tilt is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

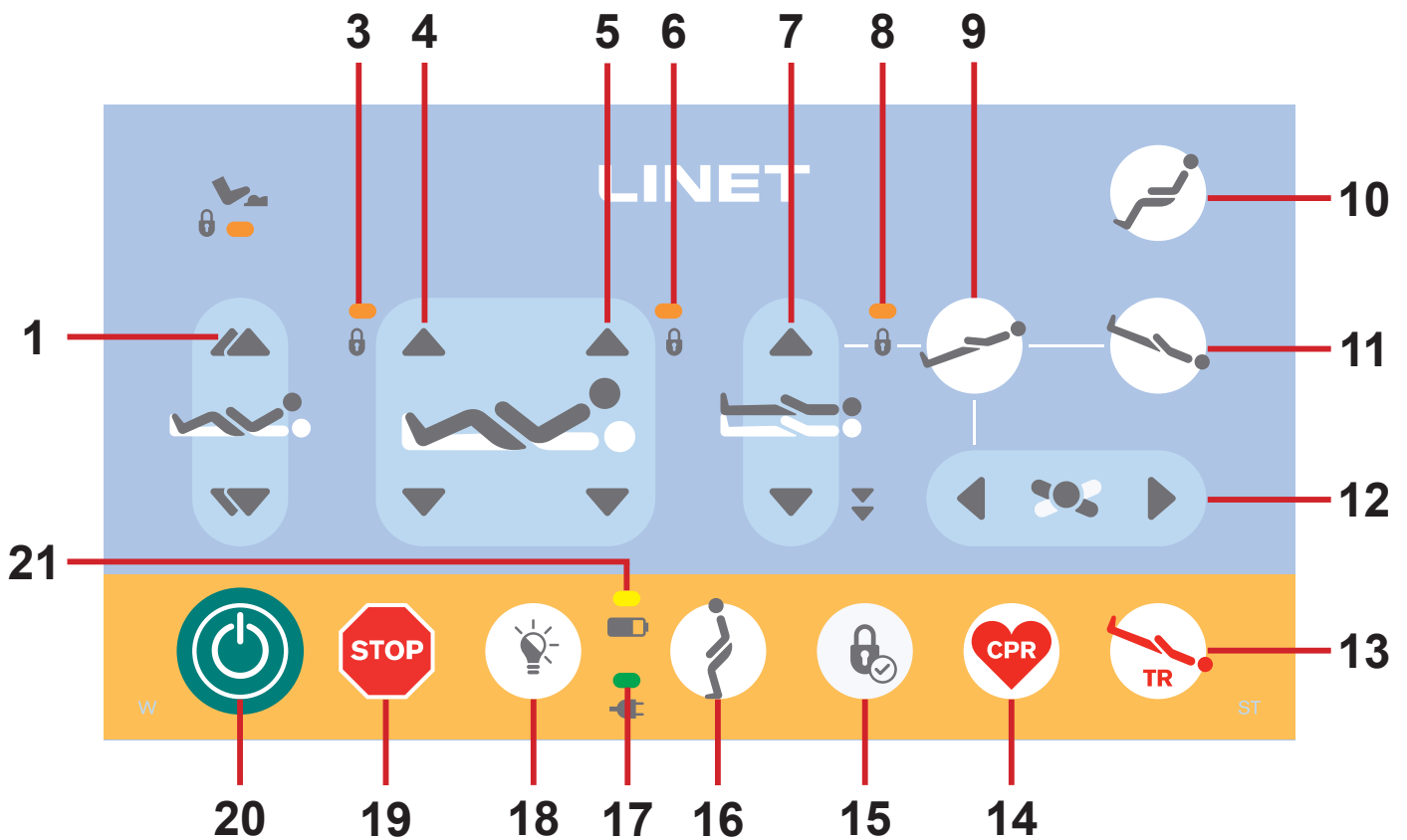


Fig. Attendant Control Panel for Essenza 300 W with Lateral Tilt

1. Autocontour Adjustment Buttons (up/down)
3. Thighrest Adjustment Lock LED Indicator (ON - locked)
4. Thighrest Adjustment Buttons (up/down)
5. Backrest Adjustment Buttons (up/down)
6. Backrest Adjustment Lock LED Indicator (ON - locked)
7. Bed Height Adjustment Buttons (up/down + Extra Low Position)
8. Bed Height Adjustment Lock LED Indicator / Lateral Tilt Lock LED Indicator / Trendelenburg Tilt Lock LED Indicator / Reverse-Trendelenburg Tilt Lock LED Indicator (ON - locked)
9. Reverse-Trendelenburg Tilt Button
10. Cardiac Chair Position Button
11. Trendelenburg Tilt Button
12. Lateral Tilt Buttons (right/left)
13. Emergency Trendelenburg Position Button
14. Emergency CPR Position Button
15. Lock Button
16. Mobilization Position Button
17. Mains Power LED Indicator
18. Bed Illumination Control Button
19. Central STOP Button
20. GO Button
21. Accumulator Charge Status LED Indicator

12.2 Attendant Control Panel for Essenza 300 W

Ensure that exclusively trained nursing staff operates the Attendant Control Panel for Essenza 300 W.

The use of the Attendant Control Panel for Essenza 300 W is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

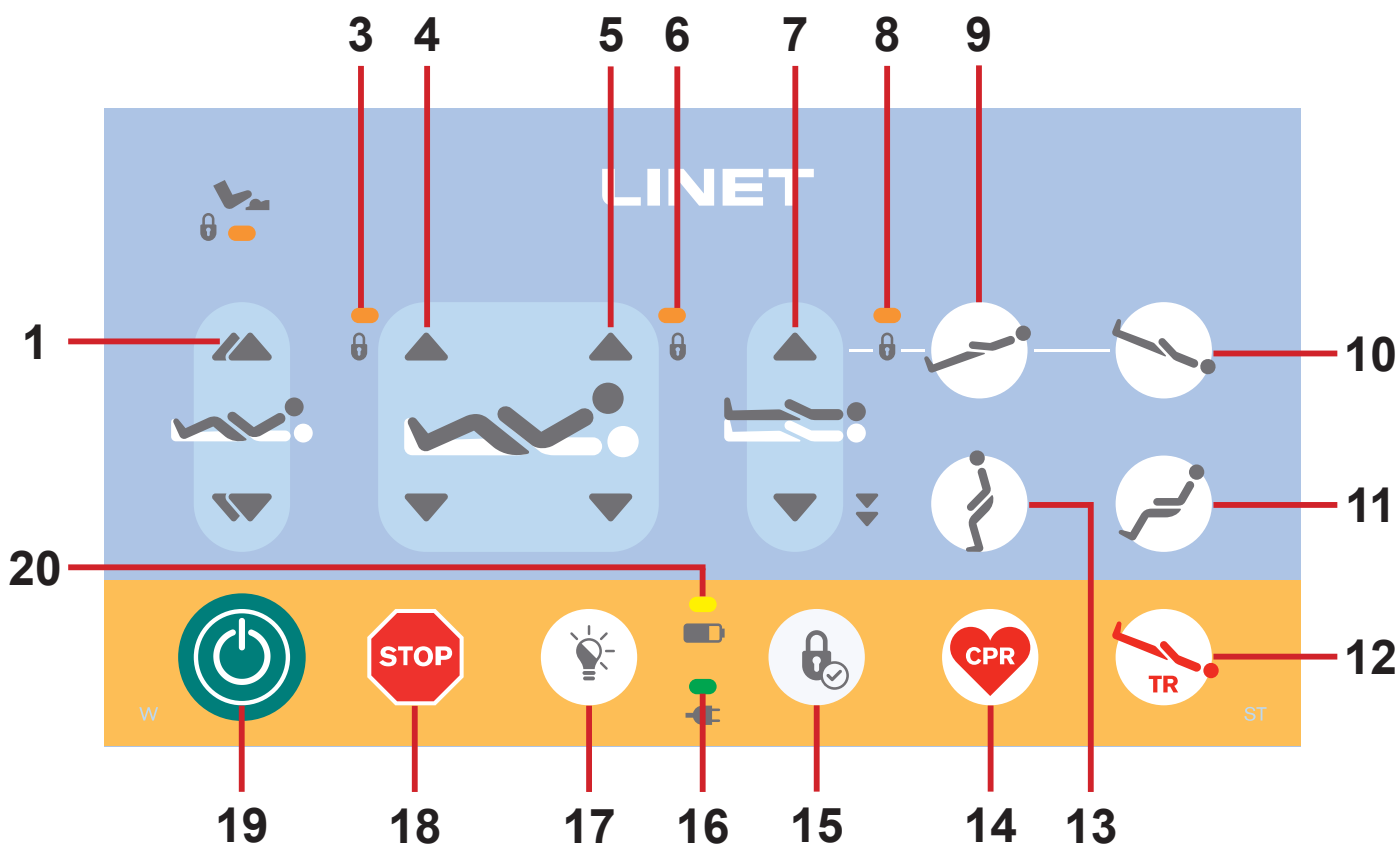


Fig. Attendant Control Panel for Essenza 300 W

1. Autocontour Adjustment Buttons (up/down)
3. Thighrest Adjustment Lock LED Indicator (ON - locked)
4. Thighrest Adjustment Buttons (up/down)
5. Backrest Adjustment Buttons (up/down)
6. Backrest Adjustment Lock LED Indicator (ON - locked)
7. Bed Height Adjustment Buttons (up/down + Extra Low Position)
8. Bed Height Adjustment Lock LED Indicator / Trendelenburg Tilt Lock LED Indicator / Reverse-Trendelenburg Tilt Lock LED Indicator (ON - locked)
9. Reverse-Trendelenburg Tilt Button
10. Trendelenburg Tilt Button
11. Cardiac Chair Position Button
12. Emergency Trendelenburg Position Button
13. Mobilization Position Button
14. Emergency CPR Position Button
15. Lock Button
16. Mains Power LED Indicator
17. Bed Illumination Control Button
18. Central STOP Button
19. GO Button
20. Accumulator Charge Status LED Indicator

12.3 Siderail Control Panel for Essenza 300 W with Lateral Tilt

Ensure that exclusively trained nursing staff operates the Siderail Control Panel for Essenza 300 W with Lateral Tilt. The use of the Siderail Control Panel for Essenza 300 W with Lateral Tilt is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

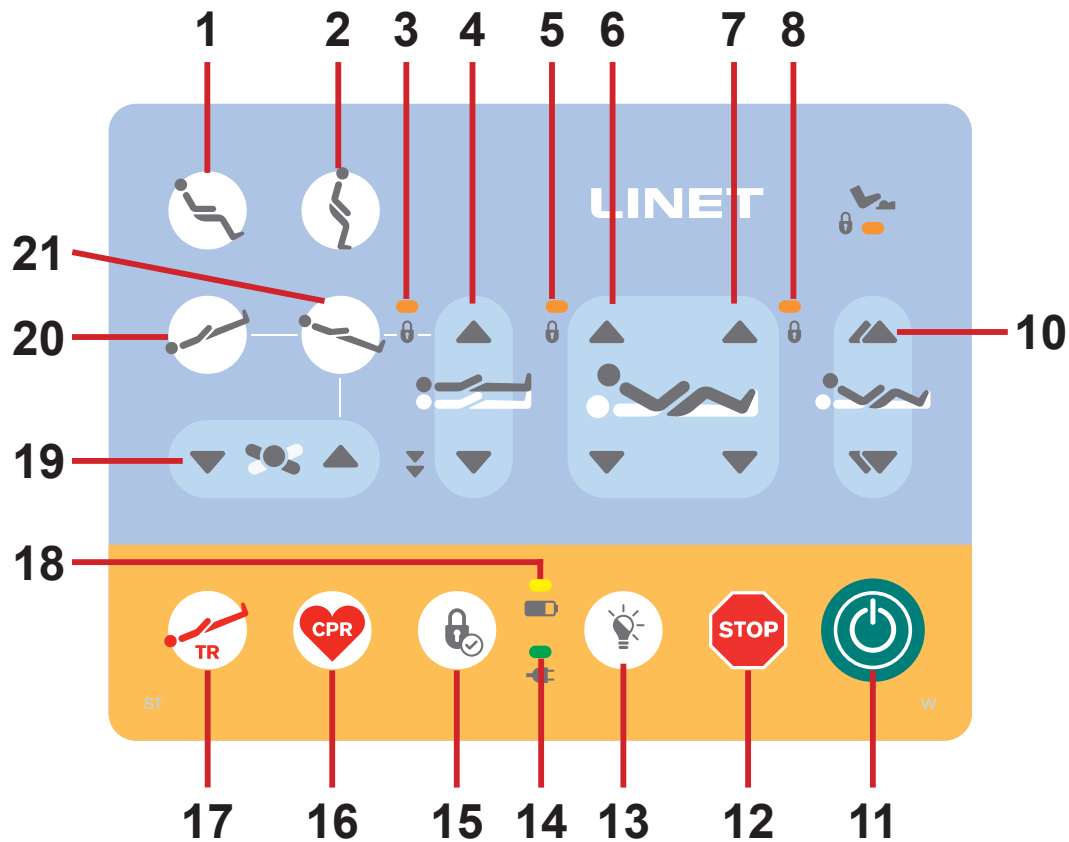


Fig. Siderail Control Panel for Essenza 300 W with Lateral Tilt

1. Cardiac Chair Position Button
2. Mobilization Position Button
3. Bed Height Adjustment Lock LED Indicator / Lateral Tilt Lock LED Indicator / Trendelenburg Tilt Lock LED Indicator / Reverse-Trendelenburg Tilt Lock LED Indicator (ON - locked)
4. Bed Height Adjustment Buttons (up/down + Extra Low Position)
5. Backrest Adjustment Lock LED Indicator (ON - locked)
6. Backrest Adjustment Buttons (up/down)
7. Thighrest Adjustment Buttons (up/down)
8. Thighrest Adjustment Lock LED Indicator (ON - locked)
10. Autocontour Adjustment Buttons (up/down)
11. GO Button
12. Central STOP Button
13. Bed Illumination Control Button
14. Mains Power LED Indicator
15. Lock Button
16. Emergency CPR Position Button
17. Emergency Trendelenburg Position Button
18. Accumulator Charge Status LED Indicator
19. Lateral Tilt Buttons (right/left)
20. Trendelenburg Tilt Button
21. Reverse-Trendelenburg Tilt Button

12.4 Siderail Control Panel for Essenza 300 W

Ensure that exclusively trained nursing staff operates the Siderail Control Panel for Essenza 300 W.

The use of the Siderail Control Panel for Essenza 300 W is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

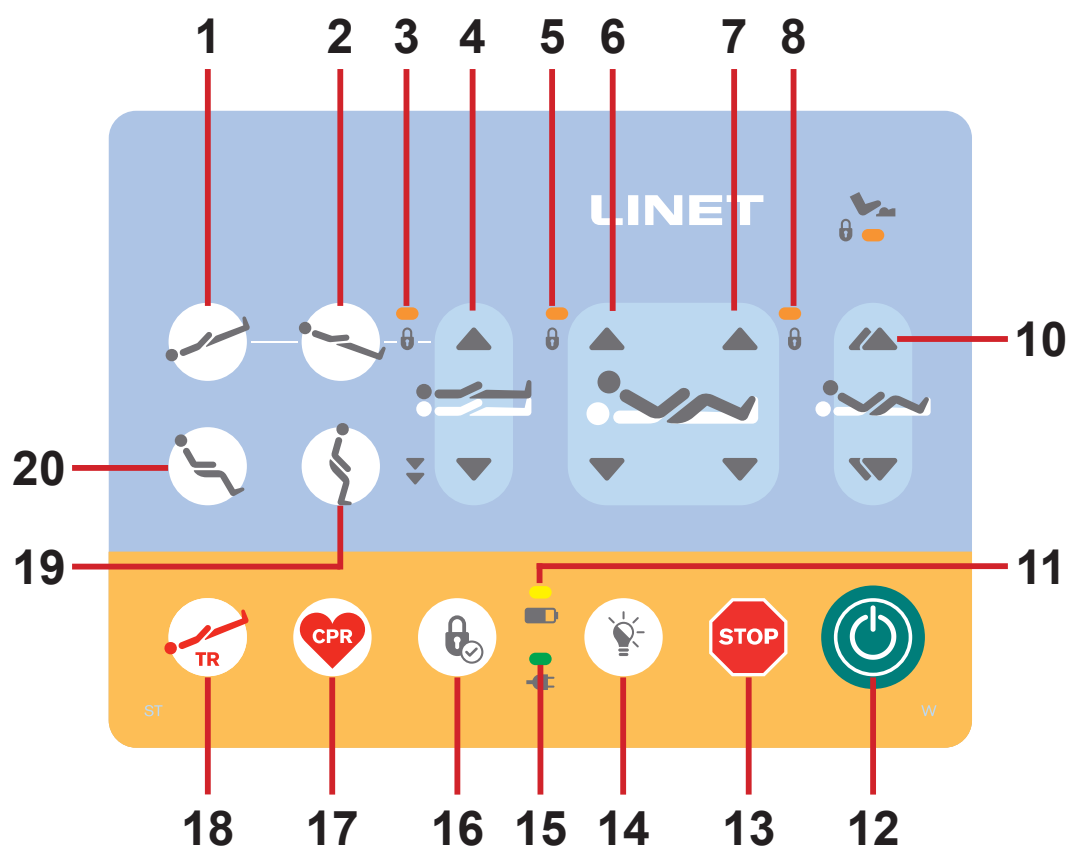


Fig. Siderail Control Panel for Essenza 300 W

1. Trendelenburg Tilt Button
2. Reverse-Trendelenburg Tilt Button
3. Bed Height Adjustment Lock LED Indicator / Lateral Tilt Lock LED Indicator / Trendelenburg Tilt Lock LED Indicator / Reverse-Trendelenburg Tilt Lock LED Indicator (ON - locked)
4. Bed Height Adjustment Buttons (up/down + Extra Low Position)
5. Backrest Adjustment Lock LED Indicator (ON - locked)
6. Backrest Adjustment Buttons (up/down)
7. Thighrest Adjustment Buttons (up/down)
8. Thighrest Adjustment Lock LED Indicator (ON - locked)
10. Autocontour Adjustment Buttons (up/down)
11. Accumulator Charge Status LED Indicator
12. GO Button
13. Central STOP Button
14. Bed Illumination Control Button
15. Mains Power LED Indicator
16. Lock Button
17. Emergency CPR Position Button
18. Emergency Trendelenburg Position Button
19. Mobilization Position Button
20. Cardiac Chair Position Button

12.5 Handset Control Panel

Patient is allowed to use the Handset Control Panel only if hospital personnel had assessed that the patient's physical and psychological state is in accordance with use of this control element and only if the hospital personnel had trained the patient in accordance with these instructions for use.

The use of the Handset Control Panel is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

Keyboard of the Handset Control Panel can optionally be illuminated.

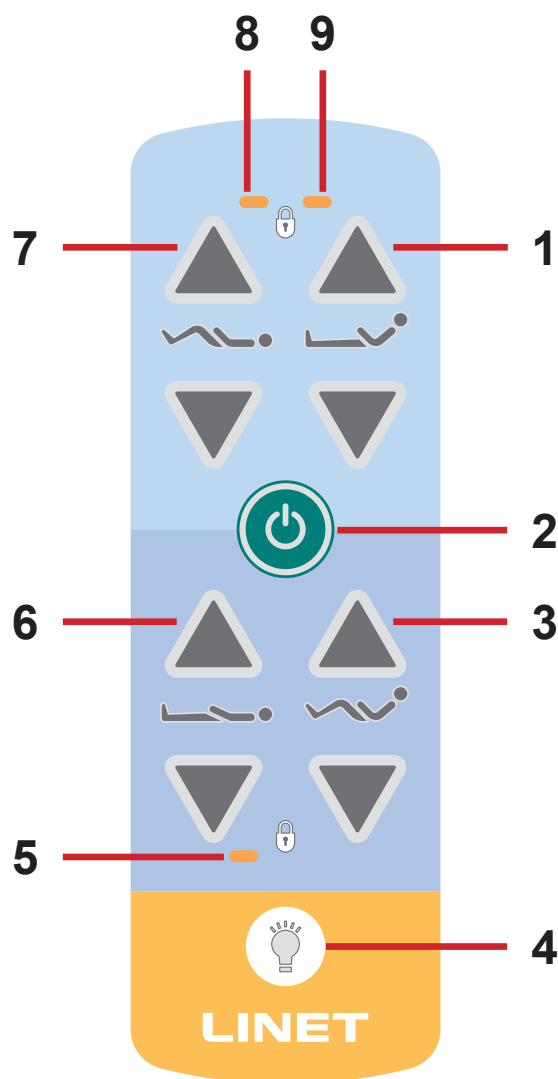


Fig. Handset Control Panel

1. Backrest Adjustment Buttons (up/down)
2. GO Button
3. Autocontour Adjustment Buttons (up/down)
4. Flashlight Button
5. Bed Height Adjustment Lock LED Indicator (ON - locked)
6. Bed Height Adjustment Buttons (up/down + Extra Low Position for Essenza 300 W with Lateral Tilt)
7. Thighrest Adjustment Buttons (up/down)
8. Thighrest Adjustment Lock LED Indicator (ON - locked)
9. Backrest Adjustment Lock LED Indicator (ON - locked)

To switch on the flashlight:

- ▶ Press and hold Flashlight button.

Flashlight emits light as long as the Flashlight button is pressed.

To disable bed positioning functions for patient, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W

12.5.1 Cable of the Handset Control Panel

Handset Control Panel Cable is equipped with a plastic lock.

The Plastic Lock is intended for fixing the Handset Control Panel on the selected bed side.

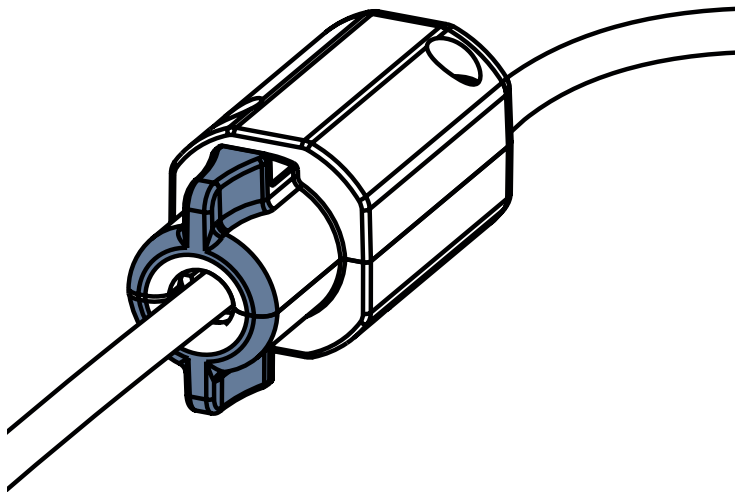


Fig. Plastic Lock of the Handset Control Panel Cable

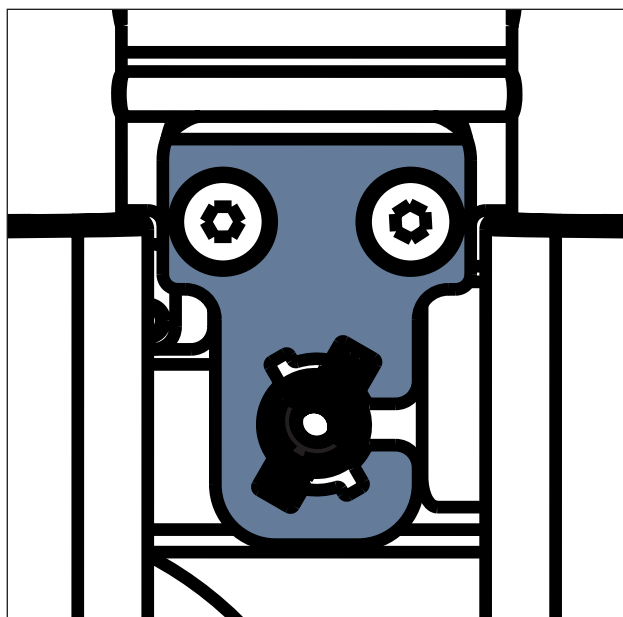


Fig. Handset Control Panel Cable Holder with unlocked plastic lock

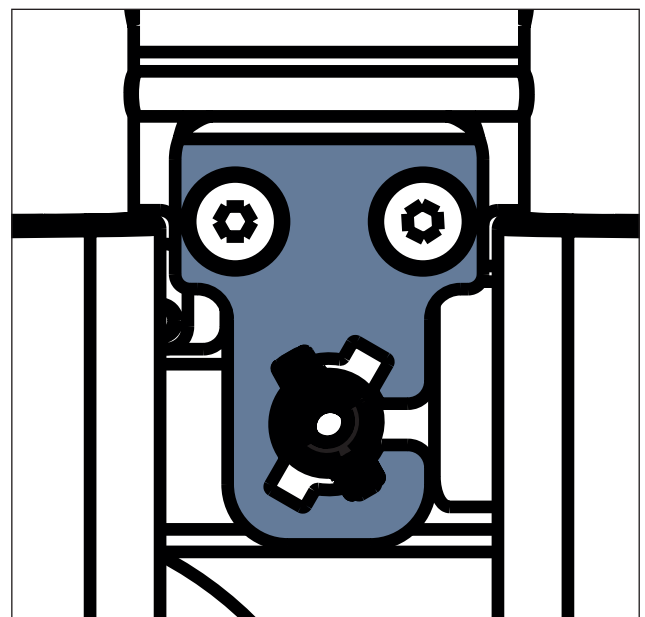


Fig. Handset Control Panel Cable Holder with locked plastic lock

12.5.2 Moving the Handset Control Panel to the opposite bed side

There are 2 types of the Handset Control Panel Cable Routing depending on the bed configuration:

- ▶ **Cable Routing through the bed Seat Section**
- ▶ **Cable Routing between the Seat Section and Thighrest**

Cable Routing through the bed Seat Section is related to the bed configurations without Bed Illumination, without Safe Height Visual Signalling, without Siderail Control Panels, without preparation for SafeSense 3 and without Backrest stop at 30 degrees.

Cable Routing between the Seat Section and Thighrest is related to the bed configurations with Bed Illumination or with Safe Height Visual Signalling or with Siderail Control Panels or with preparation for SafeSense 3 or with Backrest stop at 30 degrees.

Handset Control Panel Cable Routing through the bed Seat Section

There are 2 types of the Handset Control Panel Cable Holder position for:

- ▶ **bed with telescopic siderails or collapsible siderails**
- ▶ **bed with plastic siderails**

To fix the Handset Control Panel on the opposite side of the bed:

- 1) Remove the plastic cover on the seat section of the mattress support platform.
- 2) Unlock the Handset Control Panel Cable from the Handset Control Panel Cable Holder on the bed side.
Turn the plastic lock to the right to release the cable from the holder.
Take the cable out of the holder.
- 3) Run the Handset Control Panel upwards through the nearest opening in the seat section of the mattress support platform.
- 4) Take the Handset Control Panel Cable out of the cable clip.
- 5) Put the Handset Control Panel Cable into the cable clip on the opposite side.
- 6) Run the Handset Control Panel downwards through the farthest opening in the seat section of the mattress support platform on the opposite bed side.
- 7) Insert the Handset Control Panel Cable into the Handset Control Panel Cable Holder on the opposite bed side.
- 8) Insert the plastic lock into the holder and lock the cable in the holder by turning the plastic lock to the left.
- 9) Return the plastic cover back into the seat section of the mattress support platform.

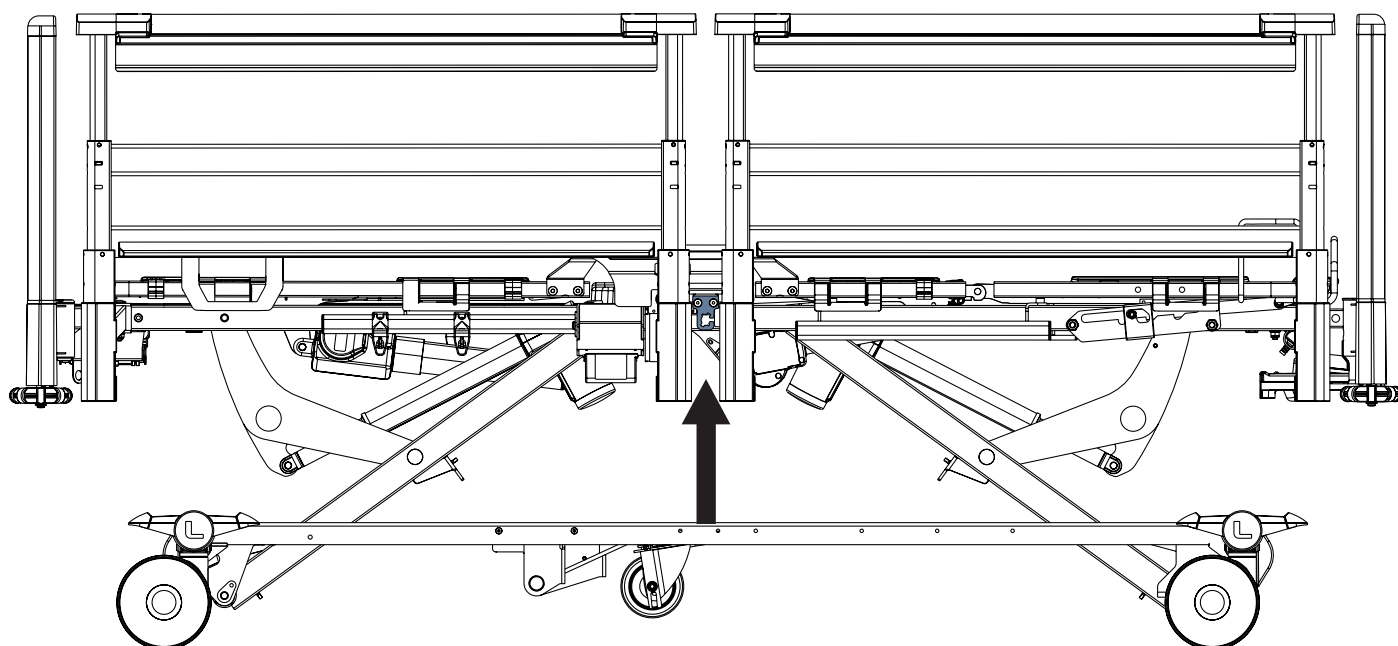


Fig. Position of the Handset Control Panel Cable Holder on the bed with telescopic siderails

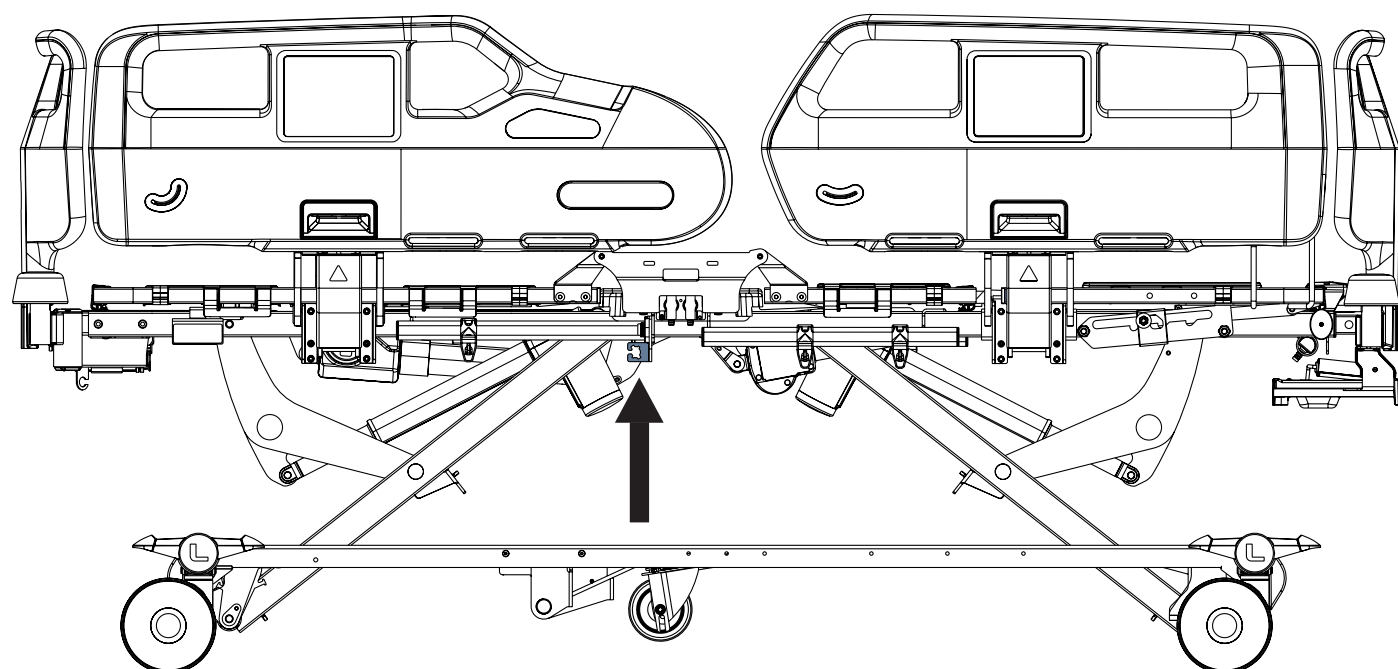


Fig. Position of the Handset Control Panel Cable Holder on the bed with plastic siderails

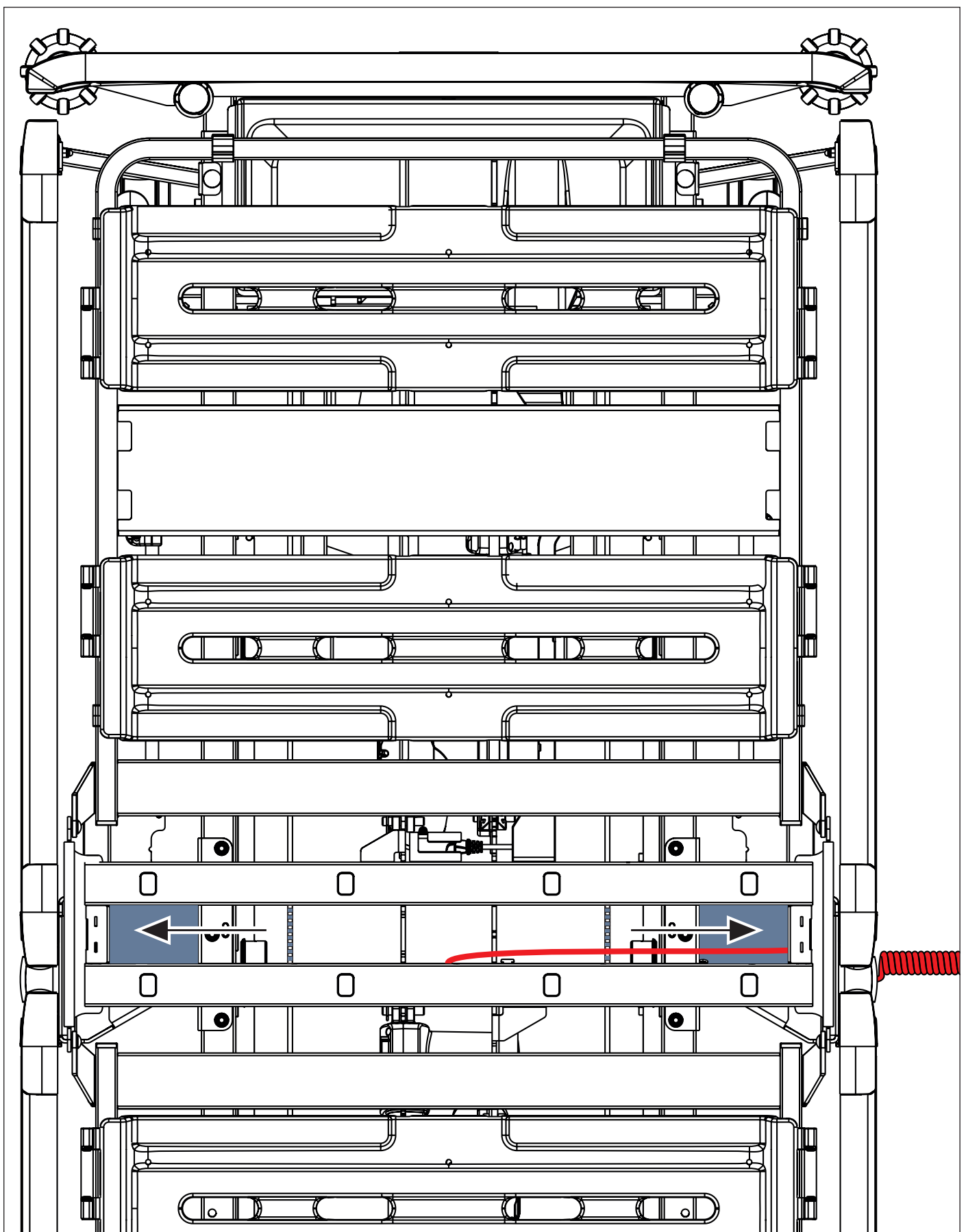


Fig. Moving the Handset Control Panel Cable to the opposite bed side through the bed seat section

Handset Control Panel Cable Routing between the Seat Section and Thighrest

There are 2 types of the Handset Control Panel Cable Holder position for:

- ▶ bed with telescopic siderails
- ▶ bed with plastic siderails or collapsible siderails

To fix the Handset Control Panel on the opposite side of the bed with telescopic siderails:

- 1) Remove the plastic cover from the seat section of the mattress support platform.
- 2) Unlock the Handset Control Panel Cable from the Handset Control Panel Cable Holder on the bed side under the thighrest.
Turn the plastic lock to the right to release the cable from the holder.
Take the cable out of the holder.
- 3) Lift the thighrest by hand to make space between the seat section and thighrest for the Handset Control Panel.
- 4) Run the Handset Control Panel upwards between the seat section and thighrest.
- 5) Take the Handset Control Panel Cable out of the cable clip on the seat section frame.
- 6) Put the Handset Control Panel Cable into the cable clip on the opposite side.
- 7) Lift the thighrest by hand to make space between the seat section and thighrest for the Handset Control Panel.
- 8) Run the Handset Control Panel downwards between the seat section and thighrest on the opposite bed side.
- 9) Insert the Handset Control Panel Cable into the Handset Control Panel Cable Holder on the opposite bed side under the thighrest.
- 10) Insert the plastic lock into the holder and lock the cable in the holder by turning the plastic lock to the left.
- 11) Return the plastic cover back into the seat section of the mattress support platform.

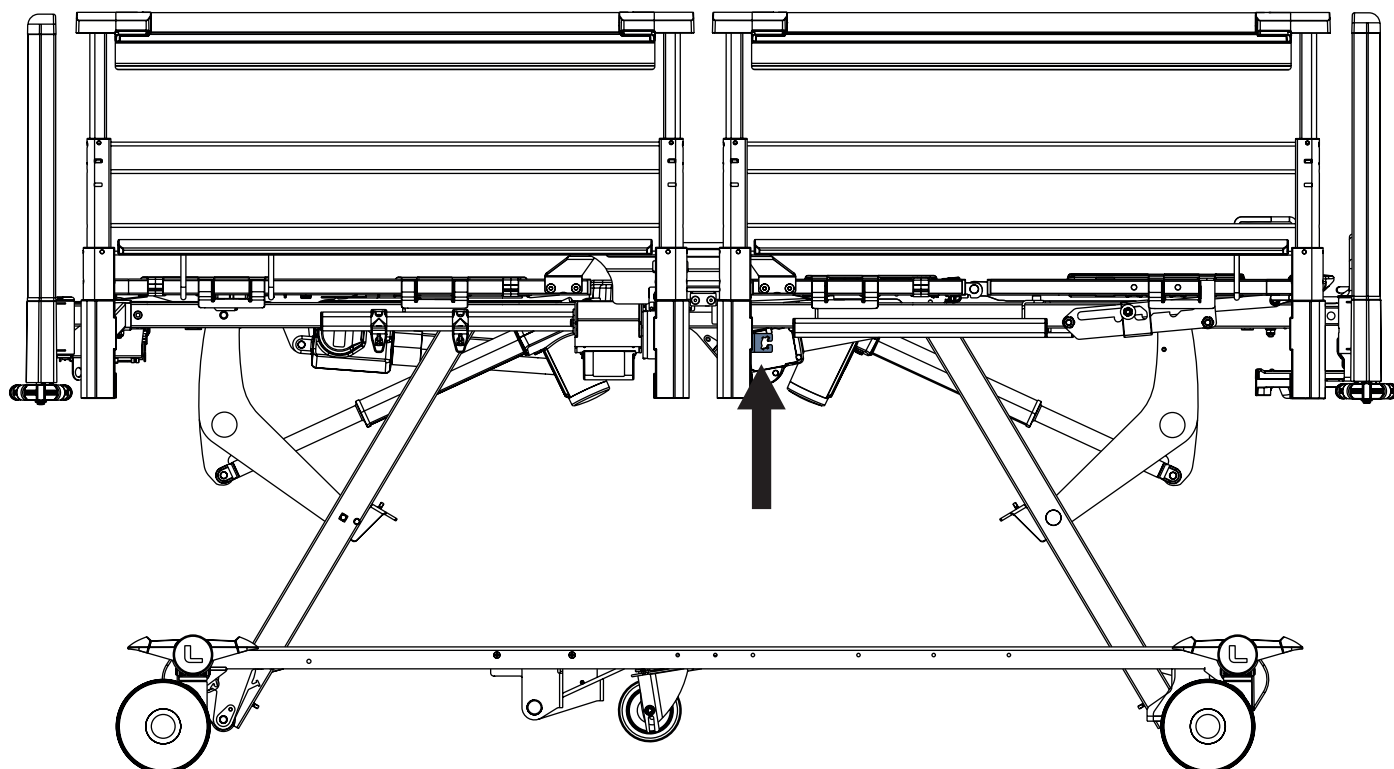


Fig. Position of the Handset Control Panel Cable Holder on the bed with telescopic siderails

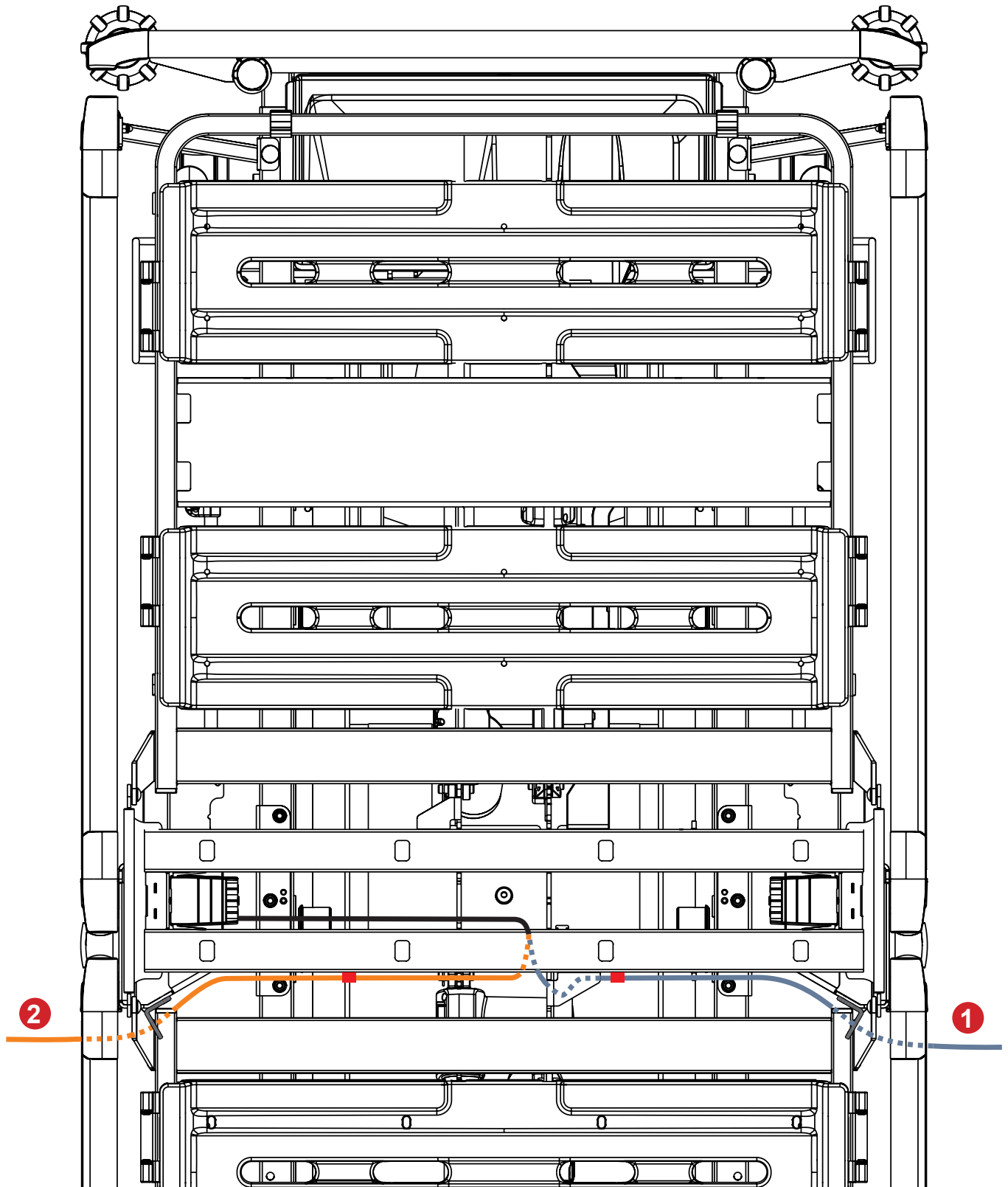


Fig. Handset Control Panel Cable Routing between the seat section and thighrest on the bed with telescopic siderails

1. Cable routing assembled in the factory
2. Optional cable routing with Handset Control Panel on the right bed side

To fix the Handset Control Panel on the opposite side of the bed with plastic siderails or collapsible siderails:

- 1) Remove the plastic cover from the seat section of the mattress support platform.
- 2) Unlock the Handset Control Panel Cable from the Handset Control Panel Cable Holder on the bed side under the seat section.
Turn the plastic lock to the right to release the cable from the holder.
Take the cable out of the holder.
- 3) Lift the thighrest by hand to make space between the seat section and thighrest for the Handset Control Panel.
- 4) Run the Handset Control Panel upwards between the seat section and thighrest.
- 5) Take the Handset Control Panel Cable out of the cable clip on the seat section frame.
- 6) Put the Handset Control Panel Cable into the cable clip on the opposite side.
- 7) Lift the thighrest by hand to make space between the seat section and thighrest for the Handset Control Panel.
- 8) Run the Handset Control Panel downwards between the seat section and thighrest on the opposite bed side.
- 9) Insert the Handset Control Panel Cable into the Handset Control Panel Cable Holder on the opposite bed side under the seat section.
- 10) Insert the plastic lock into the holder and lock the cable in the holder by turning the plastic lock to the left.
- 11) Return the plastic cover back into the seat section of the mattress support platform.

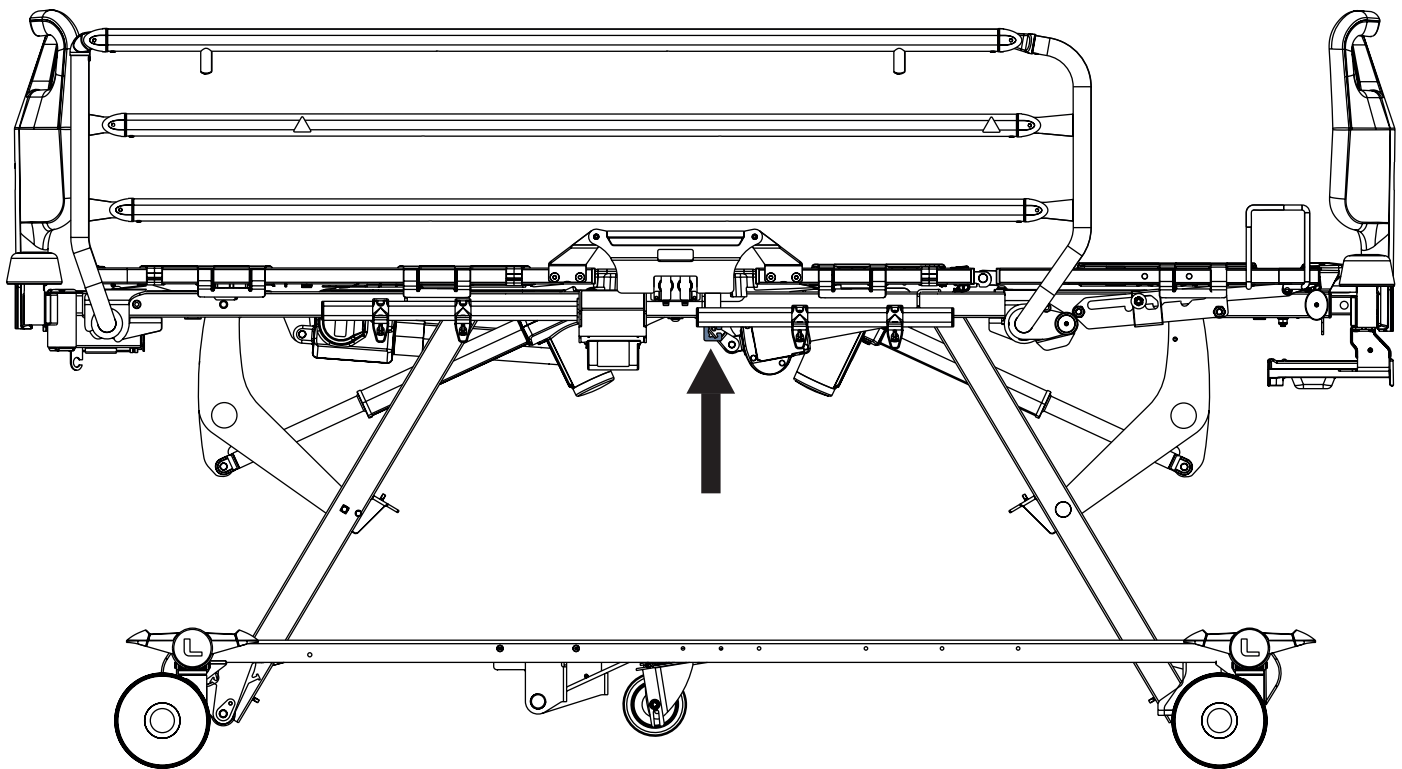


Fig. Position of the Handset Control Panel Cable Holder on the bed with collapsible siderails

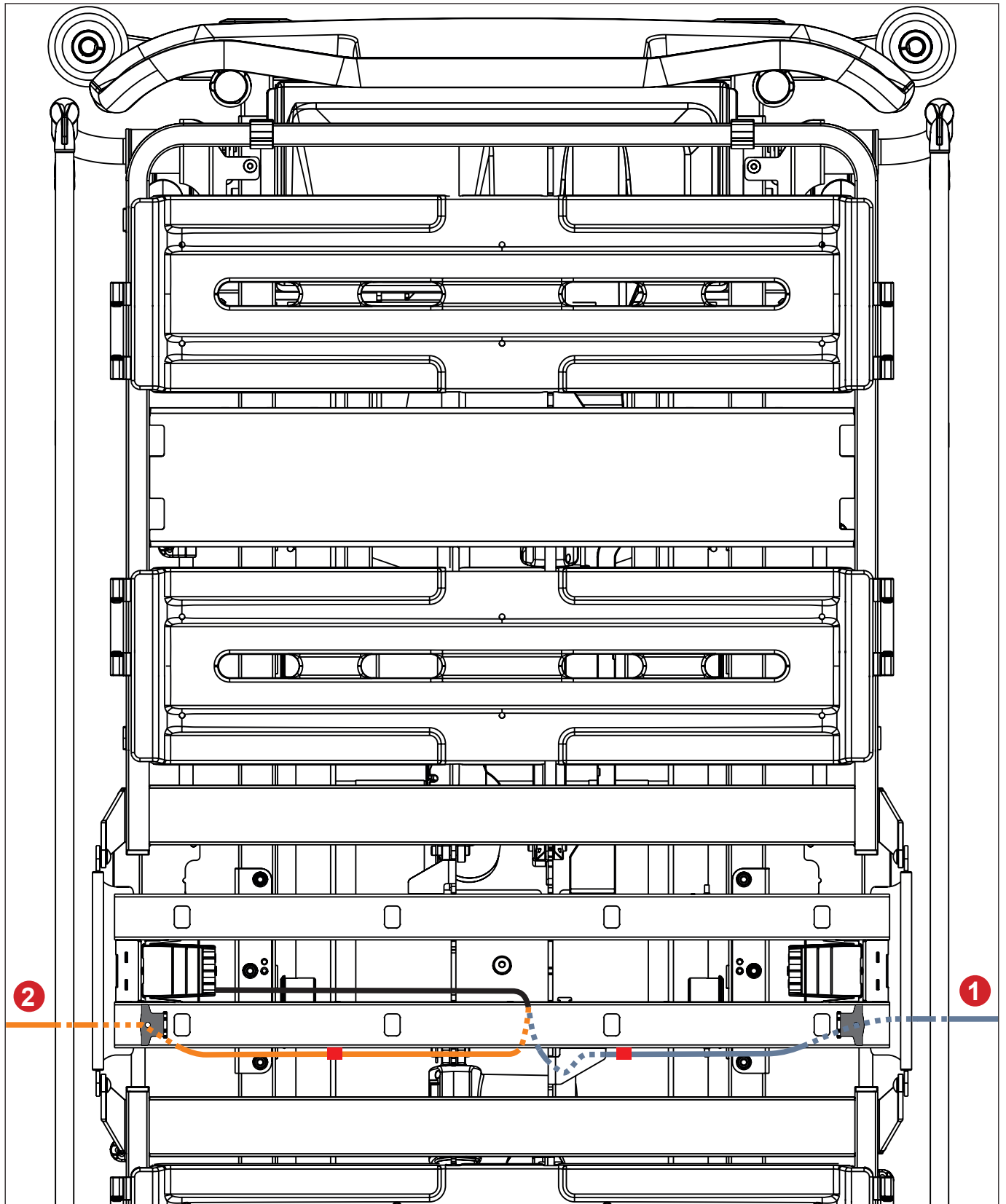


Fig. Handset Control Panel Cable Routing between the seat section and thighrest on the bed with collapsible siderails

1. Cable routing assembled in the factory
2. Optional cable routing with Handset Control Panel on the right bed side

12.6 Patient Control Panel

Patient is allowed to use the Patient Control Panel only if hospital personnel had assessed that the patient's physical and psychological state is in accordance with use of this control element and only if the hospital personnel had trained the patient in accordance with these instructions for use.

The use of the Patient Control Panel is described in the chapter Bed Positioning (see „13 Bed Positioning“ on the page 82).

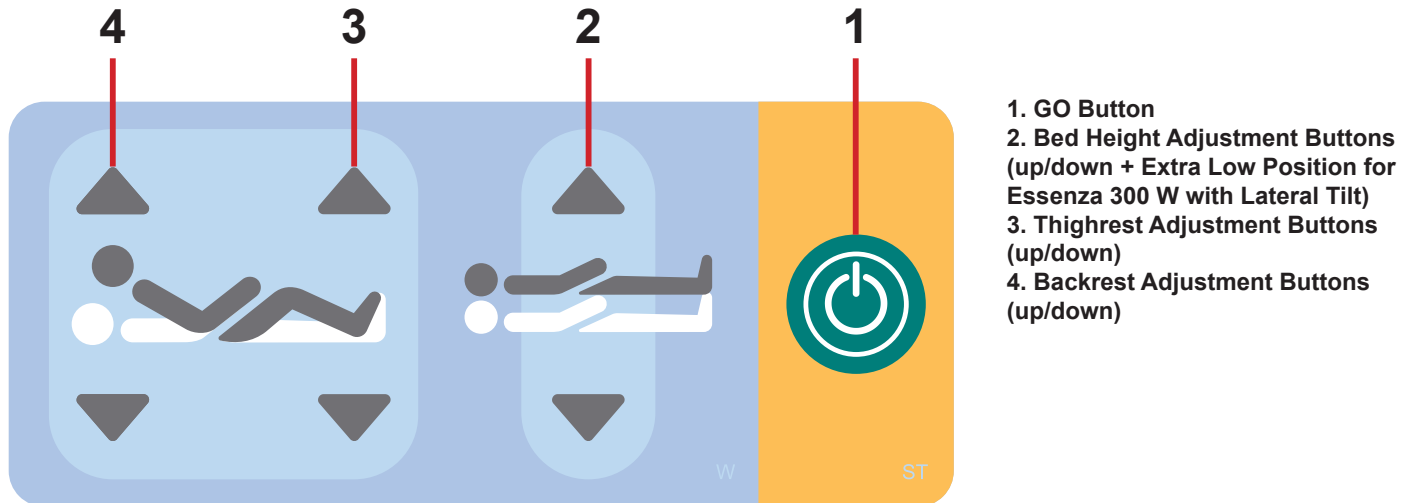


Fig. Patient Control Panel

To disable bed positioning functions for patient, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W

13 Bed Positioning



WARNING!

Risk of inappropriate reduction of the available bed positioning ways! Foot Control Pedals are not accessible in the Extra Low Position! Backrest Handles and Calfrest Handles are not ready for use when siderails are in lower position!

- ▶ The hospital personnel is responsible for adjusting the appropriate bed position for each patient!



WARNING!

Risk of injury when adjusting the bed!

- ▶ Ensure that there are no body parts between the Mattress Support Platform parts and the Mattress Support Platform frame when adjusting the bed.
- ▶ Ensure that there are no body parts below the Mattress Support Platform frame before adjusting the bed.



WARNING!

Risk of injury or material damage when adjusting the bed!

- ▶ Ensure that there are no installed accessories that may collide with moving parts of the bed.
- ▶ Extra care must be taken to avoid collisions of the accessories with bed parts during bed positioning.



WARNING!

Risk of injury due to incorrect use!

- ▶ Medical personnel is responsible for siderails adjustment and Lateral Tilt adjustment in accordance with patient's safety in bed!



WARNING!

Risk of collision due to the inappropriate combination of bed positions when bed is equipped with HPL Head Board and HPL Foot Board and is not equipped with Lateral Tilt adjustment!

- ▶ Do not adjust bed extension with maximum Trendelenburg Tilt or maximum Reverse-Trendelenburg Tilt or Cardiac Chair Position when bed is equipped with HPL Head Board and HPL Foot Board and is not equipped with Lateral Tilt adjustment in order to avoid this collision!
- ▶ Make sure that there is no risk of injury to the lower limbs of persons standing in the immediate vicinity of the bed head end and bed foot end by checking and alerting people near the bed before adjusting these bed positions!
- ▶ Do not slide out the Plastic Linen Shelf or Telescopic Linen Shelf in order to avoid even more probable collision in this combination of bed positions in this bed configuration!



CAUTION!

Risk of squeezing hand due to incorrect Backrest adjustment!

- ▶ Take extra care during Backrest adjustment to avoid squeezing your hand!



CAUTION!

Risk of squeezing hand due to incorrect Calfrest adjustment!

- ▶ Take extra care during Calfrest adjustment to avoid squeezing your hand!



CAUTION!

Risk of damaging control elements and their cables due to incorrect positions of the control elements during lowering to the Extra Low Position!

- ▶ Ensure the control elements (e. g. Attendant Control Panel or Handset Control Panel) are located on the safe positions on the bed before and during lowering to the Extra Low Position!

**CAUTION!**

Risk of damaging accessories due to incorrect positions of the accessories during lowering to the Extra Low Position!


- Ensure the accessories are located on the safe positions on the bed before and during lowering to the Extra Low Position!

To activate a selected control element :

- Press GO Button.

A positioning function can be used via the corresponding button on the activated control element.

13.1 Backrest Adjustment

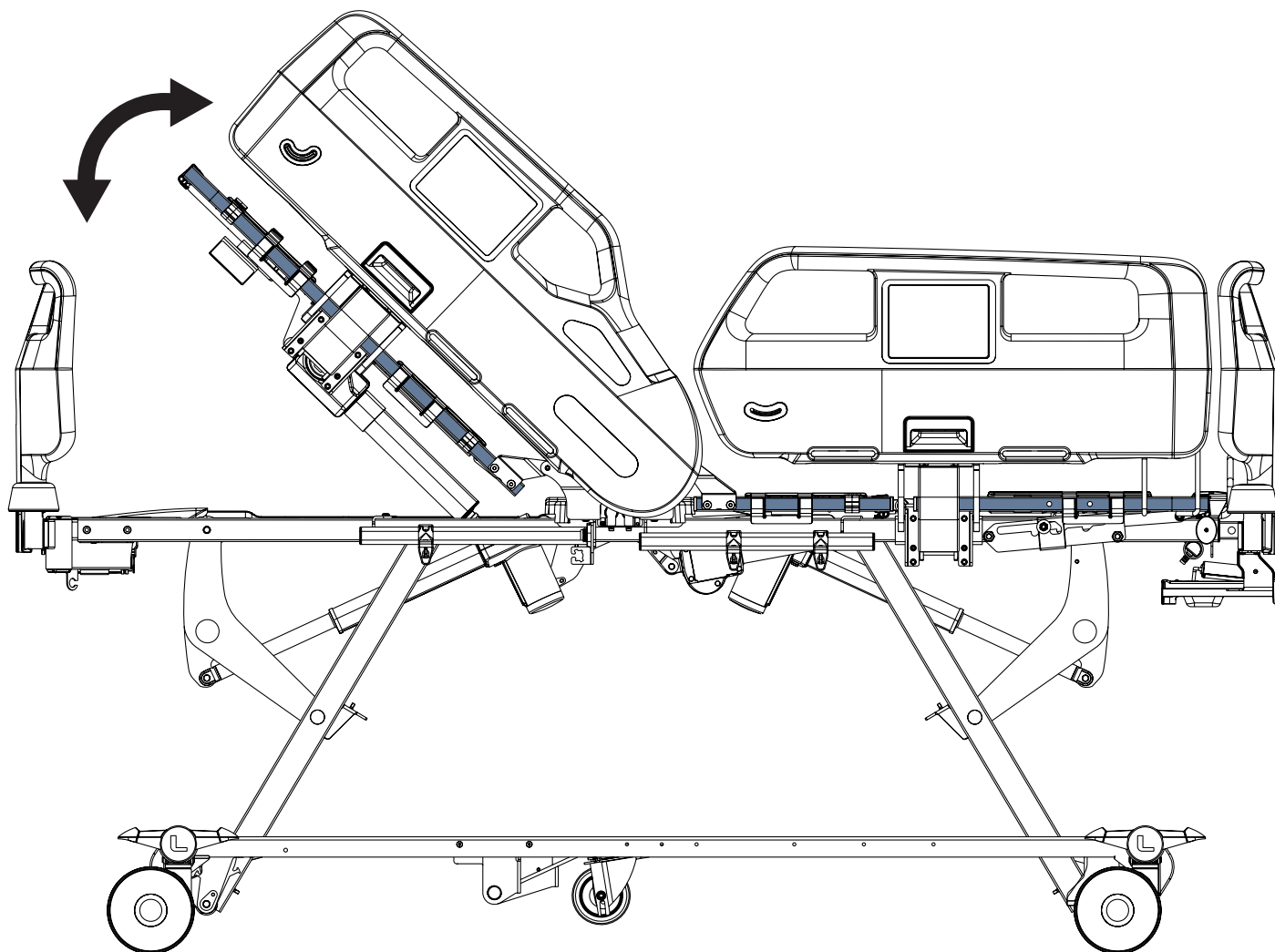



To adjust Backrest, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W
- ▶ Handset Control Panel
- ▶ Patient Control Panel



Buttons intended for Backrest adjustment:

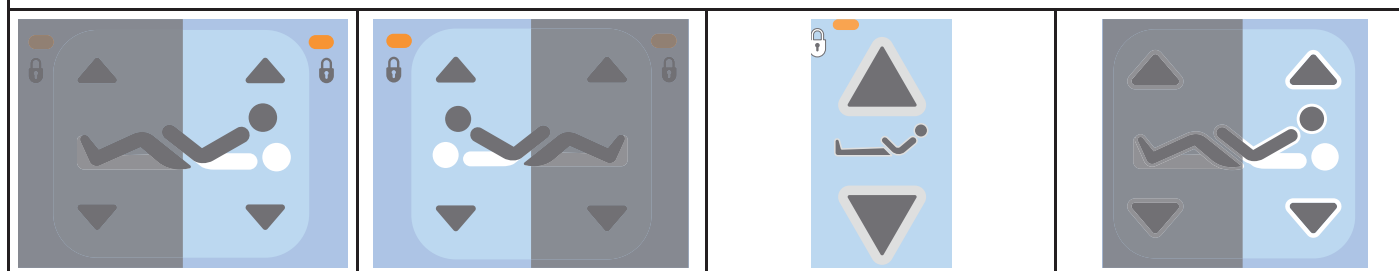

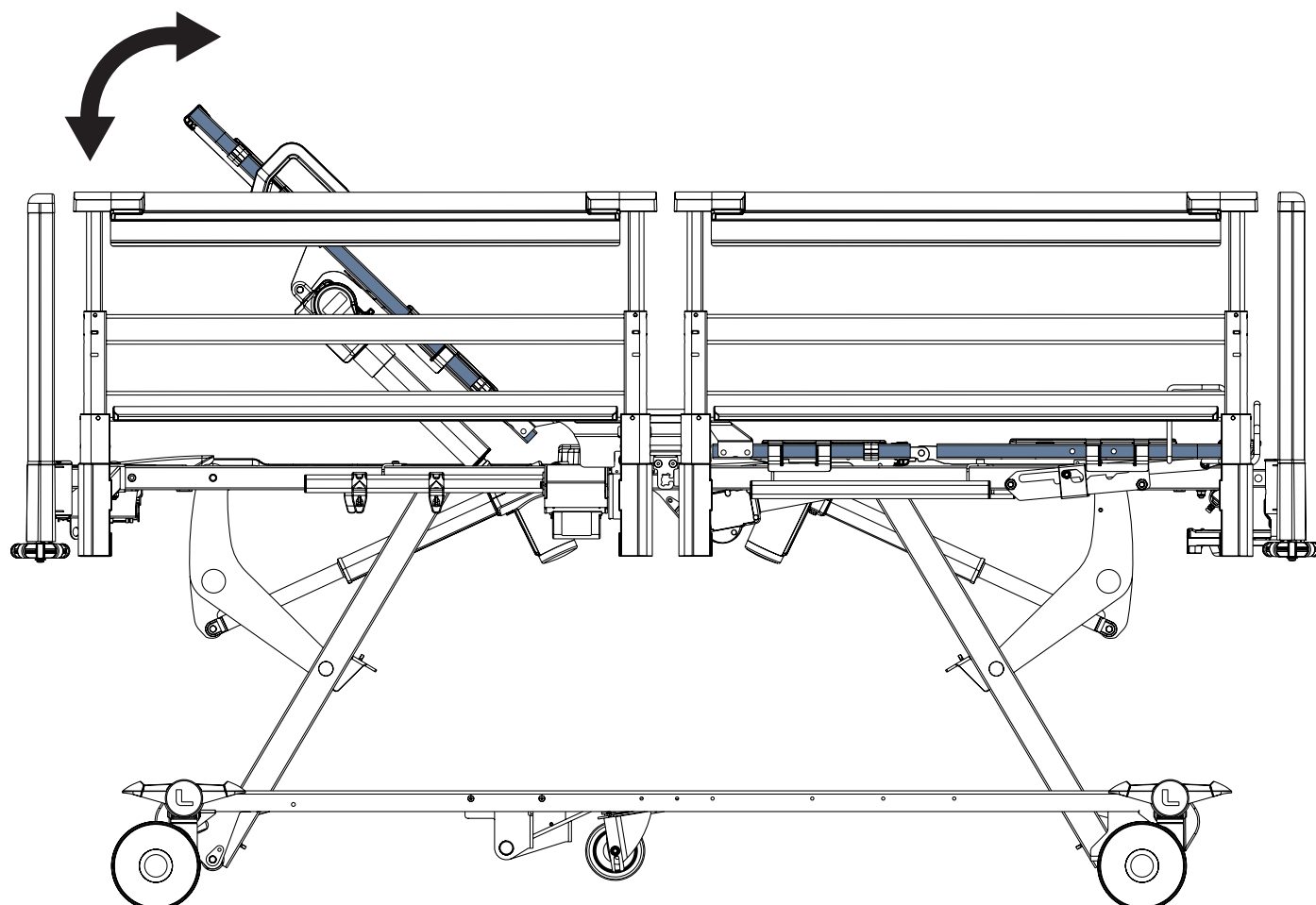


To adjust Backrest, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Handset Control Panel



Buttons intended for Backrest adjustment:

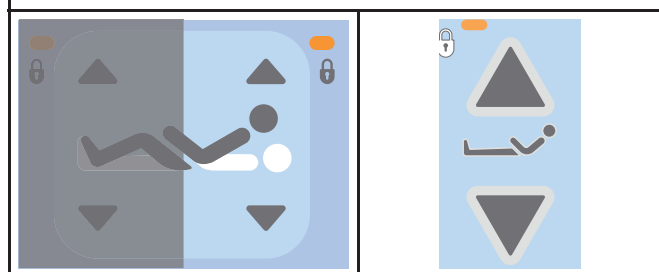



WARNING!

Risk of injury due to incorrect use!

► Ensure that there is no leg or the other body part between the Plastic Head Siderail folded down and the Plastic Foot Siderail folded down when patient sits in the middle of the bed and somebody lowers Backrest! Medical personnel are responsible for instructing the patient to avoid this risk if the patient is allowed to use Handset Control Panel or Patient Control Panel!


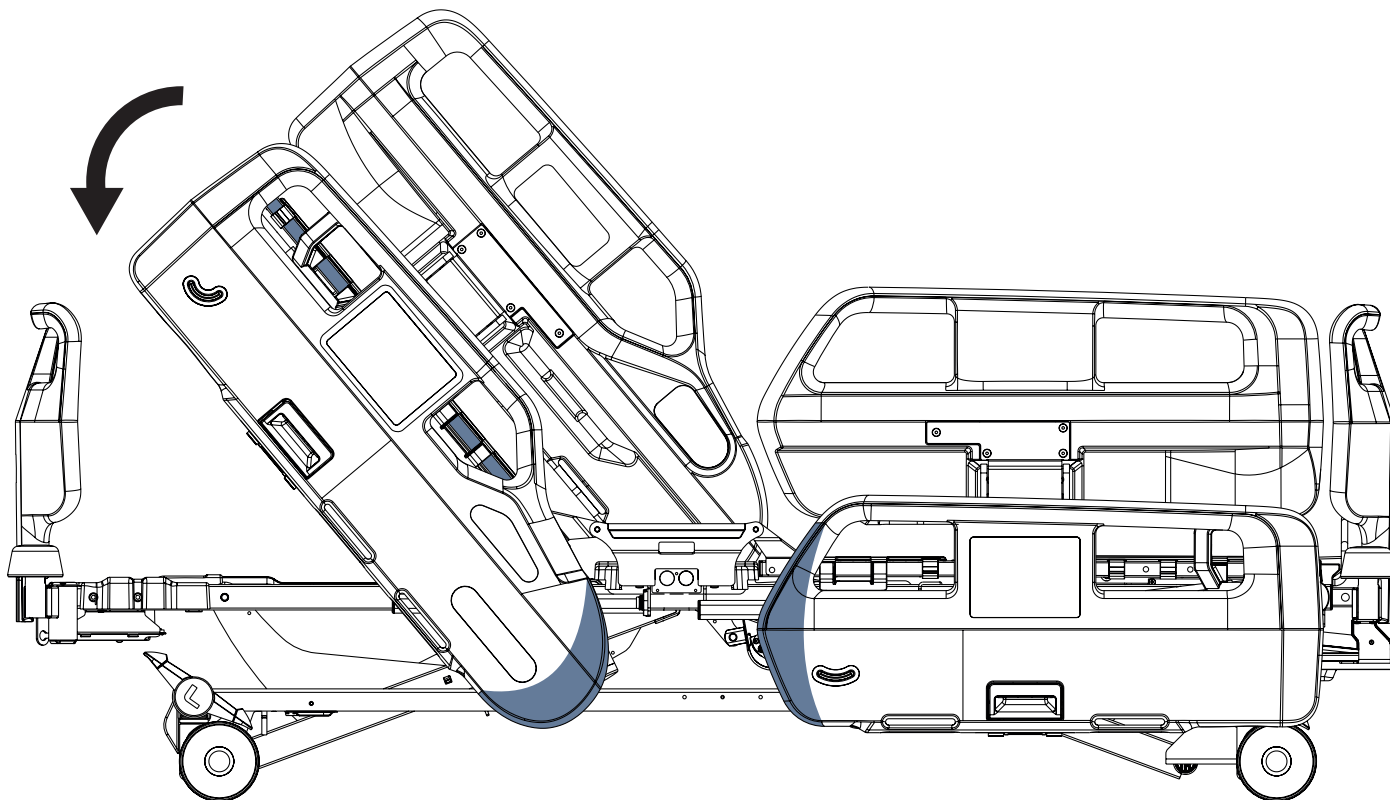


Fig. Bed adjustment unsuitable for sitting in the middle of the bed

13.2 Thighrest Adjustment

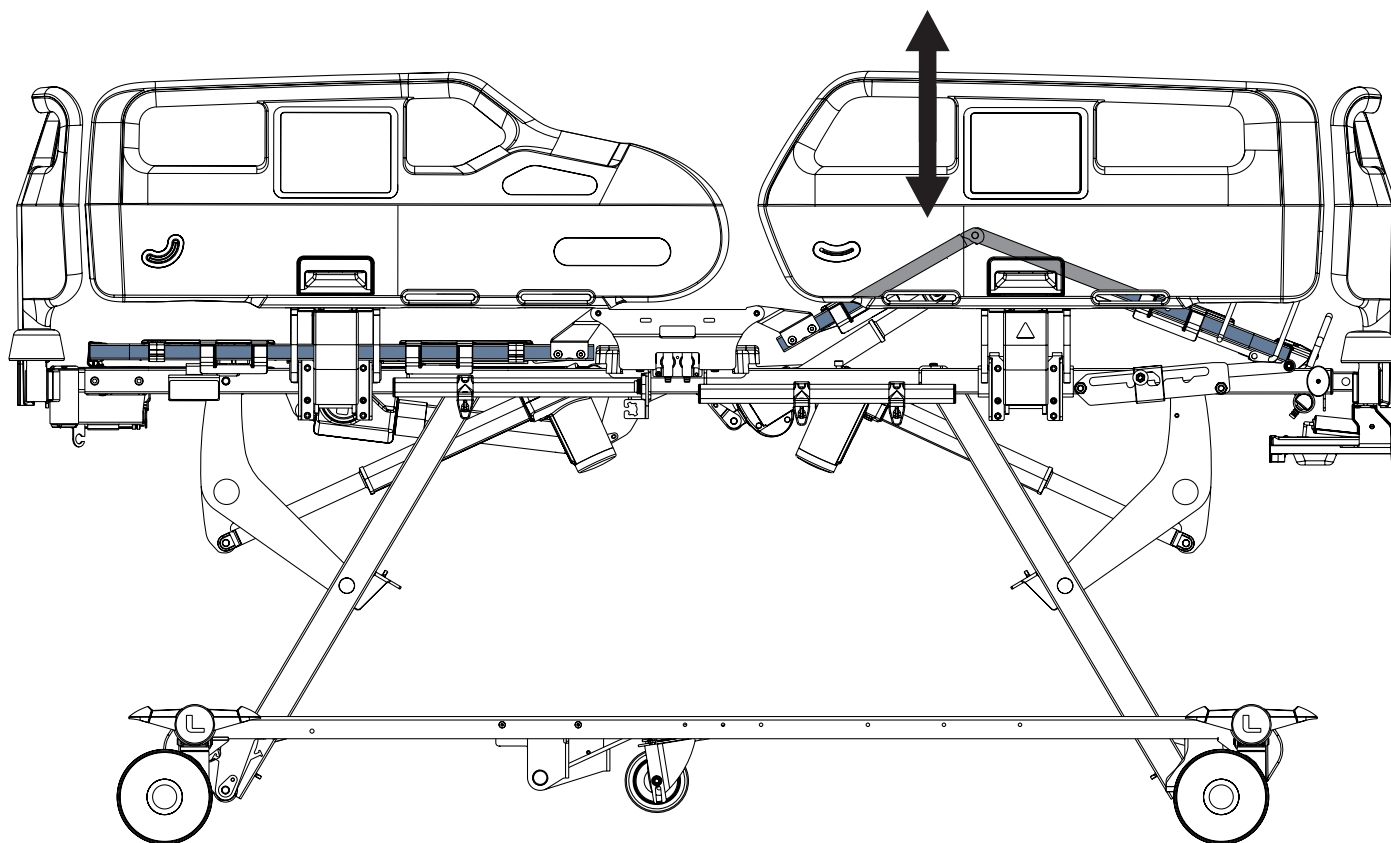



To adjust Thighrest, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W
- ▶ Handset Control Panel
- ▶ Patient Control Panel



Buttons intended for Thighrest adjustment:

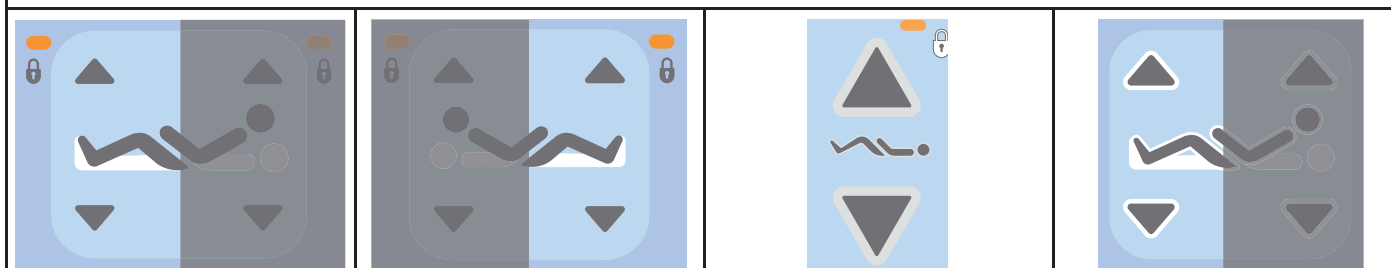

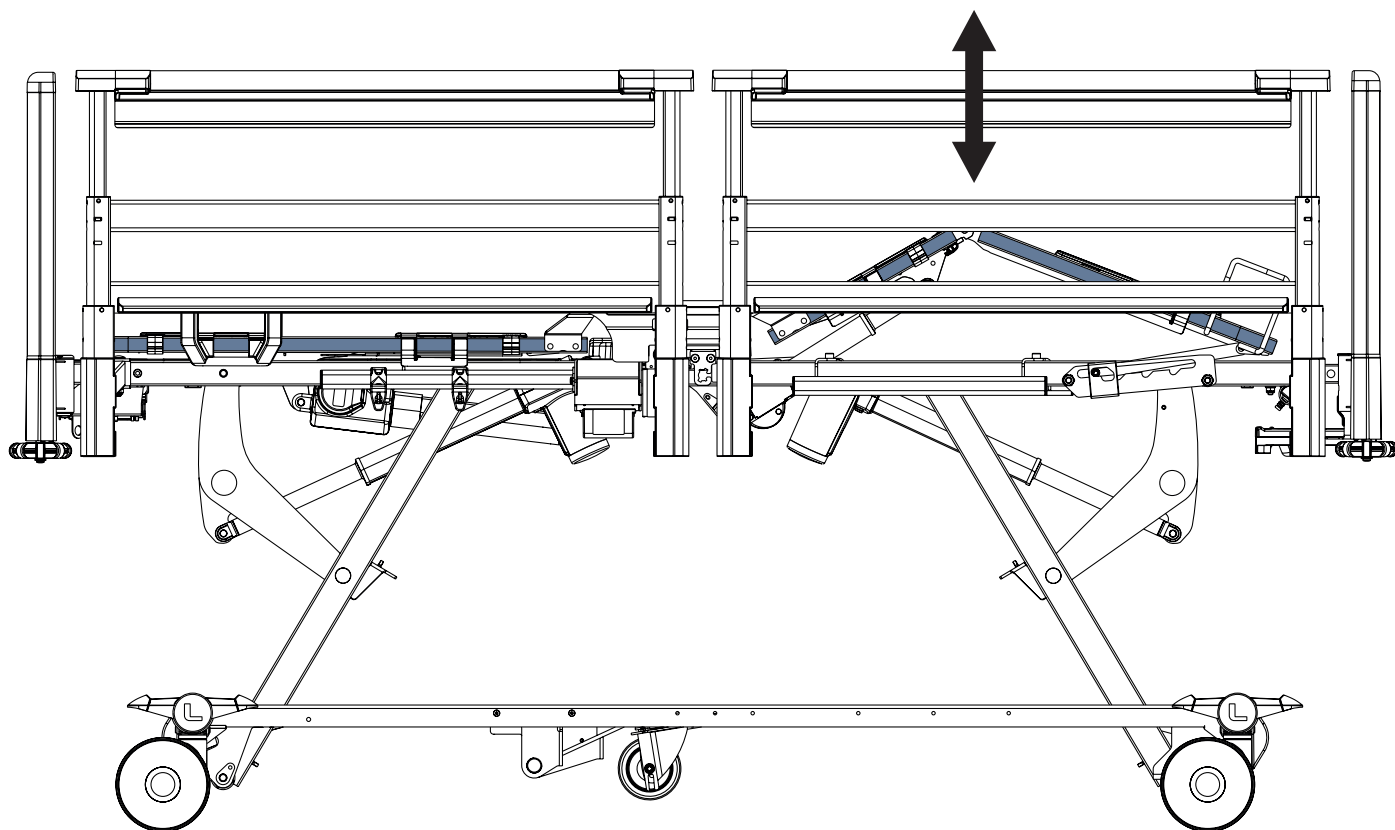


To adjust Thighrest, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Handset Control Panel



Buttons intended for Thighrest adjustment:

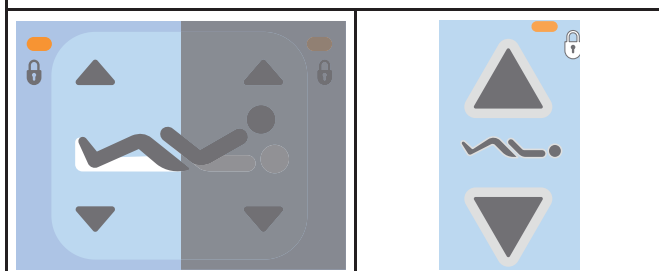

13.3 Calfrest Adjustment

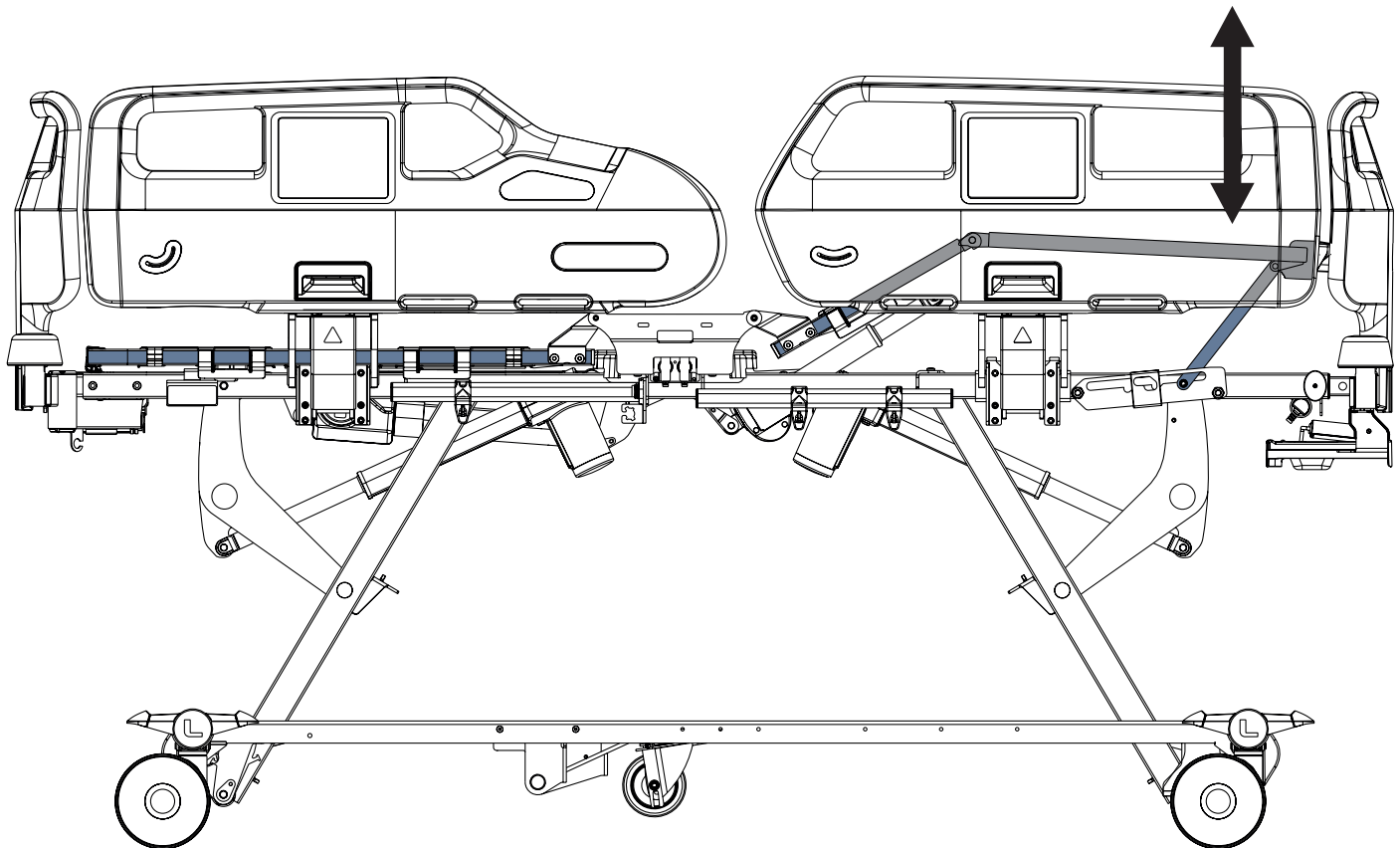

To lift Calfrest:

- ▶ Grab the bottom part of the Calfrest Handle with one hand and the Siderail Release Handle with the other hand.
- ▶ Lift Calfrest with one hand to intended position and fold the foot siderail down with the other hand.
- ▶ Lower the Calfrest so that catch fits in the ratchet-bar.

To lower Calfrest:

- ▶ Lift Calfrest slightly by the Calfrest Handle.
- ▶ Lower Calfrest to intended position.
- ▶ Ensure the catch fits in the ratchet-bar during slight lifting of the Calfrest.


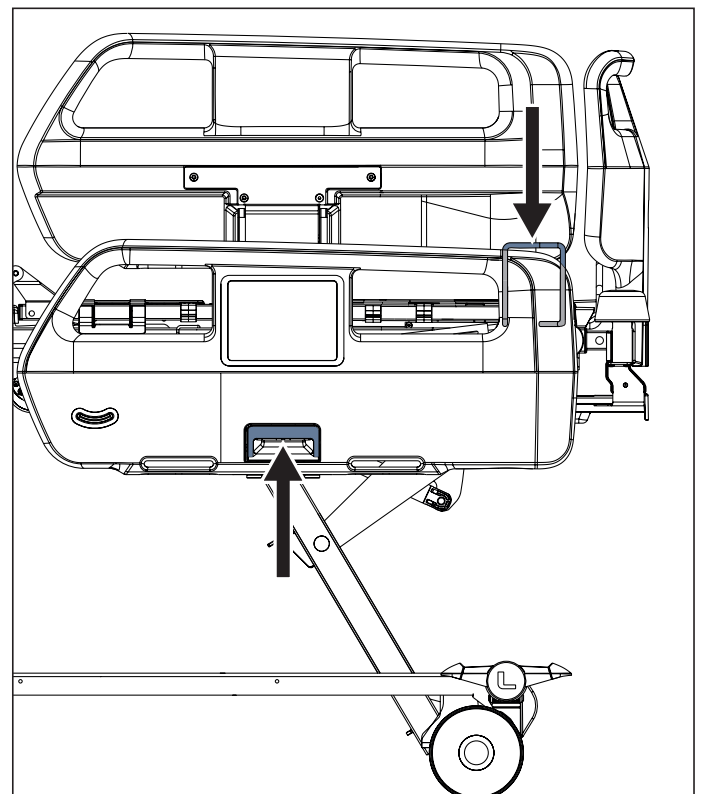


Fig. Calfrest Adjustment

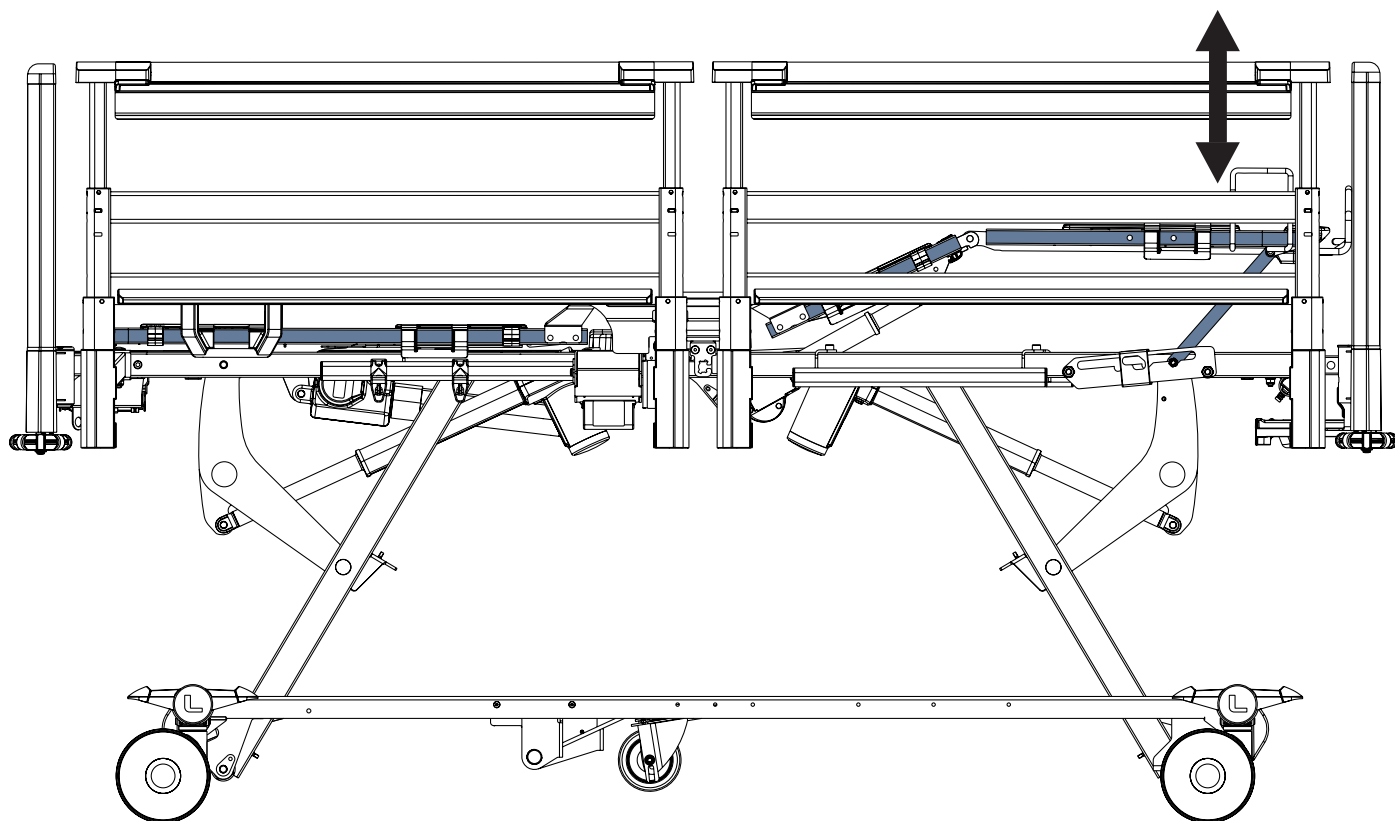

To lift Calfrest:

- ▶ Lift Calfrest by the Calfrest Handle to intended position.
- ▶ Lower the Calfrest so that catch fits in the ratchet-bar.

To lower Calfrest:

- ▶ Lift Calfrest slightly by the Calfrest Handle.
- ▶ Lower Calfrest to intended position.
- ▶ Ensure the catch fits in the ratchet-bar during slight lifting of the Calfrest.


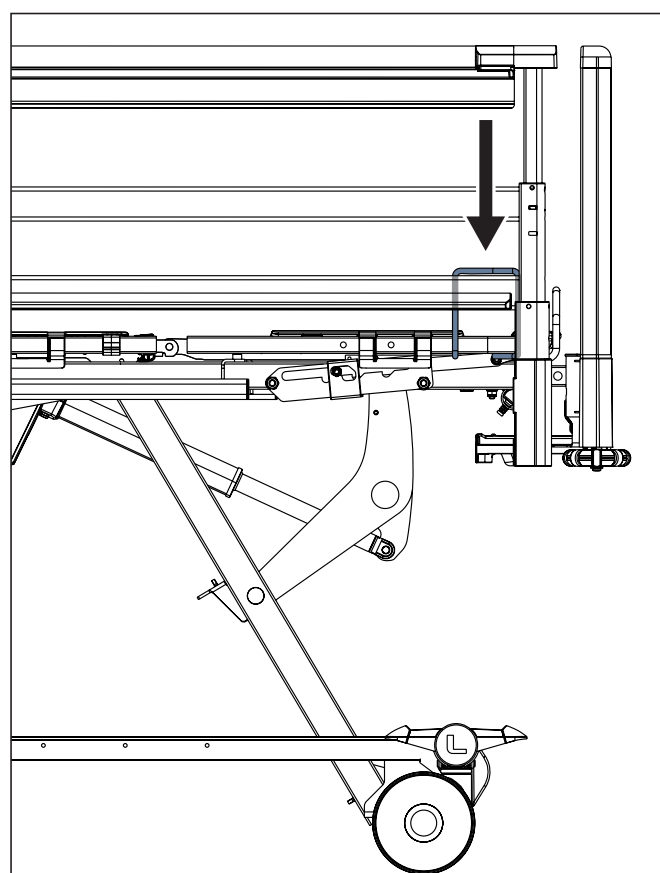


Fig. Calfrest Handle

13.4 Autocontour Adjustment

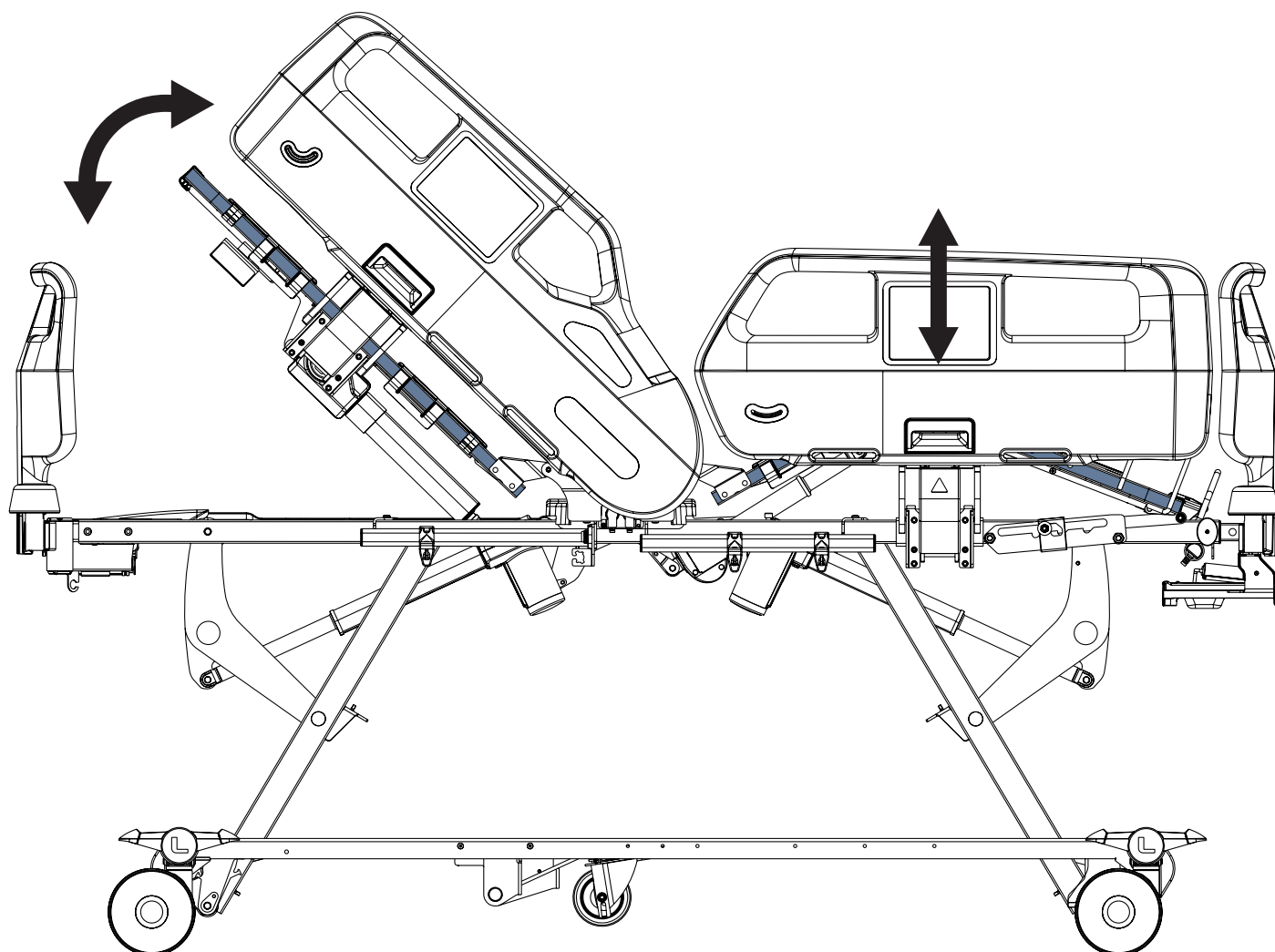



To adjust Autocontour, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W
- ▶ Handset Control Panel



Buttons intended for Autocontour adjustment:

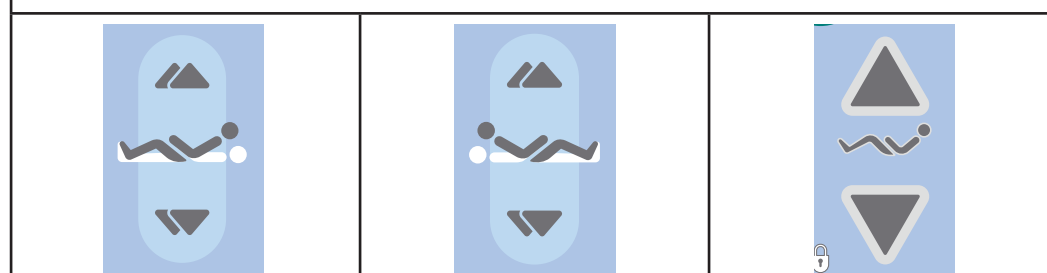

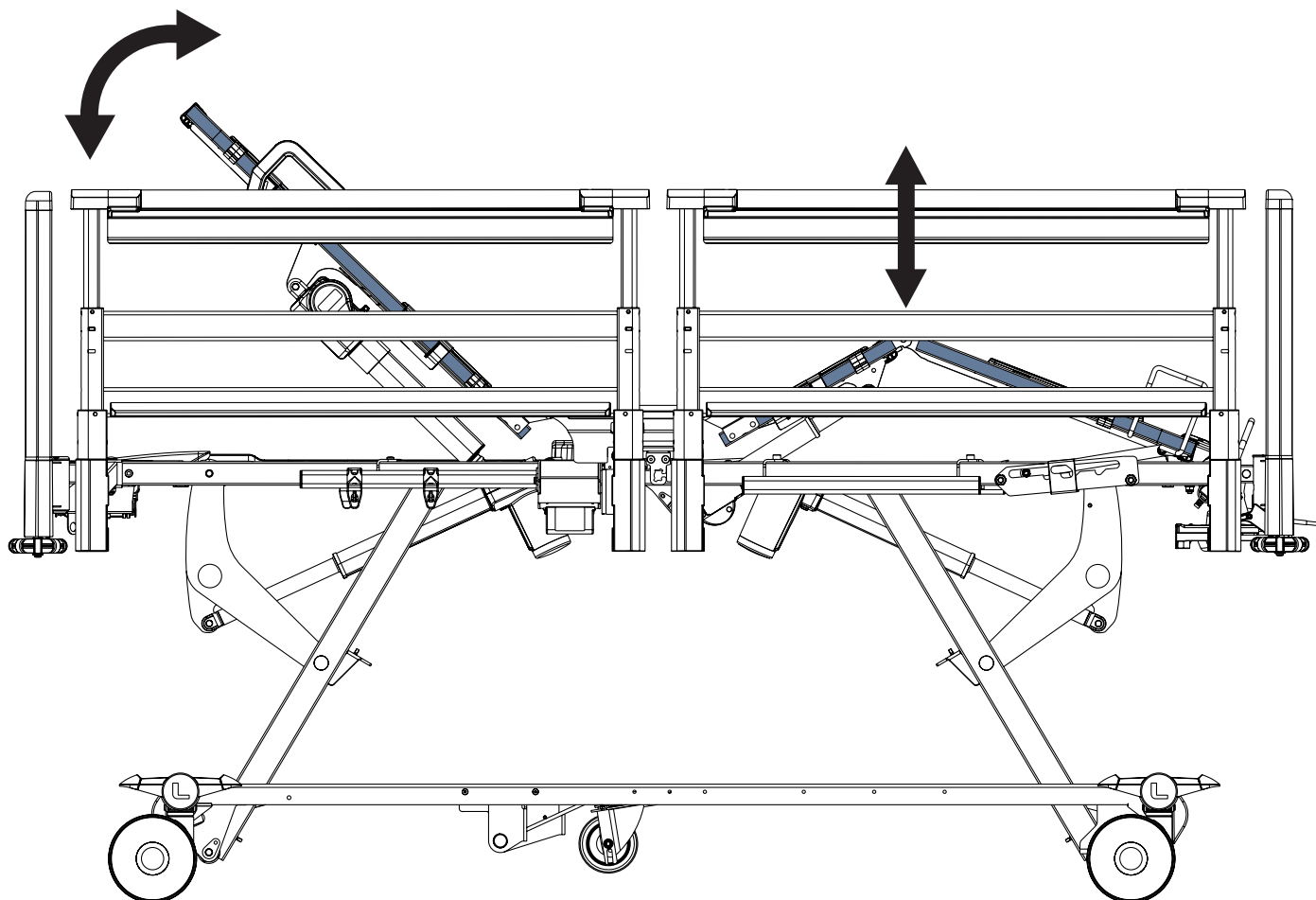


To adjust Autocontour, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Handset Control Panel



Buttons intended for Autocontour adjustment:

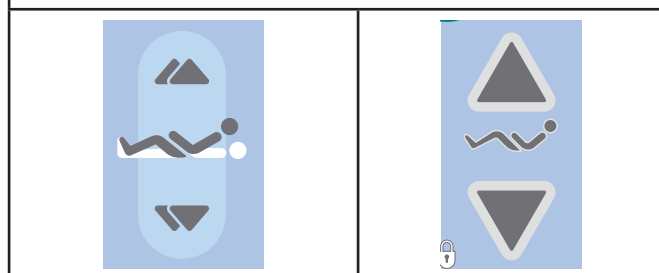

13.5 Bed Height Adjustment

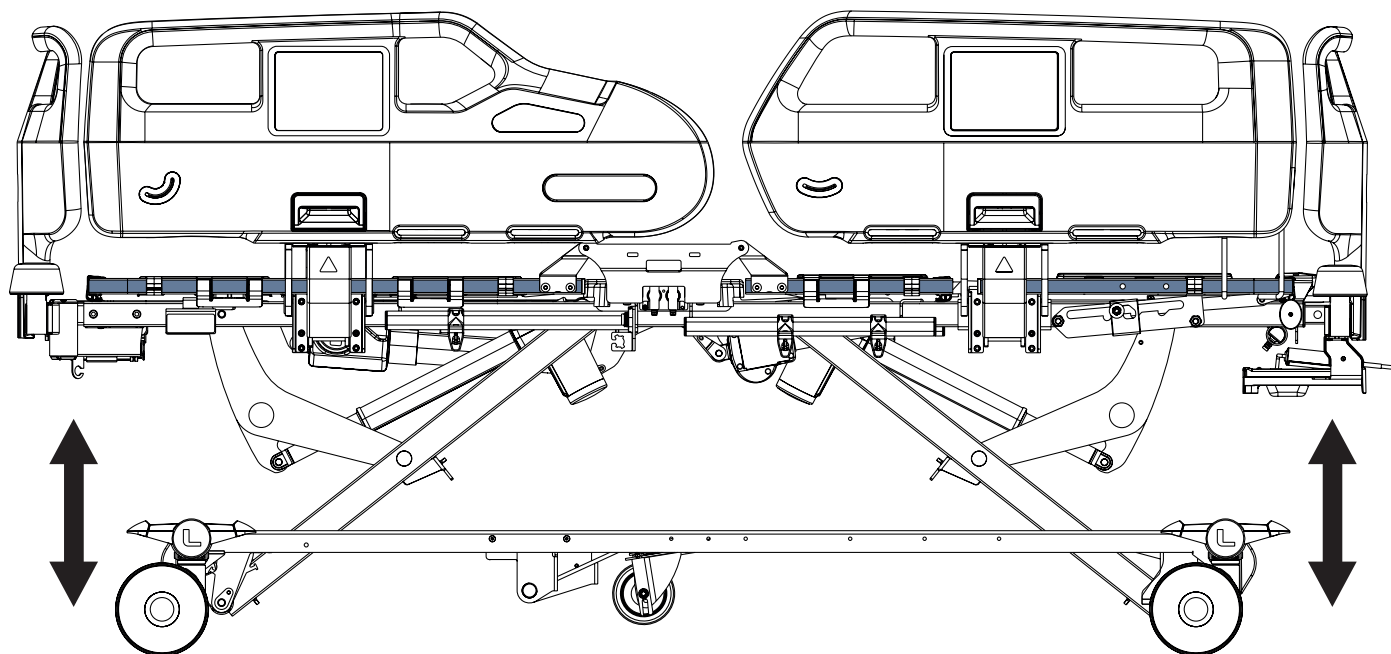



To adjust Bed Height, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W
- ▶ Handset Control Panel
- ▶ Patient Control Panel



Buttons intended for Bed Height adjustment:

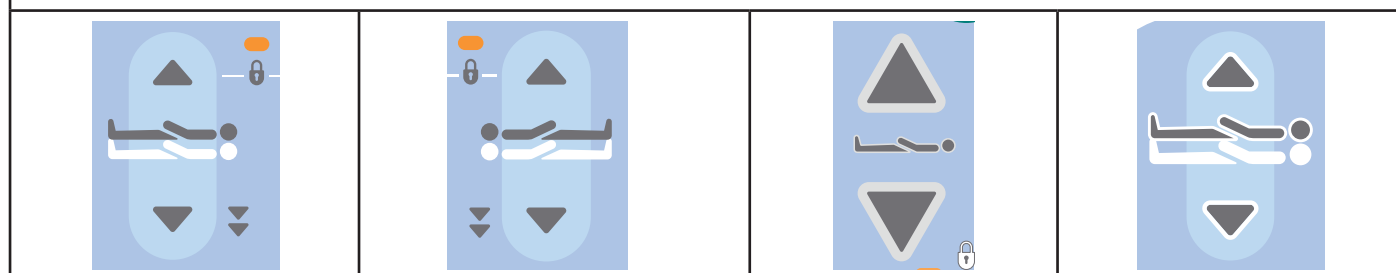

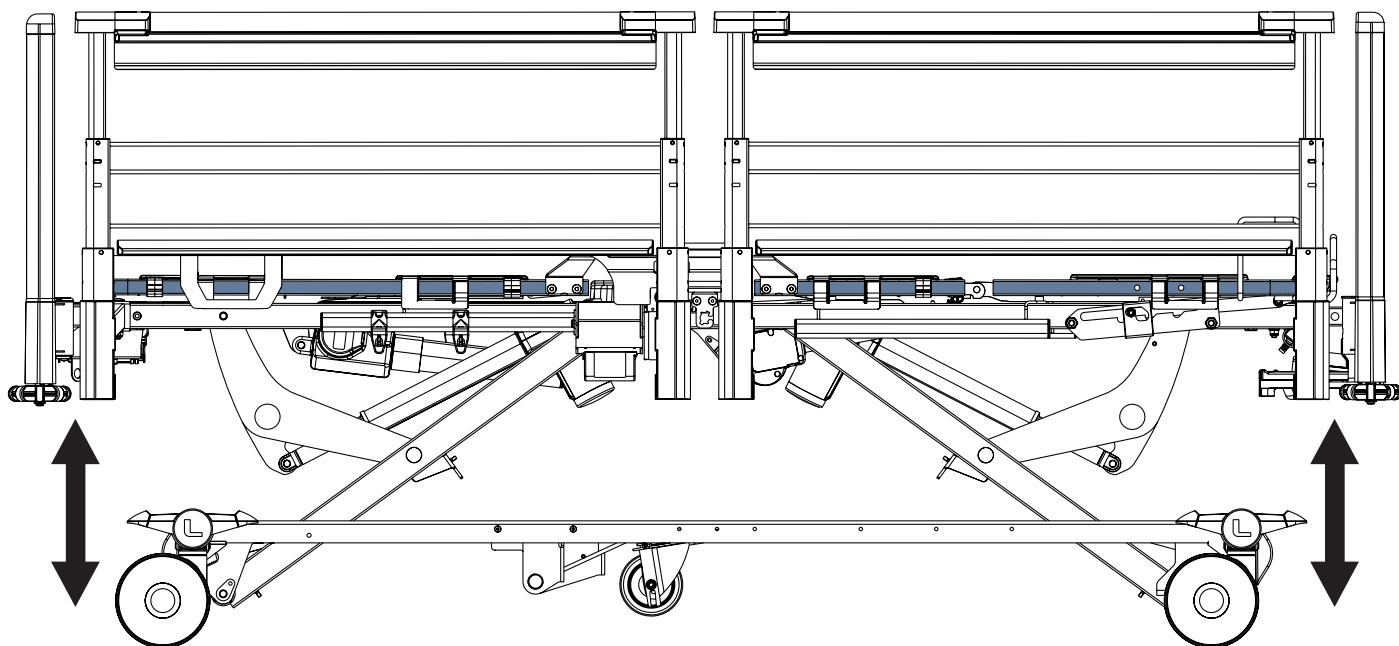


To adjust Bed Height, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Handset Control Panel



Buttons intended for Bed Height adjustment:

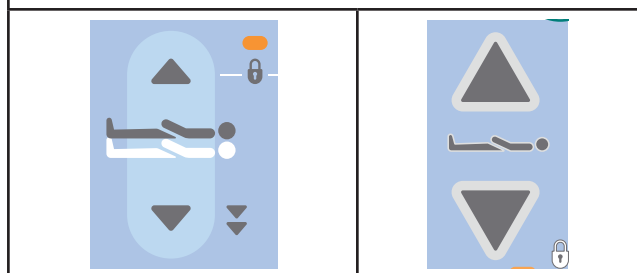

13.6 Emergency Trendelenburg Position

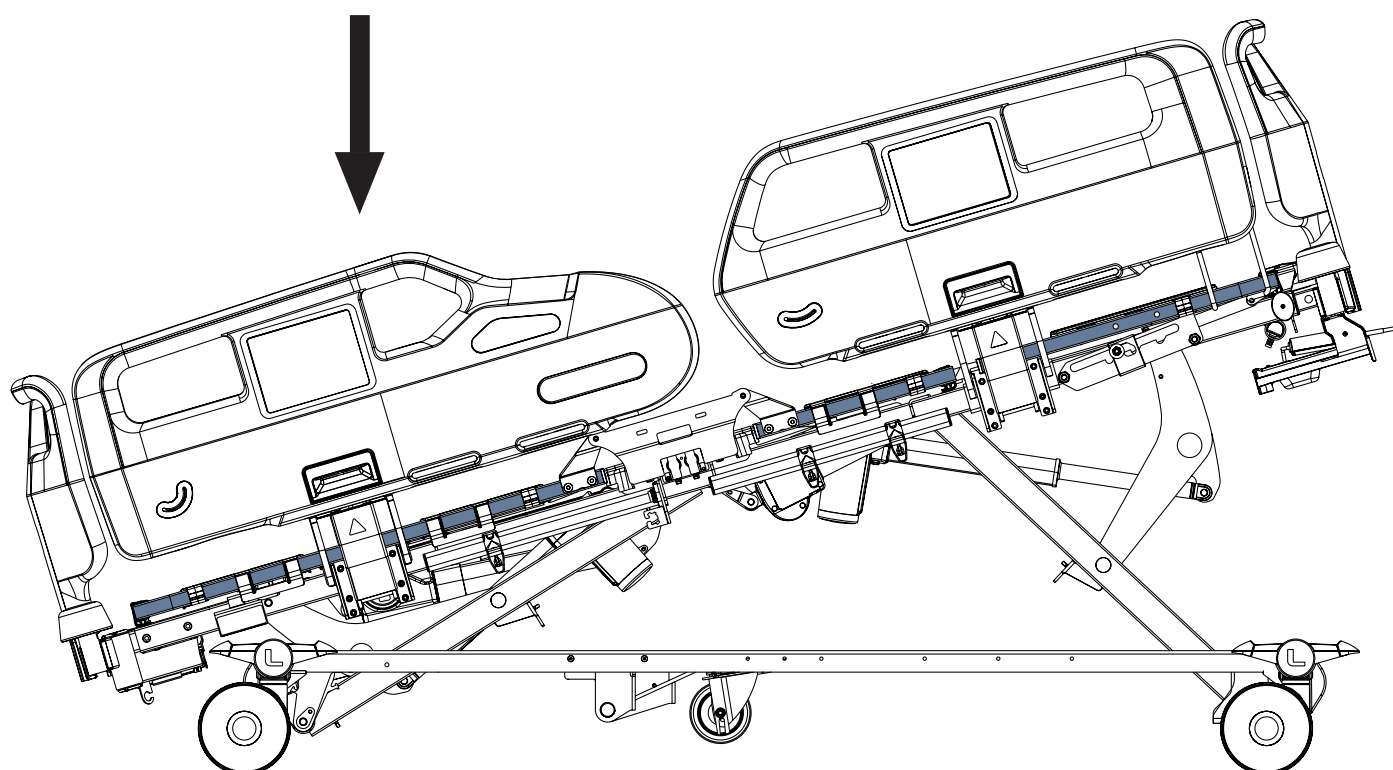

WARNING!

Risk of injury due to improper use of Trendelenburg Position!

- Hospital staff is responsible for assessing if the physical and psychological state of the patient is in accordance with use of the Trendelenburg Position.

It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.


To adjust Emergency Trendelenburg Position, use:

- Attendant Control Panel for Essenza 300 W with Lateral Tilt
- Attendant Control Panel for Essenza 300 W
- Siderail Control Panel for Essenza 300 W with Lateral Tilt
- Siderail Control Panel for Essenza 300 W



Buttons intended for Emergency Trendelenburg Position adjustment:



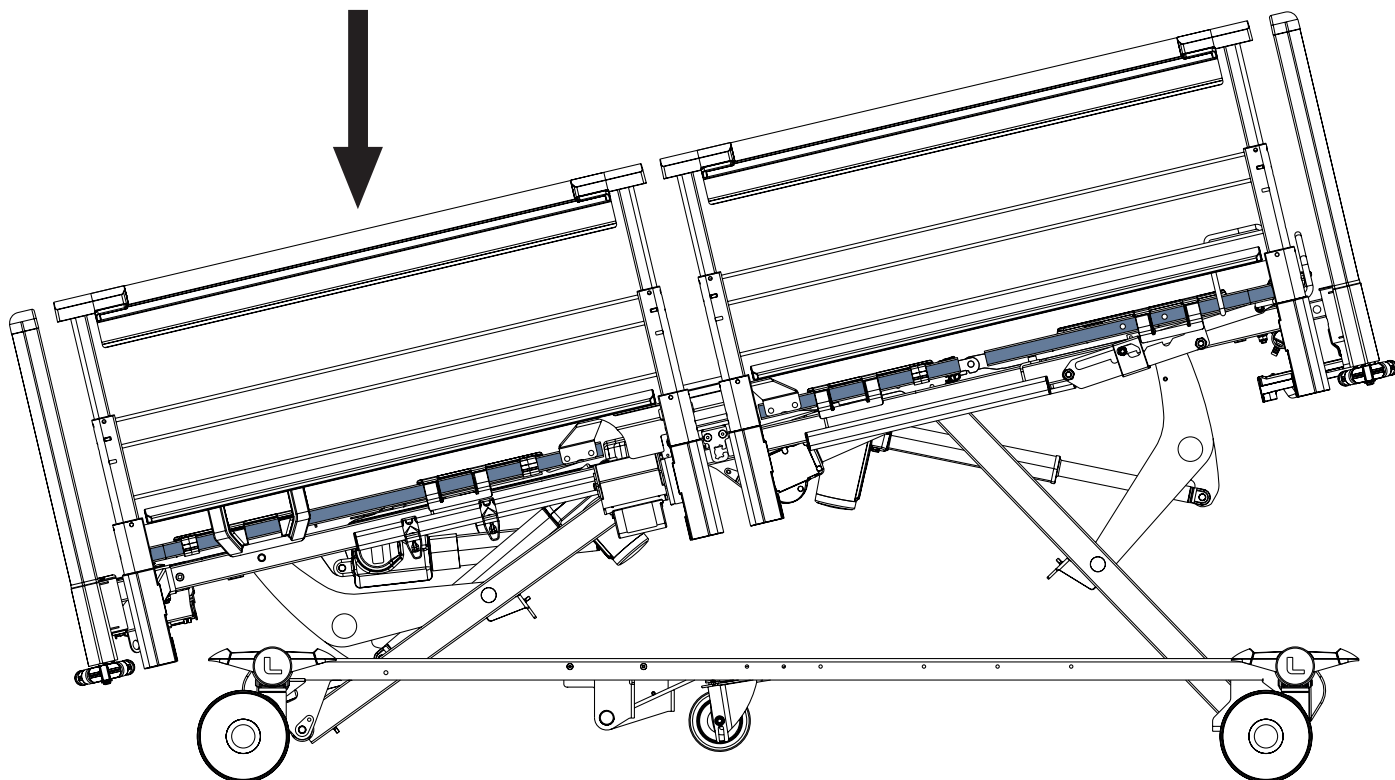

WARNING!

Risk of injury due to improper use of Trendelenburg Position!

- Hospital staff is responsible for assessing if the physical and psychological state of the patient is in accordance with use of the Trendelenburg Position.

It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.

To adjust Emergency Trendelenburg Position, use:


- Attendant Control Panel for Essenza 300 W with Lateral Tilt
- Attendant Control Panel for Essenza 300 W



Button intended for Emergency Trendelenburg Position adjustment:



13.7 Trendelenburg Tilt

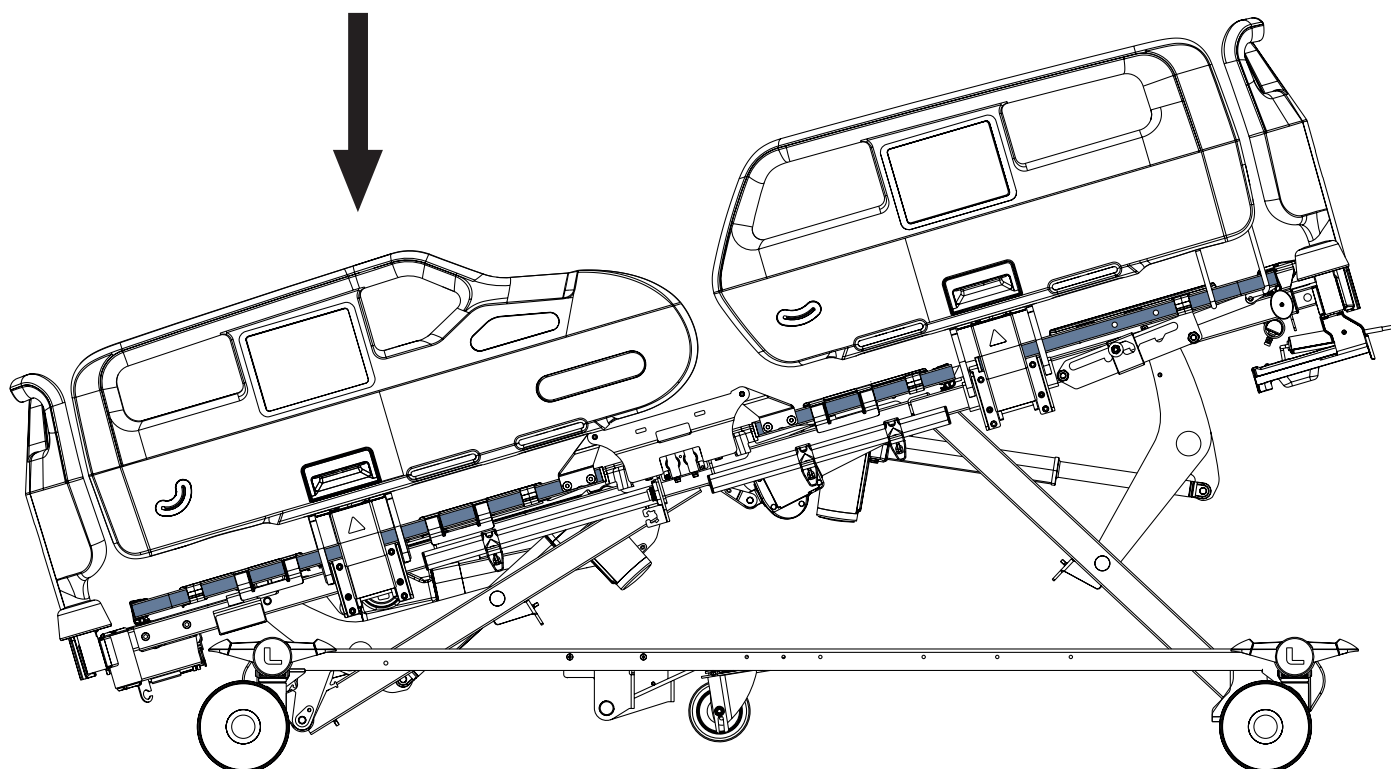

WARNING!

Risk of injury due to improper use of Trendelenburg Position!

► Hospital staff is responsible for assessing if the physical and psychological state of the patient is in accordance with use of the Trendelenburg Position.



It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.


To adjust Trendelenburg Tilt, use:

- Attendant Control Panel for Essenza 300 W with Lateral Tilt
- Attendant Control Panel for Essenza 300 W
- Siderail Control Panel for Essenza 300 W with Lateral Tilt
- Siderail Control Panel for Essenza 300 W



Buttons intended for Trendelenburg Tilt adjustment:

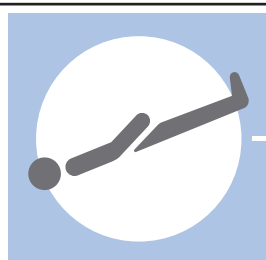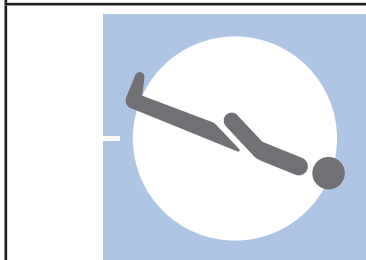

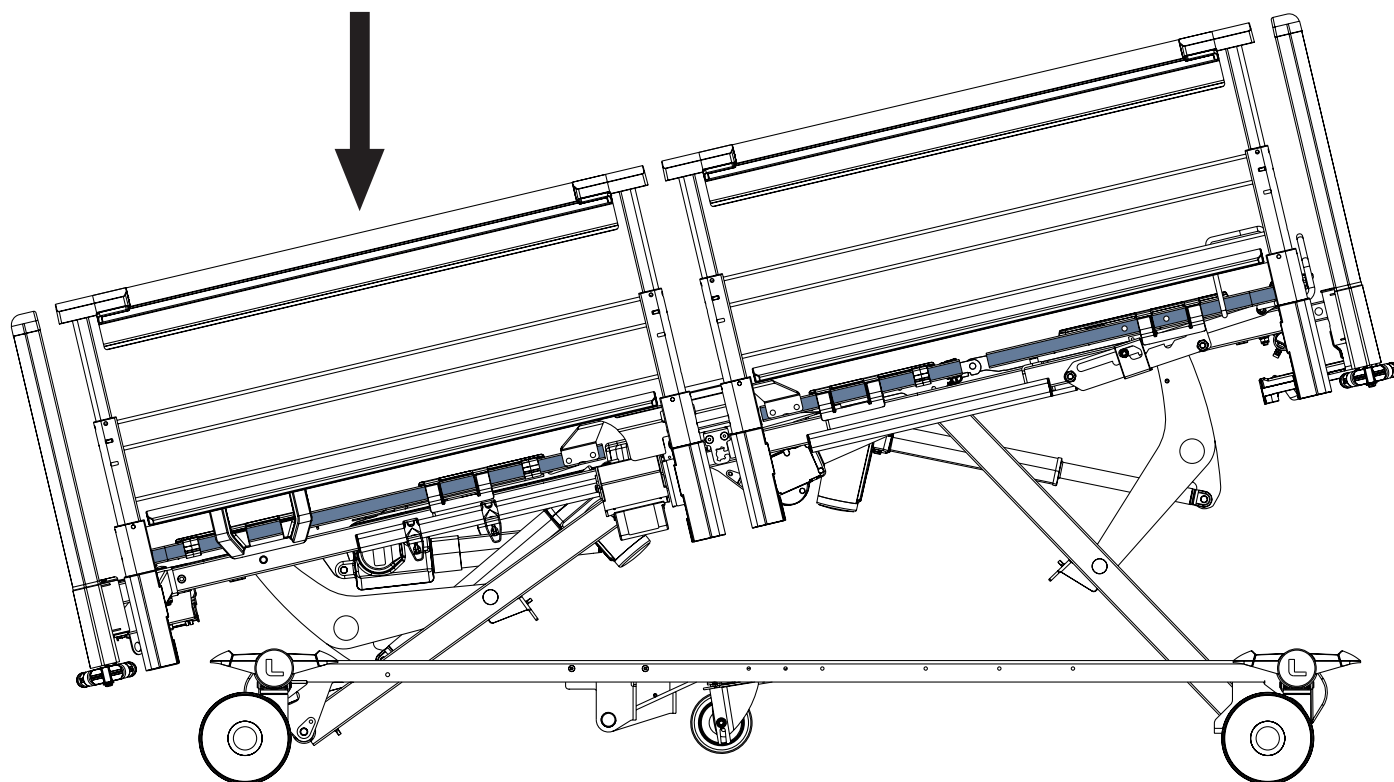

WARNING!

Risk of injury due to improper use of Trendelenburg Position!

- Hospital staff is responsible for assessing if the physical and psychological state of the patient is in accordance with use of the Trendelenburg Position.


It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.

To adjust Trendelenburg Tilt, use:


- Attendant Control Panel for Essenza 300 W with Lateral Tilt
- Attendant Control Panel for Essenza 300 W



Button intended for Trendelenburg Tilt adjustment:

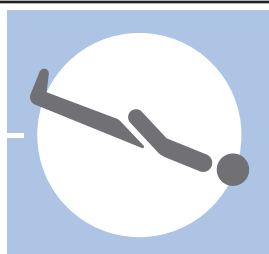

13.8 Reverse-Trendelenburg Tilt

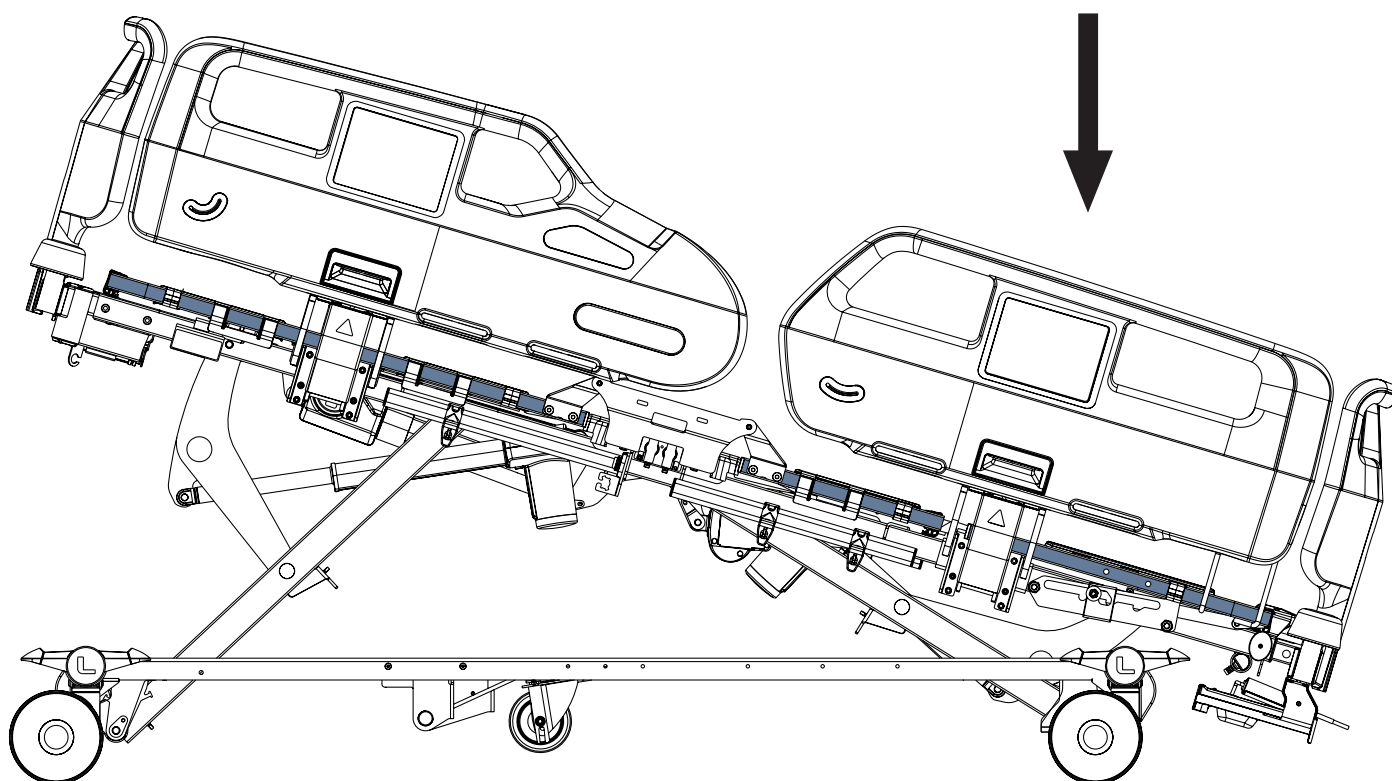



It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.


To adjust Reverse-Trendelenburg Tilt, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W



Buttons intended for Reverse-Trendelenburg Tilt adjustment:

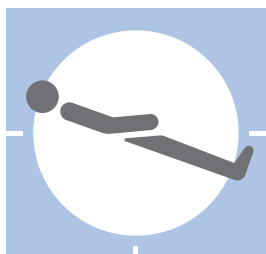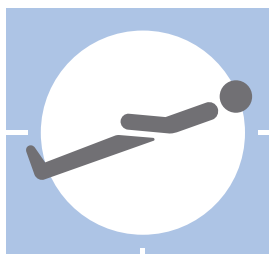

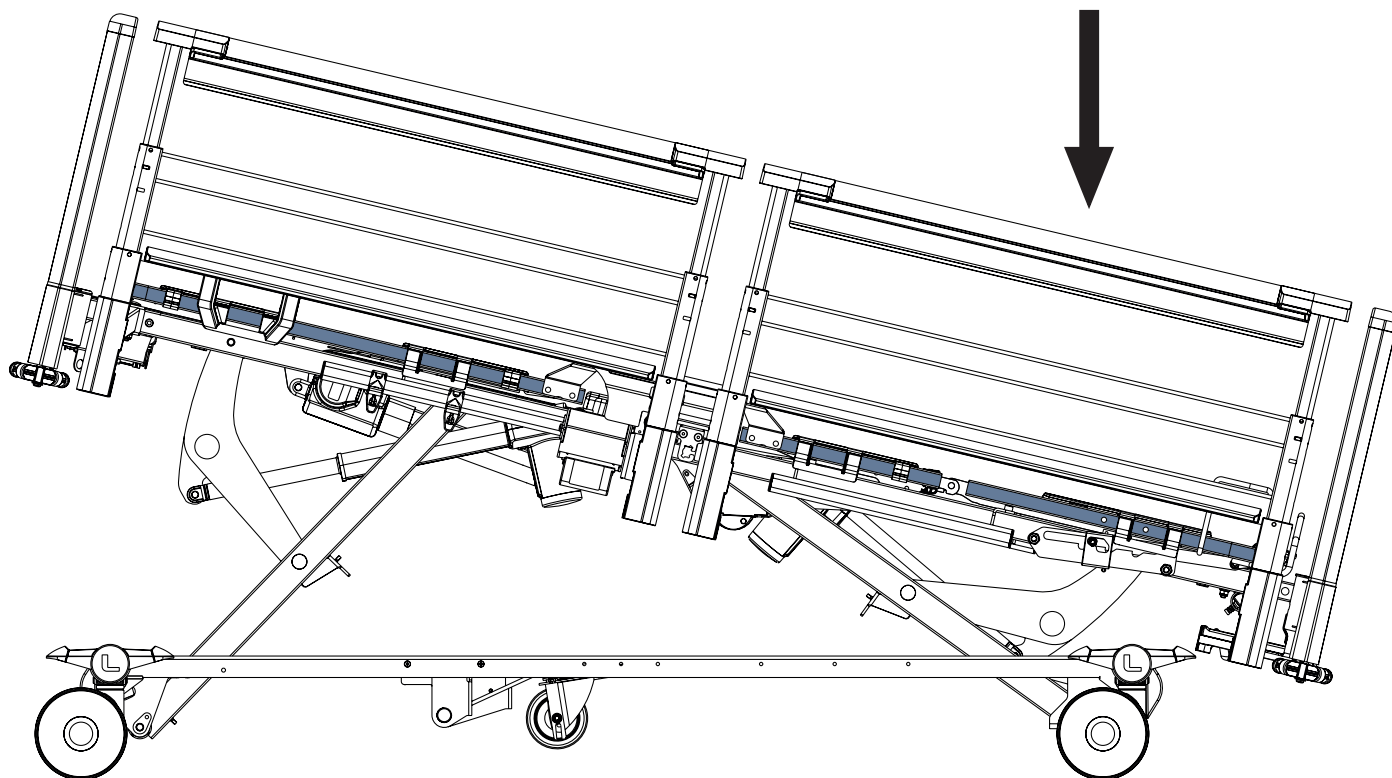


It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.

To adjust Reverse-Trendelenburg Tilt, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W



Button intended for Reverse-Trendelenburg Tilt adjustment:

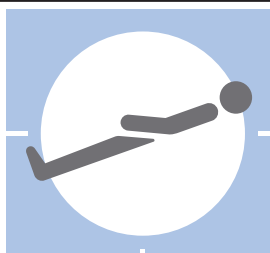

13.9 Bed Extension Adjustment

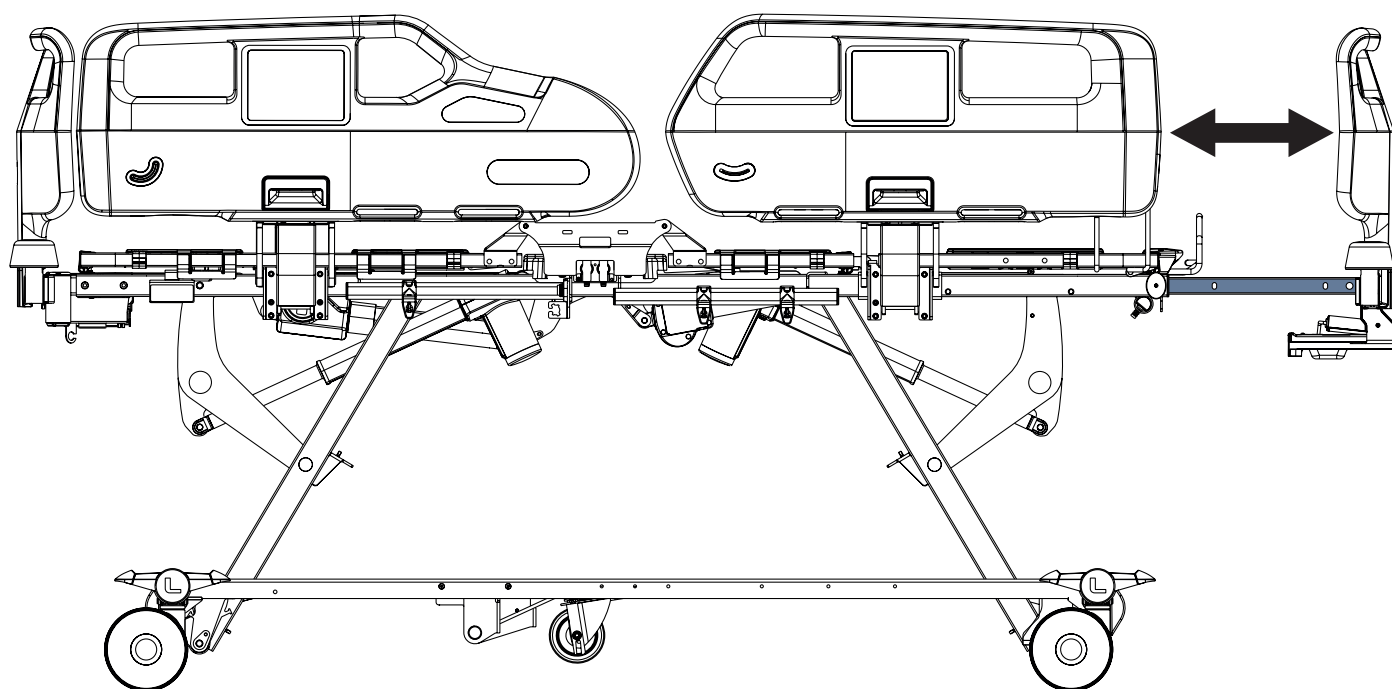

Simple Bed Extension

To extend the bed:

- ▶ Pull out the Safety Catches at Foot End on the both sides of the frame.
- ▶ Turn the Safety Catches by 90°.

The Safety Catches are released.

- ▶ Pull the Foot Board to the selected position.
- ▶ Turn the Safety Catches by 90° and push them in to fix the extension position.

To shorten the bed:

- ▶ Pull out the Safety Catches at Foot End on both sides of the frame.
- ▶ Turn the Safety Catches by 90°.

The Safety Catches are released.

- ▶ Push the Foot Board to the selected position.
- ▶ Turn the Safety Catches by 90° and push them in to fix the extension position.


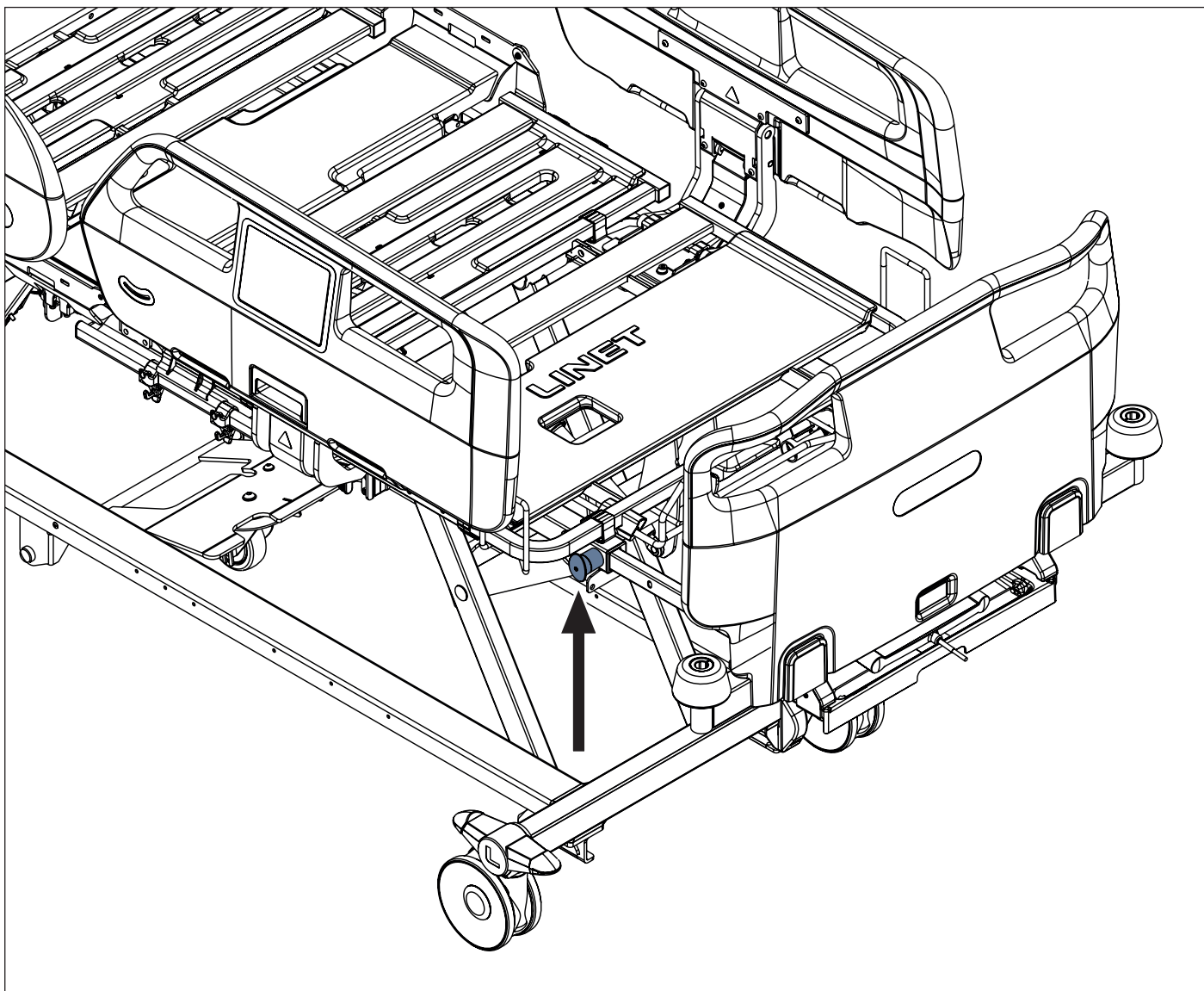


Fig. Position of Simple Bed Extension Safety Catch on the bed with plastic siderails

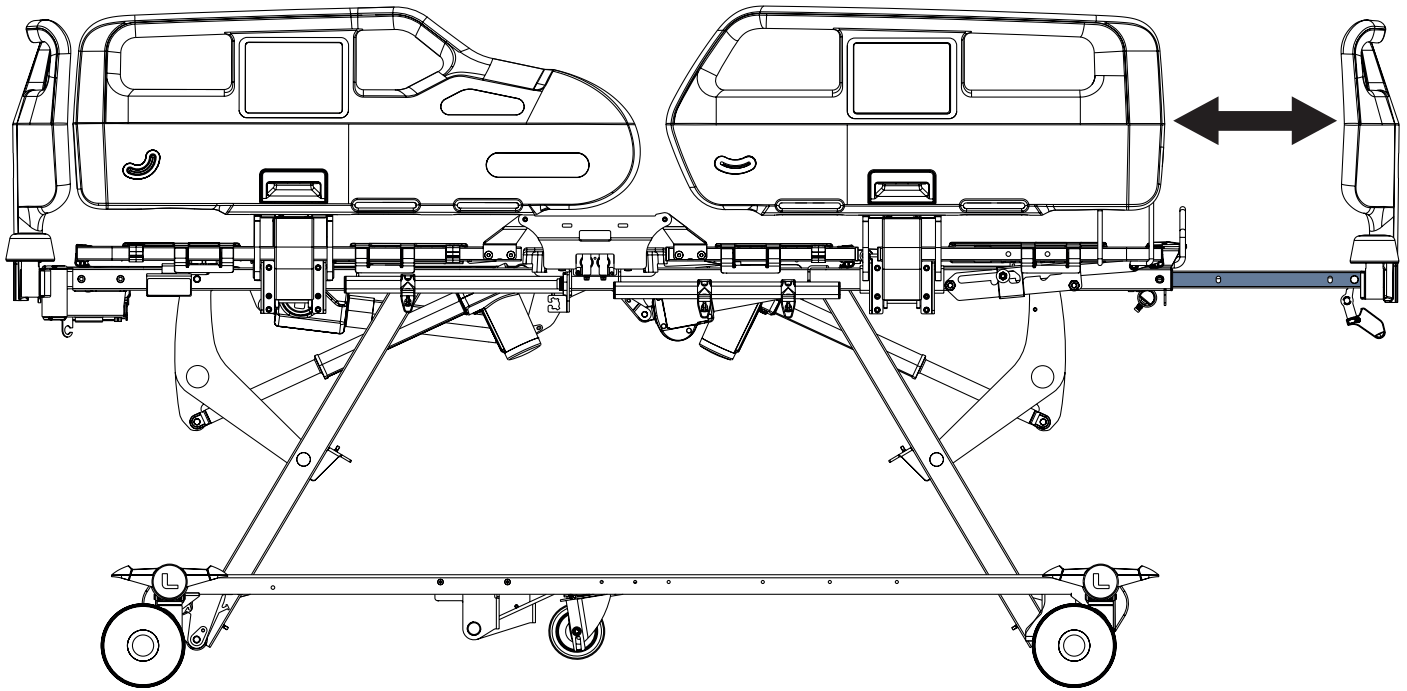

Advanced Bed Extension

To extend the bed:

- ▶ Pull Bed Extension Handle to you and keep it in this position.
- ▶ Extend the bed by pulling the foot board to you.
- ▶ Push the foot board to fix it in the required extended position.
- ▶ Move the foot board inwards and outwards to check if it is fixed.

To shorten the bed:

- ▶ Pull Bed Extension Handle to you and keep it in this position.
- ▶ Shorten the bed by pushing the foot board backwards to the required position.
- ▶ Move the foot board inwards and outwards to check if it is fixed.


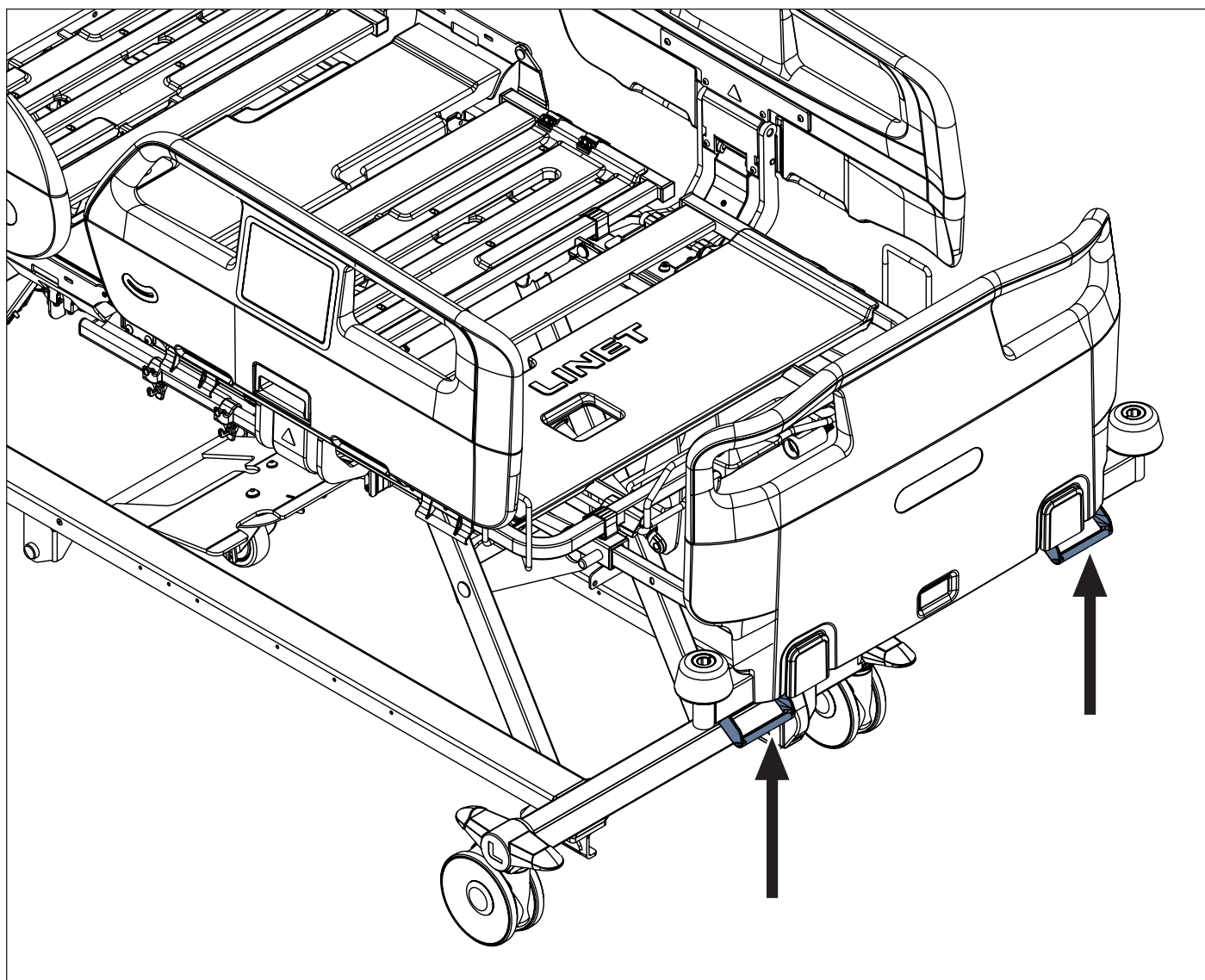


Fig. Position of Advanced Bed Extension Handles under the Foot Board on the bed with plastic siderails

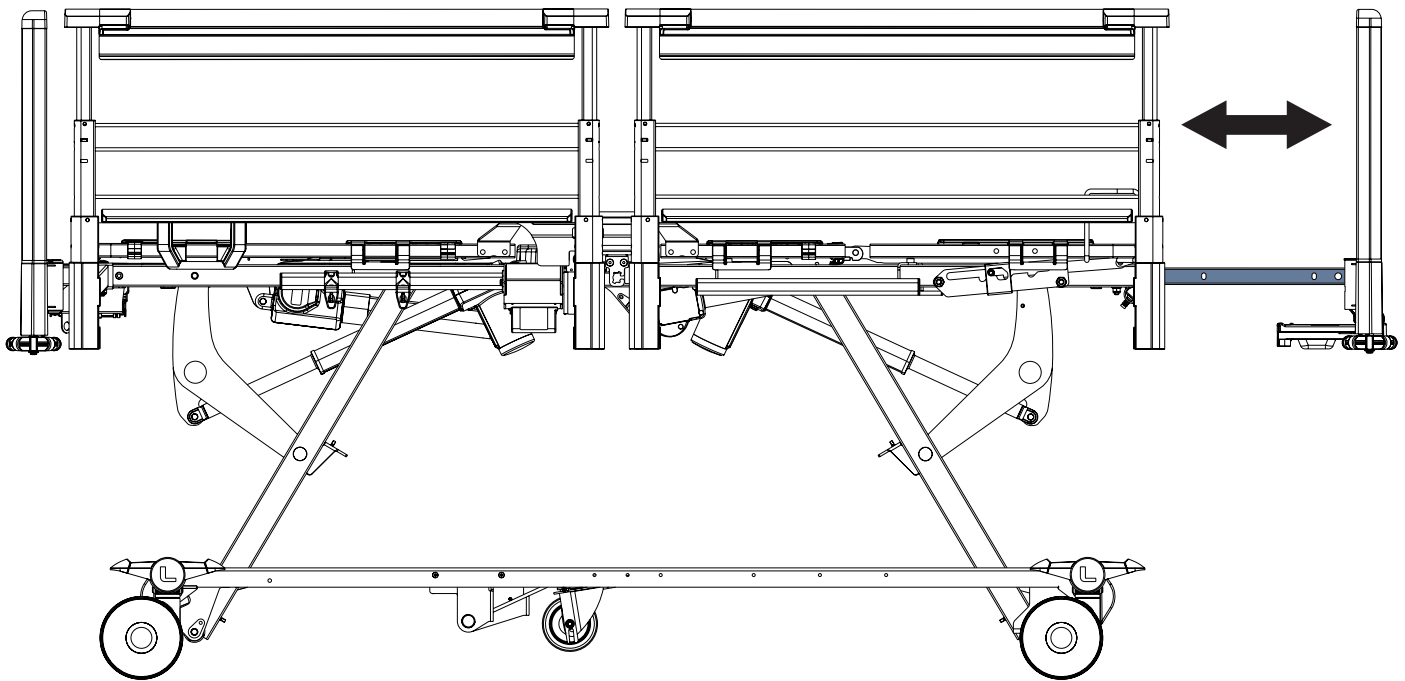

Simple Bed Extension

To extend the bed:

- ▶ Pull out the Safety Catches at Foot End on the both sides of the frame.
- ▶ Turn the Safety Catches by 90°.

The Safety Catches are released.

- ▶ Pull the Foot Board to the selected position.
- ▶ Turn the Safety Catches by 90° and push them in to fix the extension position.

To shorten the bed:

- ▶ Pull out the Safety Catches at Foot End on both sides of the frame.
- ▶ Turn the Safety Catches by 90°.

The Safety Catches are released.

- ▶ Push the Foot Board to the selected position.
- ▶ Turn the Safety Catches by 90° and push them in to fix the extension position.


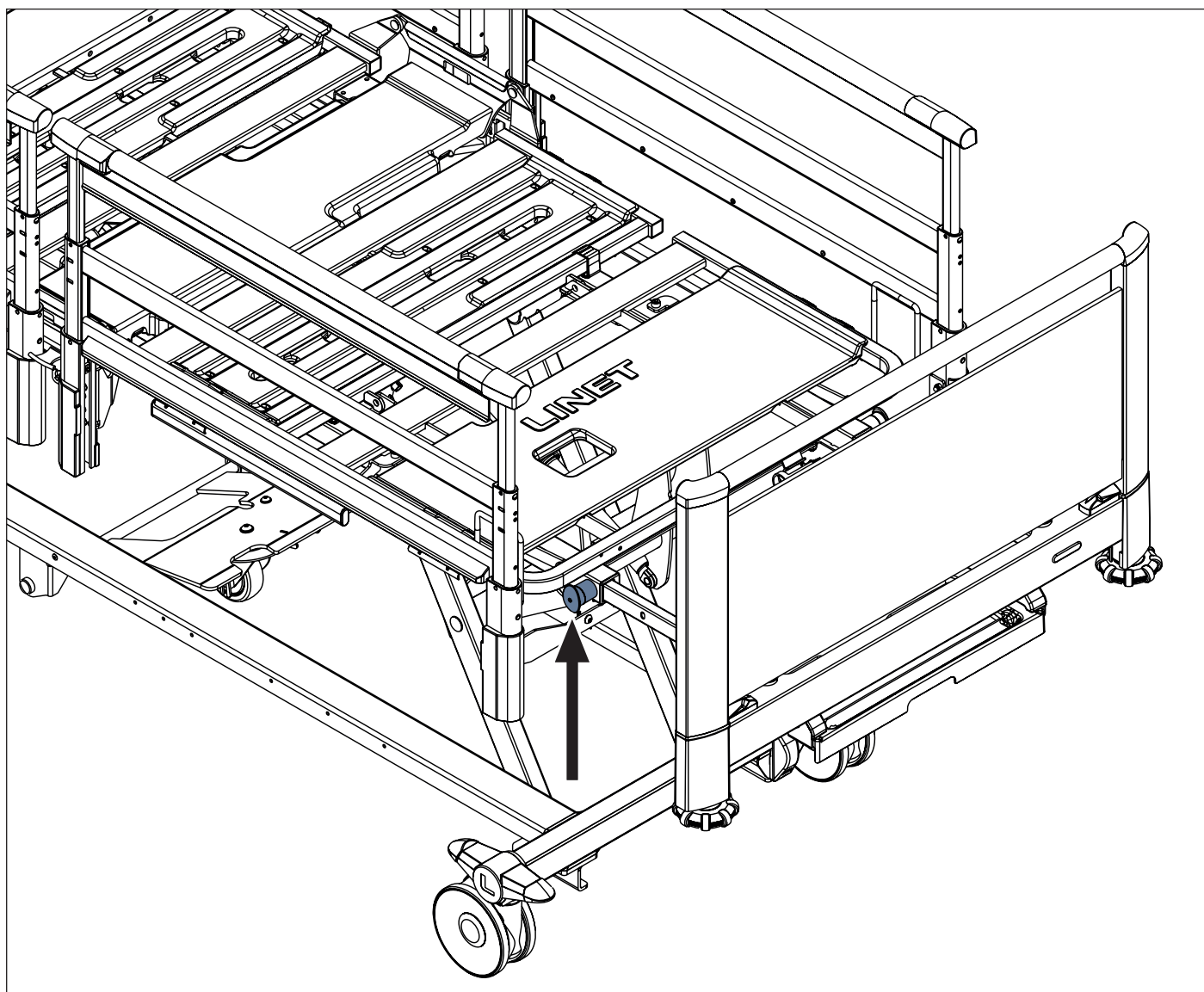


Fig. Position of Simple Bed Extension Safety Catch on the bed with telescopic siderails

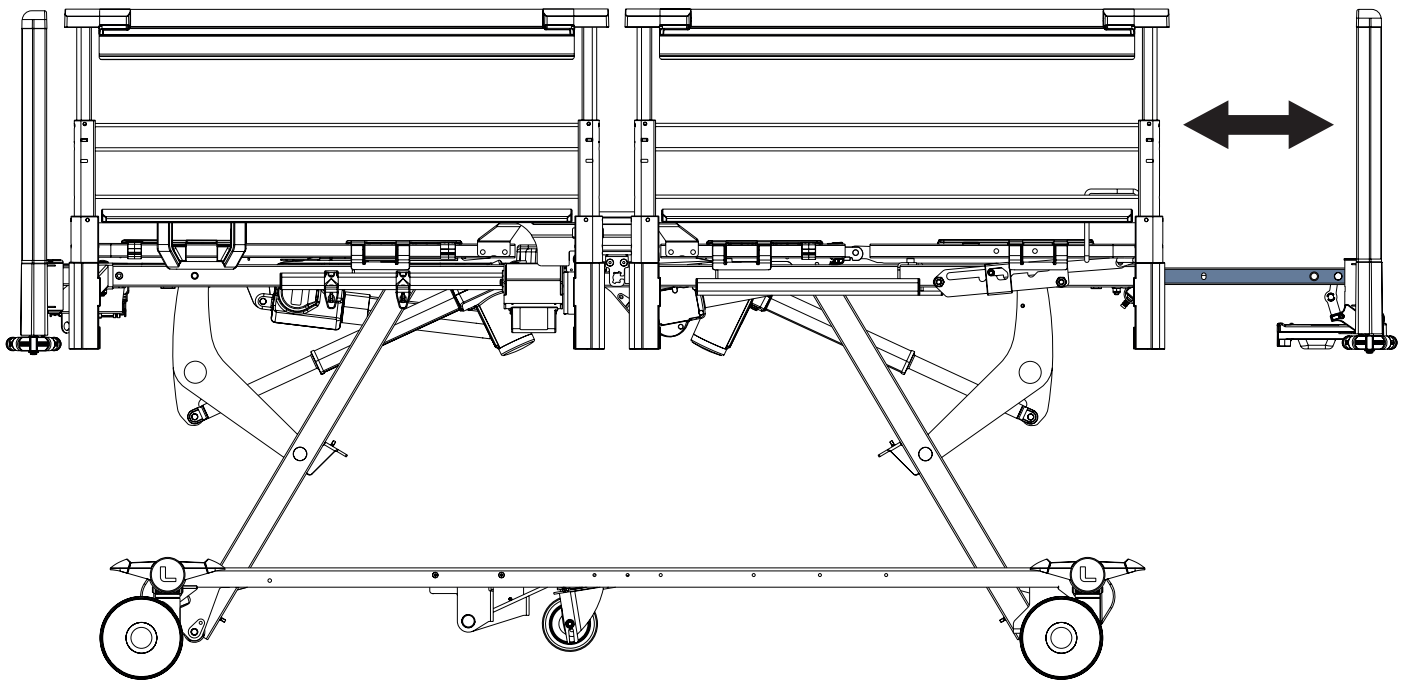

Advanced Bed Extension

To extend the bed:

- ▶ Pull Bed Extension Handle to you and keep it in this position.
- ▶ Extend the bed by pulling the foot board to you.
- ▶ Push the foot board to fix it in the required extended position.
- ▶ Move the foot board inwards and outwards to check if it is fixed.

To shorten the bed:

- ▶ Pull Bed Extension Handle to you and keep it in this position.
- ▶ Shorten the bed by pushing the foot board backwards to the required position.
- ▶ Move the foot board inwards and outwards to check if it is fixed.

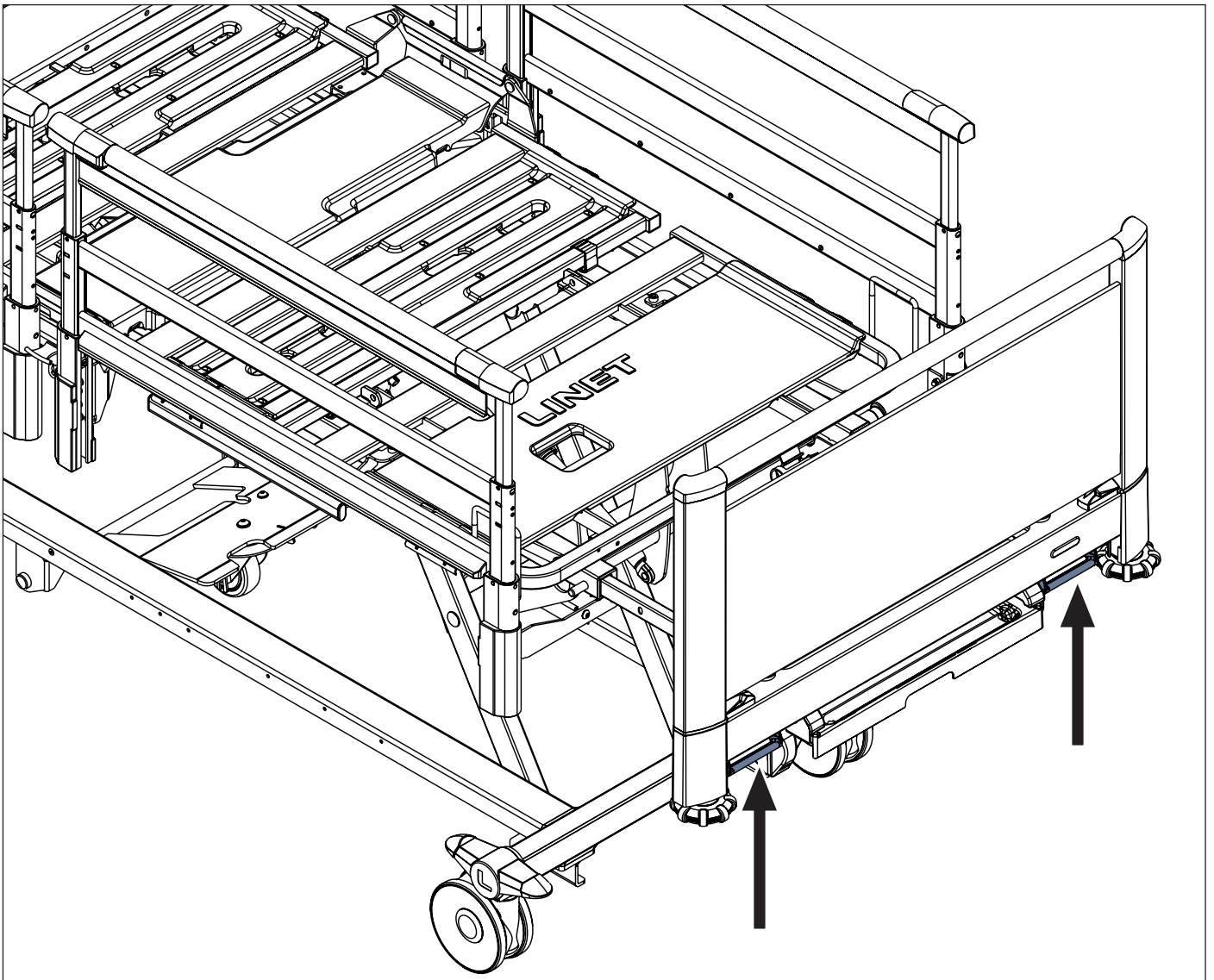


Fig. Position of Advanced Bed Extension Handles under the Foot Board on the bed with telescopic siderails

13.10 Calfrest Extension Adjustment


CAUTION!
Risk of material damage due to incorrect use!

- ▶ Avoid overloading the wire construction of Calfrest Extension and avoid shaking the wire construction of Calfrest Extension up and down violently when pulling Calfrest Extension Holders out and pushing Calfrest Extension Holders back to the original position.
- ▶ Avoid overloading the wire construction of Calfrest Extension and avoid shaking the wire construction of Calfrest Extension up and down violently when wire construction of Calfrest Extension is in extended position.

It is necessary to adjust Bed Extension before adjusting the Calfrest Extension. Calfrest Extension allows mattress support platform to be extended by 20 cm for use of the 220 cm long mattress on the Essenza 300 W bed. It is allowed to use also 200 cm long mattress in combination with 20 cm long mattress foam extension. Safe Working Load of the wire construction of Calfrest Extension in extended position is 75 kg.

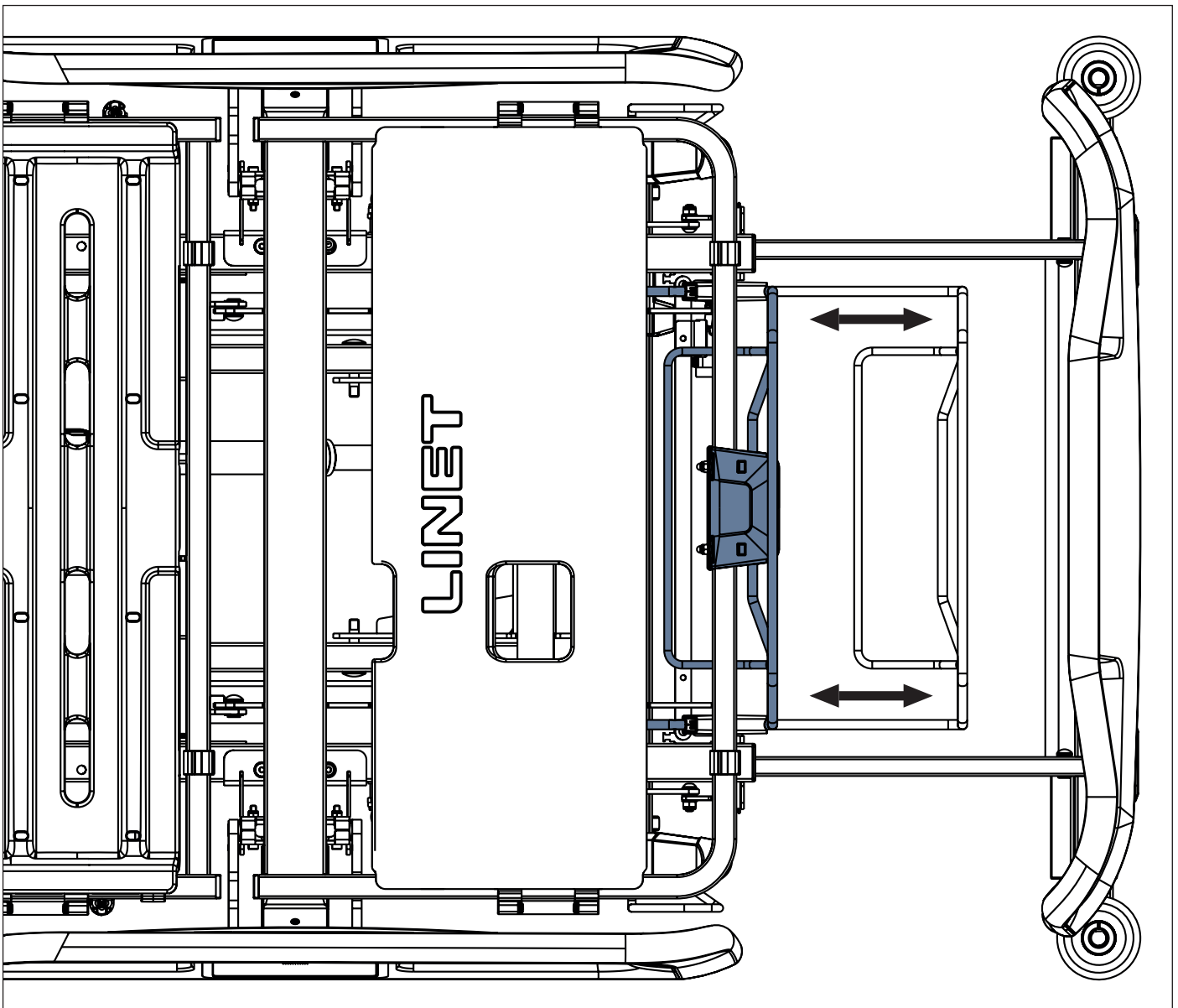


Fig. Adjustment of Calfrest Extension Holder

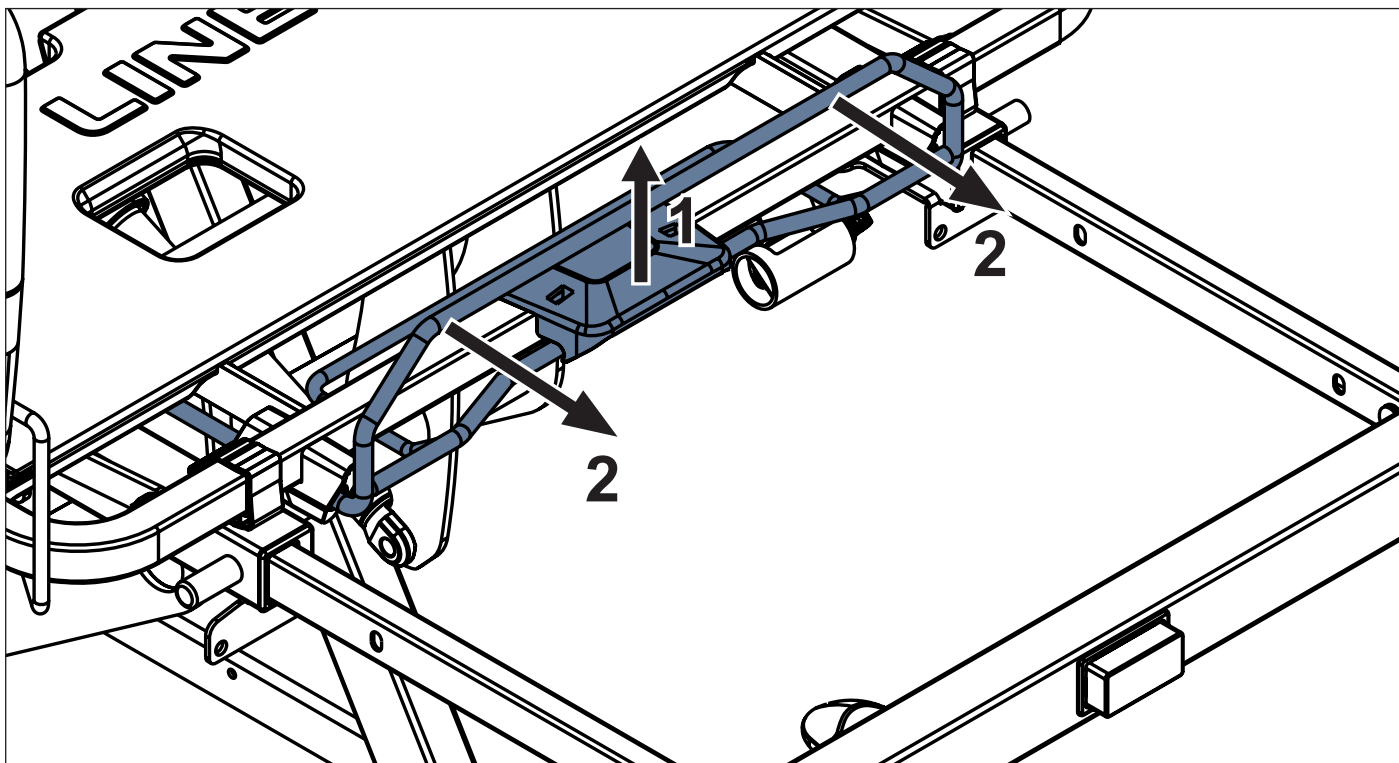


Fig. Two steps of the Calfrest Extension Adjustment

13.10.1 Calfrest Extension Adjustment

To adjust Calfrest Extension:


- ▶ Adjust Bed Extension to the second position.
- ▶ Unlock the Calfrest Extension Holder by lifting the securing element. (1)
- ▶ Pull Calfrest Extension Holder to a required position. (2)

Mattress support platform is prepared to be used with 220 cm long mattress or with 200 cm long mattress in combination with 20 cm long mattress extension.

To prepare Essenza 300 W bed for use with 200 cm long mattress again:

- ▶ Push Calfrest Extension Holder back to its original position.
- ▶ Adjust Bed Extension as required.

13.11 Emergency CPR Position

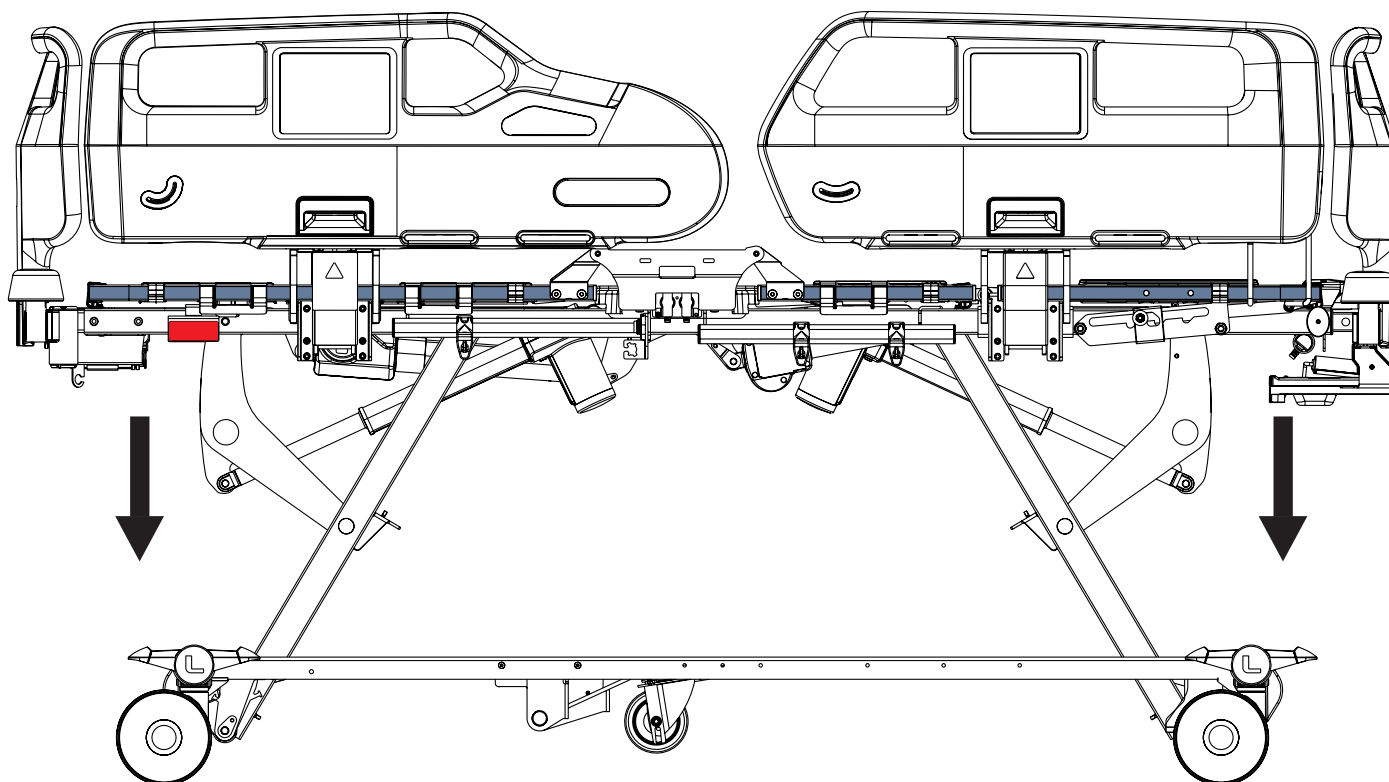

Essential performance of the Essenza 300 W medical bed consists in manual Backrest release by CPR lever.

To adjust Emergency CPR Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W

Buttons intended for Emergency CPR Position adjustment:



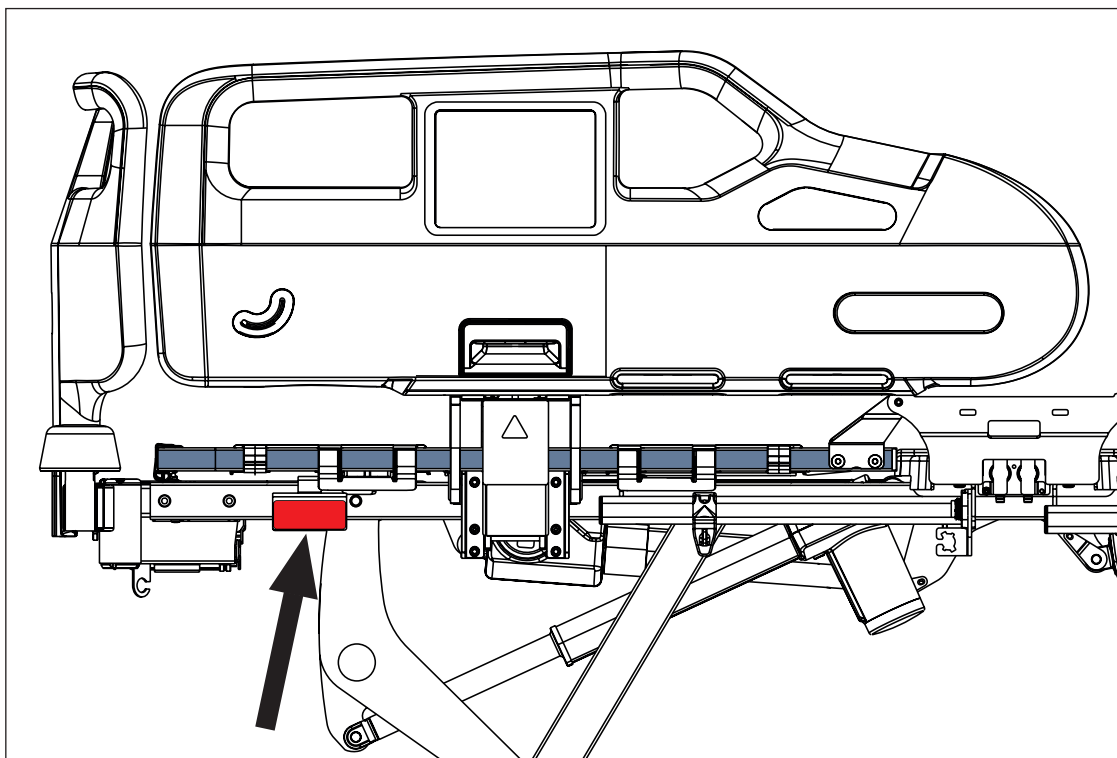


Fig. CPR lever under the Backrest of the bed with Plastic siderails



WARNING!

Risk of injury due to incorrect use!

- ▶ Ensure the CPR lever is accessible when necessary!



WARNING!

Risk of injury due to lowering the backrest too quickly!

- ▶ Ensure that there are no body parts between any movable parts of the bed.
- ▶ Push the Backrest down by the Backrest handle only.


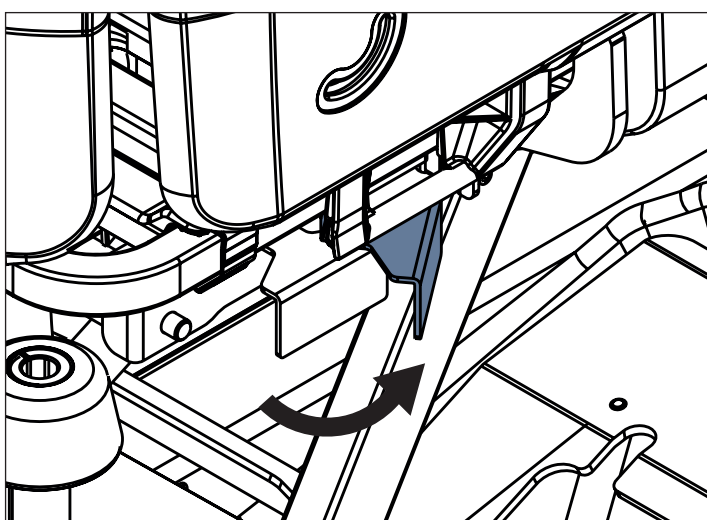


Fig. Manipulation with CPR lever in version with plastic siderails

To adjust Emergency CPR Position manually:

- ▶ Pull and hold CPR lever.
- ▶ Push Backrest down.

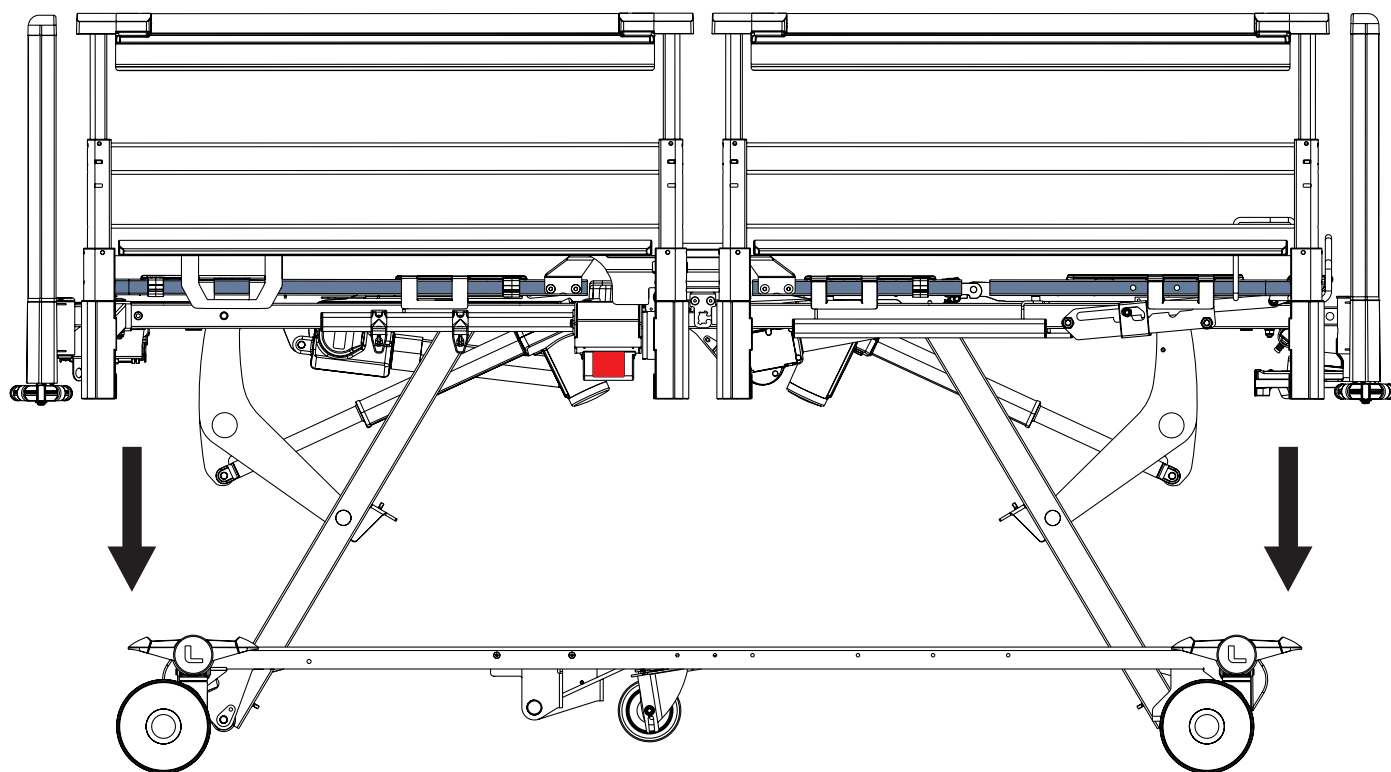

Essential performance of the Essenza 300 W medical bed consists in manual Backrest release by CPR lever.

To adjust Emergency CPR Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W

Buttons intended for Emergency CPR Position adjustment:



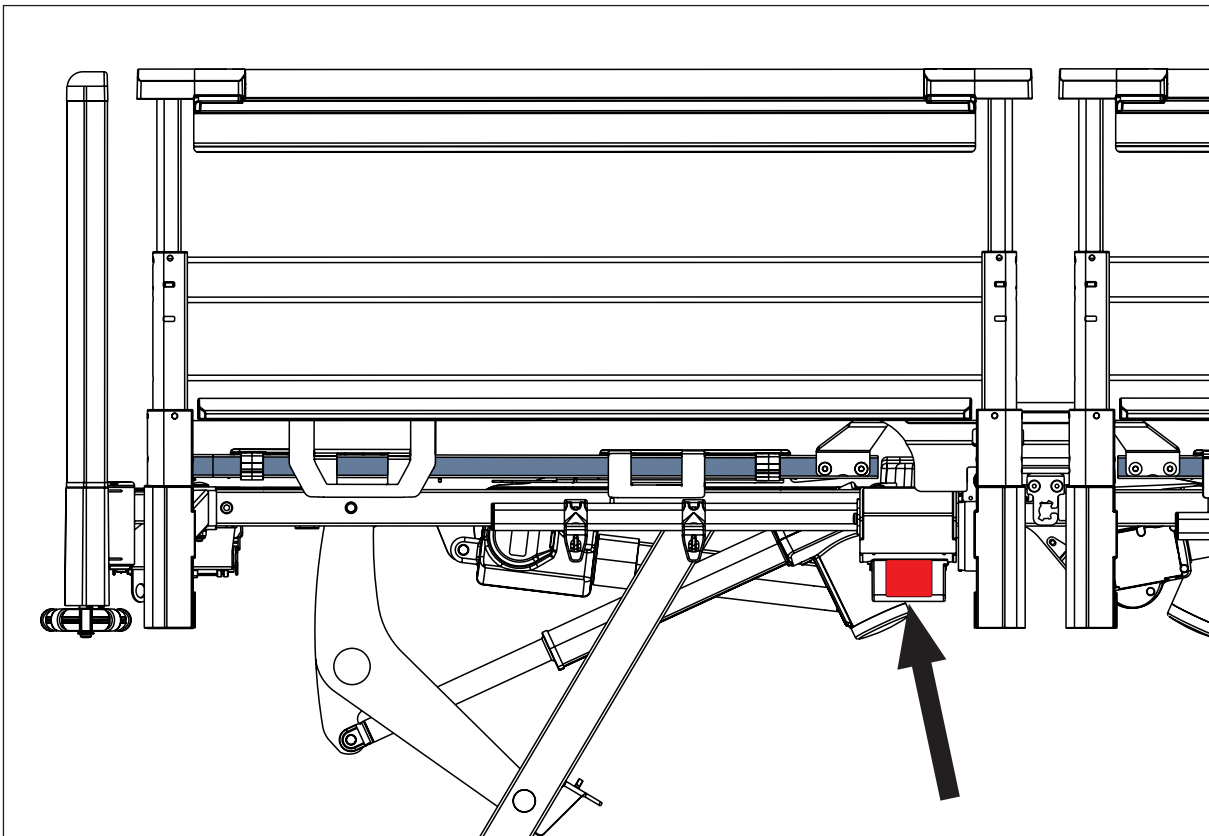


Fig. CPR lever under the Seat Section of the bed with Telescopic siderails



WARNING!

Risk of injury due to incorrect use!

- ▶ Ensure the CPR lever is accessible when necessary!



WARNING!

Risk of injury due to lowering the backrest too quickly!

- ▶ Ensure that the siderails are in the down position.
- ▶ Ensure that there are no body parts between any movable parts of the bed.
- ▶ Push the Backrest down by the mattress guard handle only.


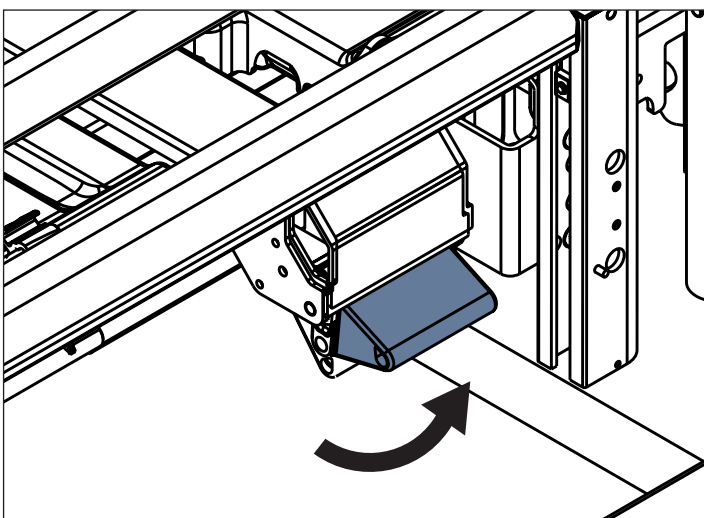


Fig. Manipulation with CPR lever in version with telescopic siderails

To adjust Emergency CPR Position manually:

- ▶ Pull and hold CPR lever.
- ▶ Push Backrest down.

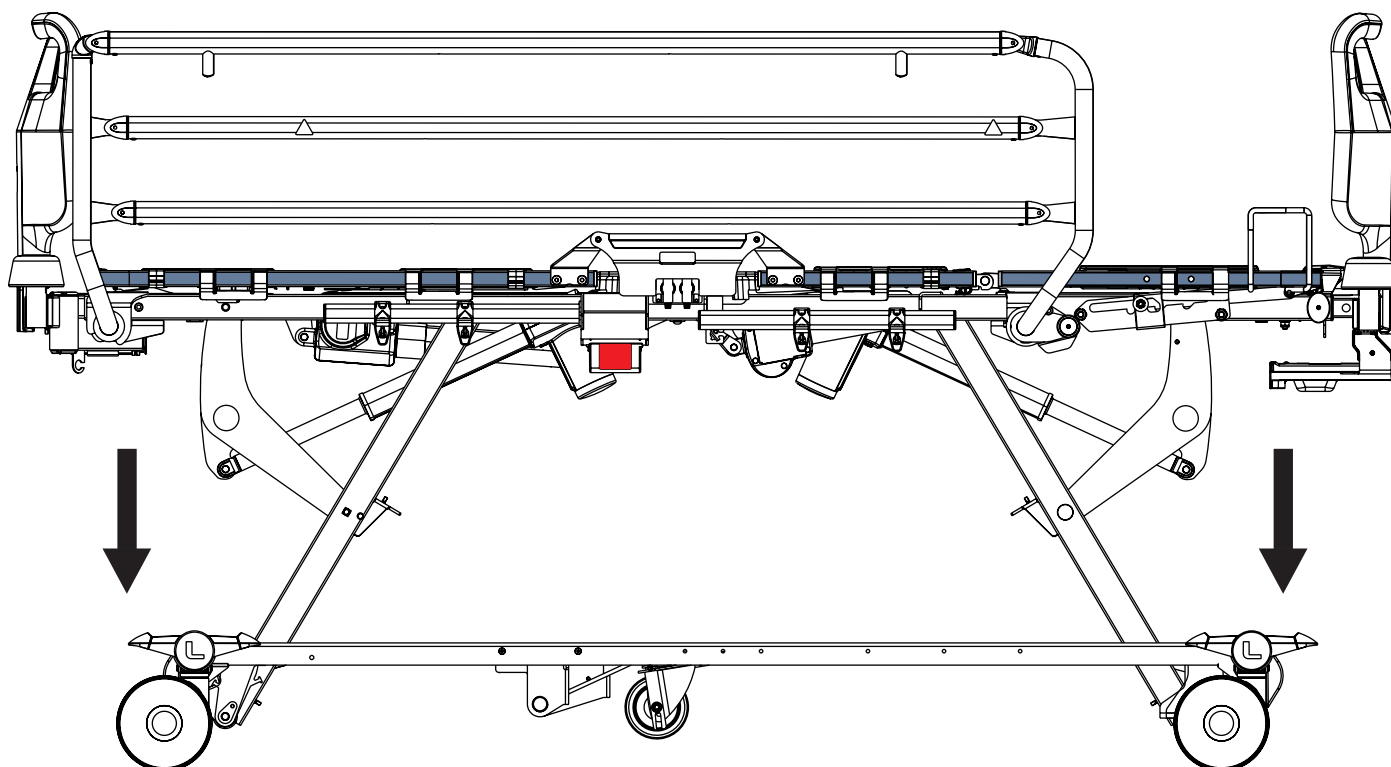

Essential performance of the Essenza 300 W medical bed consists in manual Backrest release by CPR lever.

To adjust Emergency CPR Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W

Buttons intended for Emergency CPR Position adjustment:



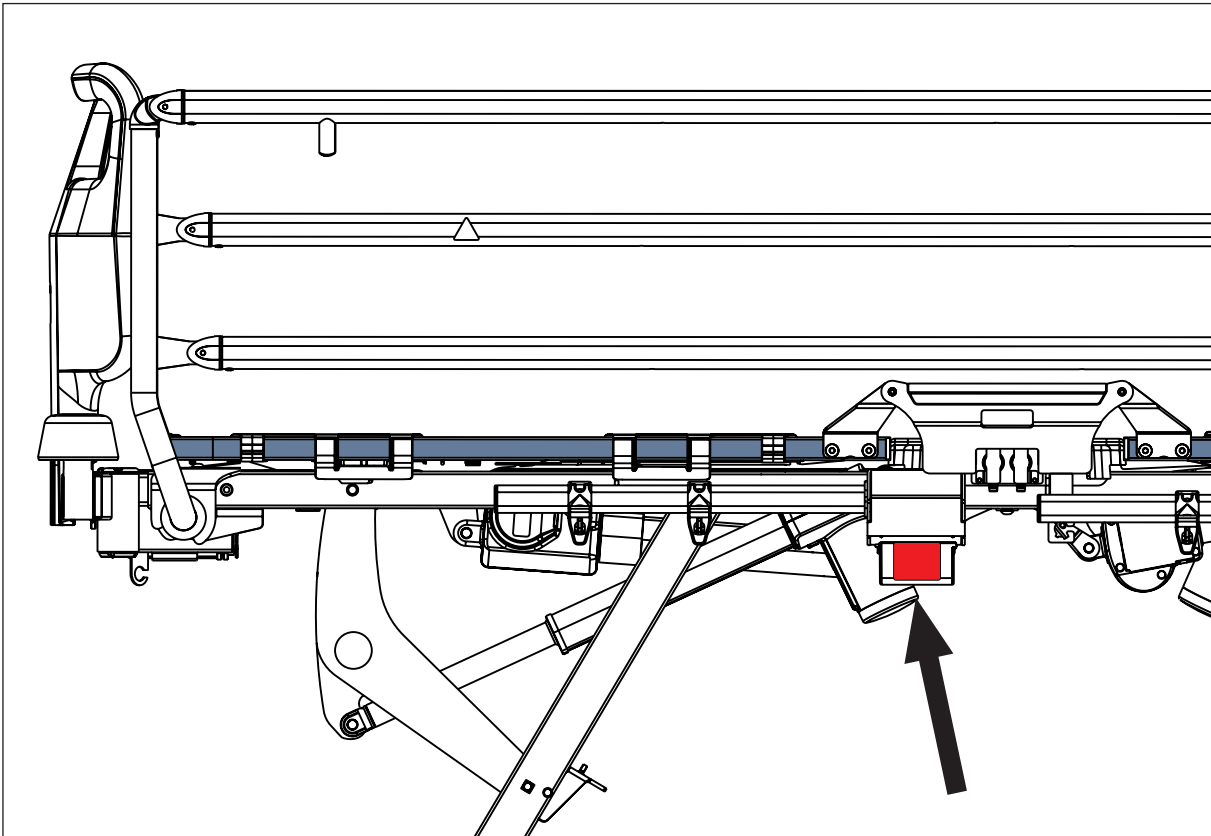


Fig. CPR lever under the seat section of the bed with collapsible siderails



WARNING!

Risk of injury due to incorrect use!

- ▶ Ensure the CPR lever is accessible when necessary!



WARNING!

Risk of injury due to lowering the backrest too quickly!

- ▶ Ensure that the siderails are in the down position.
- ▶ Ensure that there are no body parts between any movable parts of the bed.
- ▶ Push the backrest down by the mattress guard handle only.


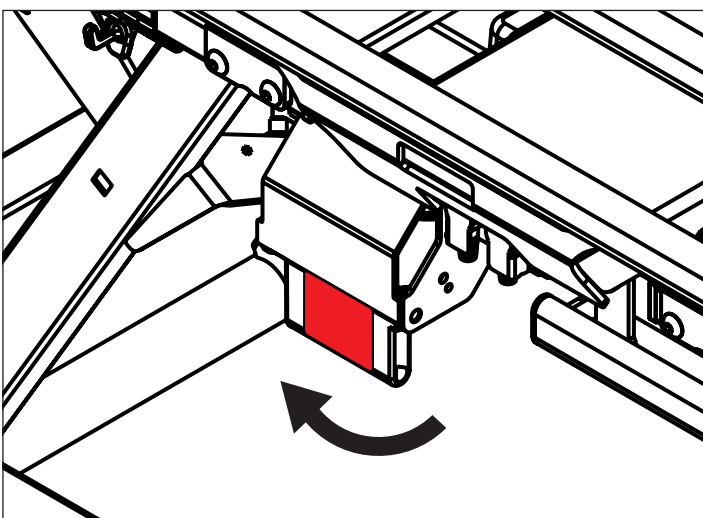


Fig. Manipulation with CPR lever in version with collapsible siderails

To adjust Emergency CPR Position manually:

- ▶ Pull and hold CPR lever.
- ▶ Push backrest down.

13.12 Cardiac Chair Position

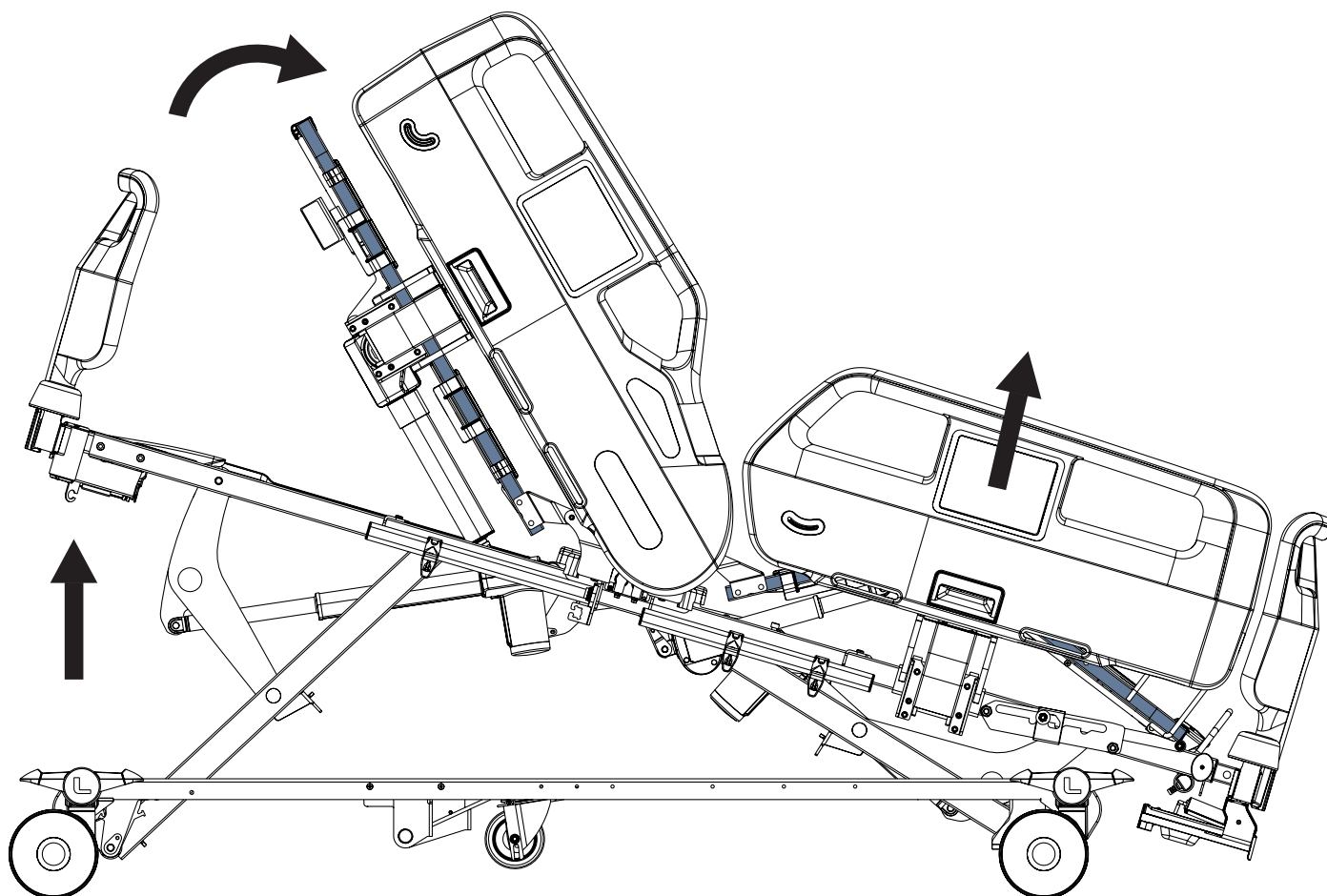


To adjust Cardiac Chair Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W



Buttons intended for Cardiac Chair Position adjustment:



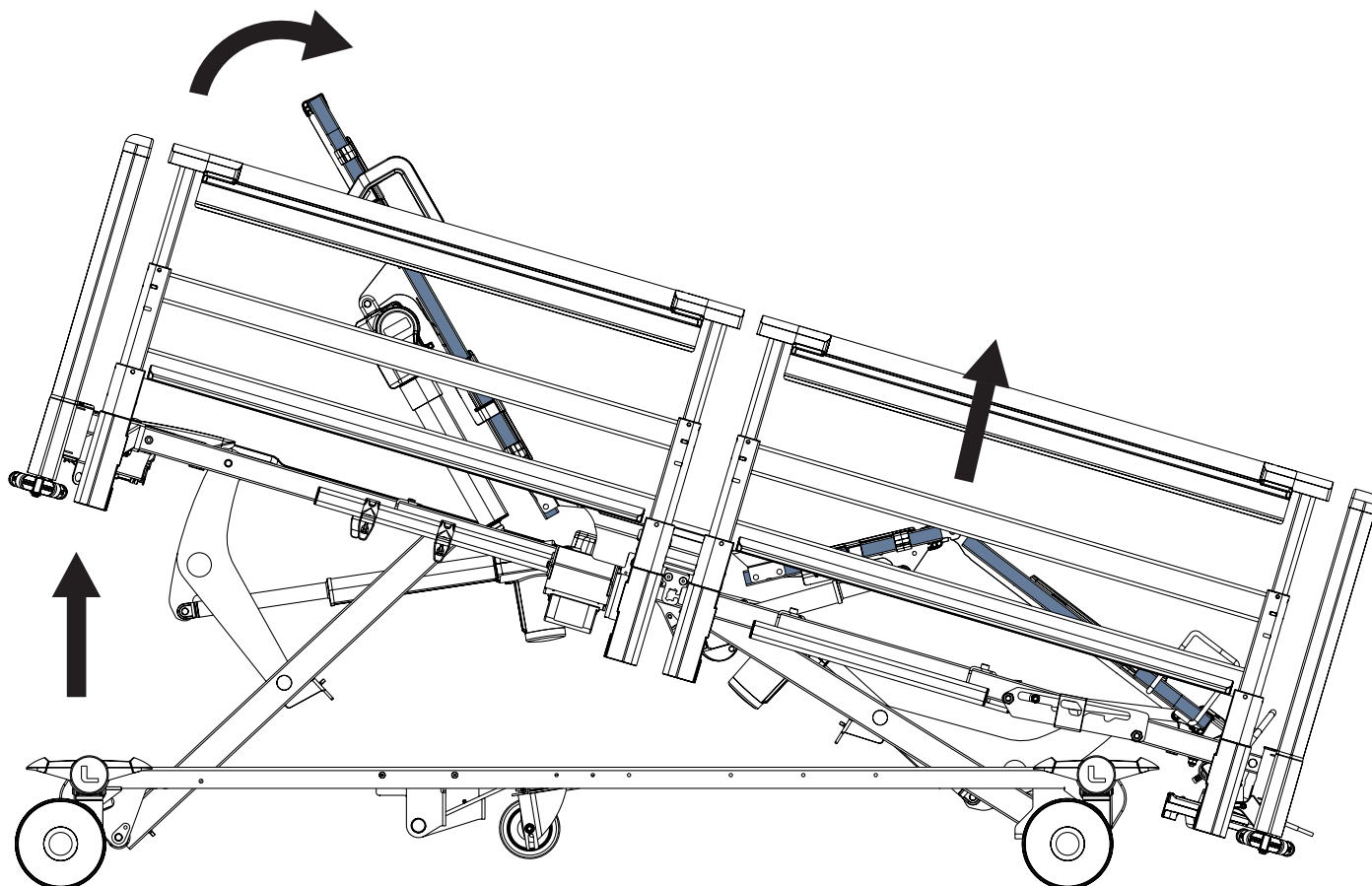

To adjust Cardiac Chair Position, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W



Button intended for Cardiac Chair
Position adjustment:



13.13 Mobilisation Position

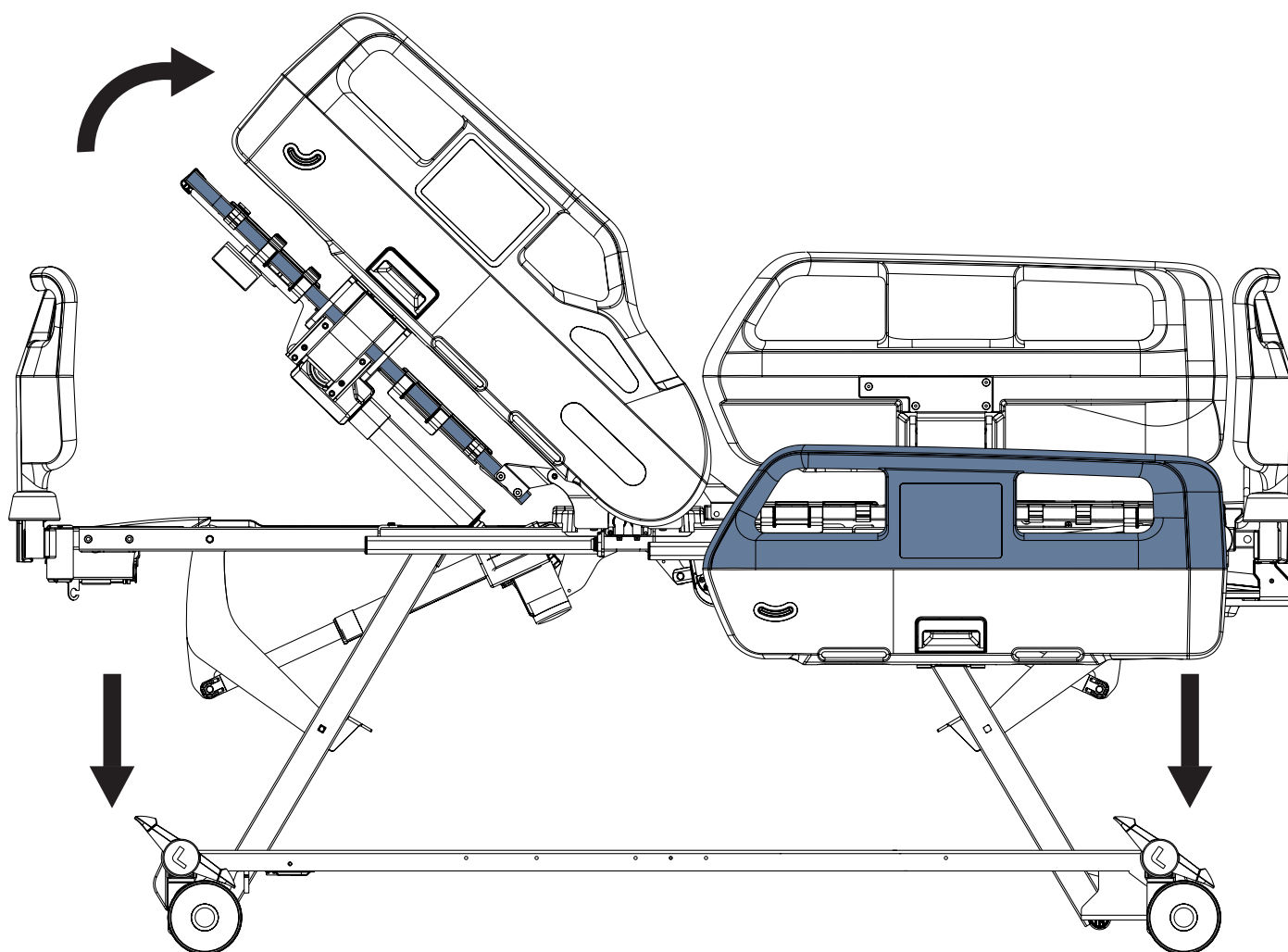


To adjust Mobilisation Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W



Buttons intended for Mobilisation Position adjustment:



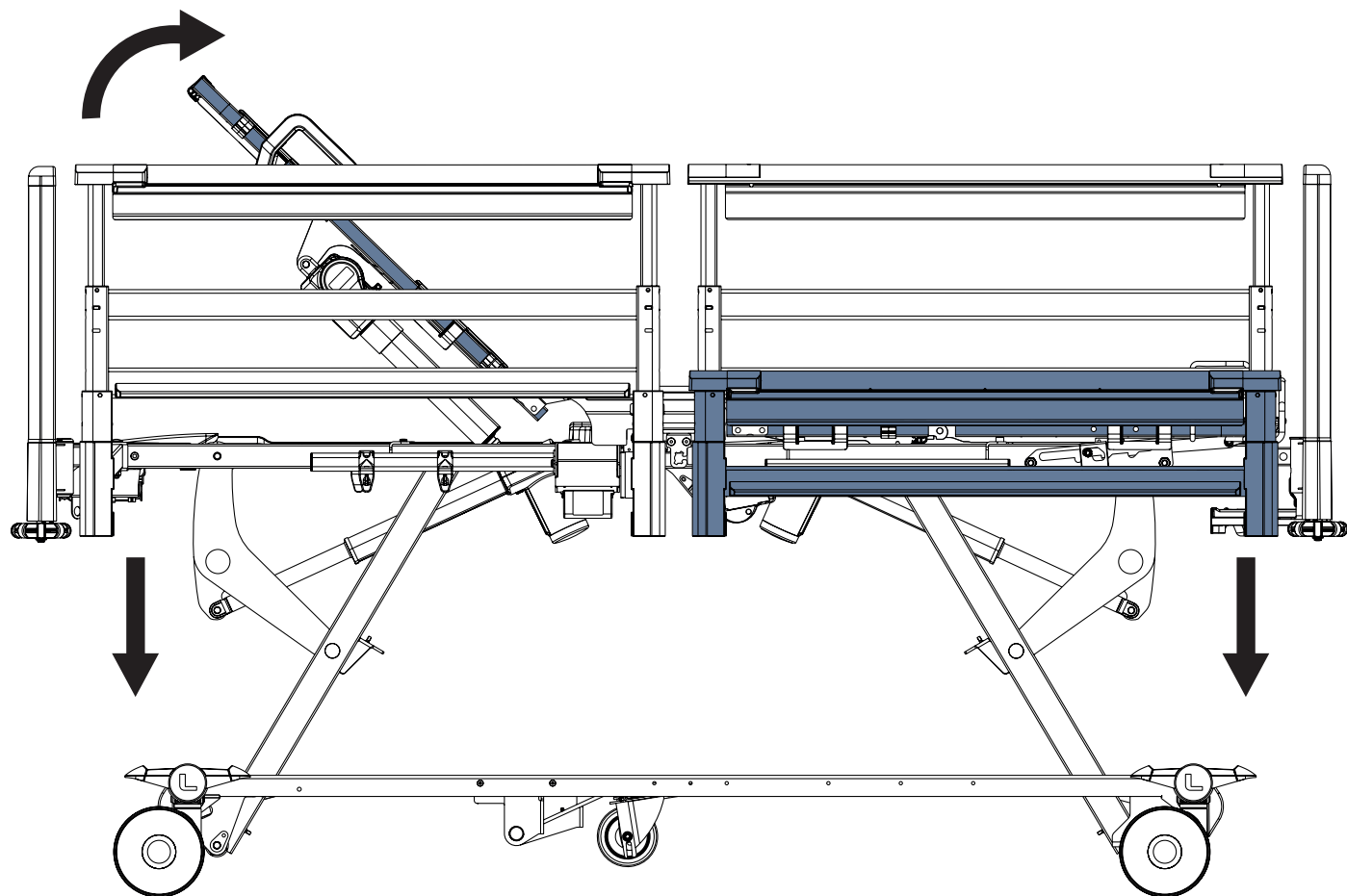

To adjust Mobilisation Position, use:


- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W



Button intended for Mobilisation
Position adjustment:



13.14 Lateral Tilt (only Essenza 300 W with Lateral Tilt)

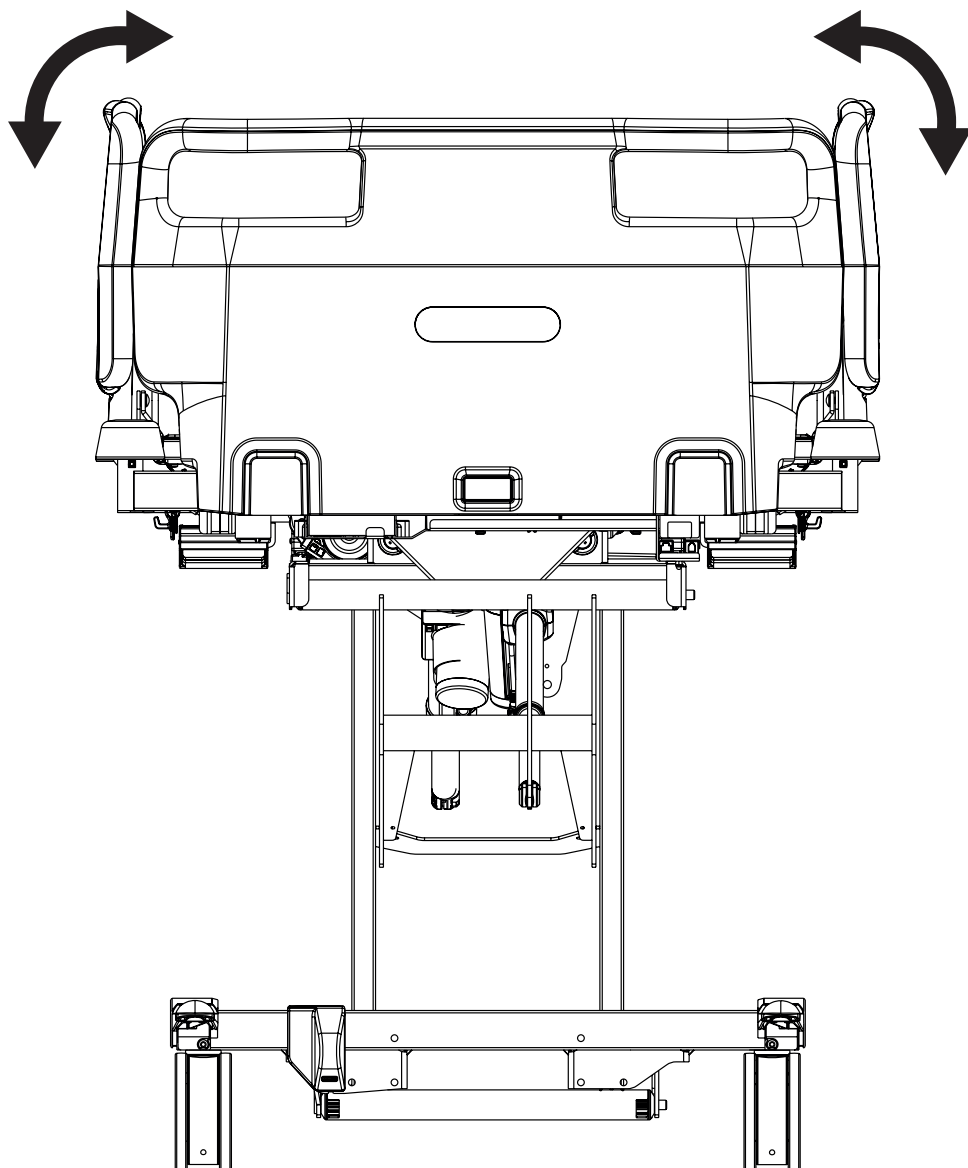



To adjust Lateral Tilt, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt



Buttons intended for Lateral Tilt adjustment:

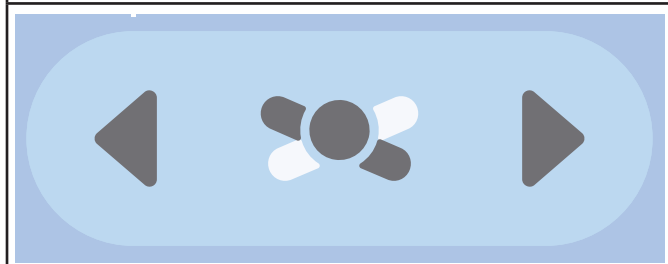

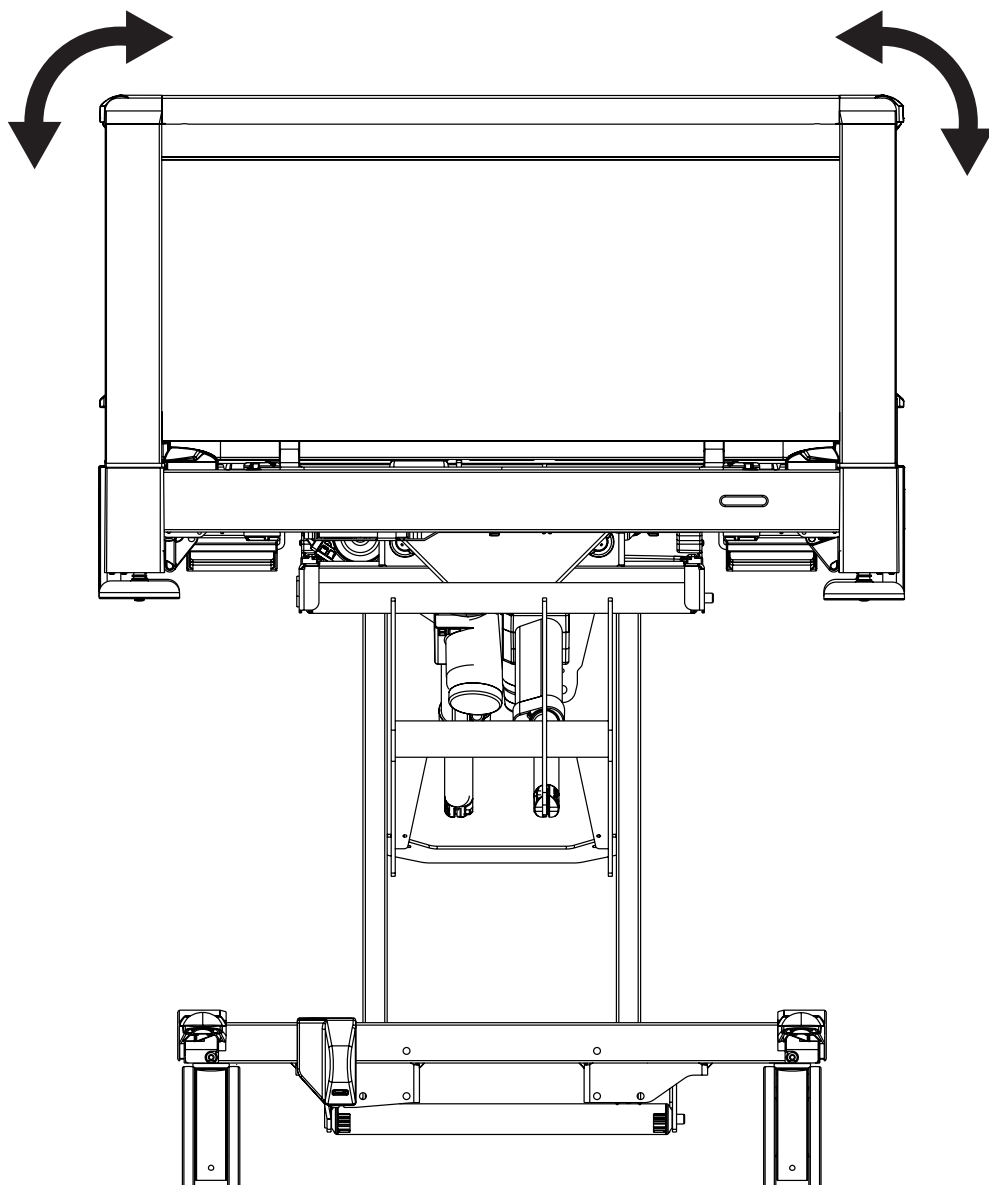


To adjust Lateral Tilt, use:

- Attendant Control Panel for Essenza 300 W with Lateral Tilt



Buttons intended for Lateral Tilt adjustment:

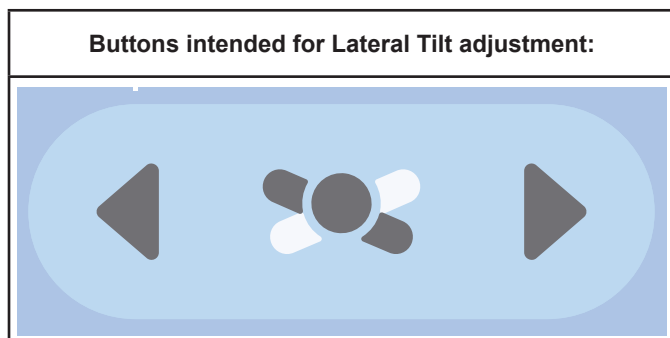

**WARNING!****Risk of injury due to incorrect use!**

- ▶ **Lateral Tilt is a positioning function intended for trained medical personnel only. The use of the Lateral Tilt, including its angle of adjustment, must be preceded by an evaluation of the patient's condition by the staff. This evaluation must take into account all risks resulting from the physical and mental state of the patient, especially the risk of fall, and loss of bed stability.**

It is possible to adjust either bed longitudinal tilt or bed lateral tilt outside the bed horizontal position, but not both at the same time in the case of the Essenza 300 W with Lateral Tilt.

The Lateral Tilt adjustment can only be performed by medical personnel who have been professionally trained by the manufacturer or by persons authorized by the manufacturer.

Lateral Tilt adjustments cannot be performed by persons without proper training and healthcare expertise, including the patient and visitors. Depending on the physical and mental state of the patient, medical personnel can lock the Lateral Tilt adjustment of the control elements.

When adjusting the Lateral Tilt, the siderails must be set in the upper position and bed must be braked. Trained personnel are responsible for braking the bed and reliably securing the siderails in the upper position.

When positioning or when turning the patient onto his side using a Lateral Tilt, proceed in such a way as not to reduce the lateral stability of the bed.

During preventive positioning the patient with the help of Lateral Tilt, place suitable positioning aids (pillows) between the siderails and parts of the patient's body that may come into contact with the siderails. In this way, it is possible to prevent pressure injuries on the patient's body parts and increase the patient's comfort.


If a patient is left unattended in a bed laterally tilted, staff must provide the following in advance:

- ▶ The siderails of the bed are reliably secured in the upper position.
- ▶ The position of the patient in the bed is stable. The patient is secured against accidental overturning with suitable positioning aids (see „17 Accessories“ on the page 149).

If the Lateral Tilt is used when getting the patient out of bed, or when sitting in bed with the siderails in the down position, medical personnel must ensure the following:

- ▶ The patient gets out of bed, or sits in bed only in the presence of staff.
- ▶ The staff ensures the patient's safety during getting out of bed, or sitting in bed to prevent the patient from falling. In addition to the staff's own assistance, suitable support aids can be used to ensure safety, e.g. siderails or support handles (Mobi-Lift).
- ▶ The height of the bed must be adapted to the height and physical and mental status of the patient.

13.15 Extra Low Position

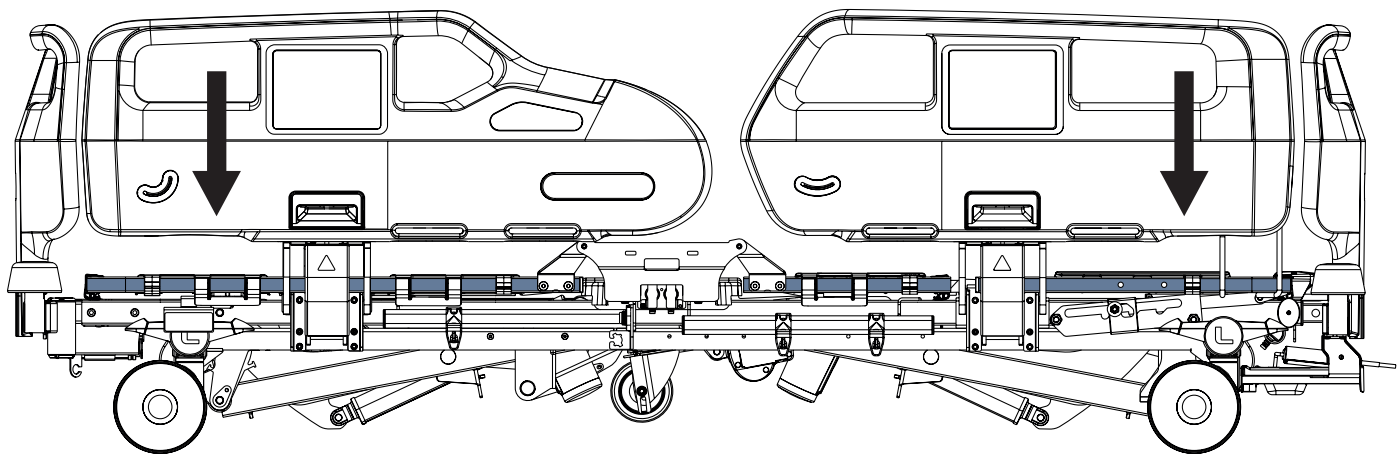


To adjust Extra Low Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W
- ▶ Handset Control Panel (only for Essenza 300 W with Lateral Tilt)
- ▶ Patient Control Panel (only for Essenza 300 W with Lateral Tilt)



Buttons intended for Extra Low Position adjustment:			

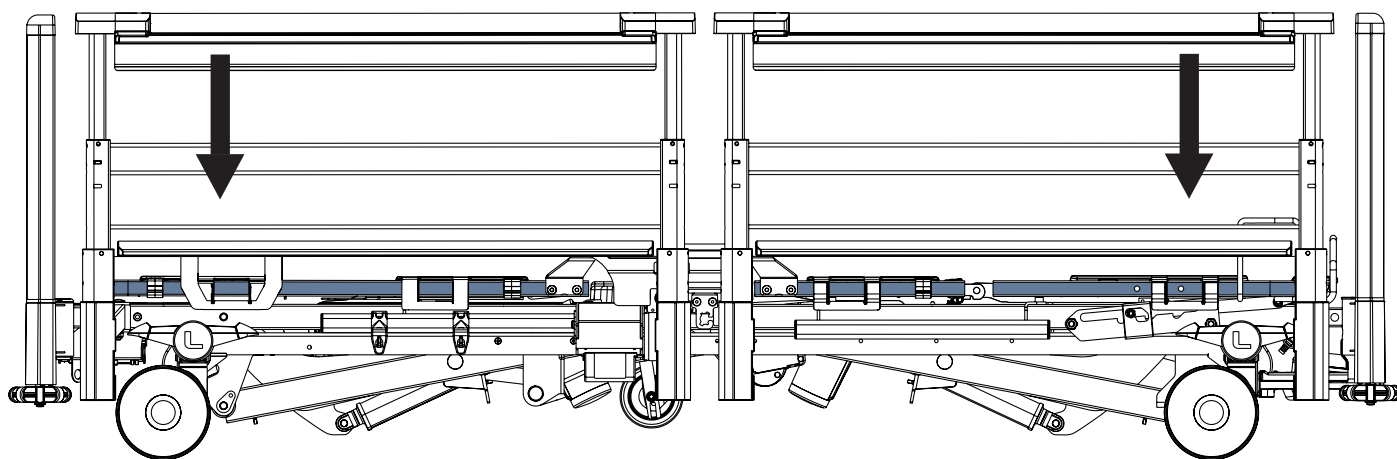


To adjust Extra Low Position, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Handset Control Panel (only for Essenza 300 W with Lateral Tilt)



Buttons intended for Extra Low Position adjustment:

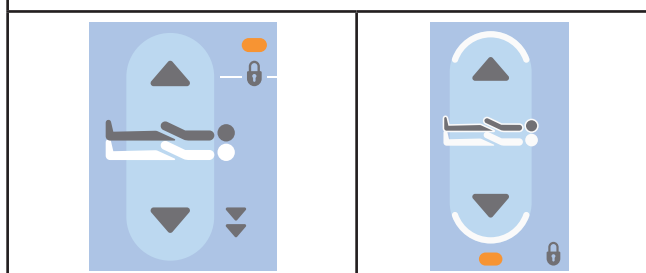



WARNING!

When setting the Extra Low Position (the lowest position of the bed) with the lowered plastic split siderails, there is a risk of injury to the lower limbs of persons standing in the immediate vicinity of the bed or the patient sitting at the bed edge (entrapment between the siderail lower edge and the floor).

- ▶ Make sure that there is no risk mentioned above by checking and alerting people near the bed before setting it to the lowest position!
- ▶ Do not adjust the lowest position of the bed if there is a risk of injury mentioned above!



WARNING!

When setting the Extra Low Position (the lowest position of the bed), once the so-called basic height is reached, downward movement will first temporarily stop and at the same time an acoustic signal (beep) will sound. If the bed is equipped with a light signal, the normal white light will change to green. An acoustic signal, a temporary stop of movement and a change in the colour of the light signal inform the personnel of the potential risk of injury if any of the siderails are in the down position when the lowest position is set. This risk only applies to plastic split siderails.

- ▶ Take extra care during Extra Low Position adjustment!

Essenza 300 W has the possibility of adjusting the bed height to the lowest position in accordance with the table of technical parameters. The lowest bed height is intended as part of measures to reduce the risk of patient falls, especially for high-risk patients.

To use the lowest position function effectively and safely, follow the following precautions:

- ▶ The decision to set the lowest position is made by the healthcare staff after evaluating the risk of falling in the specific condition of the patient (mobility, mental state, age, primary diagnosis).
- ▶ The lowest bed position shall be adjusted exclusively by trained healthcare staff using the dedicated control elements (see „12 Control Elements“ on the page 66).
- ▶ As part of measures to reduce the risk of falling, use the siderails in the upper position when setting the lowest position.
- ▶ Make sure that there is no accessory in probable collision with floor.
- ▶ Make sure that there is no cable or hose that could be pinched or damaged by bed movement to the Extra Low Position.
- ▶ Make sure that there is no need for bed transport after reaching the Extra Low Position because the castor control levers are not accessible in the Extra Low Position.
- ▶ Make sure that there is no need for quickly reached highest position of the mattress support platform after reaching the Extra Low Position because the Bed Height Foot Controls are not accessible in the Extra Low Position.

13.16 Ergoframe

Ergoframe® is the kinematic system of Backrest and Thighrest Adjustment resulting in extension of the Mattress support platform in the seat section. Ergoframe enlarges the space for pelvic area during Auto-contour. Because of increase of the space the force applied results in decrease of the pressure that can cause pressure injuries in the pelvic area. Ergoframe maintains a stable ergonomic position of the body and spine of the patient, thus limiting unwanted movement of the patient by moving down or up in beds. Unified movement eliminates the patient's shift over the mattress and thus maintains a uniform position of the patient's body that is not bound to the position of the bed parts.

14 Locking and Unlocking the Bed Positioning Functions

To disable Bed Positioning Functions for patient, use:

- ▶ Attendant Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Attendant Control Panel for Essenza 300 W
- ▶ Siderail Control Panel for Essenza 300 W with Lateral Tilt
- ▶ Siderail Control Panel for Essenza 300 W

If the bed does not react to individual bed positioning adjustment:

- ▶ Check whether the function is disabled on the Attendant Control Panel or on the Siderail Control Panel.

Disabling individual bed positioning functions on the Attendant Control Panel or on the Siderail Control Panel will affect all control elements.


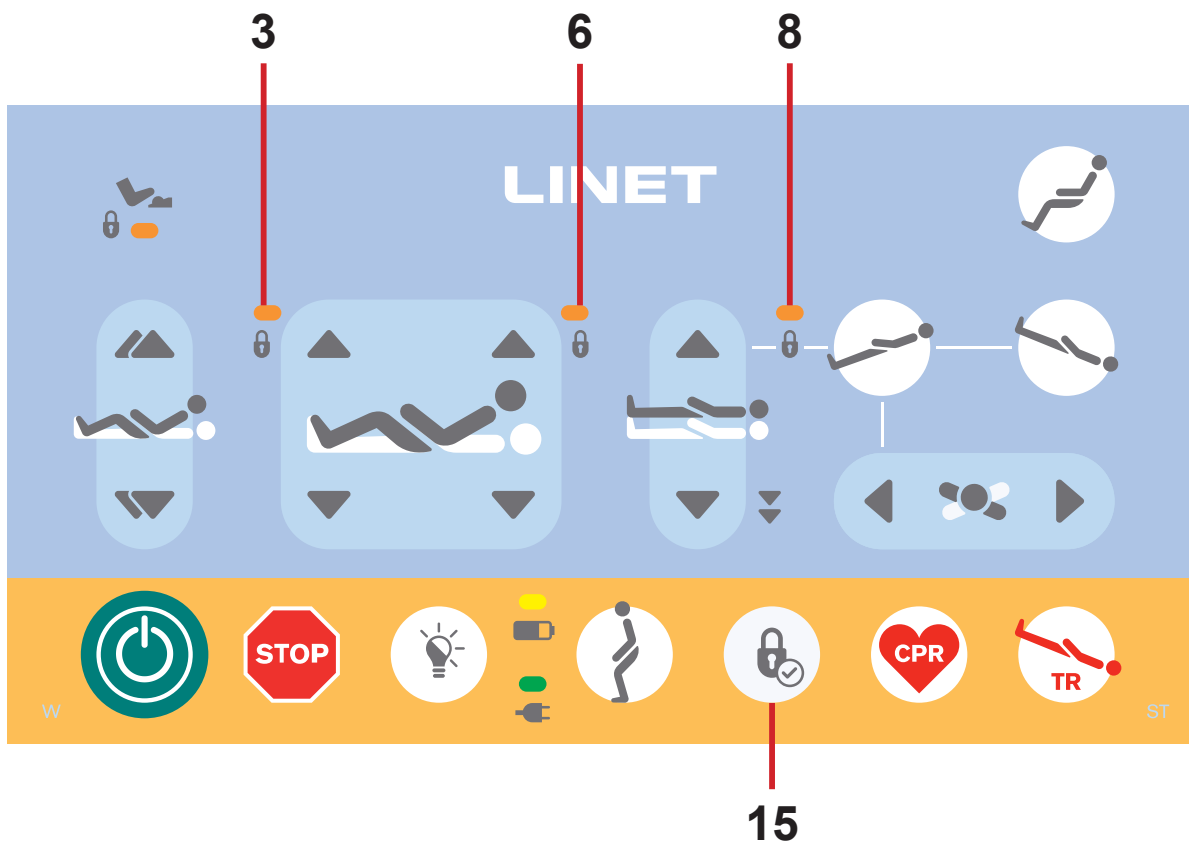


Fig. Lock Button and Lock LED Indicators on the Attendant Control Panel

- 3. Thighrest Adjustment Lock LED Indicator (ON - locked)
- 6. Backrest Adjustment Lock LED Indicator (ON - locked)
- 8. Bed Height Adjustment Lock LED Indicator / Lateral Tilt Lock LED Indicator / Trendelenburg Tilt Lock LED Indicator / Reverse-Trendelenburg Tilt Lock LED Indicator (ON - locked)
- 15. Lock Button

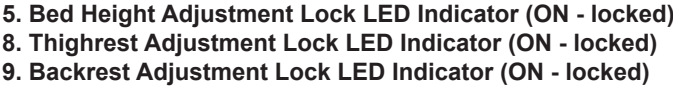

14.1 Locking

To lock **Backrest Adjustment**:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Backrest Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **6** on the Attendant Control Panel is illuminated.
 Lock LED Indicator **5** on the Siderail Control Panel is illuminated.
 Lock LED Indicator **8** on the Handset Control Panel is illuminated.
 Backrest Adjustment is disabled for the corresponding control element.

Lock Button for Locking and Unlocking:



To lock **Thighrest Adjustment**:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Thighrest Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **3** on the Attendant Control Panel is illuminated.
 Lock LED Indicator **8** on the Siderail Control Panel is illuminated.
 Lock LED Indicator **7** on the Handset Control Panel is illuminated.
 Thighrest Adjustment is disabled for the corresponding control element.

To lock **Bed Height Adjustment, Lateral Tilt Adjustment, Trendelenburg Tilt Adjustment and Reverse-Trendelenburg Tilt Adjustment**:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Bed Height Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **8** on the Attendant Control Panel is illuminated.
 Lock LED Indicator **3** on the Siderail Control Panel is illuminated.
 Lock LED Indicator **4** on the Handset Control Panel is illuminated.
 Bed Height Adjustment, Lateral Tilt Adjustment, Trendelenburg Tilt Adjustment and Reverse-Trendelenburg Tilt Adjustment are disabled for the corresponding control element.

14.2 Unlocking

To unlock Backrest Adjustment:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Backrest Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **6** on the Attendant Control Panel is not illuminated.
 Lock LED Indicator **5** on the Siderail Control Panel is not illuminated.
 Lock LED Indicator **8** on the Handset Control Panel is not illuminated.
 Backrest Adjustment is enabled again for the corresponding control element.

Lock Button for Locking and Unlocking:



To unlock Thighrest Adjustment:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Thighrest Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **3** on the Attendant Control Panel is not illuminated.
 Lock LED Indicator **8** on the Siderail Control Panel is not illuminated.
 Lock LED Indicator **7** on the Handset Control Panel is not illuminated.
 Thighrest Adjustment is enabled again for the corresponding control element.

To unlock Bed Height Adjustment, Lateral Tilt Adjustment, Trendelenburg Tilt Adjustment and Reverse-Trendelenburg Tilt Adjustment:

- ▶ Press GO Button.
 - ▶ Press Lock Button.
- Lock LED Indicators **3**, **6** and **8** on the Attendant Control Panel are flashing.
 Lock LED Indicators **3**, **5** and **8** on the Siderail Control Panel are flashing.
 Lock LED Indicators **4**, **7** and **8** on the Handset Control Panel are flashing.
- ▶ Press any **Bed Height Adjustment Button** during flashing of the LEDs.
- Lock LED Indicator **8** on the Attendant Control Panel is not illuminated.
 Lock LED Indicator **3** on the Siderail Control Panel is not illuminated.
 Lock LED Indicator **4** on the Handset Control Panel is not illuminated.
 Bed Height Adjustment, Lateral Tilt Adjustment, Trendelenburg Tilt Adjustment and Reverse-Trendelenburg Tilt Adjustment are enabled again for the corresponding control element.

15 Equipment

Product equipment depends on product configuration.

15.1 Backrest stop at 30 degrees

Backrest stops automatically in 30 degrees during continuous Backrest positioning if the Essenza 300 W bed is equipped with this optional feature.

To continue in the Backrest positioning:

- release the Backrest adjustment button, then press and hold the Backrest adjustment button until the desired position is reached.

15.2 Accessory Rail with plastic hooks

Accessory rail with plastic hooks is intended for hanging accessories. It is located on the sides of bed.

Maximum load of the Accessory Rail is 10 kg without leverage.

Maximum load of the plastic hook intended for the Accessory Rail is 2 kg.

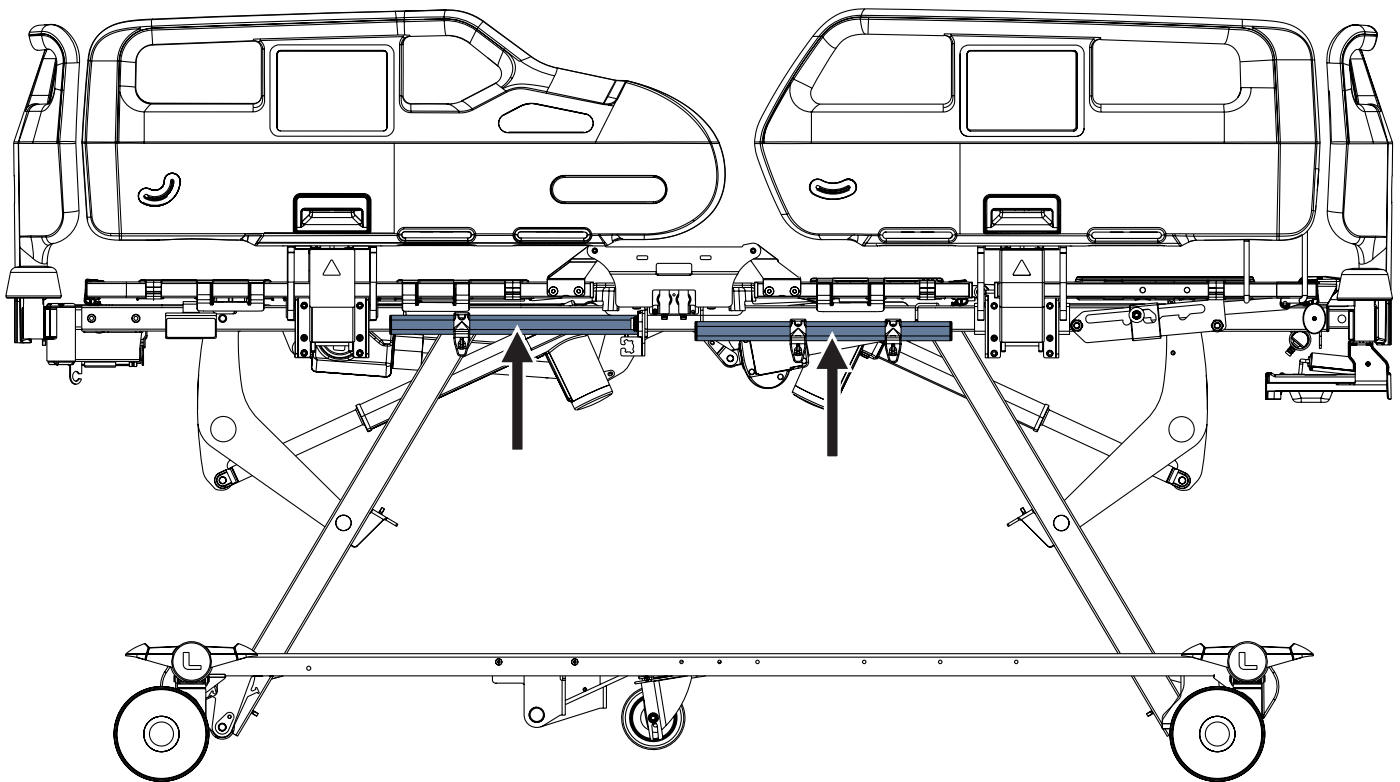


Fig. Accessory Rails with plastic hooks on the bed with Plastic Siderails

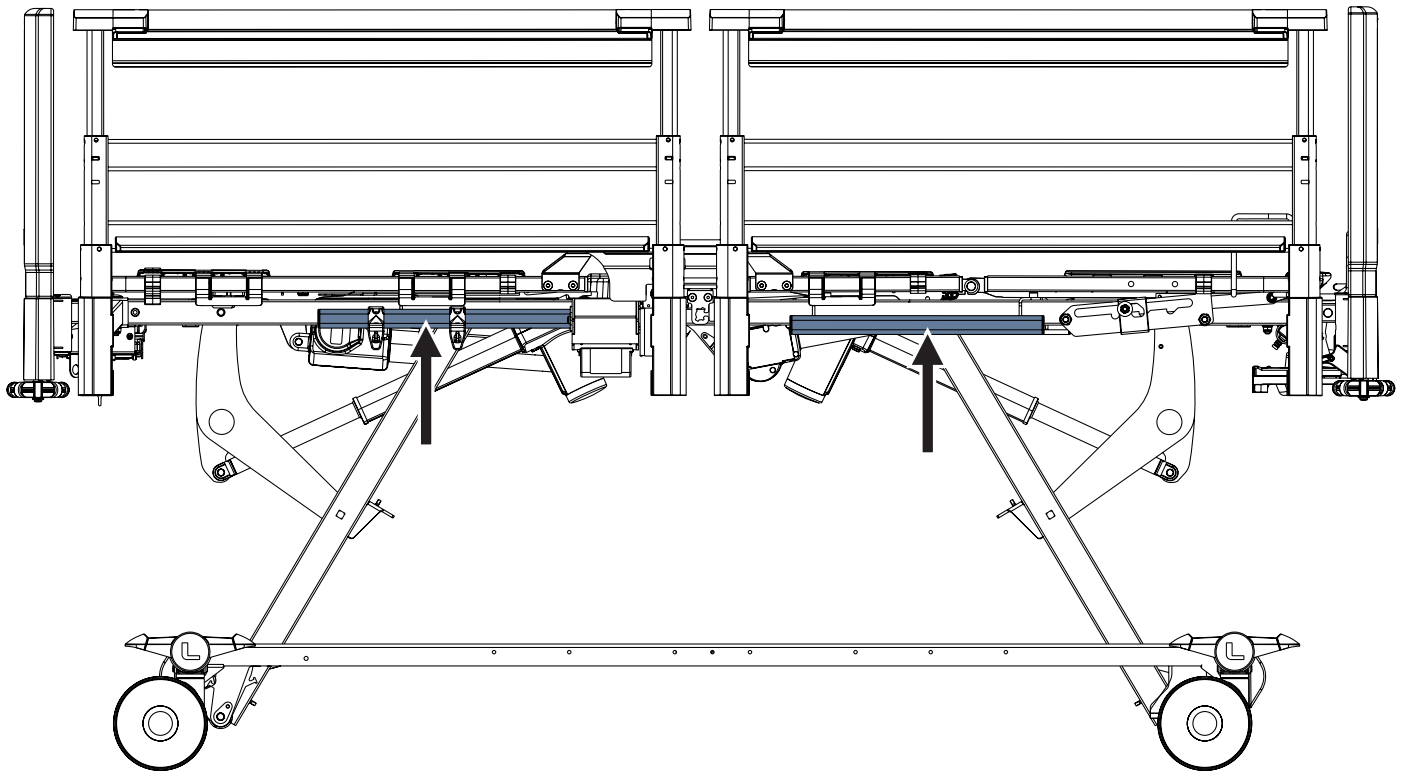


Fig. Accessory Rails with plastic hooks on the bed with Telescopic Siderails

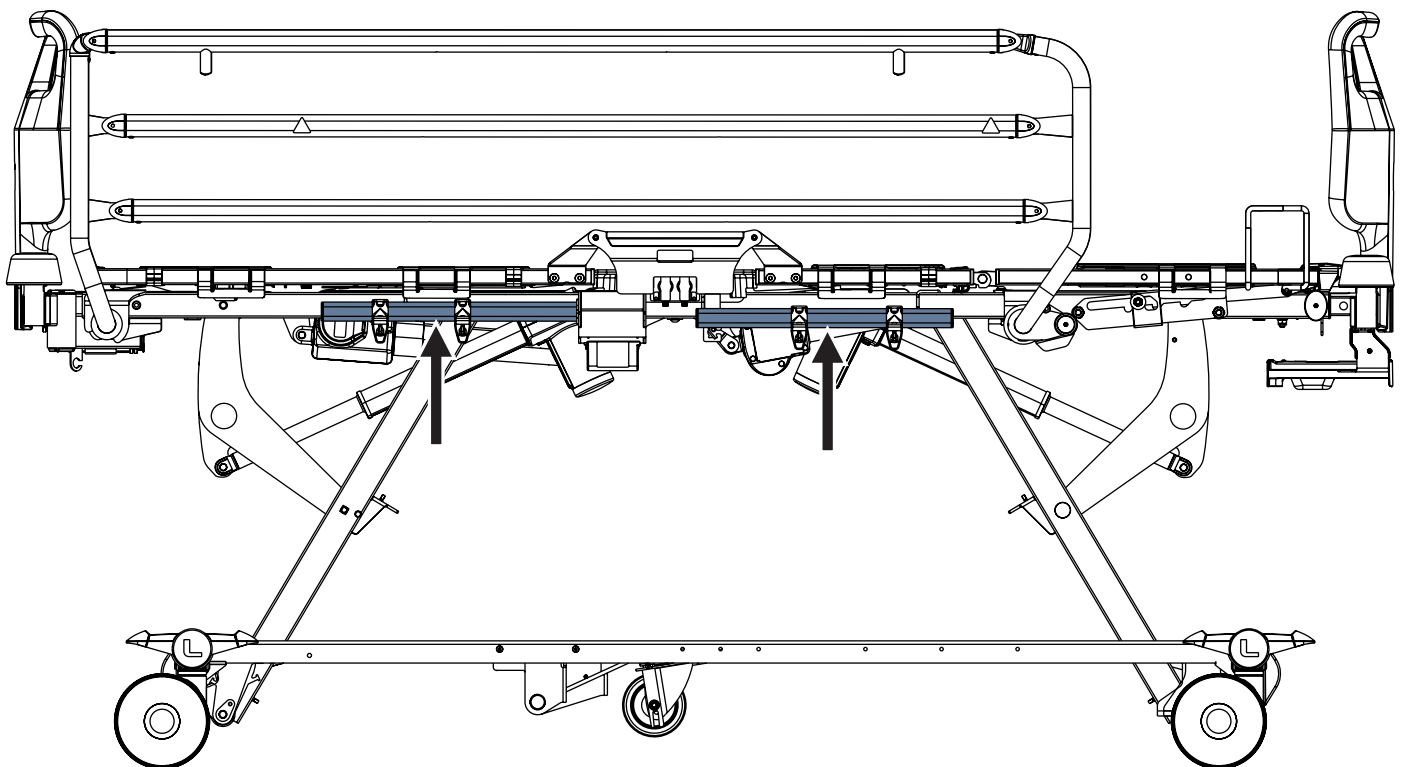


Fig. Accessory Rails with plastic hooks on the bed with Collapsible Siderails

15.3 Urinary Bag Holders in Plastic Siderails



CAUTION!

Risk of collision between Urinary Bags suspended from the Urinary Bag Holders and bed undercarriage or floor!

- Take extra care during lowering to the Extra Low Position if an Urinary Bag is suspended from an Urinary Bag Holder!

Eyelets in the Plastic Siderails are intended for hanging Urinary Bags.

Maximum Load of a compatible Urinary Bag is 3 kg.

A Plastic siderail should be lifted in up position if an Urinary Bag is suspended from an Urinary Bag Holder.

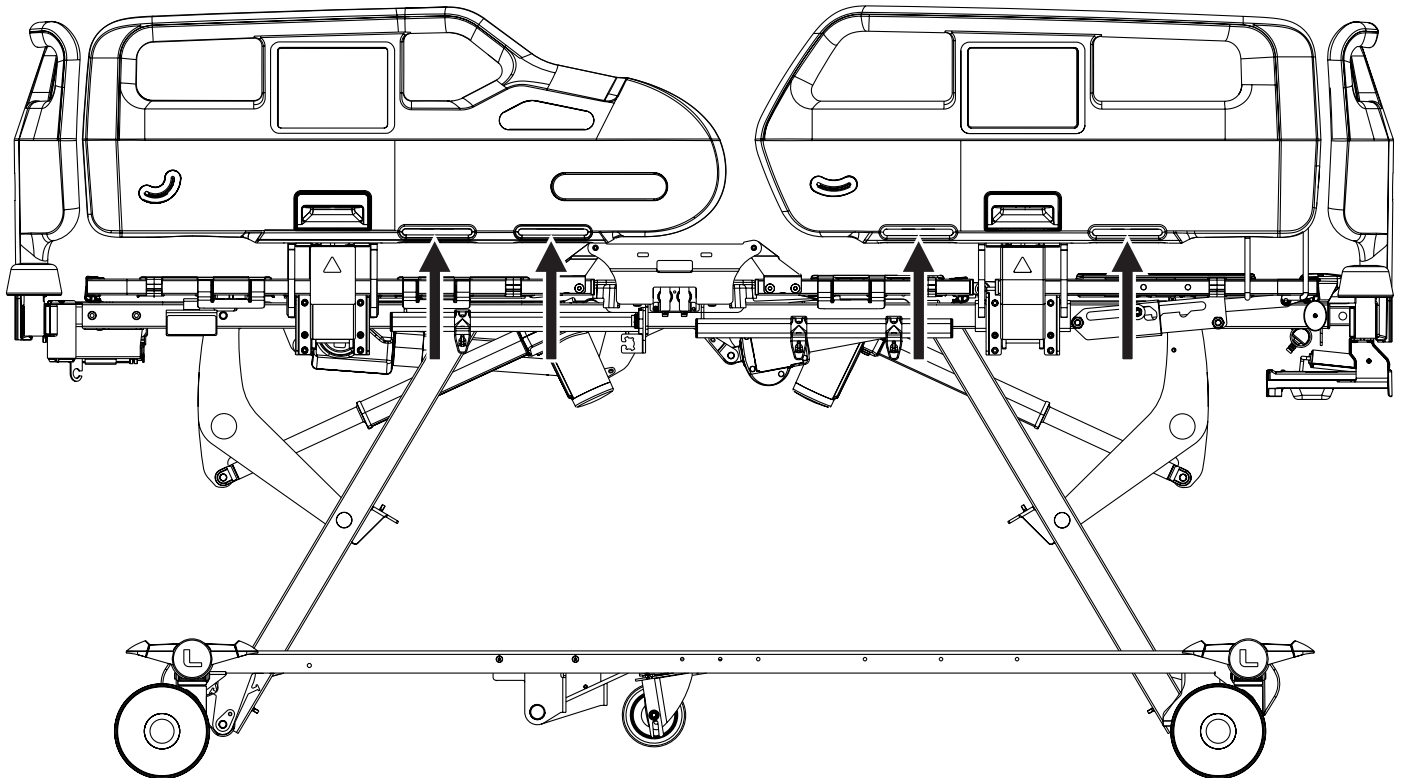


Fig. Urinary Bag Holders in Plastic Siderails

15.4 Bed Illumination and Handset Control Panel Illumination

The Essenza 300 W bed can optionally be equipped with Bed Illumination (the source of light situated under the seat section) and Handset Control Panel can optionally be equipped with illuminated keyboard. Bed Illumination and Handset Control Panel Illumination help the nursing staff as well as the patient to orientate.

Bed Illumination Control Button:

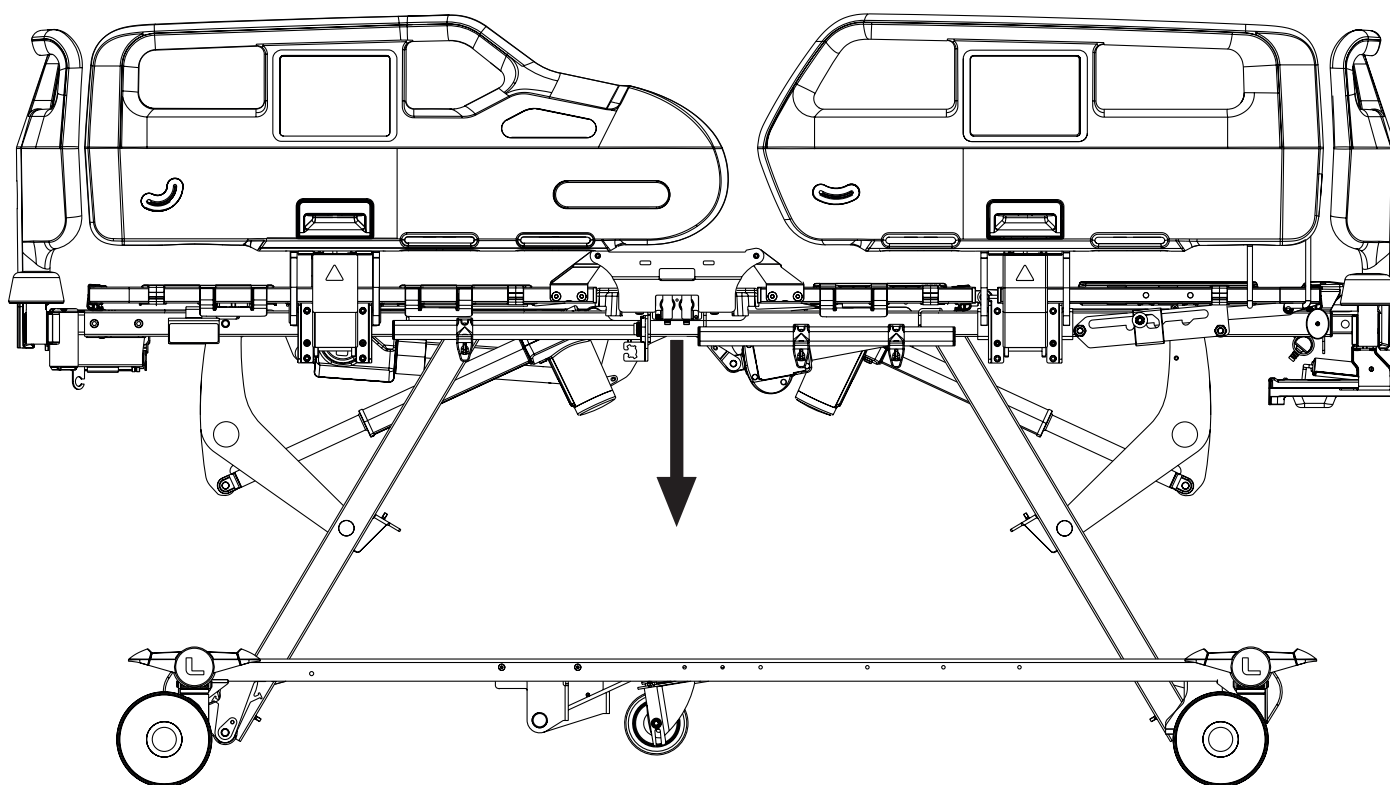


Fig. Bed Illumination source on the Essenza 300 W medical bed

The Essenza 300 W bed is optionally equipped with three-phase illumination: Full intensity of illumination, Lowered intensity of illumination and Illumination is turned off.

The lowered intensity of illumination is set up after the bed is connected to the mains. The illumination is turned off when the bed is powered by accumulator only.

After pressing any button:

- Handset Control Panel keyboard and Bed Illumination will light up at full intensity within 3 minutes.

After elapsing the full intensity period:

- The illumination intensity of the Handset Control Panel keyboard and Bed Illumination will be lowered within 7 minutes.

Turning off all illumination:

- Press **Bed Illumination Control Button** on the Attendant Control Panel or on the Siderail Control Panel and all illumination on the bed will be turned off.

15.5 Safe Height Visual Signalling

Green light emitted under the bed seat section is intended for the safe height visual signalling.
The green light cannot be activated when the bed is not connected to the mains power.
The green light is emitted in the Extra Low Position only.

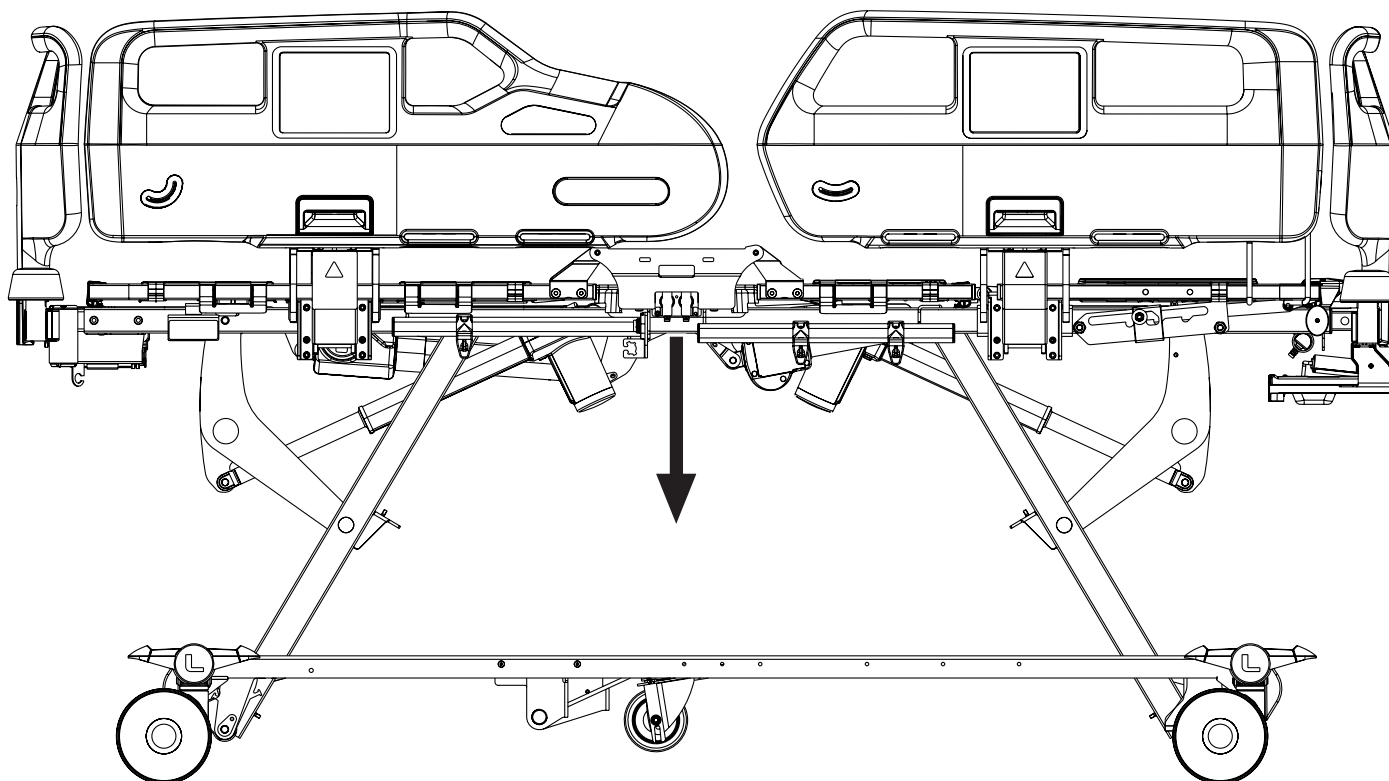


Fig. Green light source on the Essenza 300 W medical bed

15.6 Fifth Castor



WARNING!

Risk of damaging the Fifth Castor during its movement across slope greater than 15%!

- Ensure the bed does not move over the slope!



The bed can be equipped with the Fifth Castor to make it easier to manoeuvre and control the bed movement. The standard locking castor will automatically be excluded from the bed configuration if the bed is equipped with the Fifth Castor as this will replace the need for this function. The control lever for the Fifth Castor is identical to the control lever for the braking system.

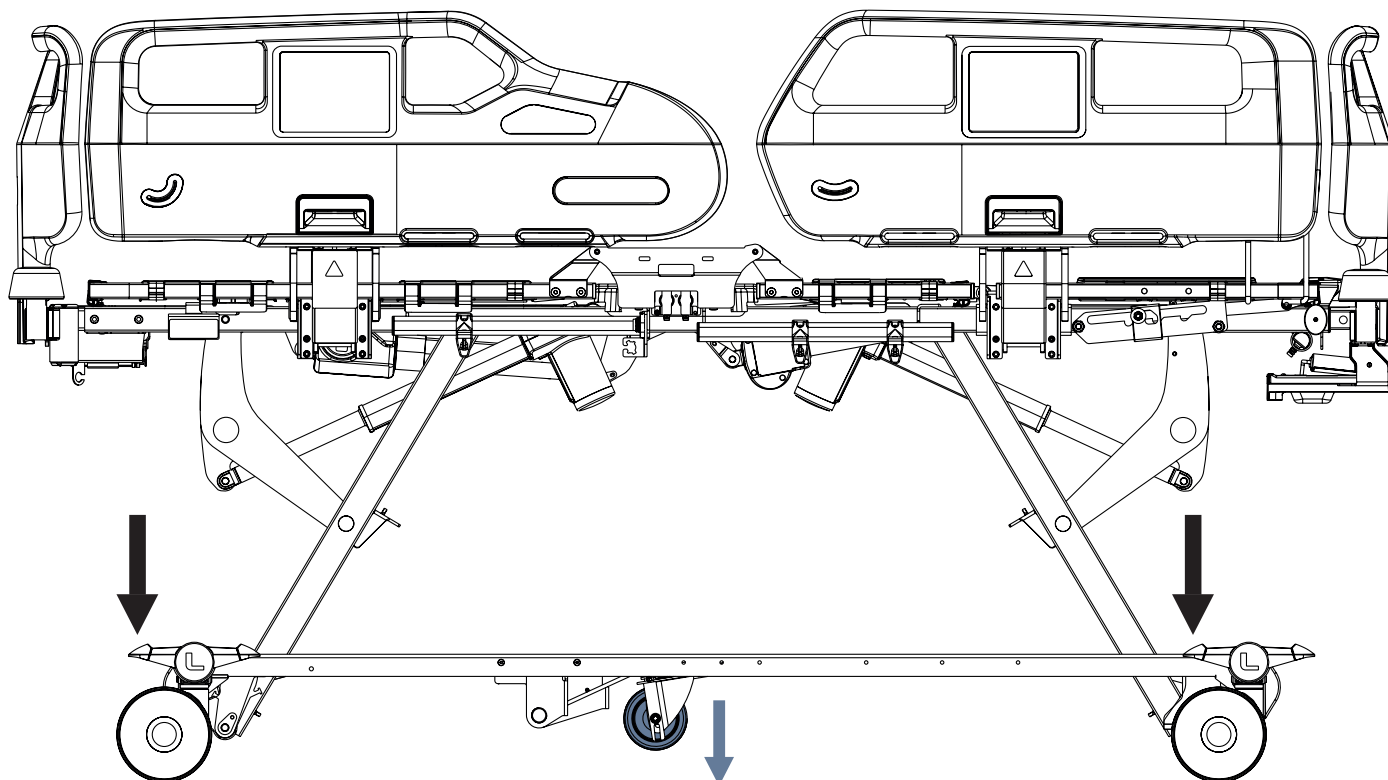


Fig. Activation of the Fifth Castor

To activate the Fifth Castor:

- Press green Drive pedal to the low position.
The Fifth Castor moves in the direction of the bed head end and in the direction of bed foot end.

To release the Fifth Castor from its limited position in the direction of the bed head end and in the direction of bed foot end :

- Leave all Brake pedals and Drive pedals unpressed in their middle position.
The Fifth Castor moves sideways and around.

15.7 Angle Indicators in Plastic Siderails

Angle Indicators are situated on both outer sides of plastic siderails. Head siderail indicators are intended for an approximate Backrest angle measurement. Foot siderail indicators are intended for an approximate measurement of Trendelenburg tilt and Reverse-Trendelenburg tilt.

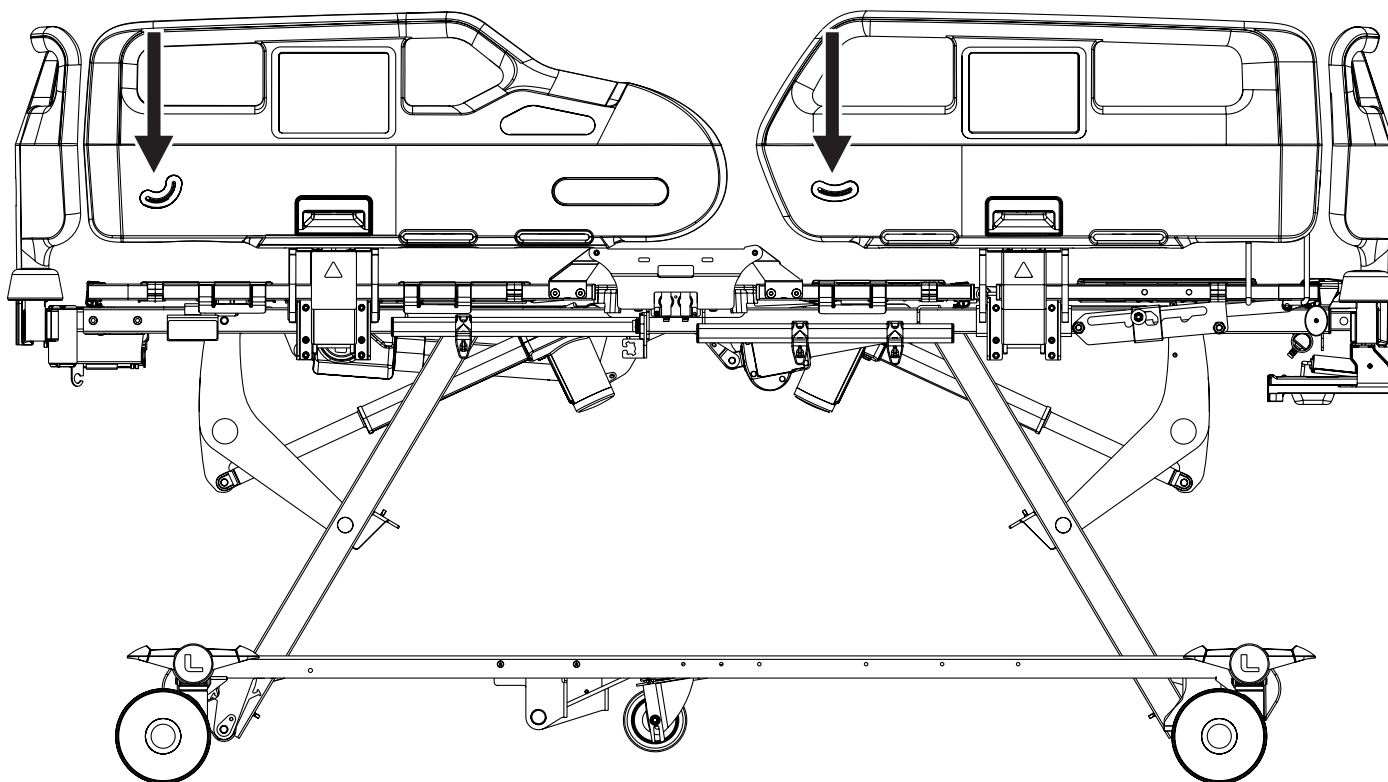


Fig. Angle Indicators in Plastic Siderails

15.8 Plastic Linen Shelf



WARNING!

- Risk of injury due to incorrect placement of the Attendant Control Panel in the Linen Shelf!
- ▶ Ensure any button of the Attendant Control Panel is not pressed in the Linen Shelf!
 - ▶ Only place the Attendant Control Panel into the Linen Shelf without any other things!

Linen Shelf is intended for storing the things corresponding to its dimensions (e.g. linens or sheets).

Wire construction assists in keeping linens and sheets on the Linen Shelf.

The Attendant Control Panel can be stored in the Linen Shelf.

The Linen Shelf is located at foot end under the Foot Board.

Safe Working Load of the Linen Shelf is 15 kg.

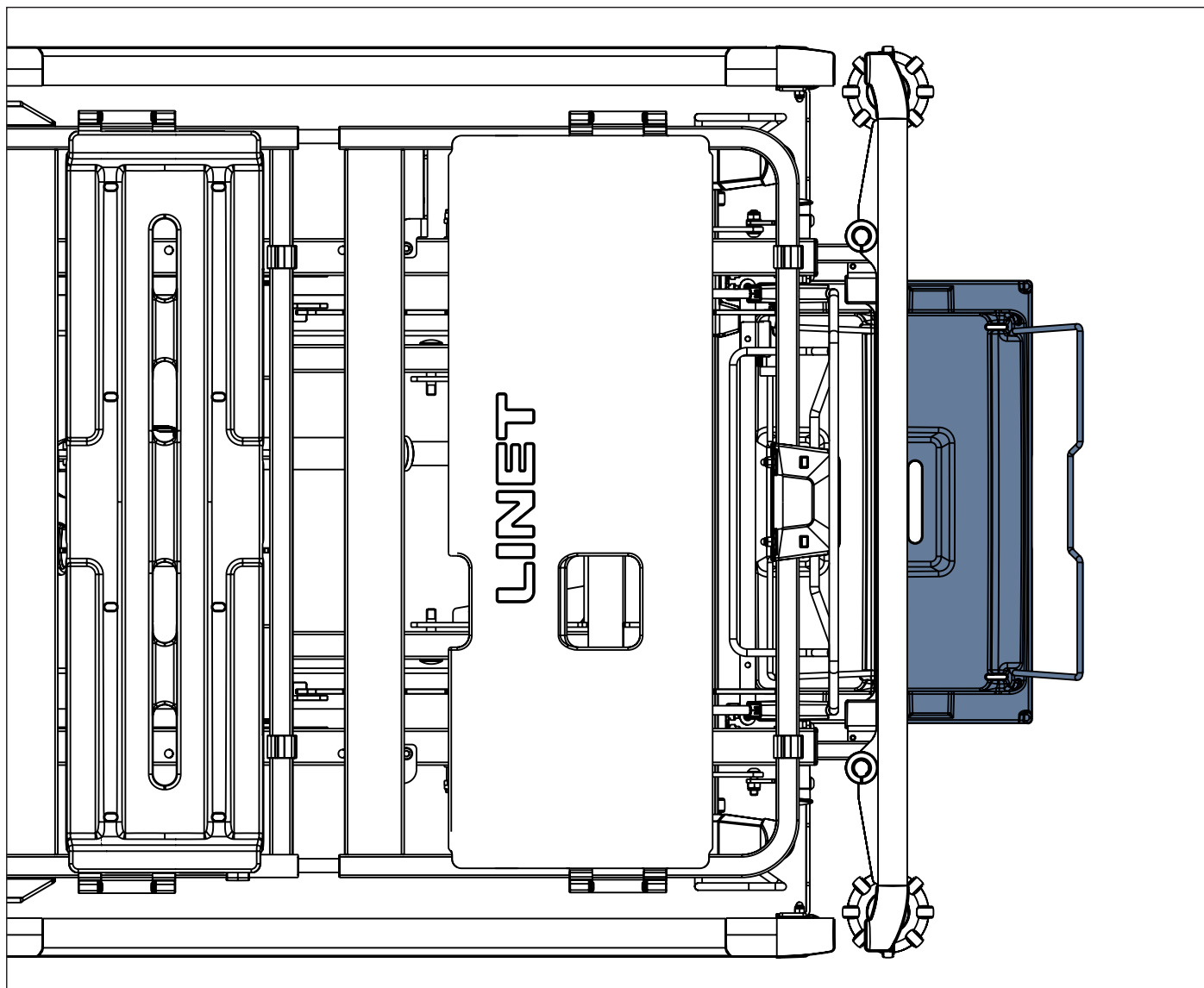


Fig. Plastic Linen Shelf at bed foot end

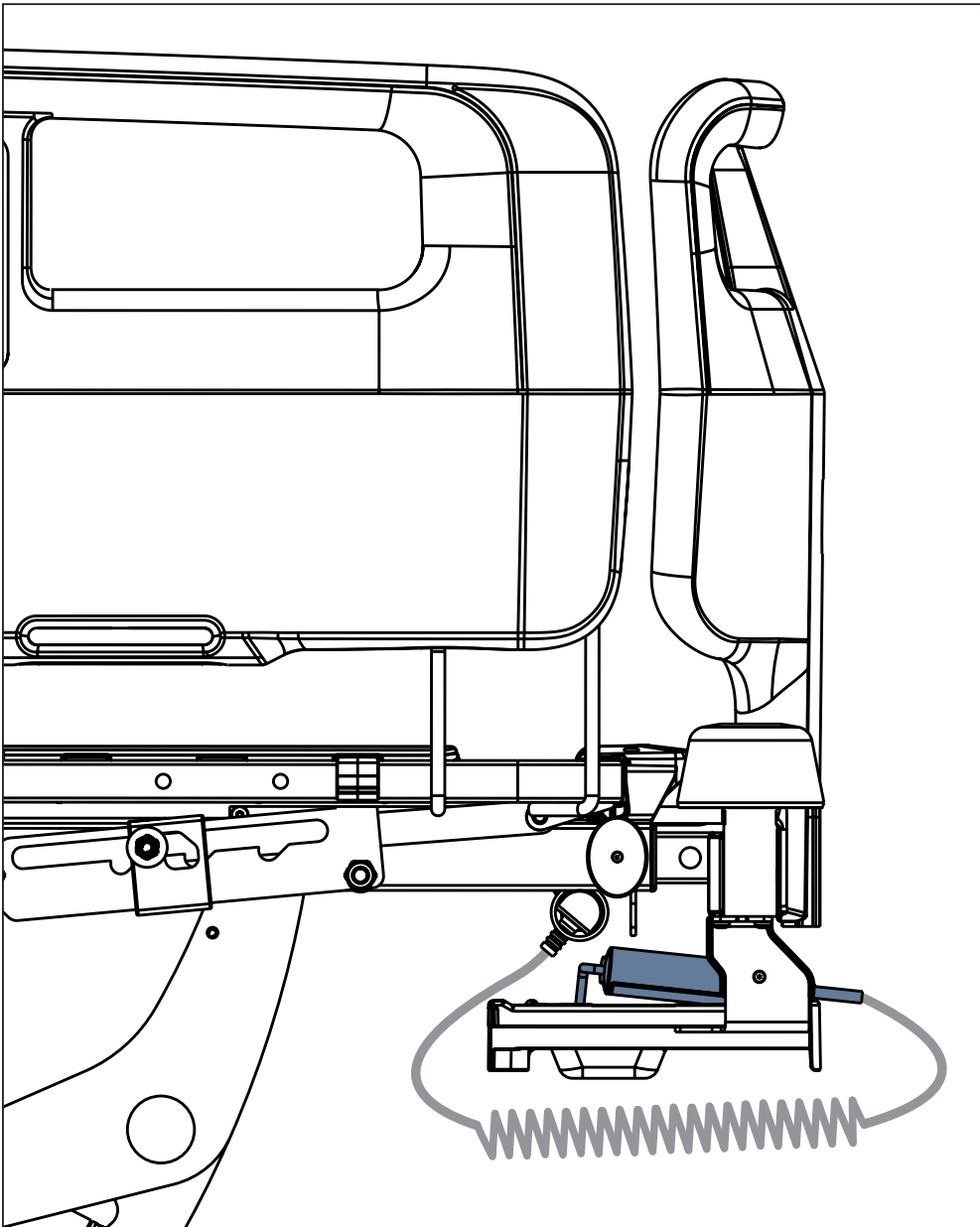


Fig. Attendant Control Panel in the Plastic Linen Shelf

15.9 Telescopic Linen Shelf



WARNING!

Risk of injury due to incorrect placement of the Attendant Control Panel in the Telescopic Linen Shelf!

- ▶ Ensure any button of the Attendant Control Panel is not pressed in the Telescopic Linen Shelf!
- ▶ Only place the Attendant Control Panel into the Telescopic Linen Shelf without any other things!



CAUTION!

Risk of material damage due to the incorrect bed position during Reverse-Trendelenburg Tilt Adjustment!

- ▶ Ensure the mattress support platform is not extended and the Telescopic Linen Shelf is not spread out during Reverse-Trendelenburg Tilt Adjustment!

Telescopic Linen Shelf is intended for storing the things corresponding to its dimensions (e.g. linens or sheets).

Wire construction assists in keeping linens and sheets on the Telescopic Linen Shelf.

The Attendant Control Panel can be stored in the Telescopic Linen Shelf.

The Telescopic Linen Shelf is located at foot end under the Foot Board.

Safe Working Load of the Telescopic Linen Shelf is 15 kg.

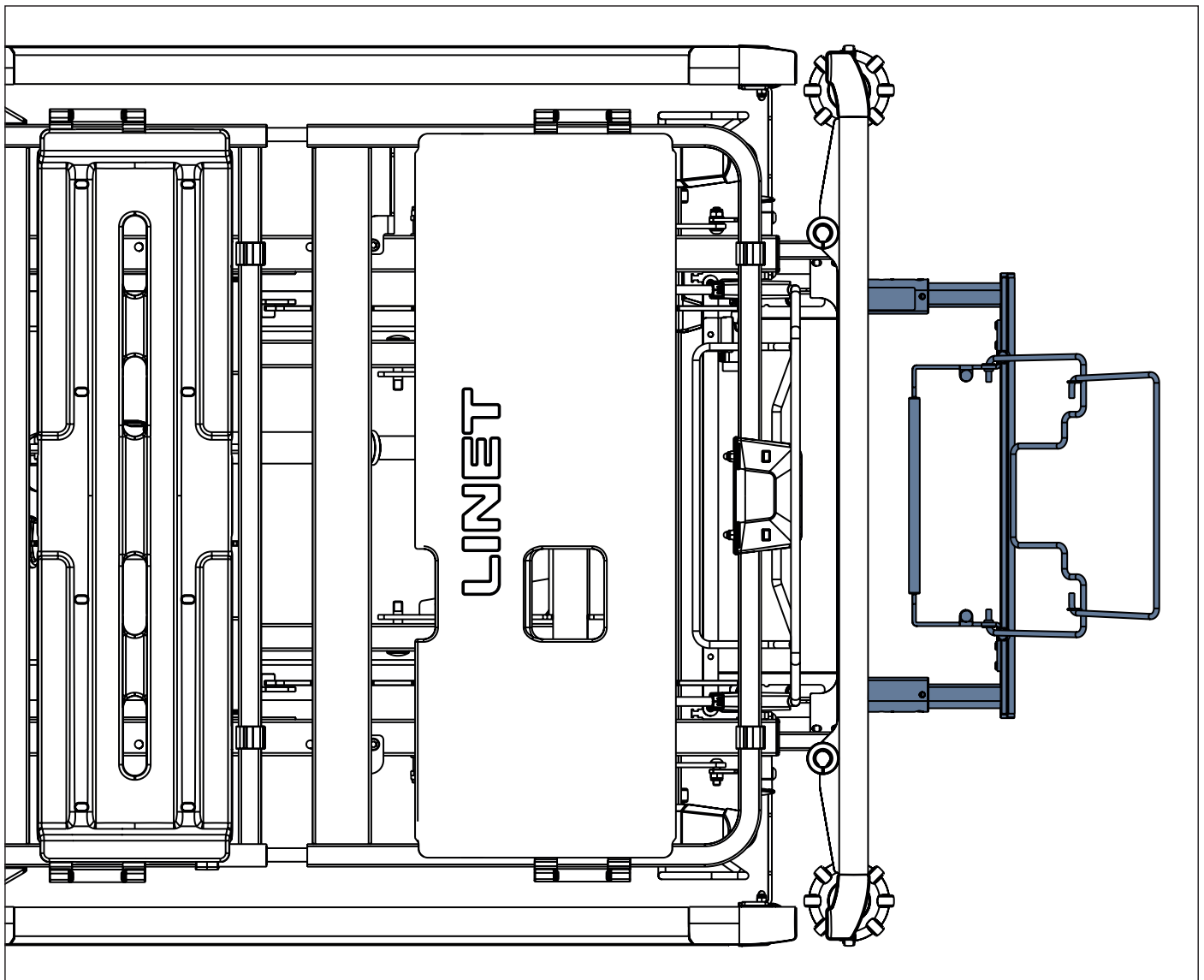


Fig. Telescopic Linen Shelf at bed foot end

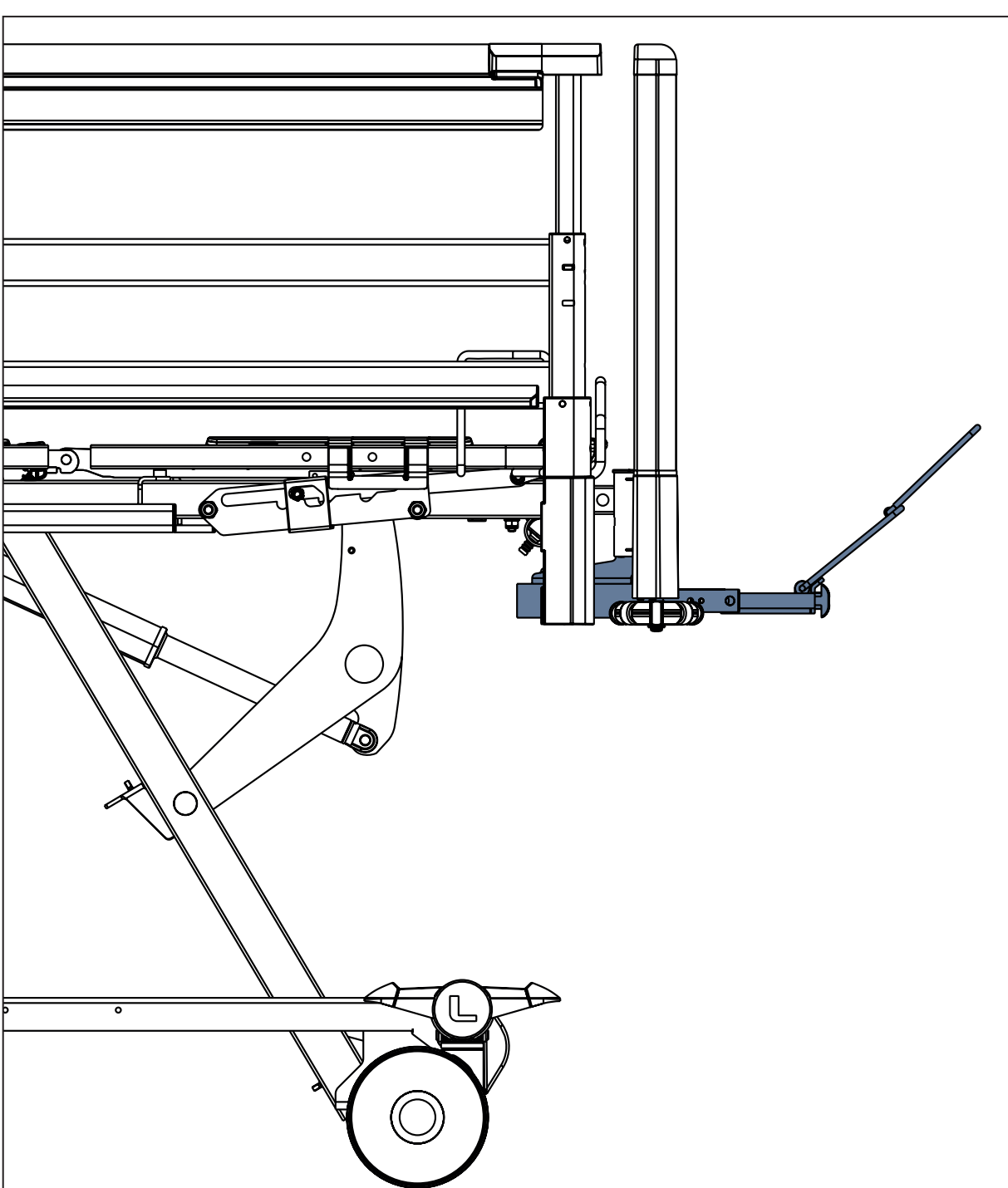


Fig. Telescopic Linen Shelf in spread out position

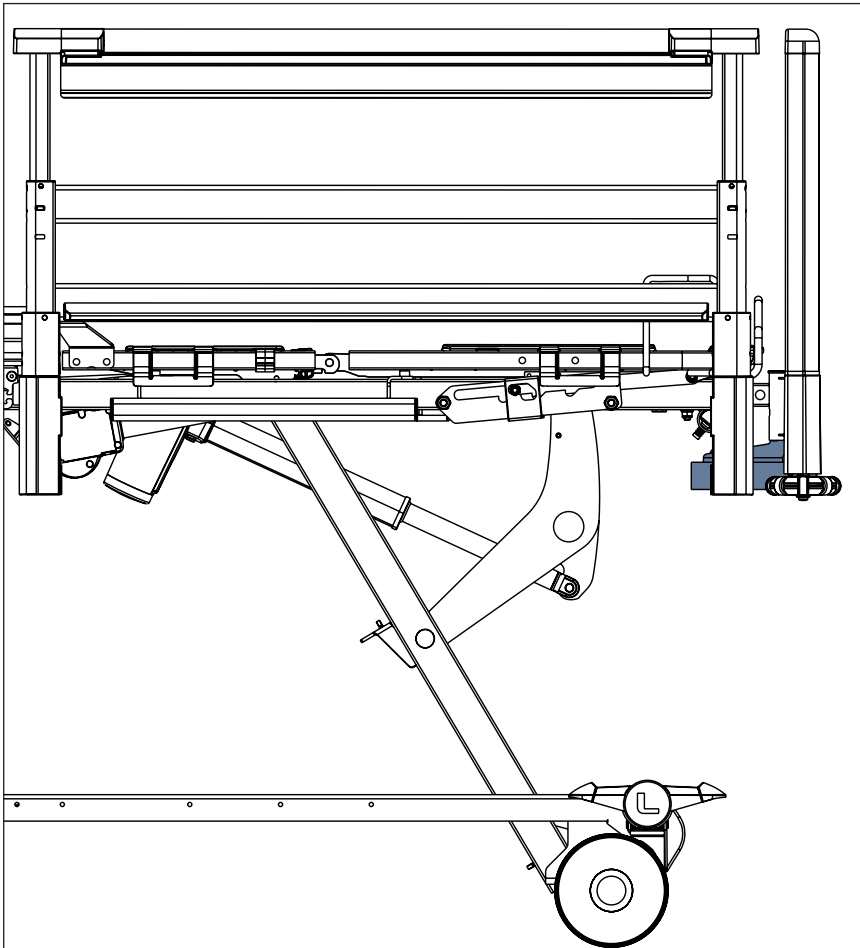


Fig. Telescopic Linen Shelf in retracted position

15.10 Brake Signal

If bed is equipped with Brake Signal and this bed is connected to the mains power, the Brake Signal sounds when the bed is not braked. For specification of the Brake Signal see „1.7 Acoustic Signalling“ on the page 14.

15.11 SafeSense® 3 ready bed


WARNING!

Before use of the SafeSense® 3 Monitoring System with Essenza 300 W medical bed read the instructions for use for SafeSense® 3 Monitoring System!

Specific configuration of the Essenza 300 W medical bed is compatible with SafeSense® 3 Monitoring system. SafeSense® 3 ready bed includes the accessory rail with Power Supply/Data cable in safety position.

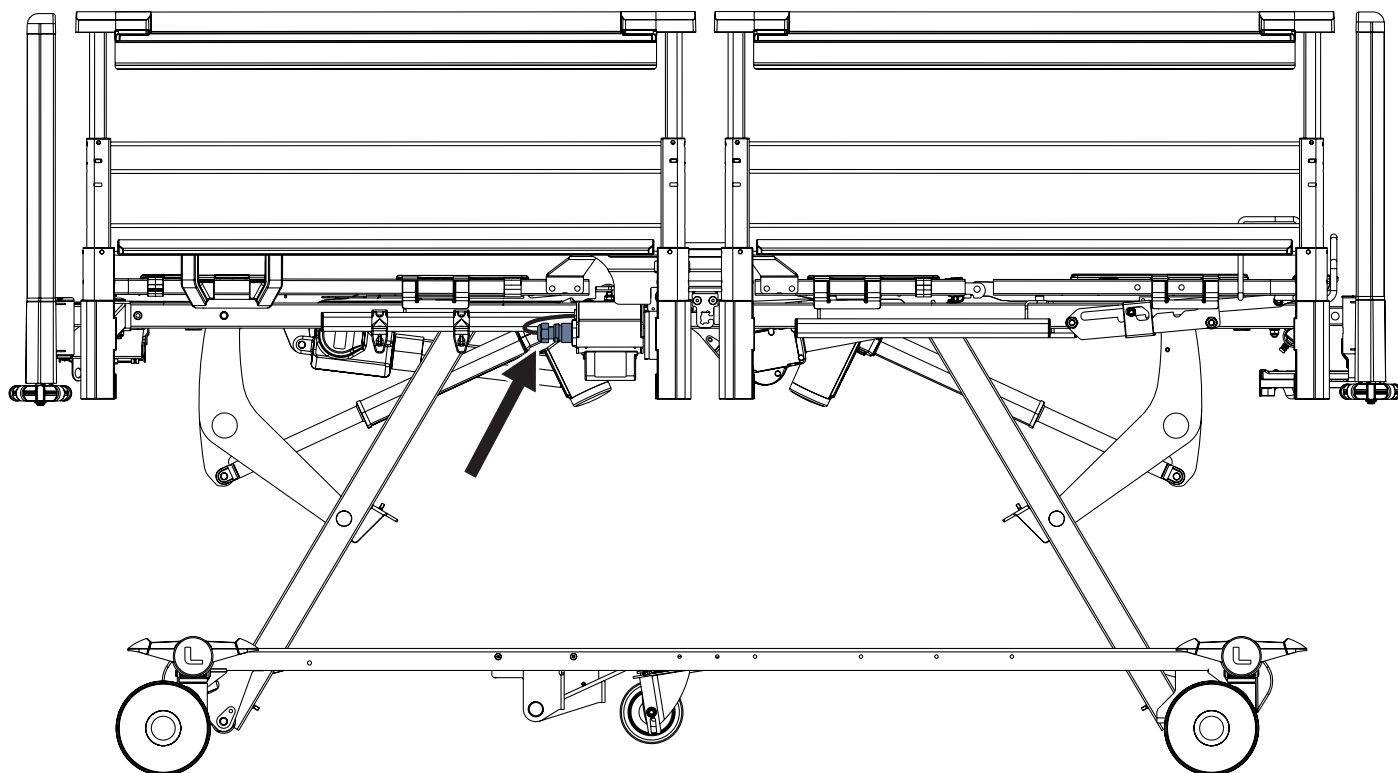


Fig. SafeSense 3 ready Essenza 300 W bed

16 Mattress



WARNING!

Risk of injury due to incorrect use!

- Ensure the compatible mattress is correctly fixed with straps on the mattress support platform after placing the compatible mattress foam extension on the Backrest in front of the mattress head end!



CAUTION!

Incompatibility with bed due to incorrect mattress dimensions!

- In the case of other mattresses check maximum approved mattress dimensions (Technical Specification chapter) and its specific shape.
- Clinical staff is responsible for assessing the risks resulting from the use of mattress with a mattress height exceeding the maximum mattress height!

Essenza 300 W bed is designed for mattresses from LINET portfolio.

Compatible Mattress Foam Extensions are intended to be placed in front of the mattress head end on the Backrest.

Compatible Mattresses	Manufacturers and Identification Numbers	Mattress Mass
EffectaCare 10	LINET: 4PCA21100AS, 4PCA51100AS	6 kg
EffectaCare 20	LINET: 4PPB21100AS, 4PPB51100AS	6,5 kg
PrimaCare 10	LINET: 4PPD21100AS, 4PPD51100AS	10 kg
MediMatt 30	WULFF MED TEC: 4PWR21100LS, 4PWR51100LS, 4PWR71100LS, WULF: 4PWRH1300LS, 4PWR21100LS, 4PWR51100LS, 4PWR71100LS	11,65 kg; 9,4 kg (4PWR51100LS)
ViskoMatt 10	WULFF MED TEC: 4PW121100LS, 4PW151100LS, 4PW171100LS, 4PW1H1300LS, 4PW121100LS, 4PW151100LS, 4PW171100LS	9,4 kg
ViskoMatt 30	WULFF MED TEC: 4PW321100LS, 4PW351100LS, 4PW371100LS, 4PW3H1300LS, 4PW321100LS, 4PW351100LS, 4PW371100LS	14,04 kg; 9,4 kg (4PW371100LS)
ViskoMatt 50	WULFF MED TEC: 4PW521100LS, 4PW551100LS, 4PW581100LS, 4PW5H1300LS, 4PW521100LS, 4PW551100LS, 4PW581100LS	17,82 kg
Air2Care 5, 6, 8	AIRFLO (XIAMEN) MEDICAL: 4VEM000340B00, 4VEM000340C00, 4VEM000440C00, 4VEM000640C00, 4VEM000740C00, 4VEM000340D00, 4VEM000540D00	6 kg (4VEM000340B00); 7,5 kg (4VEM000340C00); 8,5 kg (4VEM000440C00); 8 kg (4VEM000640C00); 8,8 kg (4VEM000740C00); 9,3 kg (4VEM000340D00); 11,2 kg (4VEM000540D00)
CliniCare 100 HF	LINET: 4VCM0001R01P0, 4VCM0001R11M0	16 kg
Virtuoso (Series III)	LINET: 1VNM000220M00 (Virtuoso Overlay), 1VNM000220L00 (Virtuoso 50), 1VNM000220F00 (Virtuoso 100), 1VNM000220G00 (Virtuoso 200), 1VNM000220H00 (Virtuoso 300), 1VNM000220K00 (Virtuoso Pro)	6,5 kg (Virtuoso Overlay); 12,5 kg (Virtuoso 50); 8,5 kg (Virtuoso 100); 10 kg (Virtuoso 200); 12 kg (Virtuoso 300); 13 kg (Virtuoso Pro)
HybriMatt 200	LINET: 4VHM201-2	15 kg
HybriMatt 100	LINET: 4VHM101-5	12 kg
ViskoMatt 30 W	WULFF MED TEC: 4PW621100LS, 4PW671100LS	

Dimensions of Compatible Foam Extensions	Manufacturers and Identification Numbers	Foam Extension Mass
86 cm x 14 cm x 20 cm	LINET: 4PRPRODMAT11	0,5 kg
86,5 cm x 14 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT25	0,5 kg
86,5 cm x 16 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT31	0,5 kg
86 cm x 20 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT35	0,5 kg
86 cm x 15 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT36	0,5 kg
90 cm x 14 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT28	0,5 kg
90 cm x 16 cm x 20 cm	WULFF MED TEC: 4PRWRODMAT34	0,5 kg

Mattress Installation

Compatible mattresses are equipped with straps intended for fixing mattress on the mattress support platform.

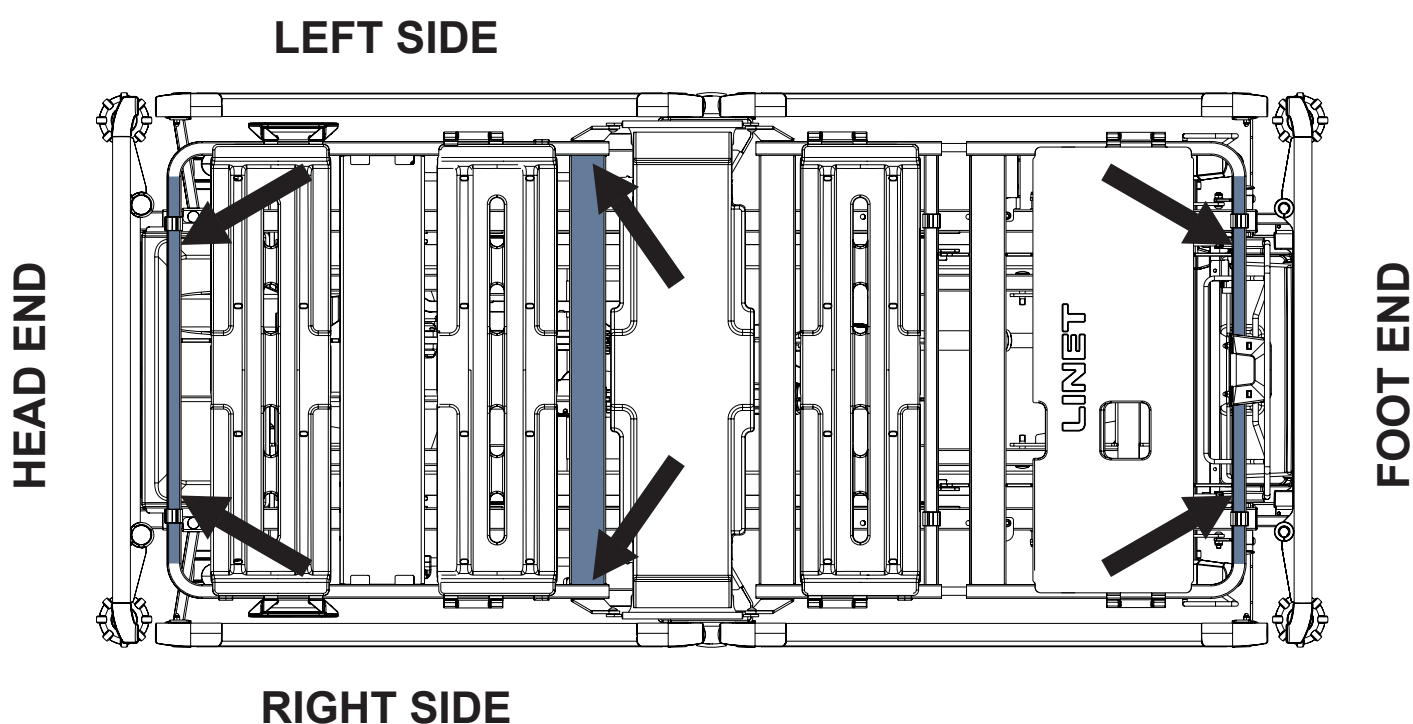


Fig. Places for fixation of the mattress on the mattress support platform of the Essenza 300 W bed



Patient on Mattress

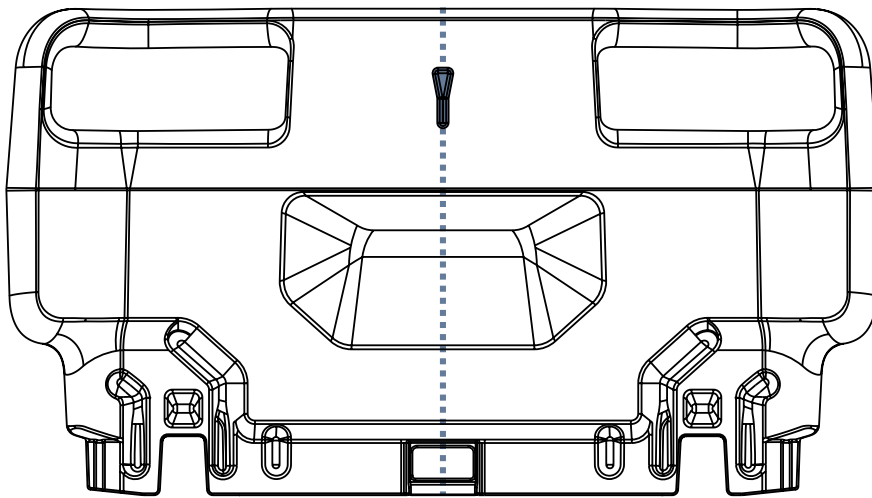


Fig. Arrow-like symbol indicating bed centre on Head Board and Foot Board

To put patient in bed with mattress:

- ▶ Head of a lying patient must point towards the bed head end.
- ▶ Put patient to the bed centre according to the arrow-like symbols on Head Board and Foot Board and according to position of the seat section on the mattress support platform.

16.1 Passive Mattress

Compatible Passive Mattresses are equipped with straps intended for fixing mattress on the Mattress support platform. Symbols on the mattress wash label serve as instructions for cleaning the mattress cover.

16.2 Active Mattress



WARNING!
Follow instructions for use of a compatible active mattress carefully!



CAUTION!
Risk of material damage due to an incorrect fixation of compatible active mattress on the mattress support platform!
▶ Adjust the bed to maximum Cardiac Chair Position before fixing all the straps of the inflated mattress to the mattress support platform!

Installation instructions:

- ▶ Remove any existing mattress.
- ▶ Observe mattress dimensions and its orientation before putting it on the Mattress support platform.
- ▶ Place SCU on the foot board of the bed or on the floor.
- ▶ Fix mattress on the Mattress support platform with straps.

16.3 Cleaning of Mattress



CAUTION!

Incorrect cleaning/disinfection can damage the mattress!

- ▶ Do not use pressure or steam cleaners.
- ▶ Follow the instructions and observe the dosages recommended by the manufacturer.
- ▶ Ensure that disinfectants are selected and applied by qualified hygiene experts only.
- ▶ The surface of the mattress should not be exposed to liquids for a long time.

16.3.1 General Guidance

For safe and gentle cleaning:

- ▶ Do not use any strong acids or alkalines, (optimum pH range 6 – 8. Do not exceed pH of 9).
- ▶ Only use detergents that are suitable for cleaning medical equipment.
- ▶ Do not use abrasive powders, steel wool, or other material and cleaning agents that might damage the mattress. Do not scrub mattress surface.
- ▶ Never use any corrosive or caustic detergents.
- ▶ Never use detergents that deposit calcium carbonate.
- ▶ Never use detergents with solvents that might affect the structure and consistency of the plastics (benzene, toluene, acetone etc.).
- ▶ Use only hospital-approved cleaners and observe local directives concerning infection control.
- ▶ Always rinse with water after cleaning and dry thoroughly before use.
- ▶ Observe local directives concerning infection control.

Mattress parts to be cleaned	Recommended Cleaning Agents (General Cleaning
Top Cover, Bottom Cover, Protect Cover	Standard hospital detergents, Alcohol or Quaternary Ammonium based disinfectants, Chlorine based disinfectants containing up to 0,1% Chlorine, followed by rinsing with water and drying thoroughly before use.
	Decontamination: Blood spills/Clostridium difficile etc Chlorine based disinfectants containing up to 0,1% Chlorine. Dwell time on surface at 0,1% of 5 minutes, followed by rinsing with water and drying thoroughly before use.
Mattress Core	Do not clean!

Due to the variety of laundry equipment, chemicals and conditions in use, customers should satisfy themselves through pretesting. It is essential that cover be thoroughly rinsed and dried after all cleaning procedures and before storage or reuse. Wet or damp PU surfaces are more prone to mechanical damage than when dry.

As stated above, after application of a suitable cleaner, the surface must be rinsed with water and dried before use. (Even if the cleaner instructions say that this is not required). This prevents a build-up of chemicals on the mattress surface which could be reactivated during use and affect biocompatibility.

NOTE Continued use of high concentration, chlorine-based disinfectants may significantly reduce the performance and the working life of a coated material.

NOTE If disinfecting is not required, cleaning with soap and water should be enough to remove dirt stains.

NOTE Cleaning and disinfecting products based on solvent, bleach, abrasives or very high (over 70%) alcohol concentrations can damage this product.

Type of Cleaning	Parts to be cleaned
Routine Cleaning and Disinfection	external of mattress cover
Full Cleaning and Disinfection	external of mattress cover

16.3.2 Routine Cleaning and Disinfection

Cleaning the mattress:

- ▶ Check mattress cover top for any signs of damage or for liquid ingress.
- ▶ Replace or repair and completely disinfect mattress cover top if damaged. Also check if the mattress core is not contaminated. In case of core contamination, do not use the mattress and dispose the core ecologically.
- ▶ Leave mattress cover on mattress.
- ▶ Clean with water with cleaning detergent.
- ▶ Rinse mattress with cold water.
- ▶ Let mattress air dry or wipe dry.
- ▶ Wipe mattress with disinfectant and rinse mattress with cold water.
- ▶ Let mattress dry or wipe dry.

16.3.3 Complete Cleaning and Disinfection

Cleaning Top/Bottom Cover:

Use standard hospital detergents, Alcohol based cleaners or Quaternary Ammonium based disinfectants. Suitable Chlorine based cleaners can be used at a concentration of 0,05%. Stronger concentrations of chlorine can be used if required, (up to 0,1%), with a maximum dwell time of five minutes followed by rinsing with water and drying thoroughly before use.

After application of a suitable cleaner, the surface should be rinsed with water and dried before use. (Even if the cleaner instructions say that this is not required). This prevents a build up of chemicals on the mattress surface that could reactivate during use and affect biocompatibility.

Cleaning Protect Cover:

Use standard hospital detergents, Alcohol based cleaners or Quaternary Ammonium based disinfectants. Suitable Chlorine based cleaners can be used at a concentration of 0,05%. Stronger concentrations of chlorine can be used if required, (up to 0,1%), with a maximum dwell time of five minutes followed by rinsing with water and drying thoroughly before use.

After application of a suitable cleaner, the surface should be rinsed with water and dried before use. (Even if the cleaner instructions say that this is not required). This prevents a build up of chemicals under the Top/Bottom Cover that could reactivate during use and affect biocompatibility.

Cleaning the mattress:

- ▶ Check mattress cover top and base for any signs of damage.
- ▶ Replace or repair and completely disinfect mattress cover top and base if damaged. Also check if the mattress core is not contaminated. In case of core contamination, do not use the mattress and dispose the core ecologically.
- ▶ Leave mattress cover on mattress.
- ▶ Clean with water with cleaning detergent.
- ▶ Rinse mattress with cold water.
- ▶ Let mattress air dry or wipe dry.
- ▶ Wipe mattress with disinfectant.
- ▶ Rinse mattress with cold water.
- ▶ Let mattress air dry or wipe dry.

Machine washing of the top/base mattress covers:

- ▶ Remove cover.
- ▶ If machine washing mattress top/base covers, the temperature should be raised during the wash cycle, to 71°C/160° F, for 3 - 10 minutes, using hospital approved detergents and rinsing agents. The temperature could be raised to 95°C in the case of Air2Care IN mattresses.
- ▶ Dry cover in tumble dryer at low temperature.

NOTE Maximum wash temperature could be 75°C for mattresses other than Air2Care IN mattresses but such a temperature reduces lifetime period of the covers.

16.3.4 Mattress Core

The entire core of the mattress does not require any major cleaning. The core does not need disinfection. Once a month it is recommended to ventilate the mattress core (remove the mattress cover and leave the mattress core on ventilated area for 12 -24 hours). The mattress core cannot be washed by water or by disinfection.

17 Accessories



WARNING!

Risk of injury due to incompatible accessories!

- Use exclusively original accessories from the manufacturer.
- The manufacturer is not responsible for the use of unapproved accessories.



WARNING!

Risk of injury due to damaged accessories!

- Use exclusively accessories in perfect condition.



WARNING!

Risk of injury or material damage due to incorrect use!

- Compatible accessories manufactured by different manufacturers have their own instructions for use. It is necessary to read instructions for use of a compatible accessory with instructions for use of the compatible LINET product to respect especially technical parameters, warning notifications, cleaning and maintenance instructions of LINET products and their compatible accessories!



WARNING!

Risk of injury or material damage when adjusting the bed!

- Ensure that there are no installed accessories that may collide with moving parts of the bed.
- Extra care must be taken to avoid collisions of the accessories with bed parts and with the other accessories during bed positioning.
- Extra care must be taken to avoid collisions of the accessories installed on the bed and things in the vicinity of the bed.

To use the compatible Essenza 300 W Accessories correctly:

- Follow the **D9U00AEZ3-01xx** Instructions for Use where „xx“ means the code of language version.

18 Cleaning and Disinfection



WARNING!

Risk of injury due to incorrect preparation!

- Ensure pedals will not be pressed accidentally during cleaning.
- Ensure the Essenza 300 W is disconnected from the mains before cleaning the Essenza 300 W.



CAUTION!

Material damage due to incorrect cleaning or disinfection!

- Do not use pressure or steam cleaners.
- Follow the instructions and observe the dosages recommended by the manufacturer.
- Ensure that disinfectants are selected and applied exclusively by qualified hygiene experts.
- Respect used materials during cleaning and disinfection! For information see the following table.
- **Check if used cleaning agents and disinfectants are compatible with materials that the product consists of! For information see the following table.**

BED COMPONENTS THAT ARE INTENDED TO BE CLEANED	MATERIALS (SURFACES OF THE MENTIONED BED COMPONENTS)		
Do not clean what is not mentioned in this column!	Competent user is responsible for check if used cleaning agents and disinfectants are compatible with mentioned materials!		
Head board and Foot board	plastic version: Polypropylene (PP)	aluminium version with HPL board: Oxidized aluminium alloy + Lacquered steel + Stainless steel + High Pressure Laminate (HPL) + Polyamide (PA)	
Plastic Head siderails and Foot siderails / Telescopic Head siderails and Foot siderails / Single Collapsible Siderails	plastic siderails: Polypropylene (PP) + Polyamide (PA) + Lacquered steel + Acrylonitrile butadiene styrene (ABS)	telescopic siderails: Oxidized aluminium alloy + Lacquered aluminium (Al) + Acrylonitrile butadiene styrene (ABS) + Polyamide (PA) + Polypropylene (PP) + Polyoxymethylene (POM) + Stainless steel	collapsible siderails: Lacquered steel + Polyamide (PA) + Polyoxymethylene (POM) + Acrylonitrile butadiene styrene (ABS) + Polyvinyl chloride (PVC)
Mattress support platform covers on Backrest, Thighrest and Calfrest	Polypropylene (PP)		
Seat section cover	Acrylonitrile butadiene styrene (ABS)		
Castors	Polyurethane (PUR) + Polyamide (PA) + Polypropylene (PP)		
Castor control levers	Polyamide (PA6)		
Frame of the mattress support platform	Polyamide (PA) + Lacquered steel		
Undercarriage	Lacquered steel		
Corner bumpers	Polypropylene (PP)		
Keyboards (Attendant Control Panel, Handset Control Panel, control elements integrated in the siderails)	Acrylonitrile butadiene styrene (ABS) + Polyoxymethylene (POM) + Polyethylene terephthalate (PET)		
CPR levers	Polyamide (PA6)		
Labels	Polyethylene terephthalate (PET)		
Accessory rails	Polyoxymethylene (POM) + Lacquered steel		
Actuator covers	Polyamide (PA6) + Aluminium (Al)		
Linen Shelf	Plastic Linen Shelf: Polypropylene (PP)	Telescopic Linen Shelf: Oxidized aluminium alloy + Zinc-coated steel + Stainless steel + Acrylonitrile butadiene styrene (ABS) + Polyamide (PA) + Polyoxymethylene (POM)	

For safe and gentle cleaning:

- ▶ Do not use any strong acids or bases (optimum pH range 6 - 8).
- ▶ Exclusively use detergents that are suitable for cleaning medical equipment.
- ▶ Do not use abrasive powders, steel wool, or other materials and cleaning agents that might damage the mattress replacement system.
- ▶ Never use any corrosive or caustic detergents.
- ▶ Never use detergents that deposit calcium carbonate.
- ▶ Never use detergents with solvents that might affect the structure and consistency of the plastics (benzene, toluene, acetone, etc.).
- ▶ Observe local directives regarding infection control.
- ▶ Make sure any cleaning agent used is approved by:
 - the facility in which the mattress replacement system is to be used.
 - by the environmental protection agency of the country in which the mattress replacement system is to be used.

18.1 Cleaning the Essenza 300 W bed

Prepare for cleaning as follows:

- ▶ Put the Mattress Support Platform in the highest position.
- ▶ Adjust the Backrest and Thighrest so that the reverse sides are accessible.
- ▶ Move the bed to the location where it will be cleaned.
- ▶ Lock the brakes on the bed.

18.1.1 Cleaning before Changing Patients

Clean the following bed parts:

- All control elements for adjusting the bed
- All handles
- Head Board and Foot Board
- Siderails (in highest position)
- Freely accessible mattress surface
- Accessory rails
- All Mattress Support Platform covers
- Mattress on all sides
- Freely accessible metal parts of Mattress Support Platform
- Infusion stand sleeve fitting
- Bumpers
- Castors
- Brakes

18.1.2 Daily Cleaning

Clean the following bed parts:

- All control elements for adjusting the bed
- All handles
- Head Board and Foot Board
- Siderails (in highest position)
- Freely accessible mattress surface
- Accessory rails

18.1.3 Complete Cleaning and Disinfection

Clean the following bed parts:

- All control elements for adjusting the bed
- All handles
- Head Board and Foot Board
- Siderails (in highest position)
- Freely accessible mattress surface

- Accessory rails
- All Mattress Support Platform covers
- Mattress on all sides
- Freely accessible metal parts of Mattress Support Platform
- Infusion stand sleeve fitting
- Bumpers
- Castors
- Brakes

18.2 Machine Washing



WARNING!

Material damage to the bed or its parts due to not following the instructions for machine washing!

- ▶ Always follow the instructions from manufacturer for machine washing due to securing maximum service life and functionality of the bed designed for machine washing. Otherwise serious damage to the bed may occur.
- ▶ Not following those instructions may lead to the void of warranty.
- ▶ Do not use any strong acids or bases. Optimum pH range 6 – 8.
- ▶ Never use detergents with solvents that might affect the structure and consistency of the plastics.
- ▶ Never exceed the dosage prescribed by the manufacturer.
- ▶ Never use high pressure cleaner for cleaning the bed.
- ▶ The cooling of the bed with cold water is not allowed.



WARNING!

Material damage to the accessories due to incorrect use in the washing machine!

- ▶ In the washing machine use accessories intended for machine washing only.
- ▶ The only allowed accessories for washing machine is lifting pole and infusion stand.
- ▶ Always check if the allowed accessory is not colliding with washing machine or if it is not preventing from its trouble-free service.

If the bed is designed for the machine washing, then the following symbol is located on the undercarriage under the foot board. The bed intended for machine washing is labeled as **1EZ3W** on the serial label which is located on the mattress support platform frame.

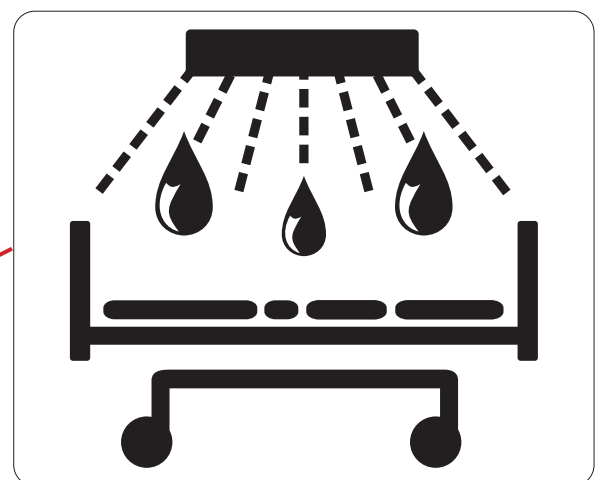
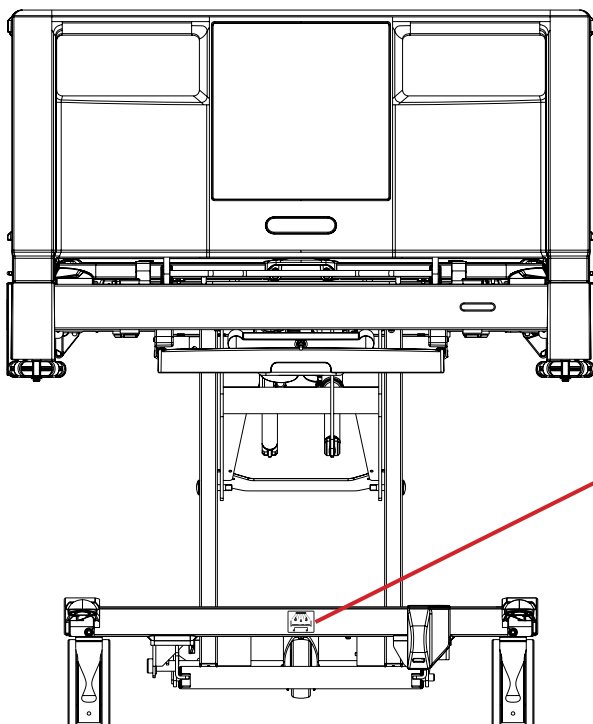


Fig. Essenza 300 W bed intended for machine washing

18.2.1 Preparing for cleaning

To prepare the **Essenza 300 W** bed for machine washing:

- ▶ Disconnect the bed from the mains.
- ▶ Ensure all connectors are correctly plugged in its place.
- ▶ Ensure the Plug&Play connector is closed correctly.
- ▶ **Allowed detergents:**

Washing Cycle	Detergent	Manufacturer
Cleaning and disinfection	Sekumatic FDR	Ecolab
	Neodisher Dekonta AF	Dr. Weigert
Rinsing	Sekumatic FKN	Ecolab
	Neodisher TN	Dr. Weigert

- ▶ **Recommended position of the bed before machine washing:**

- Height of the mattress support platform: approx. 42 cm above the ground
- Siderails: in the highest position
- Linen shelf: pulled out, as much as the washing machine allows
- Bed extension: pulled in
- Mattress support platform: flat (**If the washing tunnel has no standard tilting, clean the bed in the**

Reverse-Trendelenburg Position.)

- Attendant Control Panel: placed on the head board or foot board
- Handset: placed on the mattress support platform
- Mains cable: hanged on the head frame
- Accesories: lifting pole, infusion stand (suitable for machine washing and if compatible with size of washing tunnel)

other accessories are not designed for machine washing

- ▶ Place the bed into the washing tunnel.

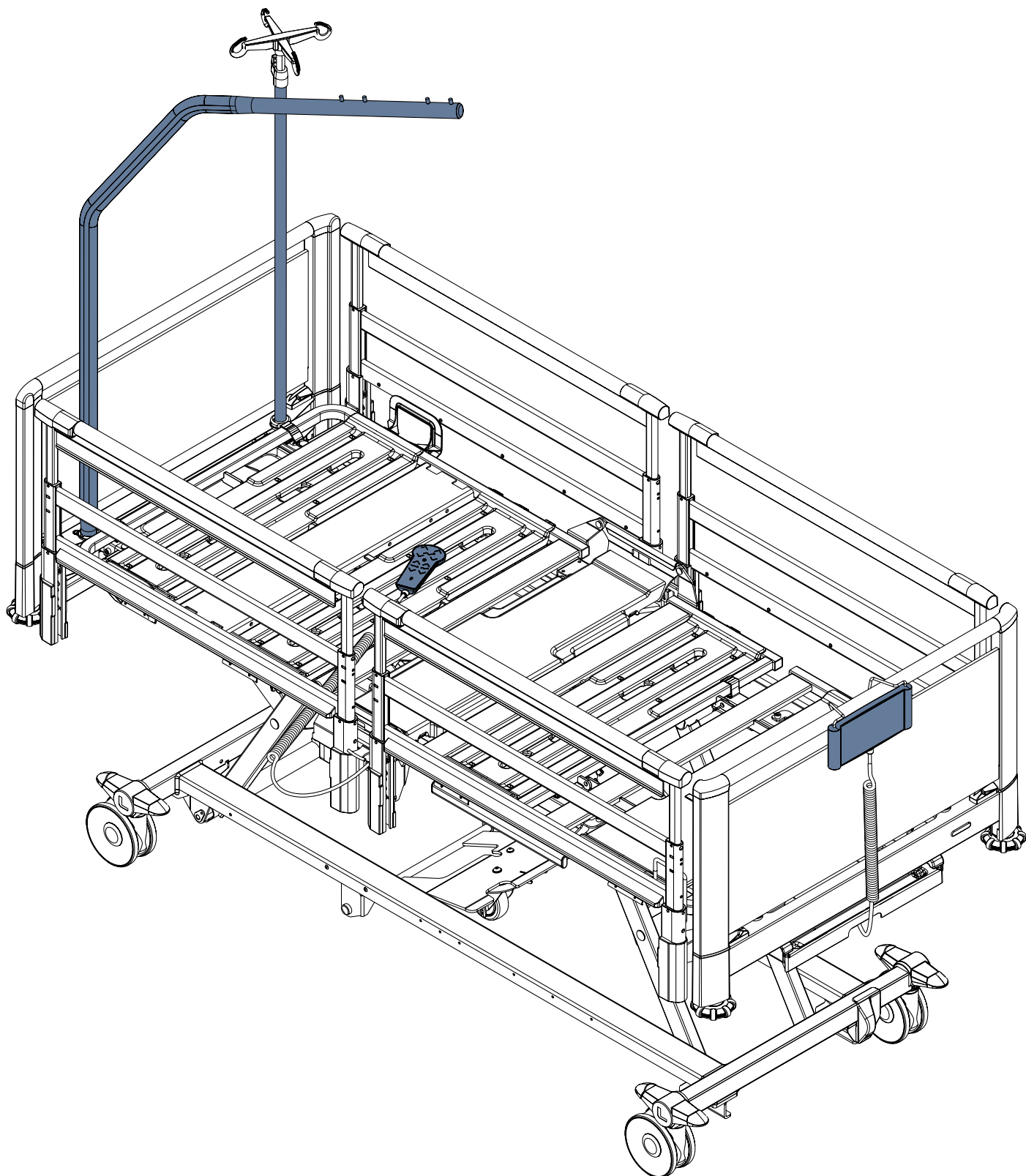


Fig. Recommended position of bed for machine washing

18.2.2 Machine washing

The following conditions for machine washing must be met:

■	Maximum length of machine washing cycle:	7 minutes
■	Maximum pressure of washing tunnel jets:	3-5 bar
■	Maximum water temperature of jets:	70 °C for 120 seconds
■	Maximum rinse water temperature:	85 °C for 20 seconds
■	Maximum drying temperature in washing tunnel:	50-60 °C for 150 seconds
■	Relative humidity:	100%
■	Maximum water hardness:	5 dH

18.2.3 After machine washing

After machine washing and before next use of the bed:

- ▶ Ensure the plug of the power supply cord is dry. Dry if needed.
- ▶ Ensure none of the electronic parts is visibly damaged.
- ▶ Wait 10 minutes with the ambient temperature of 20 °C until the bed cools down to its working temperature.
- ▶ Try some of the bed adjustment functions.
- ▶ It is recommended to visually check the bed for any signs of damage after each machine washing (especially damage on the painted parts). In case of that damage, a qualified and trained service technician authorized by the manufacturer is allowed to repair it using an appropriate paint.
- ▶ It is recommended to use Reverse-Trendelenburg Position (bed foot end is tilted down) in order to get rid of water from all parts of the bed.

19 Troubleshooting



DANGER!

Risk of mortal injury due to electric shock!

- ▶ If a fault occurs, have the electric motor, power box or other electrical parts repaired by qualified personnel exclusively.
- ▶ Do not open the protective covers of the electric motor or the control unit.

Error/Fault	Cause	Solution
Adjusting with position buttons not possible	Function disabled on Attendant Control Panel.	Enable disabled function.
	Actuators have no power. Defective actuators. Defective accumulator.	Check the mains connection. Notify the manufacturer's service department.
	Mains Plug inserted incorrectly.	Insert the Mains Plug correctly.
	Faulty Power Source.	Notify the manufacturer's service department.
	Faulty Control Element.	Notify the manufacturer's service department.
Faulty Mattress support platform Height/Tilt Adjustment	Obstacle on the undercarriage cover.	Remove the obstacle.
	Function disabled on Attendant Control Panel.	Enable disabled function.
	Actuators have no power. Defective actuators. Defective accumulator.	Check the mains connection. Notify the manufacturer's service department.
	Mains Plug inserted incorrectly.	Insert the Mains Plug correctly.
	Faulty Power Source.	Notify the manufacturer's service department.
Lowering Backrest from the upright position not possible	Obstacle under the Backrest or in the drive mechanism.	Remove the obstacle
	CPR Release Handle is defective.	Notify the manufacturer's service department.
Adjusting Siderails not possible	Obstacle in the Siderail Release Mechanism.	Remove the obstacle.
	Siderail Release Mechanism is defective.	Notify the manufacturer's service department.
Faulty brakes	Obstacle blocking brakes mechanically.	Remove the obstacle.
	The brake mechanism is defective.	Notify the manufacturer's service department.

20 Maintenance



WARNING!

Risk of injury when working on the bed!

- ▶ Ensure that the bed is disconnected from the mains power prior to installation, putting into service, maintenance and deinstallation.
- ▶ Ensure that the castors are locked prior to installation, putting into service, maintenance and deinstallation.
- ▶ No part of the Essenza 300 W ME equipment shall be serviced or maintained while in use with a patient.



WARNING!

Risk of injury due to defective bed!

- ▶ Have a defective bed repaired immediately.
- ▶ If the defect cannot be repaired, do not use the bed.



CAUTION!

Material damage due to incorrect maintenance!

- ▶ Ensure that maintenance is performed exclusively by manufacturer's customer service or by authorised service personnel certified by the manufacturer.
- ▶ If the defect cannot be repaired, do not use the bed.

LINET ® recommends attaching the maintenance plaque to the bed.

20.1 Regular maintenance

- ▶ Check regularly movable parts for wear.
- ▶ Perform regularly visual check of the product (with delivery note if necessary).
- ▶ Ask service department of the manufacturer for addition of the original spare parts if some product parts are missing.
- ▶ Ask service department of the manufacturer for replacement of any damaged product parts by the original spare parts.
- ▶ Check that the batteries are working properly. Disconnect the bed from the mains power to check signalisation of battery indicator according to the instructions for use.
- ▶ Have the batteries replaced if they are not working properly.
- ▶ Check regularly that all accessories are working properly.
- ▶ Replace damaged accessories immediately.

20.2 Spare Parts

The serial label is located on the frame of the mattress support platform. The serial label contains information for claims and ordering replacement parts.

Information about spare parts is available from:

- Manufacturer's customer service
- Sales department

20.3 Electrical Safety Check and Preventive Maintenance

20.3.1 Electrical Safety Check



WARNING!

Risk of injury due to incorrect electrical safety checks!

- ▶ Ensure that electrical safety checks are performed exclusively by manufacturer's customer service or by authorised service personnel certified by the manufacturer.
- ▶ Ensure that the electrical safety checks are recorded in accordance with local regulations and Manufacturer's recommendation.

Electrical safety check of the bed must be performed at least once every 12 months.

The procedure for performing the electrical safety check is stipulated in EN 62353:2014.

20.3.2 Preventive Maintenance

Preventive maintenance is the routine service operations of equipment to ensure safe use, prevent unexpected failures and extend life span.

To ensure the correct Preventive Maintenance:

- ▶ Ensure that preventive maintenance is performed exclusively by manufacturer's customer service or by authorised service personnel certified by the manufacturer.
- ▶ Ensure that preventive maintenance is performed in strict accordance with the Service Manual and the mandatory process steps described in the Preventive Maintenance – Checklist.
- ▶ Ensure that the electrical safety checks are recorded in accordance with local regulations and Manufacturer's recommendation.

Preventive maintenance of the bed must be performed at least once every 12 months.

NOTE On request, the manufacturer will provide service documentation (e.g. circuit diagrams, component part lists, descriptions, calibration instructions etc.) for service personnel for the repair of ME equipment designated by the manufacturer as repairable by service personnel.

20.4 Instructions related to machine washing

The manufacturer recommends performing technical check after 25 cycles of machine washing.

The technical check must check the following:

- Sufficient lubrication of the recommended bed parts
- Functionality of the adjustable functions of the bed
- Rust on the bed components
- Status of the accumulator


The technical check must grease the following:

Place on the bed	No.	Recommended grease	Manufacturer
Telescopic siderail mechanism	1	BIO FLUID	WEICON GmbH & Co. KG
Gliders of the linen shelf	2	BIO FLUID	WEICON GmbH & Co. KG
Bed height actuators in extended position	3	Klübersynth UH1 64-1302	Klüber Lubrication München GmbH & Co. KG
Mechanism of the mattress support platform lifting	4	Klübersynth UH1 64-1302	Klüber Lubrication München GmbH & Co. KG
Guideway of the undercarriage sliders	5	Klübersynth UH1 64-1302	Klüber Lubrication München GmbH & Co. KG


To complete the greasing procedure:


- position the **Essenza 300 W** bed to the highest position, to the lowest position, to the maximum Trendelenburg tilt and to the maximum Reverse-Trendelenburg tilt after greasing the five places on the bed.

20.4.1 Recommended places for greasing

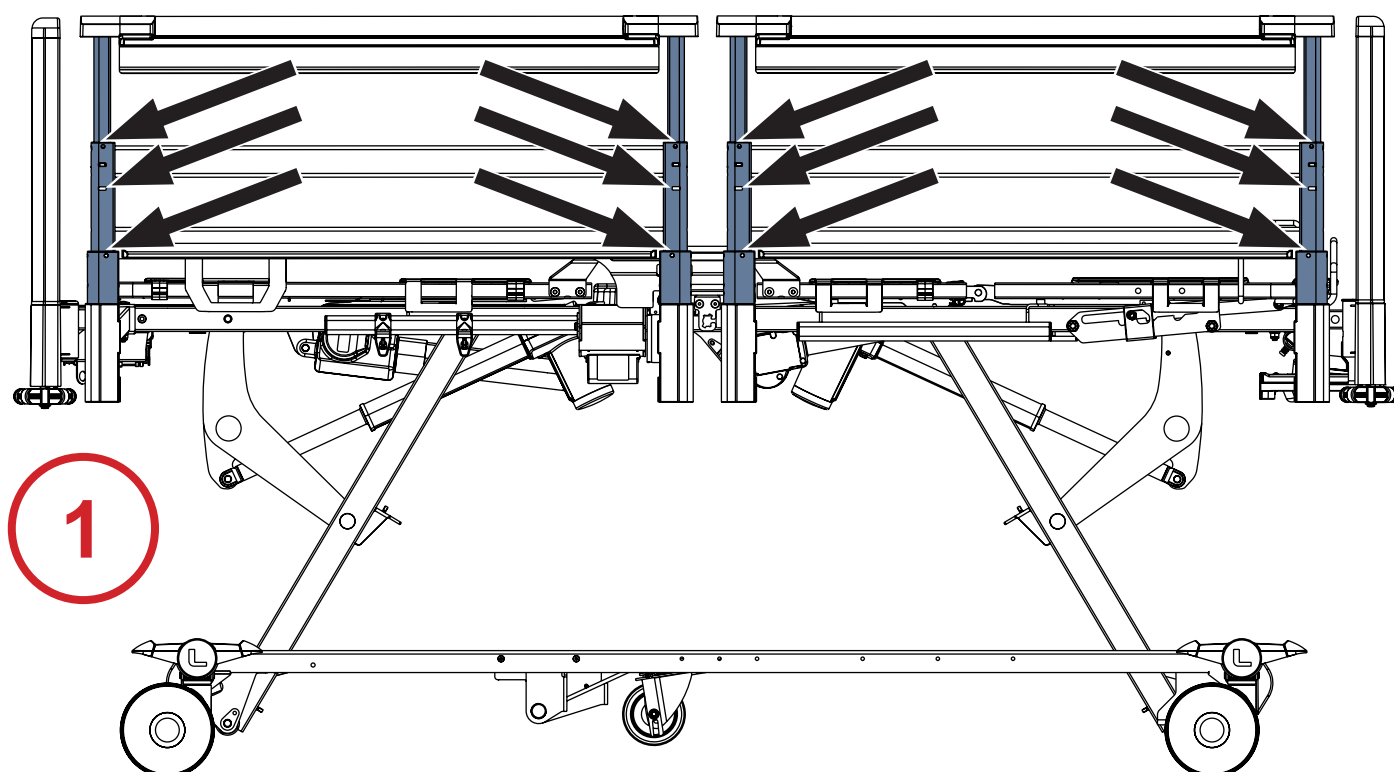



2

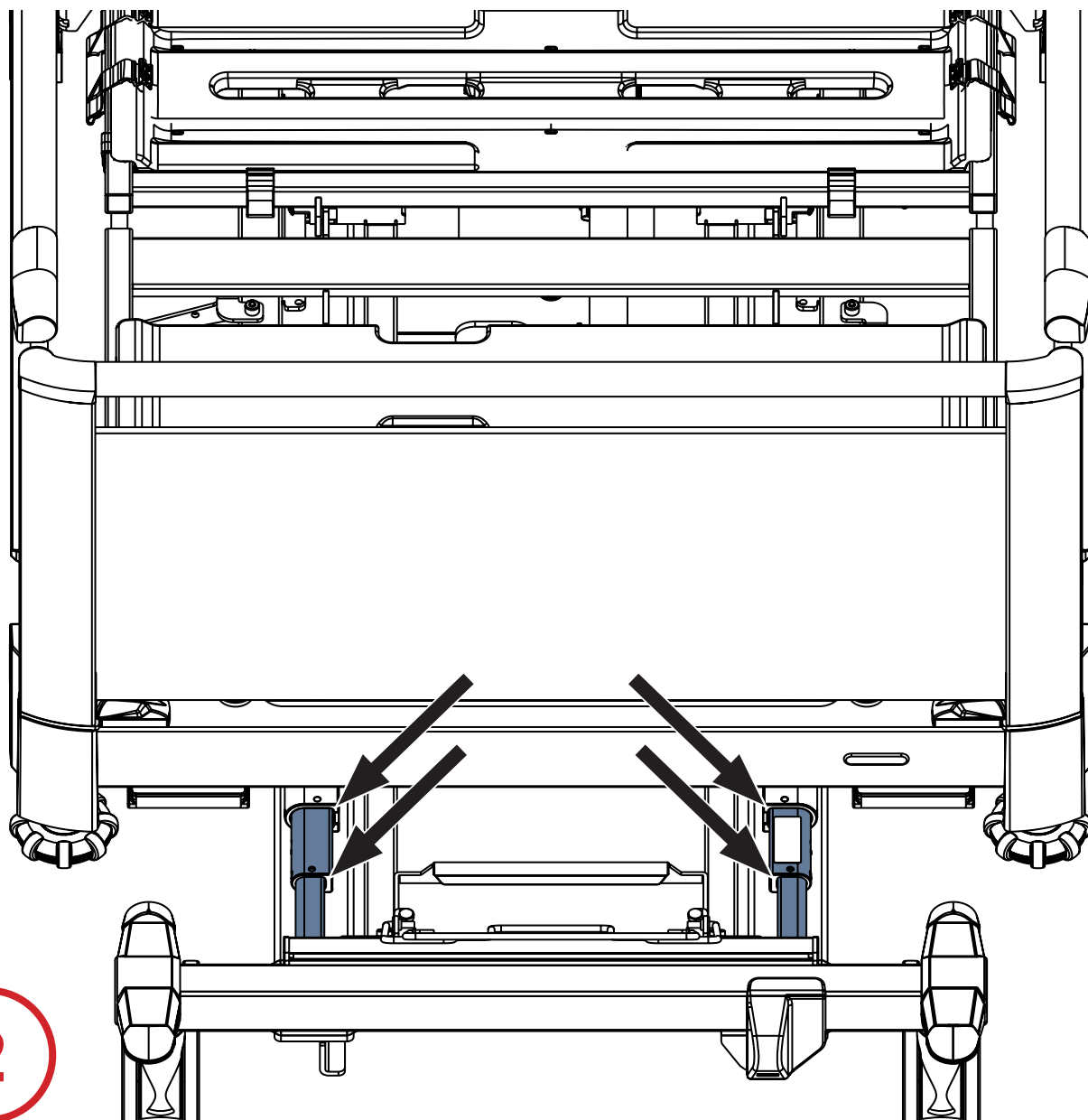

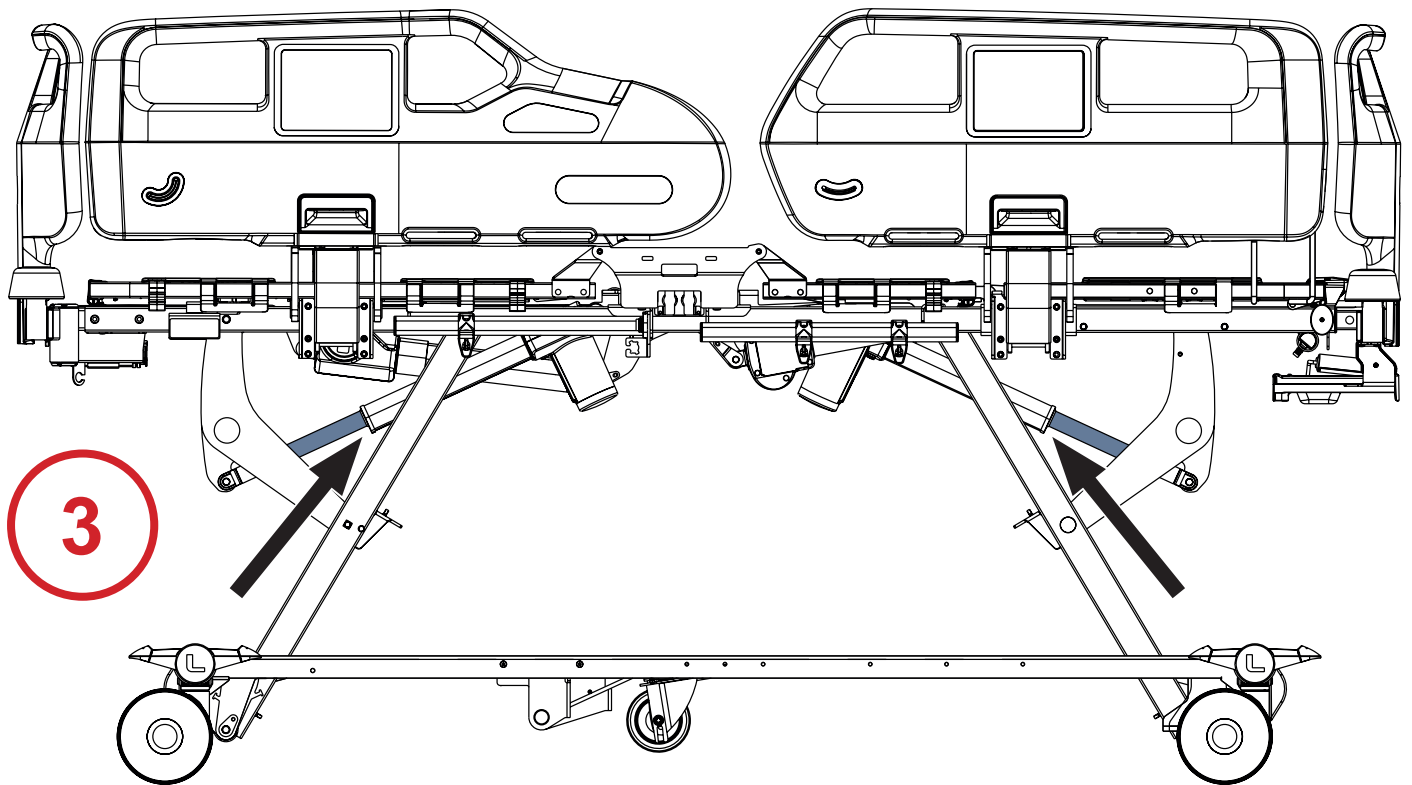

4

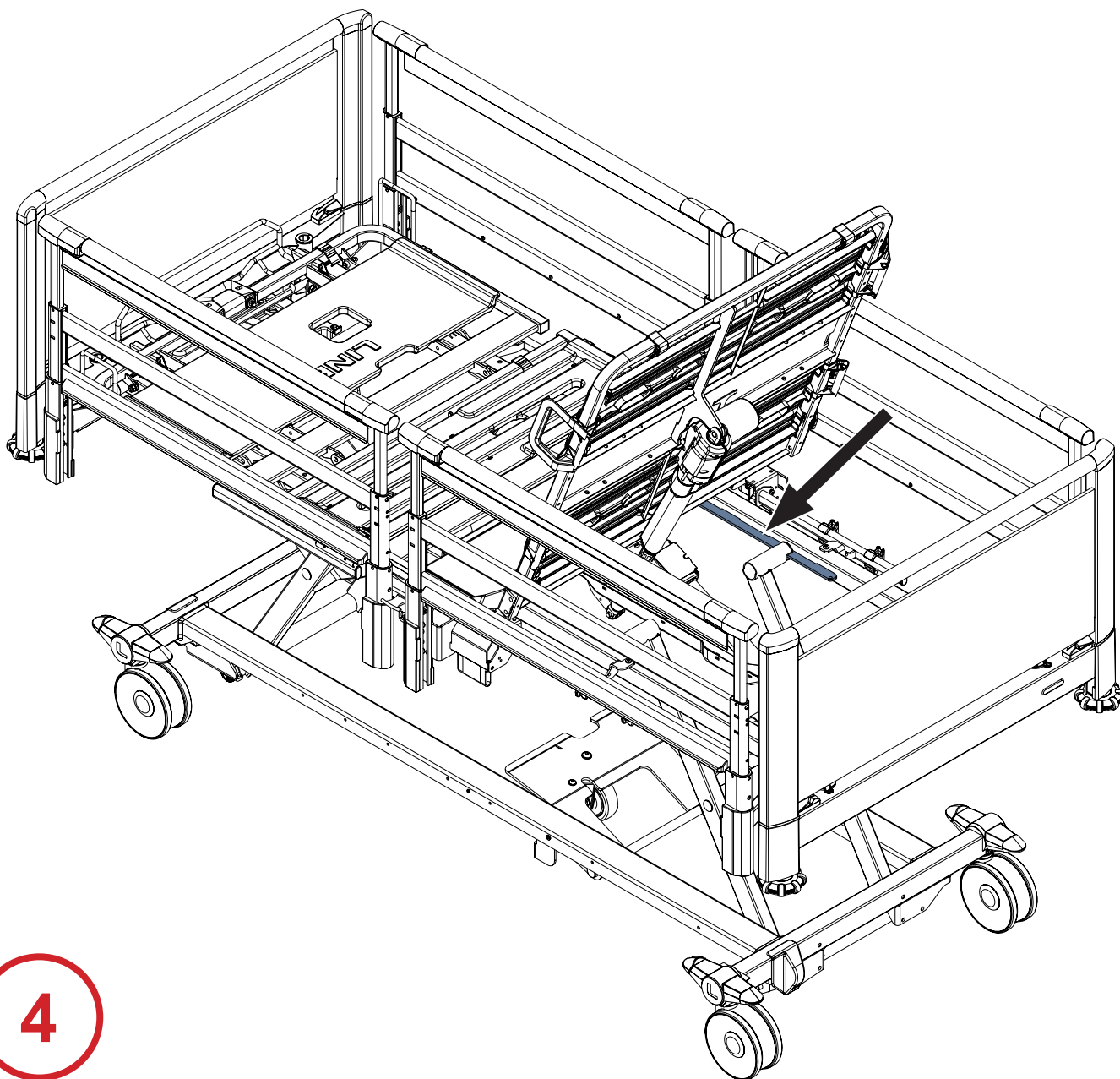

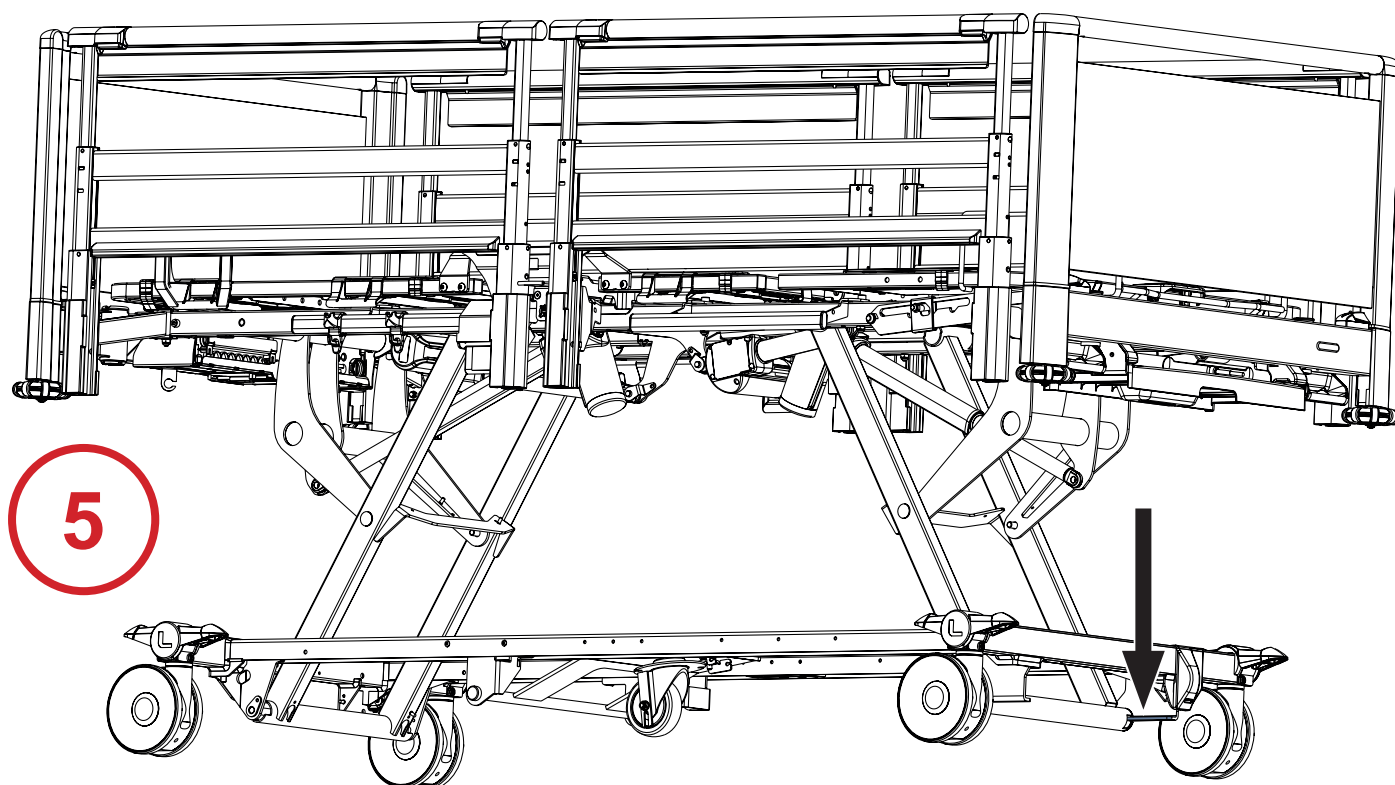

20.5 Control Unit Leak Tester (8300-9004) for the Essenza 300 W washable bed

The product has been tested and quality approved in the LINET factory.

The periodic testing is recommended also during use on the customer's side once a year or when all plugs are correctly installed after changing electric components. The easy-to-use Control Unit Leak Tester (8300-9004) is available at LINET Service Department as an equipment intended for this quality control during maintenance which must be performed by a qualified service technician who has been trained and approved by LINET.

21 Disposal

21.1 Environment Protection

The company LINET® is aware of the importance of environmental protection for future generations. Within this company the environmental management system is applied in accordance with the internationally agreed standard ISO 14001. The compliance with this standard is annually tested by the external audit executed by an authorised company. Based on the **WEEE** - Waste, Electric and Electronic Equipments Directive the company LINET, s. r. o. is registered in the List of Electric and Electronic Equipment Producers (**Seznam výrobců elektrozařízení**) on the Ministry of the Environment of the Czech Republic (Ministerstvo životního prostředí).

Materials used in this product and in the LINET® accessories are not environmentally hazardous. LINET® products and LINET® accessories meet valid requirements of national and European legislation in the areas of **RoHS** and **REACH**, so they do not contain any prohibited substances in excess quantities. None of the wooden parts is made of tropical wood (such as mahogany, rosewood, ebony, teak etc.) or made of timber from the Amazon region or from similar rainforests. Product noise (sound pressure level) meets requirements of the regulations for the protection of public health against undesirable effects of noise and vibration in protected interior spaces of buildings. Used packaging materials are in accordance with requirements of the Packaging Act (**Zákon o obalech**). For disposal of packaging materials after installation of products contact your sales representative or manufacturer's customer service about the possibility of a free take-back of packaging through an authorized company (more details on www.linet.cz).

21.2 Disposal

Materials used in this product and in the LINET® accessories burden the environment but at the same time the whole range of these materials can be very effectively reused and recycled. Mechanical disassembly of the product and material sorting into the basic waste types (plastic, metal, wooden materials) should be performed after the end of product life. The main objective of the obligations arising from the European Directive on Waste, Electric and Electronic Equipments is to increase the re-use, material recovery and recovery of electric and electronic equipment at the required level, thereby avoiding the production of waste and thereby avoiding the possible harmful effects of hazardous substances contained in electric and electronic equipment on human health and the environment. LINET® electric and electronic equipments that have a built-in battery or accumulator and LINET® accessories are designed so that the used batteries or accumulators can be safely removed by LINET® qualified service technicians. There is an information about its type on the built-in battery or accumulator.

21.2.1 Within Europe

To dispose of the electric and electronic equipment including LINET® accessories:

- ▶ The electric and electronic equipment must not be disposed of as household waste.
- ▶ Dispose of this equipment at designated collection points or take-back points.
- ▶ Dispose of the product or its components or its accessories in accordance with local laws and regulations!
- ▶ Hire an approved waste disposal company for disposal!

To dispose of the other equipment including LINET® accessories:

- ▶ The equipment must not be disposed of as household waste.
- ▶ Dispose of this equipment at designated collection points or take-back points.

LINET® participates in a collective system with take-back company REMA System (see www.remasystem.cz/sberna-mista/). By bringing electric and electronic equipment to a take-back point, you participate in recycling and you save primary raw material resources while protecting your environment from effects of unprofessional disposal.

21.2.2 Outside Europe

- ▶ Dispose of the product or its components or its accessories in accordance with local laws and regulations!
- ▶ Hire an approved waste disposal company for disposal!

22 Warranty

LINET® will only be held responsible for the safety and reliability of products that are regularly serviced, maintained and used in accordance with the safety guidelines included in the instructions for use.

Should a serious defect arise that cannot be repaired during maintenance:

- Do not continue to use the Essenza 300 W bed.

The warranty on this product and its conditions depend on the agreement between the buyer and the seller, in addition to the legal minimum guarantee provided by law.

The warranty covers all material and manufacturing-related failures and errors. Failures and errors caused by incorrect use and external effects are not covered. Justified complaints will be fixed free of charge during the warranty period. Proof of purchase, with the date of purchase, is required for all warranty service. Our standard terms and conditions apply.

23 Standards and Regulations

23.1 Essenza 300 W

Applied norms are stated on Declaration of Conformity.

23.2 Manufacturer

The manufacturer adheres to a certified quality management system in compliance with the following standards:

- ISO 9001
- ISO 14001
- ISO 13485
- MDSAP (Medical Device Single Audit Program)